



HYCD-319A

Combined Refrigerator and Freezer

Scope of Application

This product is designed to store frozen plasma, biomaterials, vaccines, electronic devices and low temperature tests of special materials, in blood stations, hospitals, epidemic prevention stations, universities, scientific research institutes, electronic chemical industry and other enterprise laboratories.

Innovative Design

Intelligent Microprocessor with Digital Display

• Audible and Visual Alert with Remote Alarm Contact

Natural Hydrocarbon Refrigerant

Stainless Steel Interior

Self-closing Door with Magnetic Gasket

Padlock to Secure Safe Storage

Qingdao Haier Biomedical Co., Ltd.

No.280 Feng Yuan Road, High-tech Zone, Qingdao, 266109, P.R. China Tel: +86-0532-88935593 E-mail: inquiry@haierbiomedical.com Website: www.haiermedical.com











Combined Refrigerator and Freezer





ন্থি

Multipurpose and Space Efficient

- One cabinet with two independent chambers, both are self-contained and independently controlled with individual doors
- Refrigerator chamber is equipped with a double-paned glass viewing window and LED interior lighting



Environment Friendly

Hydrocarbon refrigerant (R600a for refrigerator / R290 for freezer) and CFC-free polyurethane insulation for reduced global warning potential



Stainless Steel Interior

The interior of refrigerator and freezer are both made of stainless steel, which is easy to clean and corrosion-proof



Microprocessor Controller with Digital Display

- Refrigerator temperature 2~8°C with resolution of 0.1°C
- Freezer temperature -40~-10°C with resolution of 1°C
- Multiple alarm functions with remote alarm contact to ensure safe monitoring



High Performance

- \bullet Refrigerator: Forced-air cooling system with cycle defrost, temperature uniformity $\pm 3^\circ C$
- Freezer: Cold wall cooling system, manual defrost, temperature uniformity±5°C
- Standard USB port & optional printer, easy to get temperature record

Ergonomic Design 🔇



Equipped with padlocks in both refrigerator and freezer to ensure storage safety



Equipped with electrically heated glass to prevent condensation effectively



Refrigerator and freezer are equipped with self-closing hinges to avoid door ajar (less than 90° self-closing)





Two test holes for both refrigerator and freezer as standard, allow independent testing of cabinet temperature

Combined Refrigerator and Freezer







Haier Biomedical



Refrigerator Cool Down Curve

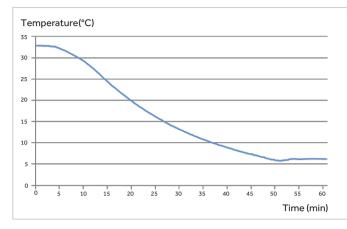


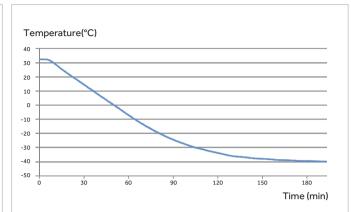


Freezer Cool Down Curve

Refrigerator Stability Curve



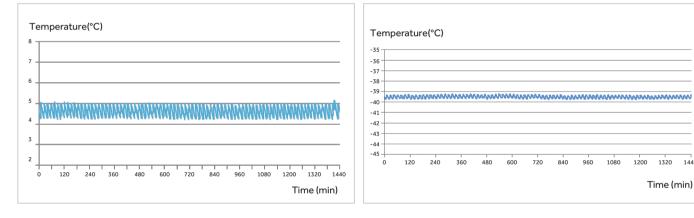



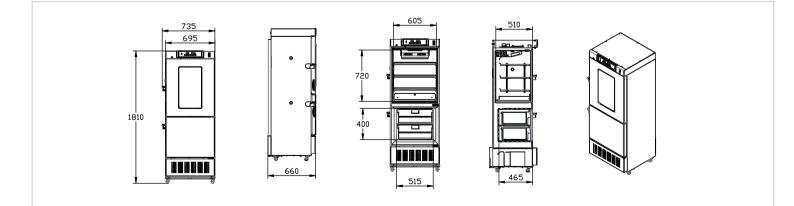

Haier Biomedical

1440

Combined Refrigerator and Freezer



Model			HYCD-319A
	Cabinet Type		Upright
Technical Data	Climate Class		Ν
	Cooling Type		Refrigerator: Forced air cooling Freezer: Direct cooling
	Defrost Mode		Refrigerator: Auto/ Freezer: Manual
	Refrigerant		Refrigerator: R600a 45g/ Freezer: R290 60g
	Sound Level (dB(A))		50
Performance	e Temperature Range (°C)		Refrigerator: 2~8 / Freezer: -10~-40
Control	Controller		Microprocessor
	Display		LED
Electrical Data	Power Supply (V/Hz)		115/60
	Electrical Current (A)		6
Dimensions	Capacity (L/Cu.Ft)		Refrigerator: 200/7.06; Freezer: 100/3.53
	Net/Gross Weight (approx)	kg	145/160
		lbs	319.7/352.7
	Interior Dimensions (W*D*H)	mm	Refrigerator: 605*510*720 Freezer: 515*465*440
		in	Refrigerator: 23.8*20*28.3 Freezer: 20.2*18.3*17.3
	Exterior Dimensions (W*D*H)	mm	735*660*1810
		in	28.9*26.0*71.3
	Packing Dimensions (W*D*H)	mm	775*720*1930
		in	30.5*28.3*76
	Container Load (20'/40'/40'H)		21/43/43
Alarms	High/Low Temperature		Υ
	Sensor Error		Y
	Power Failure		Y
	Door Ajar		Y
	High Ambient Temp		Y
	Low Battery		Y
	Remote Alarm		Y
Accessories	Caster		Y/4
	Foot		Y/2
	Porthole		Y/2
	Shelves/Drawers		Refrigerator: 3/1; Freezer:1 (For supporting the upper drawer) /2 (without sliding rack)
	Stylus printer/RS485		Optional (printer or RS 485, either-or)
	USB Interface		Optional
Others Certification			UL

Product appearance and specifications are subject to change without notice

