

**ICARE MATTRESSES & OVERLAYS** 





# M3 MEDICAL MATTRESS

The M3 has an unparalleled user comfort rating which is a crucial factor and often overlooked when. The full ActiveFloat™ layer provides edge to edge support for users eliminating ridges and increasing pressure redistribution area.


### **SPECIFICATIONS**

Feel RatingPlushMattress Height200mmRisk CategoryHigh RiskUser Weight30kg - 215kgWarranty5 Year Warranty

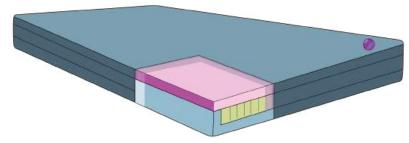

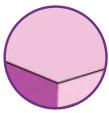

ACTIVEFLOAT™ LAYER
Our latest technology is an extremely elastic and highly breathable material designed to move in every direction even after full compression. This not only reduces shear in every direction but greatly increases user comfort.

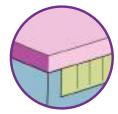

THREE LAYER CONSTRUCTION
One piece U shape firm XP
Support Foam base and side
wall support with medium feel
lateral castellated foam inlay.
Designed especially to reduce
roll-out risk and bed falls.



EDGE TO EDGE
The full ActiveFloat™ layer provides edge to edge support for users eliminating ridges and increasing pressure redistribution area.



SLIDESKIN™ ANTI-SHEAR
The unique composition of
SlideSkin™ technology
reduces shear with a
second membrane. This
enables the cover to move
or slide with the users
movements.



ZEROTEC COVER
Zerotec fabric is a medical barrier fabric that has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This cover has welded seams and waterfall flaps for fluid control.

### **BED CODES & SIZES**

LONG SINGLE	ICM3LS	2030 x 900mm
KING SINGLE	ICM3KS	2030 x 1070mm
LONG DOUBLE	ICM3D	2030 x 1350mm
QUEEN	ICM3Q	2030 x 1520mm
SPLIT QUEEN	ICM3SQ	2030 x 760mm <sub>(per side)</sub>

# **CARE & CLEANING**

- No turn design reducing manual handling
- Anticrease XP base to allow easy contour for any profiling bed
- Medical Cover BS7175: Crib 5
- Waterlow ≤20
- To launder the cover in a washing machine, set the main wash to a minimum
  of 65°C for no less than 10 minutes. The cover can be tumble dried on low
  heat and the temperature must not exceed 60°C. The foam can be
  autoclaved at 134°C.

# **PRESSURE CARE**

Risk Category High

This rating is based on our average peak pressure of <20mmHg.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. For more information refer to the Icare pressure care manual.



# **ACTIVEFLOAT™ OVERLAY**

The Icare ActiveFloat<sup>TM</sup> Overlay is the latest technology is an extremely elastic and highly breathable material designed to reduce shear and offer pressure relief.





SlideSkin™ reduces shear enables the cover to move or slide with the user movements.

ActiveFloat™ reduces shear in every direction even after full compression.

#### **SPECIFICATIONS**


Feel Rating Soft

Mattress Height 50mm

Risk Category High Risk

User Weight 30kg - 215kg

Warranty 5 Year Warranty

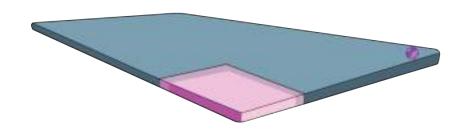



#### PRESSURE REDISTRIBUTION

The overlay combines the features the new ActiveFloat™ foam which reduces shear whilst remaining fully pliant under maximum compression. This provides a product that is responsive to temperature around inflamed areas and will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.



#### VENTILATED SIDE FABRIC

Vant fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. We use this on the side wall of this overlay cover to greatly increase air movement.



#### MEDICAL COVER

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties.

# **BED CODES & SIZES**

SINGLE	ICOFS	1900 x 900mm
LONG SINGLE	ICOFLS	2030 x 900mm
KING SINGLE	ICOFKS	2030 x 1070mm
LONG DOUBLE	ICOFD	2030 x 1350mm
QUEEN	ICOFQ	2030 x 1520mm

# **PRESSURE CARE**

**Risk Category** High

This rating is based on our average peak pressure of <15mmHg.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. For more information refer to the Icare pressure care manual.