



CM Slim

Contouring Master® HI-EMT®



Burn Fat - Build Muscle +



2 Years Warranty

Non-Invasive Contouring Treatment
High-Intensity Electromagnetic Muscle Trainer

@humanityskincare @humanityskincareaustralia humanityskincare.com.au

CONTACT US 0402909551 - ELEANOR

support@humanityskincare.com.au
100% Australian Owned & Operated
Sydney, Australia

SCAN ME!



MAKE YOUR BODY SLIM

CMSLIM®



“
Build Your Muscles

A focused electromagnetic field passes through all skin and fat to effectively stimulate the muscle


Constant Output



Maintain Strong Output

CMSlim® stimulates muscles by mixing various clinically valid frequencies depending on the treatment site or the protocol chosen.

By innovatively improving the high-voltage charging/discharging speed of the high-intensity magnetic field generator, the magnetic field energy is kept constant in the high-frequency band.

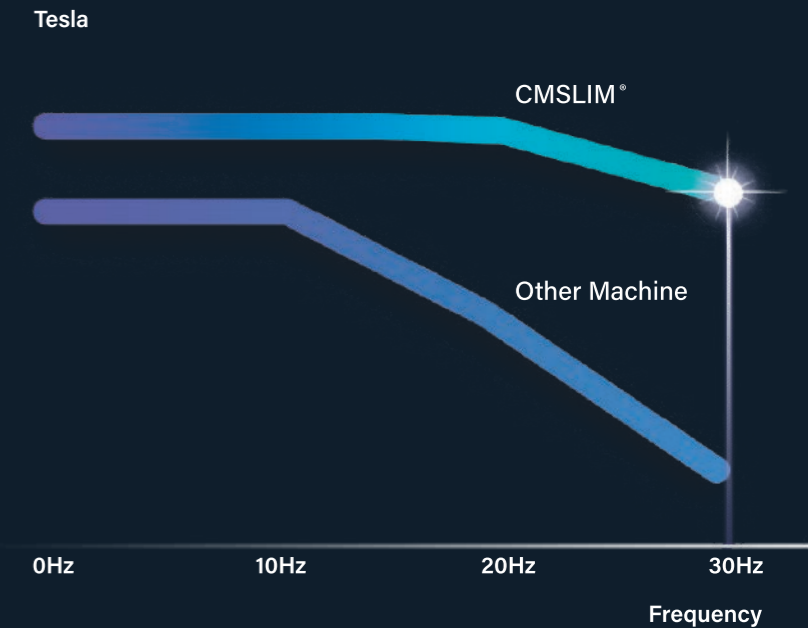

“
Tesla Output

CMSlim® maintains the strong output almost constant even as the muscle stimulation frequency increases.


HI - EMT®



High - Intensity Electromagnetic Muscle Trainer

HI-EMT® is a safe technology developed by DAEYANG MEDICAL KOREA and used in professional medical devices for rehabilitation, physical therapy, and aesthetic purposes.

Electromagnetic field passes non-invasively through the body and interacts with motor neurons which subsequently trigger supramaximal muscle contractions.

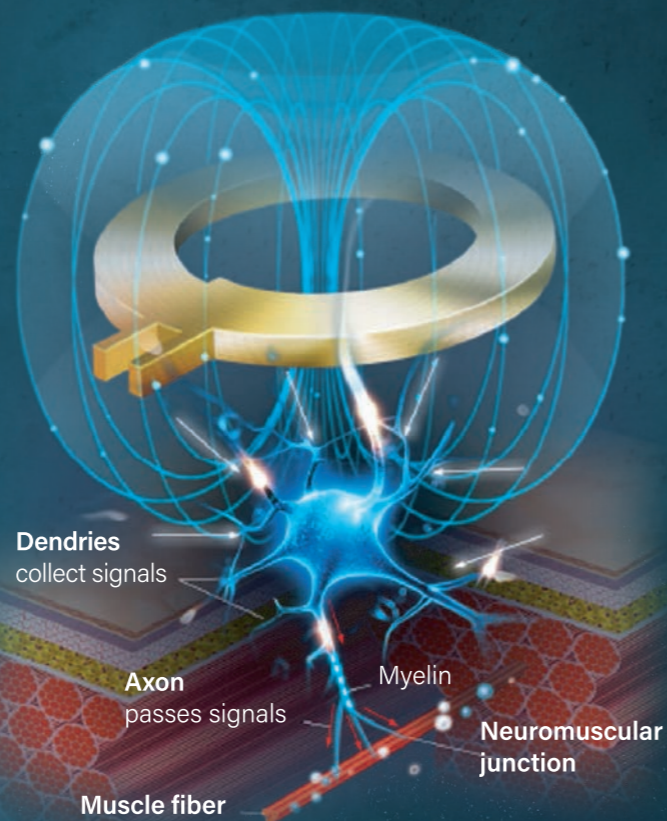

Two Applicators

Separate Control

CMSlim® is equipped with two individually controlled applicators, so it is very effective when you want to apply the same protocol simultaneously on the buttocks, thighs, shoulders or limbs.

The two applicators can be used either at the same time or separately. You can turn the left and right applicators on/off by touching the icons on the screen.

Strong Hardware + Smart Software

Optimal combination of clinically proven protocols and robust medical electromagnetic generator.



CMSlim® is a professional medical aesthetic equipment designed to produce optimal performance by combining clinically proven software protocols on durable hardware that can produce strong and stable output.

Clinical Background

In the fields of rehabilitation, physical therapy, and sports medicine, there are numerous studies on the effects of electrical stimulation such as NMES, TENS, and FES on muscle physiology.

The CMSlim® procedure protocol was developed through various papers and clinical studies on optimal frequency, duty cycle, stimulus intensity and pulse pattern.

Protocol Preset Mode

The preset mode of CMSlim® can be quickly set according to the frequently used treatment area and purpose.

The detailed protocol provides stimulation patterns optimised for gender, body part, and procedure purpose. Each stimulation pattern is designed to properly manage and offset muscle fatigue according to muscle fiber type.

Protocol Manual Mode

Depending on the patient's condition, the preset protocol may not be sufficient to achieve the intended purpose. In this case, manual mode can be the solution.

The operator can fine-tune the frequency, duty cycle, intensity and duration of stimulation. The custom protocol set in this way can be saved separately.

Workout Program

CMSlim® offers 6 workout programs based
First a stretching and warm-up phase, and after the 4-step

on muscle physiology and sports science.
main workout sessions, followed by a cool-down session.

Stretching — 1 mins

Warming Up — 5 mins

Workout 1/2/3/4 — 5 mins

Cool Down — 4 mins



This is a one-minute session for the purpose of muscle relaxation before the start of the procedure. Relieves muscle tension, increases elasticity, and makes activities smooth and comfortable.

Before starting the main program, this is a preparatory session to warm up the body and improve the condition to an adaptable state. The warm-up session lasts five minutes.

The exercise session is divided into 4 steps of 5 minutes each for a total of 20 minutes. It is a customised CMSlim® exercise program optimised according to the treatment area and purpose.


A session to relieve tight muscles after an intense, customised program. Slowly release contracted and pumped muscles to avoid muscle soreness.

Rental Mode

USB S/W Update

Cooling System

Universal Input

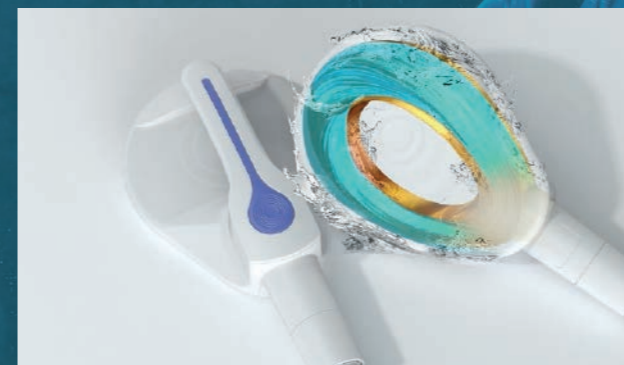

CMSlim® provides a rental mode so that you can use the product by renting it without purchasing it.

Upgrading the procedure program is easily accomplished with a single USB thumb drive.

The water cooling system optimised for CMSlim® innovatively improves heat exchange efficiency.

Stable operation is possible in 100-240V, 50-60Hz power environments around the world.

HIPEX™



HIGH INTENSITY PELVIC EXERCISE

Tighten pelvic floor muscles
Recover bladder function

Specially designed for safe pelvic floor muscle stimulation.

The HIPEX™ Chair is an add-on to CMSlim® device

Before & After



1 Month
8 Treatments



3 Weeks
5 Treatments



1 Month
2 Weeks
12 Treatments



Country of Origin : Republic of Korea
Manufacturer : Daeyang Medical Co., Ltd.
Humanity Skin Care Australia is an authorised reseller of CMSlim®.
Distributed by MM Medical Aesthetics.

Product Name	CMSLIM 1.5
Power Input	220-240V, 50/60Hz
Dimension	420(W)x660(D)x1150(H)
Weight	70kg