



Truly Easy to use.

TruZone® Peak Flow Meter (PFM)

TruZone® PFM is a true advance in helping patients monitor their condition and effectively manage their asthma. Built around the NAEPP guidelines and the ATS Standardization of Spirometry, **TruZone® PFM** delivers the precision, accuracy and repeatability health care professionals require. Designed around patient needs, the **TruZone® PFM** is reliable, easy to use and portable.

The **TruZone® PFM** is:

Truly easy to set up and use


No cumbersome calculations—just apply the unique **Color Zone®** Tape and zones are instantly and accurately established.

Truly more precise

The **TruZone® PFM** log scale provides a more easily read, precise peak flow measure in clinically critical low flow areas.

Truly versatile

One size fits all—**TruZone® PFM** measures peak flows from 60 lpm to 800 lpm. Pediatric to adult, and for lung disorders, the **TruZone® PFM** is up to the job.

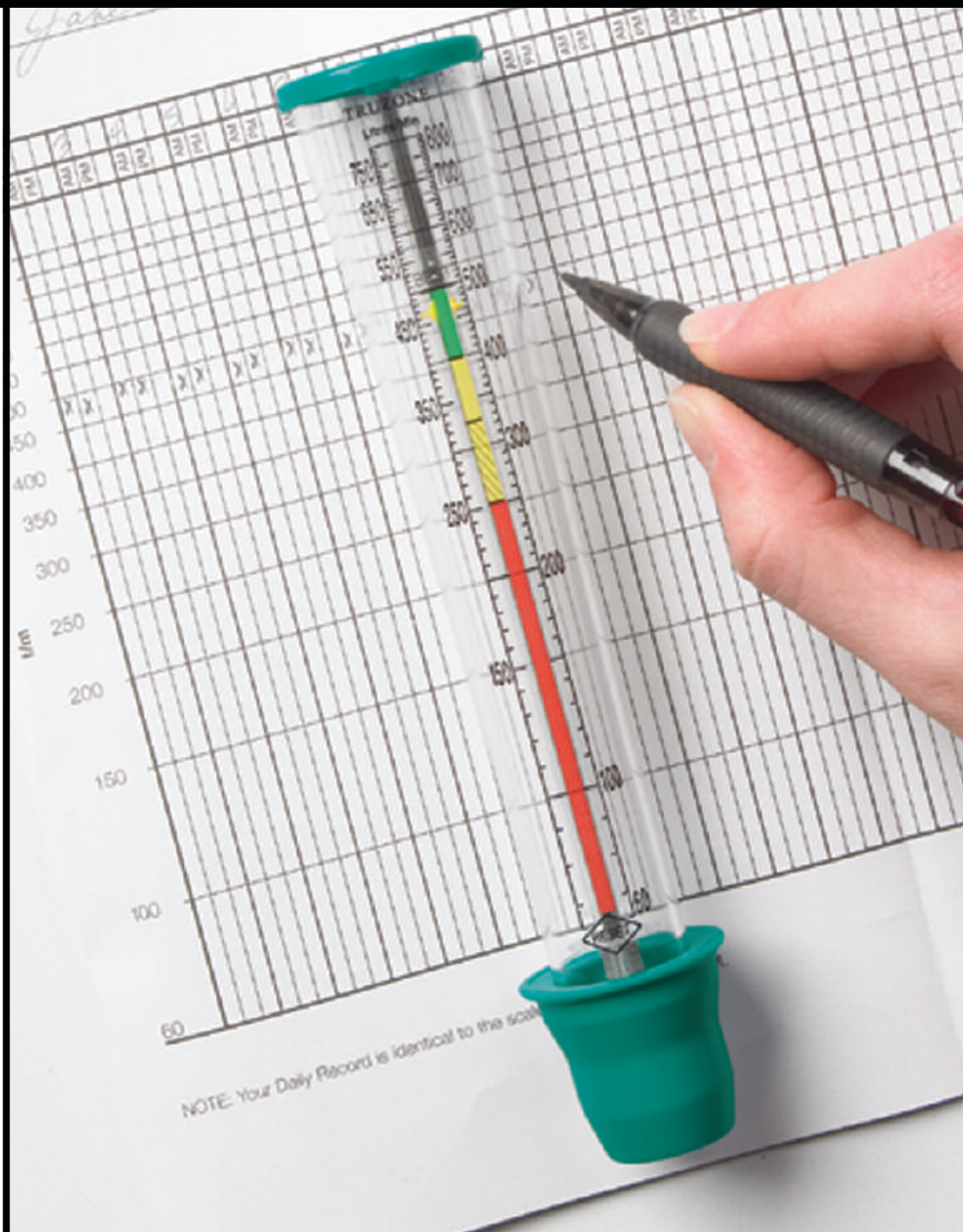

Clearly more accurate. Simple to use. Easy to read.



TruZone® PFM

Color Zone® Tape helps to make asthma management easy, identifying zones without calculations.

Color Zone® Tape allows the use of an additional fourth zone - unique to **TruZone®** PFM - to help better manage your asthma patients.

The **TruZone®** PFM comes with a one-year warranty from the date of purchase.

Internal indicator and fixed **Color Zone®** Tape deliver tamper proof PEF readings and reliable zone settings.

Clear chamber design (actual size) assures users of cleanliness.

Logarithmic scale improves precision and simplifies reading of Peak Expiratory Flow (PEF).

Unique ergonomic body shape and small size facilitate patient compliance.

Multi-patient use when used with disposable one-way valve mouthpiece.



Tested with a pulmonary waveform generator and the 25 relevant waveforms of the ATS Standardization of Spirometry, the **TruZone®** PFM meets the ATS monitoring requirements for accuracy when used with or without the disposable mouthpiece.

REGISTERED
ISO 13485
QUALITY SYSTEM

Endorsed by health care professionals



For more information or a demonstration, contact your Monaghan Medical Sales Representative.

Ordering Information

Part #	Description	
96510	TruZone® PFM	cs/10
95028	TruZone® PFM Patient Instruction Video	ea.
96200	Disposable One-way Valve Mouthpiece	cs/200
95021	TruZone® PFM Patient Instruction CD-ROM (MAC & PC)	ea.



monaghan™

Monaghan Medical Corporation
5 Latour Ave., Suite 1600 • Plattsburgh, NY 12901
Customer Service 800-833-9653
www.monaghanmed.com

™ and © are trademarks and registered trademarks of Monaghan Medical Corporation
© 2000 Monaghan Medical Corporation