

eRapid[®]
nebulizer system

MODERN INHALATION THERAPY FOR **FAST**
AND **EFFICIENT** TREATMENTS




eRapid[®] is a unique and innovative general-purpose nebulizer system developed by PARI that can significantly reduce treatment times and optimize the amount of available medicine to the lungs!

Small, lightweight, and with multiple power options, eRapid incorporates the patented PARI eFlow[®] technology - a vibrating perforated membrane system that ensures precise particle size with an incredibly fast treatment time. For more information, visit www.PARI.com/eRapid.

Modern Inhalation Therapy

The eBase® Controller is the same controller used with the drug-specific nebulizer handsets, minimizing the need to transport multiple controllers.

DRUG-SPECIFIC
NEBULIZER HANDSET†

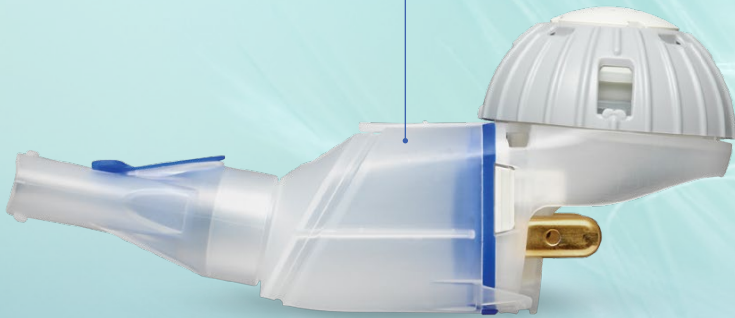

The Main Advantages at a Glance:

- Treatment time reduced approximately 50-65% compared to traditional jet nebulizers encouraging patient adherence.^{1,2}
- Delivered drug doses comparable to LC PLUS jet nebulizer.²
- Quiet operation for discrete use.
- Light, compact and portable.
- Flexible power options including off the shelf batteries, rechargeable batteries, or AC wall power.
- Two (2) complete nebulizer handsets to alternate or use while traveling.

Highest Standards of Hygiene for the Nebulizer Handset*

- Can be disinfected by boiling in distilled water.
- Autoclavable (up to 121° C).
- *easy*care cleaning aid: use to clean the eRapid aerosol head.
- *easy*care cleaning aid is designed to maintain nebulization time and aerosol performance of the eRapid Aerosol Head through a special rinsing process.

eRapid[®]

nebulizer system


UNIVERSAL CONTROLLER




GENERAL-PURPOSE ERAPID NEBULIZER HANDSET




eBase Controller provides patient with feedback including battery life, proper assembly, medication loading, and end of treatment.*

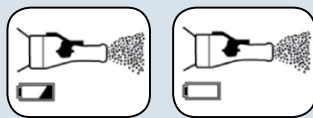

➡ **Battery charge level low**
Replace with new batteries or use AC power supply

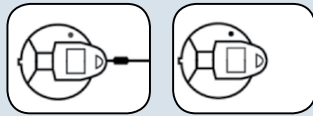

➡ **Cable connection to the control unit or on the nebulizer handset interrupted**
Check cable connections

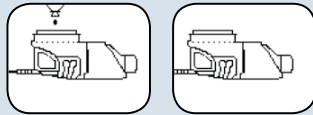

➡ **No medication added**
Add medication

1) Fischer R, Nahrig S, Kerscher M, Huber RM: Inhalation with PARI eRapid reduces treatment burden significantly in adult CF patients without altering lung function. 30th ECFC 2007.

2) Seemann S, Hug M, Lichtinghagen R, et al. Improving aerosol drug delivery in cystic fibrosis. 28th ECFC. Crete, Greece, June 2005.

*Drug-specific handsets should be used only with the drug product approved for its exclusive use.

*Please refer to the eRapid Nebulizer System Instructions for Use.

eRapid Accessories and Replacement Parts



eRapid Nebulizer System	178G1030
eRapid Nebulizer Handset with Aerosol Head*	178G8030
PARI Smartmask® Kids	078F5000
<i>easycare</i> Cleaning Aide	078G6300
Nebulizer Connection Cord	178F6009
AC Power Supply	078F0600



PARI RESPIRATORY EQUIPMENT, INC.
 TEL: 1-800-FAST-NEB (327-8632)
 PRODUCT.INFO@PARI.COM

www.PARI.com/eRapid

eFlow®, eRapid®, eBase®, and Smartmask® are registered trademarks of PARI Pharma GmbH. Pulmozyme® is a registered trademark of Genentech, Inc.

© 2015 PARI Respiratory Equipment, Inc. All rights reserved.