

“ Needless to say my eyedrops now last twice as long... in the past many drops ended up on my eye lids or cheek! ”

AutoDrop® user



Ophthalmology

AutoDrop® and AutoSqueeze™



**NORTH AMERICA**

Owen Mumford Inc.  
1755 West Oak Commons Court  
Marietta, Georgia 30062, USA  
T : 800-421-6936 F : 770-977-2866  
E : info@owenmumfordinc.com

**HEAD OFFICE**

**INTERNATIONAL SALES & MARKETING**

Owen Mumford Ltd  
Brook Hill, Woodstock  
Oxfordshire OX20 1TU, United Kingdom  
T : 01993 812021 F : 01993 813466  
E : info@owenmumford.co.uk

**UK / IRELAND SALES & MARKETING**

T : 01608 645555 F : 01608 646546  
E : info@owenmumford.co.uk

**FRANCE**

T : 02 32 51 88 70 F : 02 32 21 96 95  
E : contact@owenmumford-fr.com

**GERMANY**

T : 06026/97750 F : 06026/977529  
E : info@owen-mumford.de

[www.owenmumfordinc.com](http://www.owenmumfordinc.com)



Making eye drop application easy

**AutoDrop®**

Eye Drop Guide...



**AutoSqueeze™**

OP/6000-LTD/IN/0012/01

 OWEN MUMFORD

# AutoDrop®



## Eye Drop Guide

Eye drops can be a challenge for many people, young or old. **AutoDrop®** and **AutoSqueeze™** have been developed to make self administration as easy as possible – improving patient compliance and their reliance on others.

**AutoDrop®** and **AutoSqueeze™** are re-usable, and provide a cost effective solution to make eye drop application simple.

### Unique design prevents blinking.

A pinhole draws the patient's gaze away from the descending drop, while a small lip holds the lower eyelid down, preventing blinking and improving patient compliance.

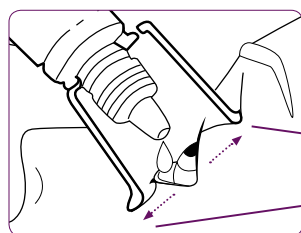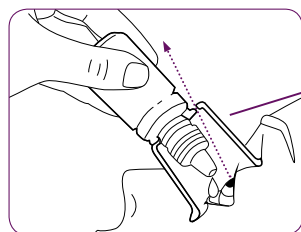

Proven **compatibility** with most eye drop bottles.

**AutoDrop®** has been tested with hundreds of eye drops to ensure compatibility.

**Prevents contamination and cross infection.**

The eye drop container is kept away from the eye at all times and the device is easy to clean.

**Delivers the drop in the right place every time.**

The nozzle of the bottle is positioned perfectly over the eye ensuring that drops are successfully delivered on the first attempt.

**AutoDrop®** and **AutoSqueeze™** have been designed to overcome common obstacles to patient compliance:

- Missing the eye with the drop.
- Blinking as the drop is delivered.
- Difficulty in squeezing the bottle to deliver the drop.

AutoDrop®	AutoSqueeze™
NCD: 08470-6000-01	NDC: 08470-6100-01



Developed together with the Royal National Institute for the Blind.



Scan the code for more information about **AutoDrop®**.

# AutoSqueeze™

Combine **AutoSqueeze™** with **AutoDrop®**, to make squeezing of the bottle easier.



Many people who use eye drops have other conditions that affect their dexterity and ability to squeeze small eye drop bottles.

**AutoSqueeze™** clips onto the bottle to provide a secure, ergonomic grip and extra leverage.

It can be used in conjunction with **AutoDrop®** or on its own.

**AutoDrop®** is compatible with most eye drop devices on the market, including those listed below:

NAME	MANUFACTURER	SIZE
Chloromycetin® (Redidrops)	Goldshield	5ml, 10ml
Ciloxan®	Alcon®	5ml
Gentamicin	Roche	10ml
Golden Eye®	Typharm	10ml
Oftaquix®	Kestral Ophthalmics	5ml
Alomide®	Alcon®	10ml
Betnesol®	UCB Pharma	10ml
Emadine®	Alcon®	5ml
Maxidex®	Alcon®	5ml, 10ml
Maxitrol®	Alcon®	5ml
Voltarol Optha®	Novartis	5ml

Tested as of November 2006. All Trademarks belong to named companies.

For a full up-to-date list of compatible eyedrops, visit [www.owenmumfordinc.com](http://www.owenmumfordinc.com)