



Breastcare





www.carbonfootprint.com

TRULIFE

Printed on
Carbon Neutral Paper

Contents

Trulife Breast Forms

- 1.** Cara
- 2.** Sublime Arís
- 3.** Impressions
- 4.** Bodicool Wave
- 5.** Silk
- 8.** Traditional
- 9.** Partial
- 11.** Swim & Activity
- 12.** Luna
- 13.** Leisure & Post Surgery
- 14.** Radiant Impressions
- 15.** Trulife Breast Form Size Guide

Trulife Bras

- 19.** Bra Fitting Guide
- 20.** Bra Measurement Guide
- 21.** Perfectly Seamless
- 22.** Stylish Fit
- 25.** Bras at a Glance
- 26.** Accessories
- 27.** Product Care





Breast Forms

After breast surgery, it can be difficult adjusting to the many changes in a woman's life and coming to terms with a new appearance. Trulife has dedicated over 60 years to the development of unique and innovative products; our focus has always been on making this transition a little easier and more comfortable for her.

Trulife breast forms are created in a wide range of shapes and sizes to complement a woman's natural beauty. Trends in surgery type and lifestyle are addressed with such products as ReCover, partial forms, leisure and swim forms and the availability of Radiant Impressions, our beautiful custom breast form.



Cara

100% Pure Silicone

Cara is customizable to provide enhanced comfort, cooling and a slight size adjustment when she needs it. This 100% pure silicone breast form can be worn several different ways to keep her comfortable and confident throughout her day.

*Cara
includes
accessories!*



Breathable Pure 100% Silicone:

Light and soft to wear against the chest.

Cooling Pad:

Unique TruCool Gel to absorb heat.

COOLMAX® Comfort Pad:

Soft absorbent foam pad to help manage perspiration.



153 Cara

Profile: Shallow

Sizes: 3-12

- Lightweight and cool
- 100% Silicone. Can be worn directly against the chest wall
- Natural skin-like feel
- Includes Comfort Pad and Cooling Pad
- Includes a COOLMAX® breast form cover



Sublime Arís

100% Pure Silicone


These revolutionary breast forms are made using 100% pure silicone and features a honeycomb core which ensures optimum breathability. The pure silicone design with tapered edges allows these breast forms to be worn directly against the skin for the most intimate fit.

Sublime offers softness and movement similar to natural breast tissue, beautiful shape and durability.

Subtle, hand engraved integrated nipple and areola

Dual layer technology with a super soft middle layer that mimics the natural movement of breast tissue and a thin front skin for an even softer and more natural feel

Full projection with smaller base for an intimate fit

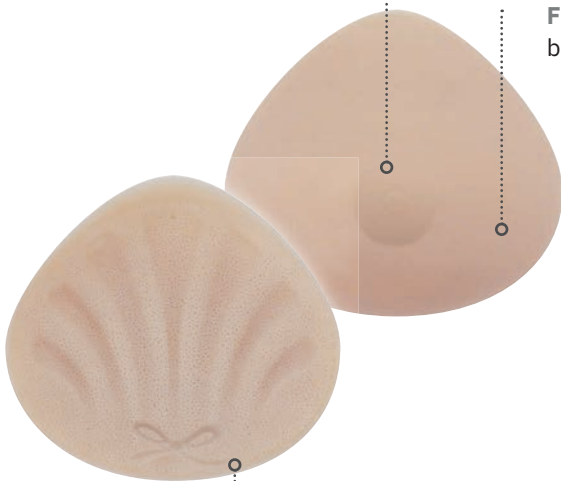

Back design features a feminine Sublime bow




151 Sublime Arís

Profile: Average

Sizes: 3-14

- Lightweight and cool
- 100% Silicone. Can be worn directly against the chest wall
- Natural skin-like feel
- Includes a COOLMAX® breast form cover

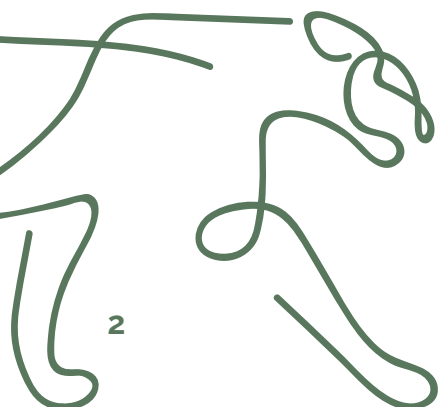

Impressions

100% Pure Silicone

Integrated nipple and areola enhances and restores femininity



Natural contours provide an intimate fit

Honeycomb core ensures optimum breathability



101 Impressions II

Profile: Shallow to Average

Sizes: 2-12

- 100% pure silicone
- Honeycomb design for enhanced breathability
- Features a shallow profile based on the world's #1 selling breast form, Trulife's 471 Silk Triangle
- Can be worn directly against the chest wall
- Includes a COOLMAX® breast form cover



102 Impressions II Encore

Profile: Shallow to Average

Sizes: 4-12

- 100% pure silicone
- Natural grip prevents shifting in bra
- Features a shallow profile based on the world's #1 selling breast form, Trulife's 471 Silk Triangle
- Includes a COOLMAX® breast form cover



Bodicool Wave

Lightweight


BodiCool Wave is made for the woman who prefers the feel of a lightweight blended silicone that offers a solution to body heat.

This breast form integrates an innovative system of TruCool Gel, a Trulife proprietary silicone blend that uses an advanced cooling technology.

BodiCool Wave provides a constant cooling effect by continuously drawing excess heat away from the body. Raised contact points of the channeling system promote enhanced air circulation allowing for optimal cooling.

TruCool Gel-filled channels provide optimal cooling

Raised channeling system promotes enhanced air circulation

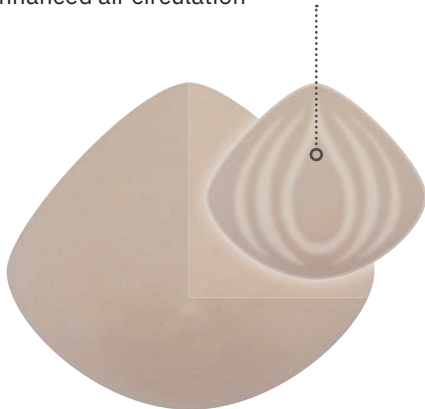

Silky-soft,
matte skin




495 Bodicool Wave

Profile: Average to Full

Sizes: 1-14

- TruCool Gel-filled channels provide optimal cooling
- Raised channeling system promotes enhanced air circulation
- Cooling properties regenerate once the breast form is removed from the chest wall
- Versatile triangle shape fits a variety of body and surgery types

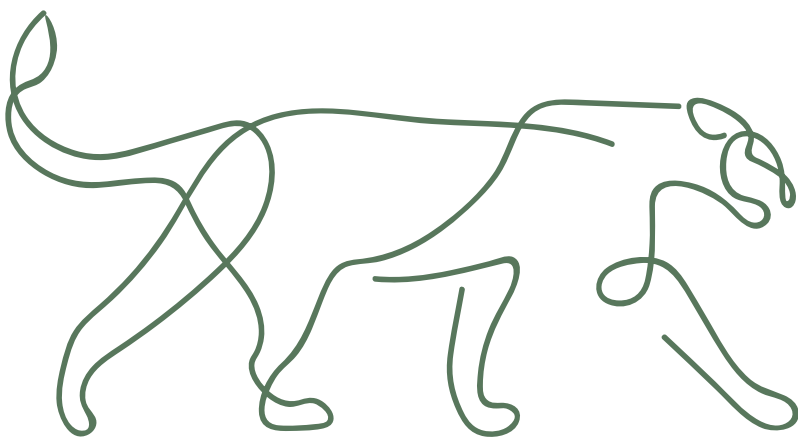

Silk

Lightweight

Silk defines what a breast form should be: soft and natural. The innovative skin is silky-soft and offers unparalleled softness and suppleness. All Silk breast forms offer a natural drape not found with other breast forms and are made from the highest quality of silicone for lasting comfort and softness.

Trulife's Silk collection offers the largest selection of shapes and sizes available on the market and are 30% lighter than traditional silicone breast forms.



471 Silk Triangle

Profile: Shallow to Average

Sizes: 1-17

- No. 1 selling breast form in the world
- Versatile triangle shape fits a variety of body and surgery types



485 Silk Curve

Profile: Shallow

Sizes: 2-14

- Elongated top transitions smoothly to the chest wall
- Soft, natural drape to match a mature breast
- Fits well in both underwire and softcup bras




473 Silk Teardrop

Profile: Shallow to Average

Sizes: 1-14

- Versatile teardrop shape tapers at top extension
- Can be positioned in any direction for optimal fit

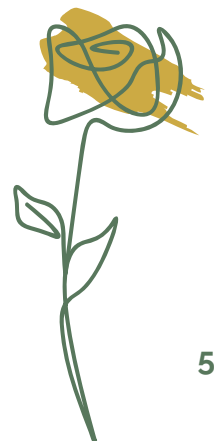

Silk

Lightweight



483 Silk Encore Triangle

Profile: Shallow to Average

Sizes: 1-17

- Triangle shape fits a variety of body and surgery types



477 Silk Flex

Profile: Average to Full

Sizes: 1-14

- Heart shape can be positioned in any direction for optimal fit
- Fits well with underwire bras



475 Silk A Supreme

Profile: Shallow to Average

Sizes: 1-14 L/R

- Asymmetrical shape replaces tissue at upper chest and underarm
- Concave back accommodates excess tissue
- Recommended for use with bilateral or extensive surgery



Silk

Lightweight



478 Silk Xtend

Profile: Average to Full

Sizes: 1-14 L/R

- Asymmetrical shape replaces tissue at upper chest and underarm
- Concave back accommodates excess tissue
- Recommended for use with bilateral or extensive surgery



472 Silk Triangle Plus

Profile: Average to Full

Sizes: 1-17

- Triangle shape fits a variety of body and surgery types
- Fits well with both underwire and softcup bras
- Comfortable back provides additional fullness and projection




481 Silk Ultima Triangle

Profile: Full

Sizes: 1-14

- Versatile triangle shape fits a variety of body and surgery types
- Fits well with both underwire and softcup bras

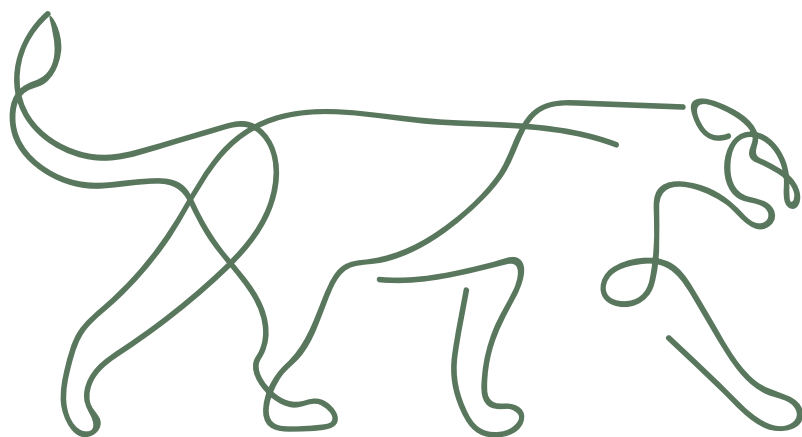

Traditional

Soft Silicone

Traditional breast forms are made from a soft silicone gel which offers more weight to mimic natural breast tissue drape and movement.



509 E Supreme

Profile: Average to Full

Sizes: 1-17

- Translucent colour
- Versatile teardrop shape tapers at top extension
- Can be positioned in any direction for optimal fit



508 Symphony

Profile: Average to Full

Sizes: 1-14

- Versatile triangle shape fits a variety of body and surgery types
- Fits well with both underwire and softcup bras



701 Duette Triangle

Profile: Average to Full

Sizes: 1-14

- Adaptable triangle shape fits a variety of body and surgery types
- Provides more fullness than other triangles
- Fits well with both underwire and softcup bras
- Flat back fits closely and securely to chest wall




571 Bella

Profile: Average to Full

Sizes: 3-14

- Translucent colour
- Symmetrical triangle shape fits a variety of body and surgery types

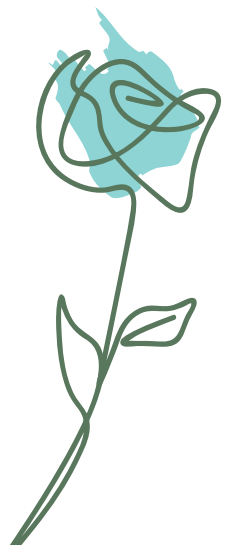

Partial

Soft Silicone

Trulife offers soft partial silicone breast forms designed to provide the lightest and most natural symmetry for post-lumpectomy, post-reconstruction and enhancement of underdeveloped breasts. All partial breast forms feature tapered edges to ensure a flawless, smooth profile under clothing.



533 Triangle Partial

Sizes: 3-9

- Triangle shape offers thin, overall coverage
- Soft silicone with an exterior silky-soft skin



531 Teardrop Partial

Sizes: 1-11


- Soft silicone with exterior silky-soft skin
- Moderate and partial coverage for a seamless look



822 ReCover Shell

Sizes: S, M, L, XL

- Designed for staged reconstruction
- Fine, thin tapered edges provide a seamless, intimate fit
- Translucent appearance to adapt to skin colour
- Three layers of silicone inserts for customized fullness and projection

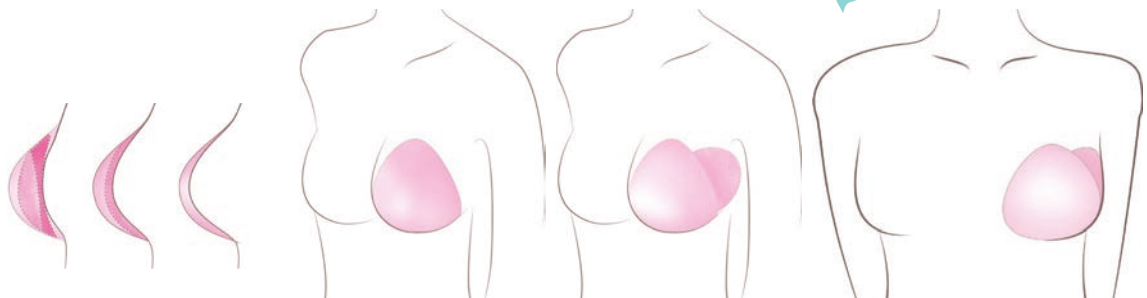

Partial

Soft Silicone



535 Triangle Partial Encore

Sizes: 3-9

- Triangle shape offers thin, overall coverage
- Soft silicone with an exterior silky-soft skin



356 Teardrop Partial Encore

Sizes: 1-11

- Soft silicone with exterior silky-soft skin
- Moderate and partial coverage for a seamless look



Swim & Activity Lifestyle

Trulife swim forms are ideal for all water activities. They retain their shape in and out of water. They can be used in chlorinated and salt water.

ActiveFlow is the newest concept in active lifestyle breast form design. The perfect combination of beads and moulded spacer fabric provide a flattering, natural profile, whether active or lying down.



630 ActiveFlow

Sizes: Original Mesh Sizes: 1-14
Black Sizes: 3-14

Colours: Original mesh, Black

- Popular for all activities including swimming and jogging
- Lightweight monoprene beads do not absorb water and dry quickly
- The beads also move and flatten when sun-bathing



641 Tropez

Sizes: 1-14

- Firm translucent triangle shape fits a variety of body and surgery types
- Concave back allows for buoyancy, water flow and shape retention
- New soft silky skin!



Luna

Soft Leisure

Trulife's newest Leisure and Post Surgery Breast Form

Our soft, lightweight foam and fibre-fill prosthesis is recommended for leisure, post-surgery and moderate activities. It is also eco-friendly! Fibre-fill is made from recycled plastic bottles. Good for the environment and good for post-surgery needs.



622 Luna

Sizes: 3-12 L/R

- Asymmetrical shape has slight extension to fit under the arm or the upper chest
- Fibre filling made from recycled plastic bottles
- Soft Bamboo material with moisture wicking properties
- A gently weighted front helps to restore balance



Leisure & Post Surgery

Soft Foam

Ideal collection of soft, lightweight foam breast forms for post-surgery, leisure and moderate activities. Each form features a removable cotton cover for ease of washing.



611 Tri-Leisure

Sizes: 3-14

- Triangle shape fits a variety of body and surgery types
- Inner weighted core provides minimal weight
- Recommended for leisure wear



616 Tri-Featherweight

Sizes: 1-14

- Triangle shape fits a variety of body and surgery types
- Soft, non-weighted foam
- Recommended for post-surgery




607 FirstFit

Sizes: 3-12

A woman's first fit should be soft and gentle, easily adjustable and lightweight for optimized comfort. Trulife's FirstFit form offers an ideal solution for post-surgery and leisure wear.

- Seamless moulded front with a breathable and moisture-wicking COOLMAX® back
- 100% hypoallergenic polyester fluff
- Easily adjusted to any fullness required
- Extra fibre-fill included for adjustment

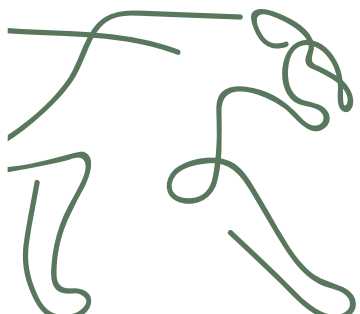

Radiant Impressions

Custom Breast Form

NON-SURGICAL BREAST & NIPPLE RESTORATION


Trulife is renowned as the quality leader in crafting a handmade custom breast prosthesis. We have over two decades of experience creating the most realistic custom breast prosthesis for women as an alternative to reconstruction, an answer to failed reconstruction, and as a beautiful prosthetic option for all women who have had breast surgery.

- Custom made to a woman's unique contour, shape and profile
- May help restore a sense of wholeness and confidence
- Enhanced softness to match natural breast tissue

With the use of a laser scanner, a 3D image is quickly created and easily sent to the Radiant Impressions Production Studio, where Design Engineers create a beautiful custom breast prosthesis that is unique to the surgery site, skin tone, colour and characteristics such as freckles, veins, nipples and tattoos.

Semi custom and custom nipples also available.

Contact Trulife about our Radiant Impressions Scan Events and/or becoming a Radiant Impressions Custom Breast Prostheses Provider.

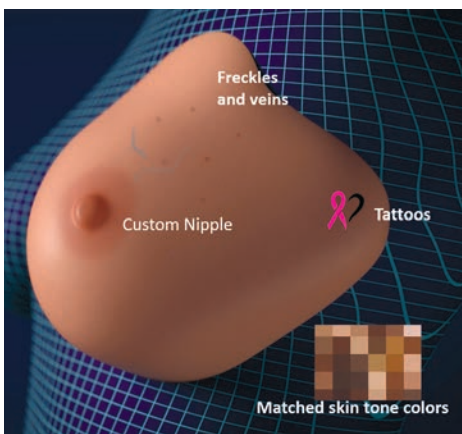

TRULIFE BREAST FORM

Size Guide

HOW TO USE THIS SIZE GUIDE

1. Select a cup size and follow the column down to where it intersects with the appropriate bra band size.
2. Follow this row to the right for suggested size by product to try on your patient.
(See example highlighted below).



								100% Pure Silicone				Lightweight					
MODEL								153	151	101	102	495	471	485	473		
NAME								Cara	Sublime Aris	Impressions II	Impressions II Encore	BodiCool Wave	Silk Triangle	Silk Curve	Silk Teardrop		
SIZES								3-12	3-14	2-12	4-12	1-14	1-17	2-14	1-14		
PAGE								1	2	3	3	4	5	5	5		
Size Guide																	
PROFILE								Shallow	Average	Shallow to Average	Shallow to Average	Shallow to Average	Shallow to Average	Shallow	Shallow to Average		
CUP SIZE																	
AA	A	B	C	D	DD	DDD											
BRA BAND SIZE (IN/CM)																	
28/60													1	1		1	
30/65	28/60									2			2	2	2	2	
32/70	30/65	28/60						3	3	3			3	3	3	3	
34/75	32/70	30/65	28/60					4	4	4	4		4	4	4	4	
36/80	34/75	32/70	30/65					5	5	5	5		5	5	5	5	
38/85	36/80	34/75	32/70	30/65				6	6	6	6	6	6	6	6	6	
40/90	38/85	36/80	34/75	32/70				7	7	7	7		7	7	7	7	
42/95	40/90	38/85	36/80	34/75	32/70			8	8	8	8		8	8	8	8	
44/100	42/95	40/90	38/85	36/80	34/75			9	9	9	9		9	9	9	9	
46/105	44/100	42/95	40/90	38/85	36/80	34/75		10	10	10	10		10	10	10	10	
48/110	46/105	44/100	42/95	40/90	38/85	36/80		11	11	11	11		11	11	11	11	
50/115	48/110	46/105	44/100	42/95	40/90	38/85		12	12	12	12		12	12	12	12	
52/120	50/115	48/110	46/105	44/100	42/95	40/90			13				13	13	13	13	
	52/120	50/115	48/110	46/105	44/100	42/95			14				14	14	14	14	
		52/120	50/115	48/110	46/105	44/100							15				
			52/120	50/115	48/110	46/105							16				
				52/120	50/115	48/110							17				

Breast Form Size Guide continued next page.

TRULIFE BREAST FORM Size Guide



HOW TO USE THIS SIZE GUIDE

1. Select a cup size and follow the column down to where it intersects with the appropriate bra band size.
2. Follow this row to the right for suggested size by product to try on your patient.

(See example highlighted below).



Lightweight						
MODEL	483	477	475	478	472	481
NAME	Silk Encore Triangle	Silk Flex	Silk A Supreme	Silk Xtend	Silk Triangle Plus	Silk Ultima Triangle
SIZES	1-17	1-14	1-14 L/R	1-14 L/R	1-17	1-14
PAGE	6	6	6	7	7	7

Size Guide												
PROFILE							Shallow to Average	Average to Full	Shallow to Average	Average to Full	Average to Full	Full
CUP SIZE												
AA	A	B	C	D	DD	DDD	▼	▼	▼	▼	▼	▼
BRA BAND SIZE (IN/CM)												
28/60							1	1	1	1	1	1
30/65	28/60						2	2	2	2	2	2
32/70	30/65	28/60					3	3	3	3	3	3
34/75	32/70	30/65	28/60				4	4	4	4	4	4
36/80	34/75	32/70	30/65				5	5	5	5	5	5
38/85	36/80	34/75	32/70	30/65			6	6	6	6	6	6
40/90	38/85	36/80	34/75	32/70			7	7	7	7	7	7
42/95	40/90	38/85	36/80	34/75	32/70		8	8	8	8	8	8
44/100	42/95	40/90	38/85	36/80	34/75		9	9	9	9	9	9
46/105	44/100	42/95	40/90	38/85	36/80	34/75	10	10	10	10	10	10
48/110	46/105	44/100	42/95	40/90	38/85	36/80	11	11	11	11	11	11
50/115	48/110	46/105	44/100	42/95	40/90	38/85	12	12	12	12	12	12
52/120	50/115	48/110	46/105	44/100	42/95	40/90	13	13	13	13	13	13
	52/120	50/115	48/110	46/105	44/100	42/95	14	14	14	14	14	14
		52/120	50/115	48/110	46/105	44/100	15				15	
			52/120	50/115	48/110	46/105	16				16	
				52/120	50/115	48/110	17				17	

TRULIFE BREAST FORM

Size Guide



Traditional				Swim & Activity		Leisure & Post-Surgery			
509	508	701	571	630	641	622	611	616	607
E Supreme	Symphony	Duette Triangle	Bella	ActiveFlow	Tropez	Luna	Tri-Leisure	Tri-Feather-weight	First Fit
1-17	1-14	1-14	3-14	1-14	1-14	3-12 L/R	3-14	1-14	3-12
8	8	8	8	11	11	12	13	13	13
Average to Full	Average to Full	Average to Full	Average to Full	N/A	N/A		N/A	N/A	N/A
▼	▼	▼	▼	▼	▼	▼	▼	▼	▼
1		1		1				1	
2	1	2		2	1			2	
3	2	3	3	3	2	3	3	3	3
4	3	4	4	4	3	4	4	4	4
5	4	5	5	5	4	5	5	5	5
6	5	6	6	6	5	6	6	6	6
7	6	7	7	7	6	7	7	7	7
8	7	8	8	8	7	8	8	8	8
9	8	9	9	9	8	9	9	9	9
10	9	10	10	10	9	10	10	10	10
11	10	11	11	11	10	11	11	11	11
12	11	12	12	12	11	12	12	12	12
13	12	13	13	13	12		13	13	
14	13	14	14	14	13		14	14	
15	14				14				
16									
17									






Bras

Trulife bras offer women everyday confidence and comfort, fashionable choices and supportive styles for a variety of lifestyles.

Our bras are hand crafted in North America using quality fabrics and components that are inspected to meet Trulife's rigorous specifications. Each bra is carefully sewn by highly skilled workers, followed by a thorough inspection before being packaged for delivery.

All Trulife bras feature pockets on both sides made from breathable soft mesh or COOLMAX® fabric.

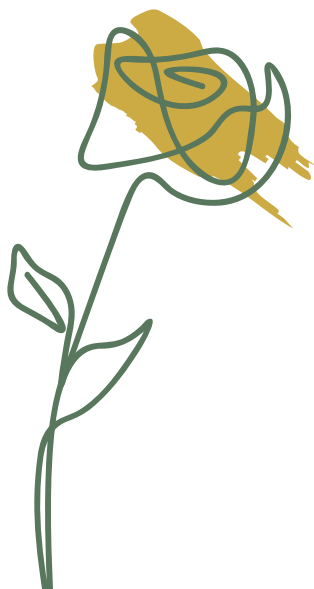


Bra Fitting Guide

ACHIEVING THE PERFECT FIT

When fitting a bra, there are three distinct shapes that describe the wearer's profile: shallow, average and full.


The type of profile is not dependent on cup capacity. For example a woman who wears a 'D' cup can have a shallow profile, just as a woman who wears an 'A' cup can have a full profile.

Below is a guide to help demonstrate the different profile types:

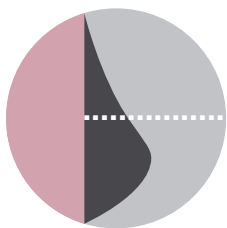


Shallow

The breast tissue drapes in the lower cup of the bra and is often seen on the mature woman with soft breast tissue. The shallow profile woman requires a bra with a large bottom cup and small upper cup.

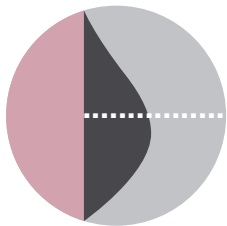

Average

There is an equal amount of breast tissue on the top and bottom of the breast, with a more round appearance in the upper cup area. The average profile woman usually has firm breast tissue and a round profile.

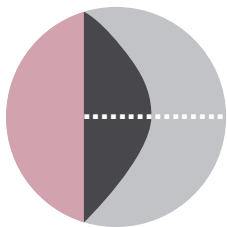


Full

Producing the most round appearance, the full profile woman requires a bra that provides not only bottom cup support but also support and coverage for the upper part of the breast. The full profile can also be broad without much projection, needing a bra that provides containment such as M-frame styles that separate as well as support.

The bras appearing on the following pages come with a suggested profile fit. Knowing your customer's profile and the best Trulife bras for each profile can provide a personalized, efficient and rewarding fitting experience.

Additionally, each Trulife bra featured in this catalog specifies innovative design features and benefits, for example COOLMAX®, and other performance enhancing materials

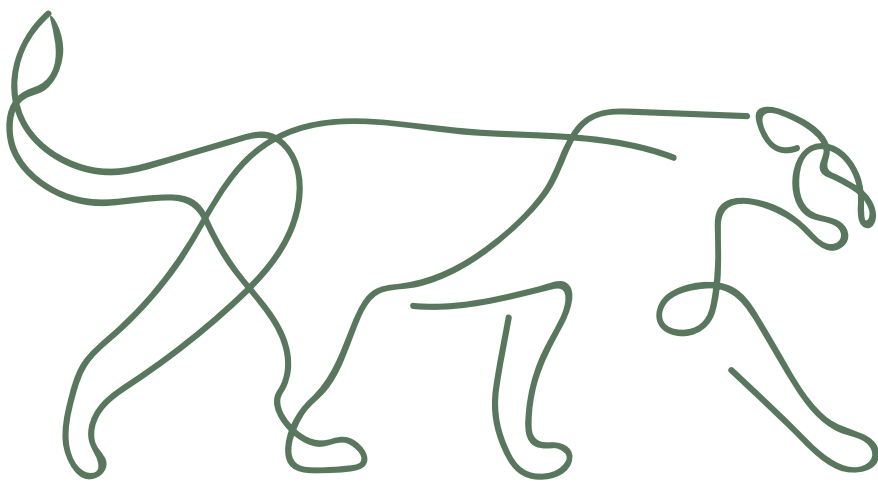

BRA MEASUREMENT GUIDE

Finding the perfect bra is essential for a woman to look and feel her best. It's important to choose a bra style that will adequately support the weight of the breast form and hold it securely against the body.

The following guidelines on measuring for bra size are only a suggested technique and starting point. Taking time to find the correct size will ensure a comfortable fit. Keep in mind that body changes due to weight gain or loss, can also change the bra size. Measurements should be checked again within six months to a year.

Bra Band Measurement

To find your bra size, measure all around the body underneath the bust. If the measurement is an EVEN number add 4 inches, e.g. 30 + 4 = 34 Band Size. If the measurement is an ODD number add 5 inches. e.g. 29 + 5 = 34 Band Size.

Cup Measurement

Measure over the most projected part of the remaining breast from center of sternum to center of spine; double the number.

Bra Size

Determine the difference between the cup measurement and band measurement and refer to chart below:

Difference between Bra Band and Cup Measurements							
CUP SIZE	AA	A	B	C	D	DD	DDD
DIFFERENCE (INCHES)	-1	0	1	2	3	4	5
DIFFERENCE (CM)	8/10	10/12	12/15	15/18	18/20	20/23	23/25

Trulife Bra Features

All Trulife bras provide features essential for a beautiful fit:

1. Adjustable shoulder straps
2. Elastic binding around neckline
3. Full cup coverage
4. High center front fabric
5. Generous underarm support
6. Good separation between cups
7. Soft breathable pockets
8. Supportive underband
9. Two or more hook & eye fasteners (not shown)



Lisa is wearing 210 Barbara in Nude (page 23)

Perfectly Seamless



4030 Emily

Style: Seamless Smooth Soft Cup

Profile: Shallow to Average

Colours: White, Latte, Black

CUP	A	B	C	D
IN	34-40	34-42	34-42	34-42
CM	75-90	75-95	75-95	75-95

Emily seamless softcup brings together all the features that make a perfect bra: secure fit, stylish and comfortable.

- Moisture-wicking extended COOLMAX® pockets with underarm opening for comfort and breathability



Stylish Fit



4019 Jessica

Style: Softcup

Profile: Shallow to Average

Colours: White, Nude, Black

CUP	A	B	C	D	DD	DDD
IN	34-42	34-44	34-48	34-48	36-48	36-46
CM	75-95	75-100	75-110	75-110	80-110	80-105

This fashionable and preferred style features a built-in lace camisole panel tailored to give you all-day security and comfort. Seamed soft cup design offers excellent support. An absolute favorite!

- Lace accents and lace camisole panel for an exceptionally feminine style
- Back-adjusting elastic straps and hook closure
- Pockets feature COOLMAX® fabric that manages moisture to help keep you cool and dry



Stylish Fit

210 Barbara

Style: Lace Accent Softcup

Profile: Shallow to Average

Core Colours: White, Nude, Black

CUP	A	B	C	D	DD	DDD
IN	32-42	32-52	34-52	34-52	34-48	34-48
CM	70-95	70-120	75-120	75-120	75-110	75-110

Fashion Colour: Amethyst

CUP	A	B	C	D	DD	DDD
IN	34-42	34-44	34-48	34-48	36-48	36-46
CM	75-95	75-100	75-110	75-110	80-110	80-105

Our best-selling style! Three core and Amethyst fashion colour make this a favorite. This seamed style offers unparalleled support, and is highlighted with a lace trim that accentuates its loveliness.

- Extremely versatile fit
- Back-adjusting elastic straps and hook closure
- Soft and breathable moisture-wicking COOLMAX® pockets for coolness and comfort



Stylish Fit

S210 Barbara (Soft)

Style: Lace Accent Softcup


Profile: Shallow to Average

Colour: Exotic Black

CUP	A	B	C	D	DD	DDD
IN	34-42	34-44	34-48	34-48	36-48	36-46
CM	75-95	75-100	75-110	75-110	80-110	85-105

The Soft 210 Barbara offers the same great fit as our Original style. Offered in a fun exotic Black fashion colour! The soft touch material provides a luxurious feel that is lace-free. A new wardrobe must-have!

- Extremely versatile fit
- Back-adjusting elastic straps and back closure
- Soft and breathable moisture-wicking COOLMAX® pockets for coolness and comfort

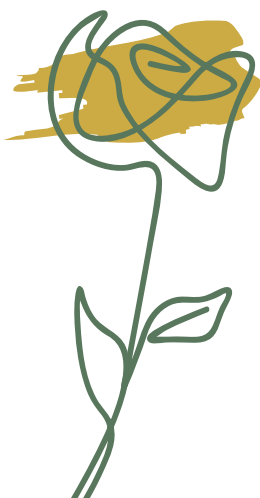

Trulife

Bras at a Glance

Key	
○	Softcup
◡	Underwire
S	Shallow to Avg
○	White
●	Nude
●	Latte
●	Black
●	Amethyst
⊗	Exotic Black

	Perfectly Seamless	STYLISH FIT			
MODEL	4030	4019	210*	210*	S210*
NAME	Emily	Jessica	Barbara	Barbara	Barbara (Soft)
STYLE	○	○	○	○	○
PROFILE	S	S	S	S	S
COLOURS	○ ● ●	○ ● ●	○ ● ●	●	⊗
SIZES					
A	34-40 in 75-90 cm	34-42 in 75-95 cm	32-42 in 70-95 cm	34-42 in 75-95 cm	34-42 in 75-95 cm
B	34-42 in 75-95 cm	34-44 in 75-100 cm	32-52 in 70-120 cm	34-44 in 75-100 cm	34-44 in 75-100 cm
C	34-42 in 75-95 cm	34-48 in 75-110 cm	34-52 in 75-120 cm	34-48 in 75-110 cm	34-48 in 75-110 cm
D	34-42 in 75-95 cm	34-48 in 75-110 cm	34-52 in 75-120 cm	34-48 in 75-110 cm	34-48 in 75-110 cm
DD		36-48 in 80-110 cm	34-48 in 75-110 cm	36-48 in 80-100 cm	36-48 in 80-110 cm
DDD		36-46 in 80-105 cm	34-48 in 75-110 cm	36-46 in 80-105 cm	36-46 in 80-105 cm
PAGE	21	22	23	23	24



Accessories



Breast Form Covers

- Most available in light and dark tone
- Replacement covers
- Sized to match breast form
- One unit per package



Product Care

Breast Forms

Trulife silicone breast forms are designed to be worn by women of all ages, with different lifestyles and surgery types. Our breast forms are available in an extensive range of shapes, sizes, textures and weights to ensure every woman can find the breast form that best suits her needs. A little care and attention will ensure long-lasting comfort and satisfaction.

To clean a silicone breast form, hand wash gently in lukewarm water with a mild soap. Dry with a soft towel, do not wring or rub. Avoid using powders, perfumes or hot water directly on the breast form. Care should be taken to avoid sharp objects that may damage the breast form.

Routine Cleaning of 100% Pure Silicone products

Wash the prosthesis as often as required by following the steps below:

1. Apply a pea size amount of soap to a wet cloth and use this to wash the breast form, front and back. Wipe with a dry cloth.
2. Alternatively use a wipe.
3. As often as needed, bathe the breast form by submerging in warm water and soap.
4. The honeycomb core of your 100% silicone product will retain water. To extract this water, place the back of the breast form on an absorbent towel and press firmly downwards. Repeat until most of the water is extracted, then leave to dry overnight, with the aerated back facing up.
5. To further enhance the silky soft feel, you may wish to apply talc evenly over the front of the breast form.

It is important to wear the correct size bra that fits properly and does not force the breast form into a different shape. This may cause the breast form to crease or permanently damage the breast form.

It is important to handle the breast form gently when inserting and removing it from the cover or pocket of the bra. Repeated squeezing or pulling may cause damage over time.

When not wearing the breast form, store it in the original packaging and cradle provided, away from heat. This will allow the breast form to maintain its shape and prevent accidental harm.

Bras

Trulife bras feature quality fabrics, components and construction. We recognize that every woman is an individual with unique needs and desires.

Proper care will prolong the life of Trulife bras and provide lasting comfort, security and an overall natural look and feel. We recommend hand washing in warm water with mild soap. Hang dry for best results. Avoid use of fabric softeners when bra contains COOLMAX® fabric. Do not bleach.

With care, your Trulife bra will beautifully support and secure your breast form, for a natural, feminine look.



www.trulife.com