



SAY HELLO
TO EXPERTISE

iStir MS10
(Magnetic Stirrer)

iStir MS10 (Magnetic Stirrer)

When looking for magnetic stirrers, users look for a device which delivers reliable performance and years of durability, but most manufacturers tend to oversell and under deliver which results in a frustrating experience for the user. Magnetic stirrers come with promises of exceptional performance, quality components and high durability, but in just a short time, the reality becomes apparent. Users face grave issues like the top getting corroded or the motor is not powerful enough to stir large quantities of solutions, its speed remains low etc.


With detailed insights about the requirements and needs of users and their feedback of existing stirrers, the R&D department at Neuation came up with the design of iStir MS10 magnetic stirrer to develop a product that delivers high performance without compromising on flexibility or durability. With its stainless steel top surface, the device is extremely resistant to corrosion while its long lasting DC motor runs at top speed without needing much of a maintenance. The same motor is capable enough to stir large quantities of solutions at once for higher productivity and less wastage of time and effort.



Product Features



DC Motor For Quite &
Long Life Performance

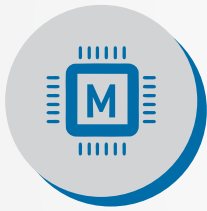

Microprocessor Controlled
For Efficient Performance
Even Under High Loads.



Variable Speed Setting
From 200-2200 RPM With
Increments Of 10 Rpm



Variable Timer Setting
From 1 To 999 Minutes
Along With An Infinite Mode.

Stirring Capacity
Up To 10 Litres



Pulse Mode Function With
Run-pause-run Function For
Thorough Mixing Of Solutions.



Closed Design Allows
Easy Cleaning And
Protects The Equipment
From Accidental Spills.



Product Features



Non-slip Strength Rubber Feet For Maximum Stability Even At High-speeds.



IP43 Class Compliant Rating For Unmatched Durability.



No corrosive SS top increases life of instrument

Compact Footprint Saves Valuable Bench Space And Also Makes It Suitable For Use In Hoods And Cold Room.



Digital Display Allows The User To Gather All The Relevant Information In A Single Glance.



Specifications

Specifications	iStir MS10
Capacity	10 Litres
Variable Speed	200-2200 RPM
Run Time	1-999 mins & infinite mode
Pulse Mode	Yes, (Run - pause - Run Mode)
Digital Display	Yes
Weight	1.4 kg
Dimension	156 x 248 x 104 mm
Protection class	IP 21

Ordering Information

iStir MS10		
Product Name	Cat #	Details
iStir MS 10	HS20-2314 EU/US/AU/UK	Magnetic Stirrer

Standard Accessories	
Part Name	Cat #
Stir Bar	B1



AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

* OEM request can be entertained. Specifications can be changed without notice for quality improvement.