



IT CAN'T GET
BETTER THAN THIS

iSwix VT
Digital Vortex Mixer

iSwix VT

Digital Vortex Mixer

Vortex Mixers play a crucial part in any laboratory setup hence the need for reliable performance and durability, but most mixers have issues with lackluster performance owing to sub-par quality components and raw-materials used. Hence, these devices leave a lot to be desired and end-up being a waste of space and investment.

Developed after years of R&D, Neuation's best-in-class iSwix VT vortex mixer provides the best of performance due to its top-of-the-line internal components such as long lasting maintenance free BLDC motor. With class-leading features like inbuilt counter balance and 3-way switch, the device makes working even easier.



Product Features




Brushless DC motor for quiet performance and a maintenance-free long life

Inbuilt counter balance ensures stable operation for complete peace-of-mind



Adjustable speed setting from 300 to 4200 RPM

Microprocessor based closed loop control for efficient performance even under high loads

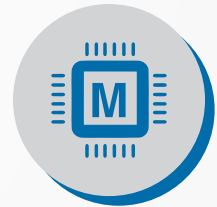

Digital Display shows the speed and time toggles after every 5 seconds

3-way switch for Touch, Continuous & Stand-by mode for unmatched convenience




Pulse Mode function with Run-Pause-Run function for thorough mixing of solution

Product Features



Orbital diameter
of 4 mm

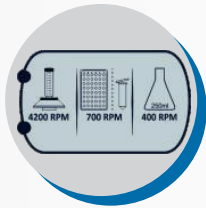


Visual indication
marker for max speed
with attachments



IP 42 Provides An
Assurance Of
Durability And Quality.


Small footprint saves
valuable bench space and
can also be used inside
hoods and cold rooms

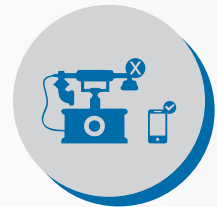

Variable timer setting
from 1 to 999 minutes and
continuous mode



Standard Attachment

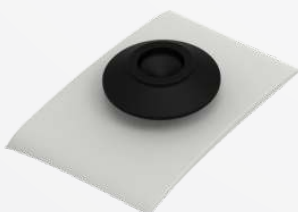


Cup Attachment


Optional Attachments



Micro Centrifuge Tube
20 x 1.5/2ml

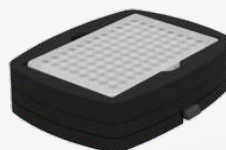

96 Well Micro Plate

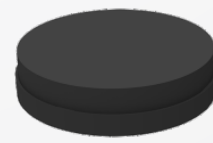

Round Disc



Flask





Specifications

iSwix VT	
Max Speed	4200 RPM
Orbital diameter	4 mm
Motor type	Brushless DC
Digital Display	Yes
Dimensions (W X D X H)	205 x 136 x 138.5 mm
Speed Setting	300-4200 RPM
Timer Setting	0-999 mins (min/sec display)
Pulse mode	Yes, Programmable Pulse
Max Load Capacity	500 gm
Weight	2.8 kg
Optional Attachments (not included)	Universal platform for attaching either - microplate / microtube / flask Attachments

Ordering Information

iSwix VT		
Product Name	Cat #	Details
iSwix VT	VM-2110 EU/US/UK/AU	Digital Vortex Mixer
Attachment	VM-A1	Cup Attachment

Accessories

 <p>Microplate CAT. NO. # VM-A2</p>	 <p>Flask CAT. NO. # VM-A3</p>	 <p>Micro Centrifuge Tubes CAT. NO. # VM-A4</p>	 <p>Round Disc CAT. NO. # VM-A6</p>
---	--	--	---



AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

* OEM request can be entertained. Specifications can be changed without notice for quality improvement.