


Pleated Style Face Mask Technical Data Sheet

The MM100 is a non sterile disposable face mask manufactured in Korea and is best suited for general public use. This mask has a multilayer design and is intended to cover the wearer’s nose and mouth.

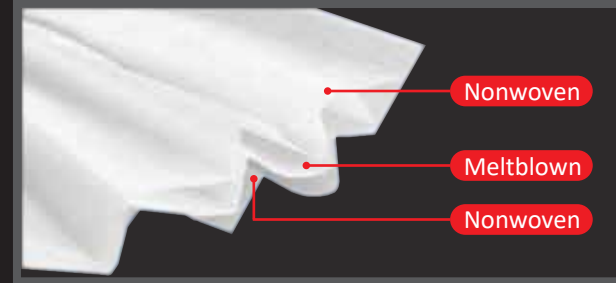

Model Number: MM100

“CDC recommends that people wear a face covering to cover their nose and mouth in the community setting. This is an additional public health measure people should take to reduce the spread of Covid-19 in addition to (not instead of) social distancing, frequent hand cleaning and other everyday preventive actions.”

	Face Mask	Surgical Mask
User	General Public	Health Care Workers and patients in healthcare settings
Risk Factor	Low Risk Situations	Moderate Risk Situations
Design	Designed to contain coughs and sneezes and prevent disease transmission to others – NOT to protect the wearer from breathing in aerosols.	Designed for health care workers to stop droplets being spread by the wearer – NOT to protect the wearer from breathing in aerosols.
When To Wear	Recommended for use when a person can't perform social distancing. • When coughing or sneezing. • Using public transportation. • Shopping and working at essential businesses like grocery stores and pharmacies.	Recommended for health care providers during single or multiple patient interactions or routine health procedures. Surgical masks are also recommended when N95s are not available.
Fit Testing Required	No Scarves and bandanas can be used if necessary.	No
Use Limitations	Homemade/Cloth Masks: • Should be washed after each use. • Should not be worn damp or when wet from spit or mucus.	Ideally should be discarded after each patient encounter. Extended use is preferable to reuse. Both are important and viable options when supplies are limited. Extended use is preferable. Follow CDC guidance on Strategies to Optimize the Supply of PPE Equipment



Features

- In process of approval to be used as a surgical mask alternative
- Flat, easy on/off mask with standard pleats
- 3 filter protection
- Adjustable nose clip
- Comfortable earloop design
- Easy to breathe through
- Made in Korea

Intended Use

Using a face mask may limit the spread of User’s respiratory droplets. For single use.



MADE IN KOREA


Materials	
Outer Layer	Polypropylene
Filter Layer	Polypropylene
Inner Layer	Polypropylene
Nose Clip	Stainless/Polypropylene
Elastic Band	Polypropylene/Polyester

Packaging	
Qty/ Box	Boxes/ Case
50	20


Fitting Instructions



1 Sanitize your hands thoroughly before touching the mask

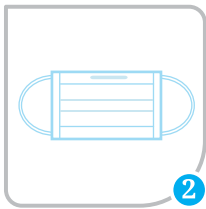

2 Inspect the overall condition of the mask to ensure there are no defective parts

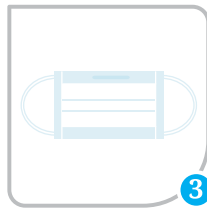

3 Face the inside of the mask and turn upright



4 Place the mask over your mouth and nose and place straps behind ears



5 Using both hands, press down the top metal edge over the bridge of your nose to create the best possible seal



6 Adjust the mask as needed and fit check one last time

Practice Social Distancing



Keep at least 6 feet apart



Warnings and Limitations

- Always confirm the product is suitable for the application
- Keep the mask on during all periods of exposure
- Fit testing should be performed before entering exposed areas
- Inspect all parts of the mask to ensure it is in good condition before use
- Dispose according to applicable regulations
- Do not alter the product in any way
- Masks should be stored in a cool and dry area