


# GINGIVAL CUFF LINK

## PRESERVING THE EMERGENCE PROFILE



GCL EMERGENCE PROFILE

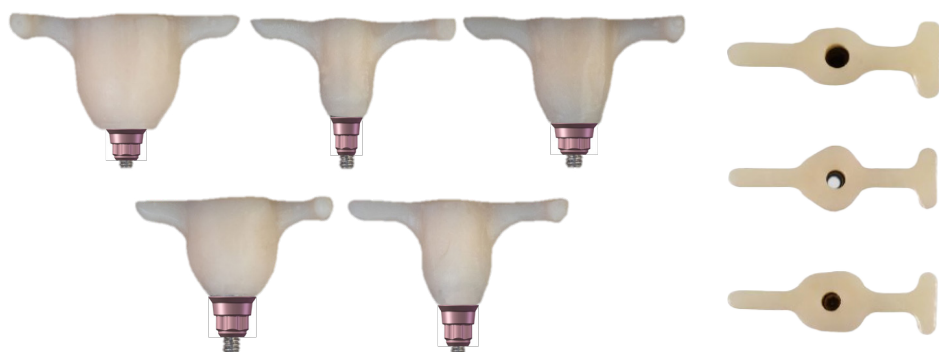

CONVENTIONAL ABUTMENT

### WHAT'S IN YOUR TOOLBOX?

A unique provisional healing cuff that allows the dentist to quickly shape and form tissue into an anatomically designed profile.

- A completely customizable chairside solution
- Tooth specific shapes to maximize optimal gingival healing
- Compatible with flowable composite and other restorative materials
- Use at any stage of the implant process
- Use in analog and digital workflows



GENESIS

PRIMA PLUS

PRIMA CONNEX

# POWER OF ONE

## OVER A DECADE OF RESTORATIVE PREDICTABILITY

The TiLobe® connection, with over a decade of restorative predictability, provides unlimited analog and digital flexibility. The self-sealing conical taper minimizes the microgap between implant and abutment, maintaining bone and soft tissue health. The rounded six-lobe design provides even-load distribution and a solid anti-rotation mechanism while the built-in platform shifts across all diameters, stabilizing bone and soft tissue.



# GINGIVAL CUFF LINK

## CASE STUDY



Socket of tooth number 8.



Seated Gingival Cuff Link.



The "T-handle", used for manipulating and orienting the Gingival Cuff Link, has been removed.



Try-in of the provisional crown prior to gingival grafting.



Healing with provisional crown healing at four months.



Provisional crown removed. Soft tissue healing at four months with enhanced buccal soft tissue thickness.

### Gingival Cuff Link, TiLobe Connection

DESCRIPTION	TOOTH #	REF #
Ø 3.8 Upper Central Incisor	8, 9	GCL-45370K
Ø 3.5/3.8 Upper Lateral Incisor	7, 10	GCL-45371K
Ø 3.8 Upper Canine	6, 11	GCL-45372K
Ø 3.5/3.8 Lower Incisor	23, 24, 25, 26	GCL-45376K
Ø 3.8 Lower Canine	22, 27	GCL-45373K
Ø 4.5 Upper Central Incisor	8, 9	GCL-45374K
Ø 4.5 Upper Canine	6, 11	GCL-45375K
Ø 4.5 Lower Canine	22, 27	GCL-45377K



Keystone Dental Group is an authorized distributor of GCL Systems.