

Dynamic CONICAL MC

Platform Switching


The platform switch maintains crestal bone and soft tissue volume around the implant platform.^{1,2,3}

Thread Design

Double-lead thread with a reverse buttress profile with 0.8 mm thread pitch.

Cutting Apex

Designed for controlled implant placement in challenging sites.

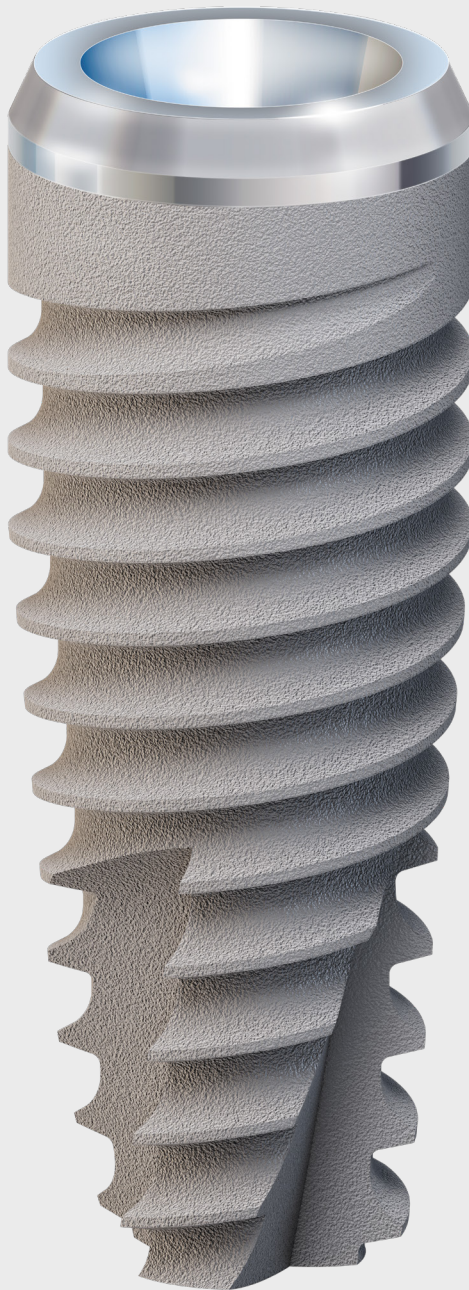

Conical Connection

Tightly sealed conical interface with a deep internal hex connection offers high mechanical strength and platform switching to optimize bone preservation and soft tissue.^{4,5,6}

Machined Collar

0.8 mm machined collar for different clinical applications and surgical flexibility.

KLEAN Surface

Sandblasted acid-etched surface, with an extensive multi-stage cleaning process, utilizes ultra-pure water (UPW) which removes undesired residues, providing a clean surface and maintaining an intact oxide layer.^{7,8,9}


Ø 4.2 mm Dynamic Conical MC Implant Shown

*Disclaimer: Some products are not available for sale in all markets including the USA.

References

1. Al-Nsour MM, Hsun-Liang C, Wang HL. Effect of the Platform-Switching Technique on Preservation of Peri-implant Marginal Bone: A Systematic Review. The International Journal of Oral & Maxillofacial Implants Volume 27, Number 1, 2012.
2. Cappiello M, Luongo R, Di Iorio D, Bugea C, Cocchetto R, Celletti R. Evaluation of peri-implant bone loss around platform-switched implants. Int J Periodontics Restorative Dent. 2008 Aug;28(4):347-55.
3. Prosper L, Redaelli S, Pasi M, Zarone F, Radaelli G, Gherlone EF. A randomized prospective multicenter trial evaluating the platform-switching technique for the prevention of post-restorative crestal bone loss. Int J Oral Maxillofac Implants. 2009 Mar-Apr;24(2):299-308.
4. M. Schmitt C, Nogueira Filho G, C. Tenenbaum H, Yuan Lai J, Brito C, Döring H, Nonhoff J. Performance of conical abutment (morse Taper) connection implants: A Systematic review. J Biomed Mater Res Part A: 102A: 552–574, 2014.
5. Sang-Woon L, Min-Sang C, Ji-Hye L, Lee-Ra C, Chan-Jin P. Joint stability of internal conical connection abutments with or without hexagon indexes: an in vitro study. Journal of Dental Rehabilitation and Applied Science, 10.14368/jdras.2020.36.2.95, 36, 2, 95-103, 2020.
6. Tsuruta K, Ayukawa Y, Matsuzaki T, Kihara M, Koyano K. The influence of implant-abutment connection on the screw loosening and microleakage. Int. J. Implant Dent. 2018, 4, 1–6.
7. Klein M, Tarnow D, Lehrfeld L. Marginal Bone Changes on Ultraclean, Micro-Threaded Platform-Switched Implants Following Restoration: 1- to 4-Year Data. Compendium of Continuin Education in Dentistry. Article Vol. 41 No. 4. 2020.
8. Singh Dhaliwal J, Rani Nakka David S, Ramizah Zulhilmil N, Kaur Sodhi Dhaliwal S, Knights J, Ferreira de Albuquerque Junior R. Contamination of titanium dental implants: a narrative review. SN Applied Sciences volume 2, Article number: 1011. 2020.
9. Duddeck D, Maghaireh H, Faber FJ and Neugebauer J. SEM surface analyses of 120 sterile-packed implants. EDI Journal. 2014; 64-75.

Conical Platform

 Conical Connection	Ø 3.25 mm	Ø 3.75 mm	Ø 4.2 mm	Ø 5.0 mm
	8.0 mm	-	24-70017MC	24-70006MC
10.0 mm	24-70018MC	24-70001MC	24-70007MC	24-70013MC
11.5 mm	24-70019MC	24-70002MC	24-70008MC	24-70014MC
13.0 mm	24-70020MC	24-70003MC	24-70009MC	24-70015MC
16.0 mm	24-70021MC	24-70004MC	24-70010MC	24-70016MC