



THE
**CANADIAN
SHIELD**

Procedural Mask with Comfort Earloops

Level 1

TCS-PMC-0001

LEVEL 1



Made-in-Canada:

All of our Medical

Procedural Masks are

100% designed and manufactured
in Ontario, Canada by our
Canadian workforce.



Full Coverage:

Our pleated design

offers full coverage of


the nose and mouth and added
nose stiffeners provide a tight
seal to minimize gaps.

The Canadian Shield's ASTM
Level 1 Procedural Masks offer
triple-layer protection from fluid
exposure.

Level 1 Procedural Masks feature
3 layers of 100% polypropylene:

- Moisture repellent, Blue
spunbond outer layer
- Moisture absorbent, White
spunbond inner layer
- Meltblown filtration layer

- Adjustable nosebridge
ensures a tight seal to
minimize gaps
- Earloops are made of 70%
polyester and 30% spandex
- ASTM F2100 Compliant

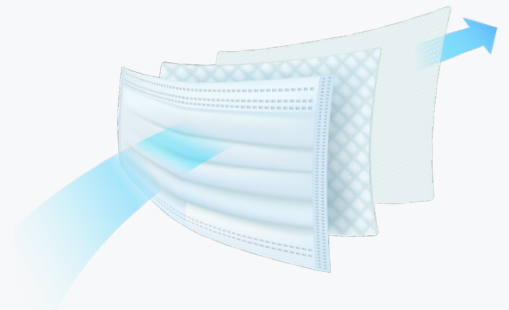

CHARACTERISTIC	DETAILS
Bacterial filtration efficiency, %	≥95
Differential pressure, mm H ₂ O/cm ²	<6.0
Sub-micron particulate filtration efficiency at 0.1 micron, %	≥95
Splash Resistance/Synthetic Blood Resistance, mmHg	80
Microbial cleanliness	<30
Flame spread	Class 1



The Canadian Shield Contact Info

Info@canadianshield.ca | 1-844-465-7684

canadianshieldppe.ca



Procedural Masks with Comfort Earloops - Level 1 - V02.22.21