



All-day ECG Monitors

Hearts Deserve
the Best Holters



CardiUP!3 and 12 ECG Holter family

- stay connected for up to 7 days for better cardiovascular analysis
- CardiUP!3: 3 recording channels
- CardiUP!12: 3/12 recording channels
- large LCD for real-time ECG display
- menu operation buttons for easy setup
- both manual & PC-based programming
- new, sensitive PM detection
- mini USB connection
- event button
- free, all-inclusive software

Warranty information

- 3-year warranty period for the recorder unit (additional one-year warranty period optional)
- 1-year warranty period for non-disposable accessories



Default kit contents


- 1 × Holter recorder
- 1 × easy to clean pouch
- 1 × ECG patient cable*
- 1 × package of disposable ECG electrodes
- 1 × high capacity AA battery
- 1 × mini USB cable
- 1 × data carrier with the latest software and user documentation
- 1 × quality certificate



Packed in black protecting case. The User Guide is also included with each kit.

*(CardiUPI12 standard kit contains a 12-channel cable by default, but 3-channel cable can also be supplied upon request; CardiUPI3 is compatible with the 3-channel cable only, therefore the standard kit comes with a 3-channel cable)

Easy to clean pouch

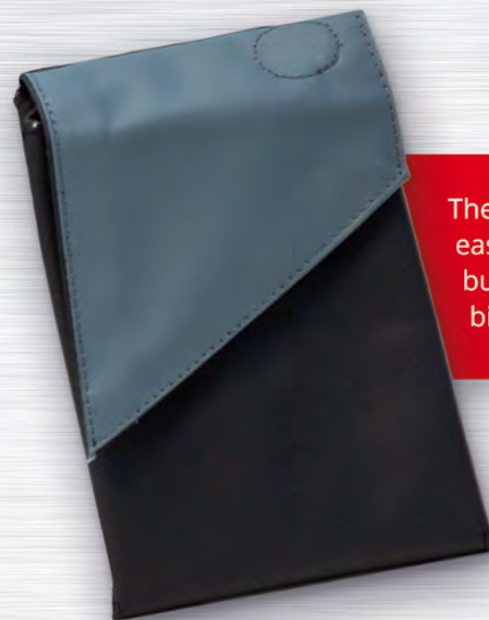


The suitable pouch is not only easy to clean and comfortable to wear, but its fabric also meets the related biocompatibility requirements.

CardioMera ECG Holter

- 1-, 2-, 3- or 5-channel recording
- up to 96 hours of continuous ECG
- individual, easily replaceable patient cables
- small size & light weight
- event button
- free, all-inclusive software

Warranty information

- 3-year warranty period for the recorder unit (additional one-year warranty period optional)
- 1-year warranty period for non-disposable accessories

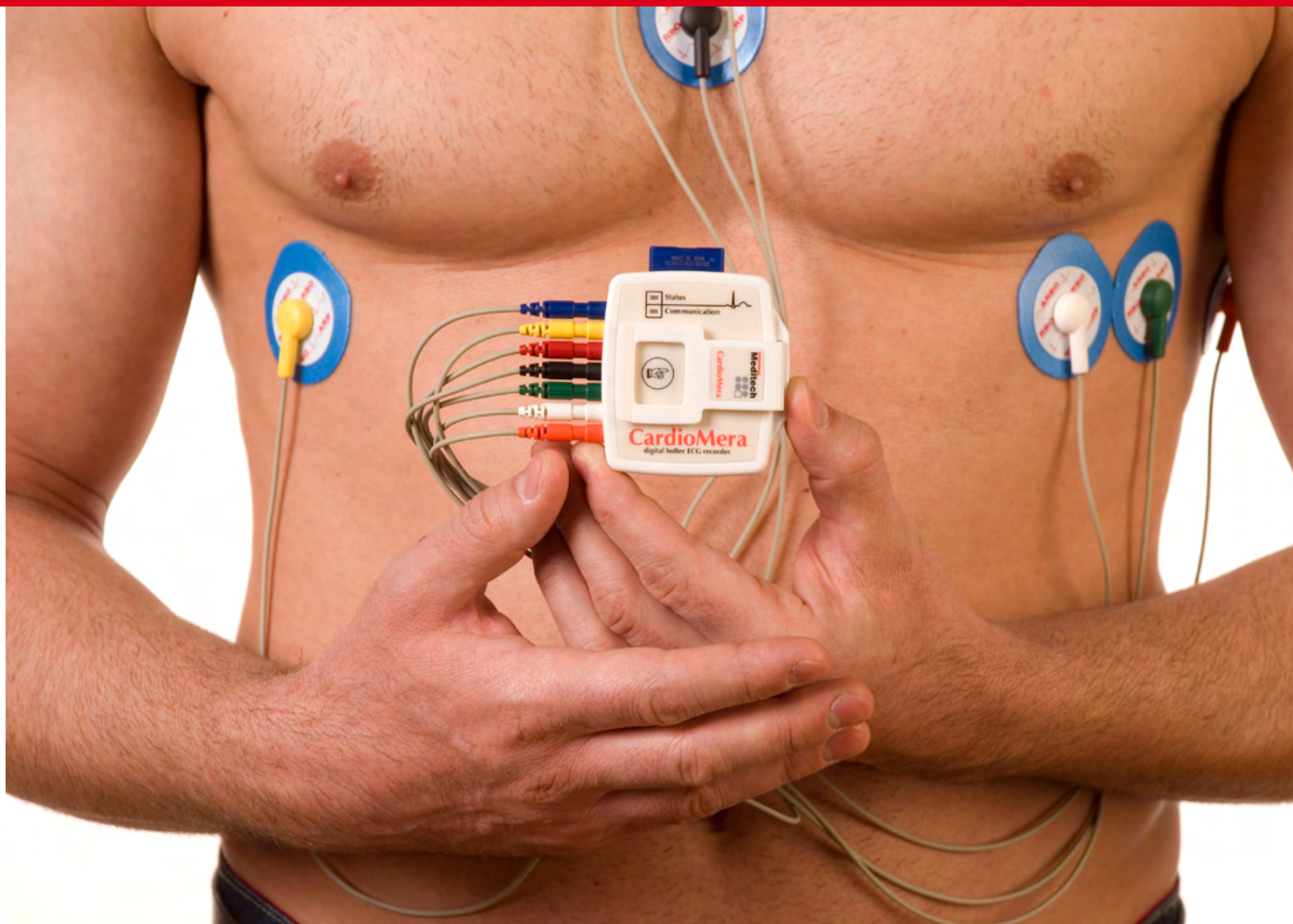

Default kit contents

- 1 × Holter recorder
- 1 × set of 7 lead ECG patient cable
- 1 × package of disposable ECG electrodes
- 1 × memory card
- 1 × high capacity AAA battery
- 1 × USB optical interface
- 1 × data carrier with the latest software and user documentation
- 1 × quality certificate



Packed in black protecting case. The User Guide is also included with each kit.

Technical specifications

Specifications	CardiUP!12	CardiUP!3	CardioMera
dimensions	67×93×21 mm		66×59×17 mm
weight (excl. batteries)	92g		50g
memory	internal		external (a card reader unit is needed for retrieving data)
data transfer	USB 2.0 (mini USB)		USB optical cable
power supply	1 AA battery		1 AAA battery
LCD display	48×28 mm		-
number of keys	4 (event, up, down, select)		1 (event)
number of channels	3/12 independent channels (7 or 10 leads)	3 independent channels (7 leads)	1, 2, 3, 5 channels (7 leads)
operation	up to 168 hours of continuous ECG		up to 96 hours of continuous ECG
sampling, A/D resolution	300, 600 or 1200 Hz, 12 bit		600 or 1200 Hz, 12 bit
storage quality	300, 600 or 1200 Hz		150, 300 or 600 Hz
conformity certificates	CE		

Meditech does not warrant the accuracy, completeness or fitness for a particular purpose of this brochure, and it cannot be considered an official proposal. For more information on the intended use of the products please check the relevant User Guide.

The soul in the body

CardioVisions free software for ECG analysis

Specifications	CardiUP!12	CardiUP!3	CardioMera
accuracy		99,90%	
beat-by-beat ECG annotation		✓	
event distribution		✓	
ST analysis		✓	
QT analysis		✓	
PQ analysis		✓	
PM analysis		✓	
time domain HRV analysis		✓	
frequency domain HRV analysis		✓	
arrhythmia analysis		✓	
graphical and tabular display of data		✓	
customizable screens		✓	
customizable & printable report		✓	
personal & network use		✓	





Meditech Kft. • H-1184 Budapest, Mikszath K. u. 24.
phone: +36 1 280 8232, +36 1 280 8233
email: meditech@meditech.eu
www.meditech.eu

MARAA0018