


Steam Sterilization Chemical Indicator Label

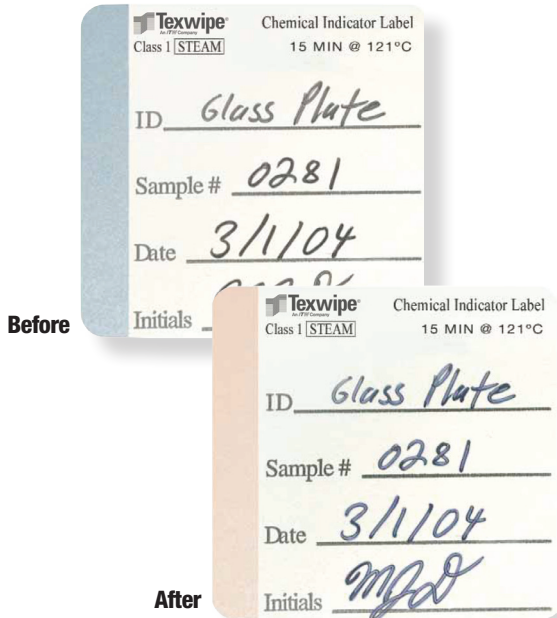

Texwipe's Steam Sterilization Chemical Indicator Label is designed and engineered for easy, accurate record keeping of materials subject to the steam sterilization process.

Autoclave

Texwipe's unique label material prevents removal of written data by the steam sterilization process. Our adhesive allows the label to be removed after the autoclave process and reapplied to a tracking form for permanent recordkeeping.

Meets AAMI ST60 and ISO11140-1

The Texwipe Steam Sterilization Chemical Indicator Label meets the performance specifications outlined in AAMI ST60 and ISO11140-1 "Sterilization of Health Care Products – Chemical Indicators – Part 1: General Requirements," Section 6 for Class 1 Process indicators.

Features

- Designed and formatted for data recording and sample tracking
- Indicator ink signals autoclave completion
- Removable from glass, metal and sterilization bag surfaces
- No residue

Benefits

- Supports accurate data recording and sample tracking
- Easy-to-use format for recording information
- Recorded data remains intact through steam sterilization process
- Signals completion of autoclave process (ink color change – blue to pink)
- Eliminates post steam sterilization cleaning (residue-free design)

Applications

- Metal trays
- Glassware
- Instruments and tools
- Autoclave bags

Material Information

- Substrate: Teslin®
- Adhesive: Autoclavable; pressure-sensitive adhesive
- Release liner: Polyester
- Category ink: Water-based
- Indicator ink: Water-based steam sterilization ink

For more information, contact:

- Your ITW Texwipe Regional Sales Manager
- Customer Service at 800.839.9472

Product

Number	Description	Packaging
TX522A	Autoclave Label 2"x2" Pre-printed, Perforated	500 labels/roll



Americas
1210 South Park Drive
Kernersville, NC 27284
Tel 336.996.7046
www.texwipe.com
info@texwipe.com

Europe/Middle East
Saffierlaan 5
2132 VZ Hoofddorp
The Netherlands
Tel +31 88 1307 410

Asia/Pacific
3A International Business Park
Icon @IBP #06-16
Singapore 609935
Tel +65 9115 8496

Steam Sterilization Chemical Indicator Label

Property	Value
Color Change	Blue to Pink
Critical Parameters	Saturated Steam; 121°, 15 minutes
Classification	STEAM
Storage Conditions (Before Use)	Store in cool, dry place
Product Code	TX522A
Storage Conditions (After Use)	Store at room temperature