


# BIODEGRADABLE NITRILE EXAM GLOVES

CHEMO DRUG TESTED | TEXTURED FINGERTIPS | POWDER-FREE | PURPLE

REF PM6-202X

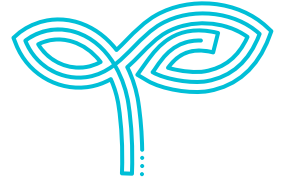

Sustain™

## ABOUT THE GLOVE

PRIMED Sustain™ nitrile gloves are specially formulated to attract microbes found in landfills to enhance the enzymatic biodegradation of the material. Sustain™ achieves 90% biodegradability in 490 days versus 0.8% in 626 days for standard nitrile formulations. Enjoy the same reliable performance and protection you expect from a nitrile glove while supporting environmental sustainability.

## PROTECTION STANDARDS

**Exceeds the current ASTM D6319:** Standard Specification for Nitrile Examination Gloves for Medical Application.

**Tested against chemotherapy drugs to the current ASTM D6978:** Standard Practice for Assessment of Resistance of Medical Gloves to Permeation by Chemotherapy Drugs.  
(List of tested chemotherapy drugs available upon request.)

**Tested to ASTM F1671 for Viral Penetration.**

## FEATURES & BENEFITS

### LIGHTWEIGHT NITRILE TECHNOLOGY

Combines strength with tactile sensitivity, allowing intricate and precise movement.

### FLEXIBLE AND COMFORTABLE FIT


Moves with the wearer to facilitate natural dexterity.

### FULLY MICRO-TEXTURED WITH ENHANCED FINGERTIP TEXTURING

Provides reliable grip when handling small instruments and objects.

### BIODEGRADABLE NITRILE FORMULATION

Specially formulated to attract microbes found in landfills to enhance the enzymatic biodegradation of the material.

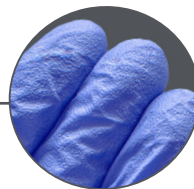

## AVAILABLE SIZES



## QUANTITIES

GLOVES/BOX	BOXES/CASE	GLOVES/CASE
200	10	2000



PRIMED Medical Products Inc.  
200, 2003-91 Street SW, Edmonton, AB, Canada, T6X 0W8  
TF +1.877.877.4633 | F +1.780.497.7670 | W PRIMED.com | E info@PRIMED.com

NA-REVA1