



QUICKIE®

QM-710



EXPERIENCE
it to **BELIEVE** it.

The QUICKIE® QM-7 Series power wheelchairs combine function and luxury for a unique driving experience. Make your day-to-day a thrill with the QM-7 Series' SpiderTrac™ Suspension, reliable 4-pole motors, PG Drives electronics, and our most inspired portfolio of power seating options ever.



QM-710

DESIGNED FOR
PERFORMANCE



POWER BASE PERFORMANCE

- Intuitive mid-wheel drive handling
- Linix 10kph, 4-pole motor packages
- SpiderTrac™ Suspension for superior traction and handling across multiple terrains
 - Gas dampening system for smooth, stable transitions
 - 3" curb climbing capability
 - 22° transition angle helps prevent high centering on curb cuts



Switch-It Dual Pro™
Proportional Head Array



PM120 Controller

ELECTRONICS OPTIONS FOR INTUITIVE DRIVING

R-net Expandable and Non-Expandable Controls

- Standard PM120 Controller offers 30% more power and improved drive algorithm
- Full programmable to ensure safe and comfortable mobility over even rough terrain.
- Compatible with iOS and Android technology - making sure the digital world is at your finger tips
- Expandable controls feature assignable buttons for control of common functions and 8 programmable profiles for transitions between environments

Alternative Drive Controls

- Alternative drive controls available on-chair and will ship complete with QM-7 Series wheelchairs
- 28 Switch-It and ASL® proportional and non-proportional (switch) input control options available
- NEW Switch-It Dual Pro™ proportional head array offers fine-tuned adjustments and customization to provide an easy-to-learn, precise, and intuitive driving experience

NEW
OPTIONS FOR
ENHANCED USER
EXPERIENCE




Power Centermount Footrest

- Length of articulation increased to 7.5"
- 3" to 19" of independent lower leg length adjustment
- Footplate extends to the floor for easier transfers



Push Button

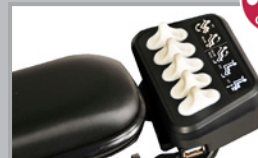

Toggle

Ctrl+5 Buttons and Toggles

- Convenient access to 5 seating functions
- USB port charges mobile phones and other devices
- May be mounted in various locations



ASAP II SEAT FRAME

- In-field width, depth, armrest, legrest, and back height adjustability. No parts required.
- Wide range of seat sizes
- Manual reclining backrest for increased function and sitting tolerance
- Optional fold down backrest for easy transport
- 16.7" Seat Height

RECLINE SEAT FRAME

- Ultra low seat-to-floor heights – as low as 16" when combined with tilt and elevate!
- Integrated JAY® J3™ back for seamless functionality with power recline and shear reduction system
- Low profile electronics enclosure shell for a clean look and compatibility with accessories

POWER SEATING OPTIONS

- Power Recline with Power Shear Reduction
- 50° of Power Tilt to optimize positioning
- 12" Power Seat Elevate for increased access to cabinets and eye-to-eye interaction with peers
- Power Elevating Legrests with Articulation
 - Swing-Away
 - NEW centermount with 7.5" of articulation



PG Drives Colour Joystick

- 3.5" colour display
- Paddle switches conveniently control on/off function and chair speed
- Speaker for increased horn volume
- Programmable shortcuts for custom control of other functions



Retractable Joystick Mount

- Adjusts from -135° to 135°
- Provides comfortable position for users whose hands turn inward/outward
- Retracts inward to clear tight spaces
- Retracts outward to access desks/tables

Comfort Reclining Armpads

- Softer armpad
- Cleaner style
- 2" width option
- 4" width option



QM-710 SPECIFICATIONS

Specifications	QM-710 BC	QM-710 SC	QM-710 MPC	QM-715 HD BC	QM-715 HD SC	QM-715 HD MPC
Drive Wheel Position	Mid	Mid	Mid	Mid	Mid	Mid
Speed	10 kph	10 kph	10 kph	10 kph	10 kph	10 kph
Overall Weight	141kgs	159kgs	170.5kgs	141kgs	159kgs	170.5kgs
Overall Width	25"	25"	25"	25"	25"	25"
Overall Length	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests	36" w/o footrests
User Weight Capacity	136kgs	136kgs	136kgs	181kgs	181kgs	181kgs Recline: 158kgs
Seat Widths	12" to 22"	12" to 22"	14" to 22"	16" to 24"	16" to 24"	16" to 22"
Seat Depths	12" to 22"	12" to 22"	16" to 22"	16" to 22"	16" to 22"	16" to 22"
Seat To Floor Height	16.3" to 21"	16.1" to 20.6"	16" to 19.3"	16.3" to 20.3"	16.1" to 20.8"	16.6" to 19.6"
Seat Frame Options	ASAP II	ASAP II or Recline	Recline	ASAP II	ASAP II or Recline	Recline
Power Seating Options	Power ELRs, Elevate	Tilt or Recline, Power ELRs, Elevate	Tilt or Recline, Power ELRs, Elevate	Power ELRs, Elevate	Tilt or Recline	Tilt & Recline

CONVENIENCE OPTIONS

- RAM® Phone and Tablet Mounts allow mobile phones and/or tablets to be connected and disconnected for easy viewing.
- RAM® Mounts Self-Leveling Cup Holder keeps drinks level during travel or when stopped on inclines.
- Active Controls ReCharge™ with Socket eliminates the challenges of connecting and disconnecting the battery charger by providing an easy-to-reach adapter.
- Bluetooth for controlling Apple and PC-Android devices.



Phone Mount

Tablet Mount

ReCharge™

JAY® Seating Selection



Union™



J2™/ J2 DC



Fusion™



J3™



J3™ PA



J3™ PL



J3™ PD, PDC, PDL