




X-Ray Grid




Radiography System

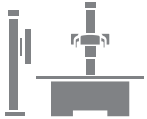

Software



Computed Radiography

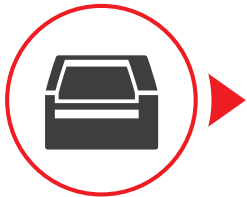

Processor



X-Ray Accessories



RF Treatment



Diagnostics



JP-33

Automatic X-ray Film Processor



L 635mm (25.0 inches)
H 385mm (15.0 inches)
W 525mm (20.6 inches)

Compact size



**Easy to use,
maintain and service**



**Retractable
feed-tray for
small spaces**

Features

New retractable feed tray provides space savings for small dark rooms

Exceptionally quiet operation

Microprocessor controlled

High-quality developing

The control panel digitally displays

Water level sensor prevents overflowing of tank

Stand-by mode for energy saving

One year warranty on all parts

Specifications

Processor type	Continuous roller transport system	
Processing time	Standard: 90 seconds Option: 90 to 230 seconds by variable speed controller	
Film size	From 4" x 4" (10cm x 10cm) up to 14" x 17" (35cm x 43cm)	
Processing capacity	100 sheets/hr (mixed sizes) - at 95 seconds standard processing time	
Developer temperature	Automatic microprocessor controlled - Temperature control range: 28-38°C (82-100°F) - Temperature can be easily changed and reset	
Dryer temperature	Automatic microprocessor controlled - Temperature control range: 45-55°C (113-131°F) - Temperature can be easily changed and reset	
Circulation mode	Anti-crystallization by automatic circulation pump	
Replenishment mode	Automatic microprocessor and manual-controlled	
Dryer type	Blower-heater	
Film feed sensor	Micro-switch sensor	
Tank capacities	Developer: 5.0 liters (1.32 gallons) Fixed: 5.0 liters (1.32 gallons) Washer: 6.0 liters (1.58 gallons)	
Water consumption	1.5 liters (0.4 gallon) / min	
Power consumption	Operation mode: 1.8 kw Stand-by mode: 0.8 kw	
Electrical connections	AC 220-240 Volt, 50Hz/60Hz, Single phase AC 110-120 Volt, 50Hz/60Hz, Single phase	
Dimensions	Closed feed tray	L 63.5cm (25.0 inches) H 38.5cm (15.0 inches) W 52.5cm (20.6 inches)
	Opened feed tray	L 86.5cm (34.9 inches) H 38.5cm (15.0 inches) W 52.5cm (20.6 inches)
Weight	Net: 35 kg (77.1 lbs) Gross: 77 kg (169.8 lbs)	
Standard accessories	Replenishment tanks: 25 liters (7.0 gallons) / 2 pieces Hoses / Stand / Manual / Spare parts kit	

Optional Features

- Variable speed controller
- Closed loop system



The Control Panel

Displays working status, and provides easy control of all functions and settings.



- Film size
- from 4" x 4" up to 14" x 17"
- Water consumption
- 1.5 liters (0.4 gallon) / min
- Processing capacity
- 100 sheets/hr

★ Specifications are subject to change without notice

JPI Healthcare Co., Ltd.

608-ho, 28, Digital-ro 33-gil, Guro-gu, 08377 Seoul, Korea

82-2-2108-2580 82-2-2108-1180

www.jpi-korea.com

