

Technical data sheet


Level 3 Surgical Mask



Face masks can be an effective control measure as they prevent droplets of bodily fluids from escaping via the nose and mouth. They also protect the wearer against splashes and sprays (sneeze and coughs) from other people.

AusMed Health Level 3 Surgical Face Mask

A single-use face mask designed to meet Australian Standard AS/NZS 4381:2015. The lightweight mask is designed to promote comfort and excellent breathability, with adjustable nose wire and soft ear loops to hold the mask in place.

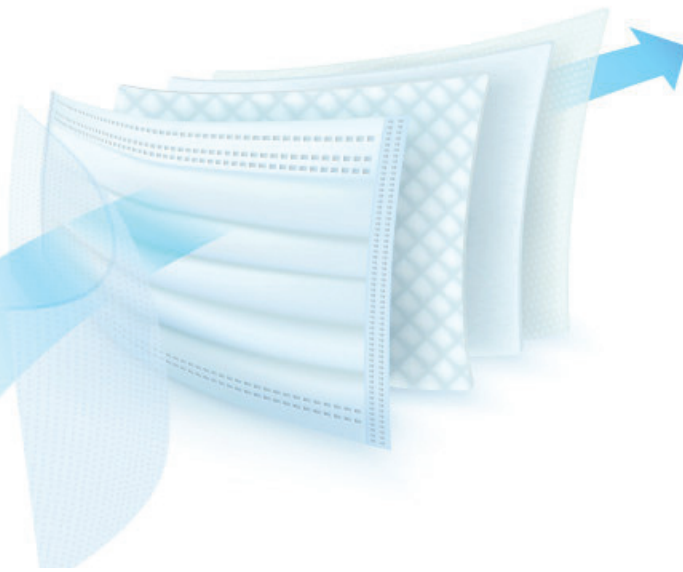

The AusMed Health Surgical Face Mask (non-sterilised) is constructed of four layers of material to ensure compliance to the highest level of protection.

The composition is two outer layers of spunbond non-woven material, a middle layer of meltblown and one inner layer of hypoallergenic spunbond nonwoven material for the wearers comfort.

Key Features

- Exceeds Australian Standard AS4381-2015
- Bacterial filtration efficiency (BFE) $\geq 98\%$
- Differential pressure, $\text{mmH}_2\text{O}/\text{cm}^2 < 5.0$
- Resistance to penetration by synthetic blood to 160mmHg
- 17.5cm x 9.5cm in size

Listed on the Australian Register of Therapeutic Goods (ARTG ID 341864).

Suitable for: general community, public transport, personal use, aged care, GPs, food service, industrial use, general hospital use.

Sold in cartons of 500 surgical face masks. Each carton contains ten boxes of 50 masks.

The AusMed Health Level 3 Surgical Face Masks have been manufactured in Melbourne, Australia.

Caution: The mask is designed for personnel who are properly trained in mask use and limitations. Misuse may result in sickness or death. DO NOT use in any manner not indicated in the Instructions for Use. Always read the label. Discard the mask if it becomes damp, soiled, damaged, or difficult to breathe through. Do not use if the packaging is damaged.



ausmedhealth.com.au
AusMedHealth PTY LTD | ABN: 71643 436 177
6 Precision Lane, Notting Hill, Victoria 3168

