


Dräger DCD 5000 Oral fluid collection system

Simple, fast and secure detection of intoxicants and drugs:
The Dräger DCD 5000 does not only enable an easy collection of saliva (oral fluid), you can also store the samples and send them to the laboratory for analysis – without the sample being falsified along the way.

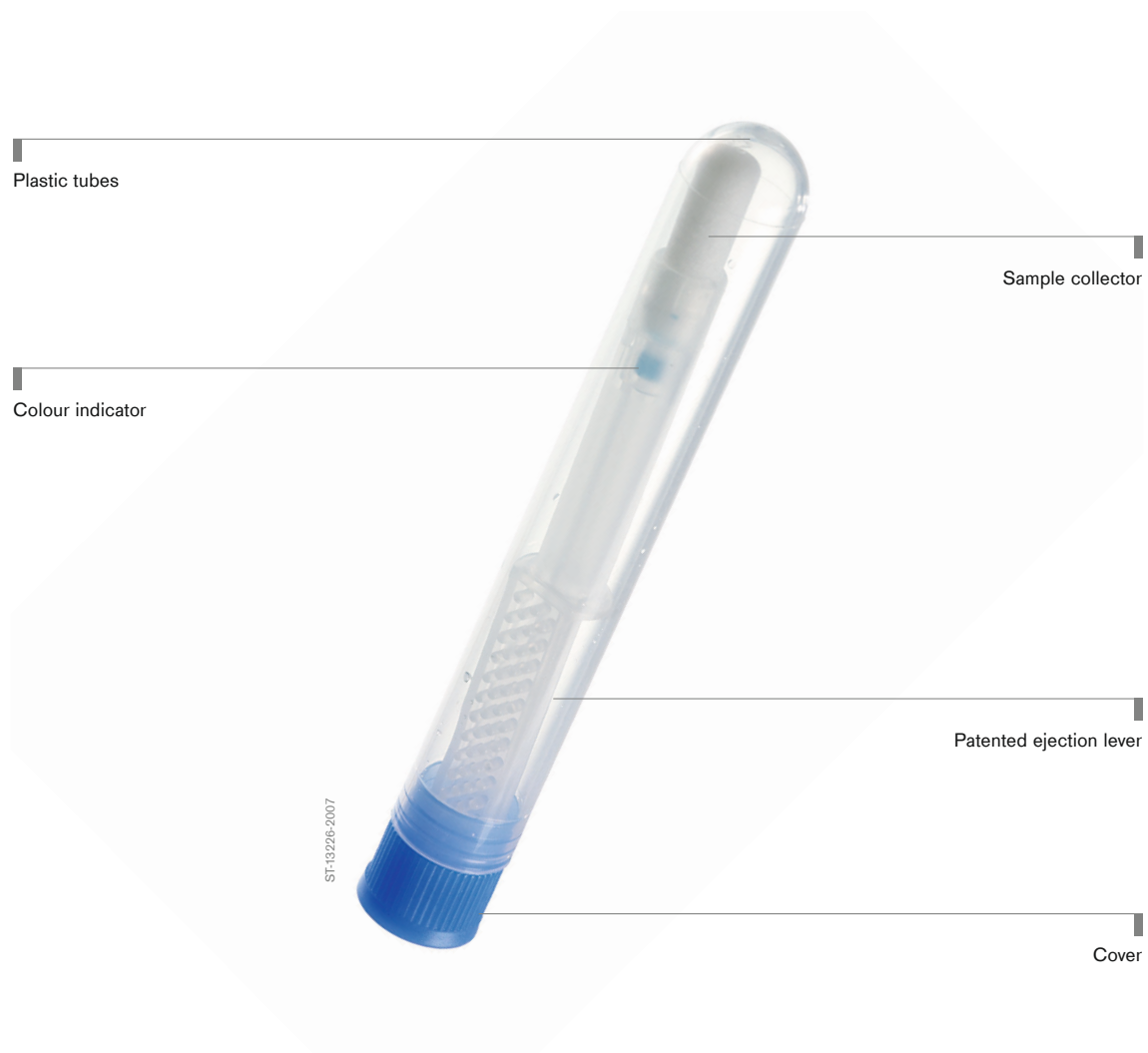

Benefits

Ultra filtrate of blood

Saliva is an established ultra filtrate in clinical chemical diagnostics. Saliva samples are approved as sample material for the detection and diagnostic analysis of substance abuse.

Check without interruptions

The continuous check of sampling is facilitated with the method based on saliva, compared to urine sampling. Thereby, manipulation by the test person is eliminated.

Easy, clean and safe

The results of a pretest are confirmed reliably via a laboratory analysis. Storage, transport and recovery of the sample must not impair this result. The DCD 5000 sampling kit allows for this security.

Double sample

Combining two DCD 5000 samplers, it is very easy to collect two separate samples simultaneously. Afterwards, the samples can be separated and serviced separately, e.g. for transport and storage.

Details



Remove the sampler from the tube and give it to the test subject.



To collect the saliva sample, swab from one inner cheek to the other.



When the colour indicator turns blue, the sampling can be ended.



The sample is secured inside the tube for further use.



Patented ejection mechanism to release the mouthpiece from the handle.



Twin sampler to collect two saliva samples at once.

Accessories



Dräger SCSK

In the SCSK confirmation set, the DCD 5000 is enhanced with a stabilising component (isopropanol) to ensure high-accuracy detection of drugs in the sample, even after lengthy periods between sampling and laboratory testing.

Accessories



D-15278-2017

Dräger SCK

The SCK confirmation set contains two samplers and stabilising components (isopropanol) as well as documentation forms and shipping materials, ensuring a chain of custody through clear documentation and secure shipment of samples to the lab.

Related Products



D-42702-2015

Dräger DrugCheck® 3000

Use the Dräger DrugCheck® 3000 to find out within minutes if a person recently consumed certain drugs. The compact and quick saliva (oral fluid) based drug test yields reliable results affordably and easily. The non-invasive drug test kit does not require electricity and can therefore be used anywhere.



D-54720-2012

Dräger DrugTest® 5000

No pipetting, no drips, no timing: With the Dräger DrugTest® 5000, a drug test is carried out simply and quickly. The collected saliva sample can be analyzed immediately – for accurate results on the spot.

Technical Data

Dimensions of the individual sampler (H x W x D)	120 mm x 15 mm x 15 mm
Working range	Storing an unused Dräger DCD 5000: 4 to 30 °C
Sample volume	approx. 0.4 ml of saliva (oral fluid)
Licenses/Standards	Approved in the EU as a medical product in accordance with Regulation 98/79/EC for in-vitro diagnostics (IVD), as well as in partner states of the EU and in countries in which a local permit exists. Outside the EU for non-medical purposes (non-IVD) or forensic applications (USA only)

Ordering Information

	Packaging unit	Part no.
Dräger DCD 5000 Sampler Kit	20 saliva samplers, packaged individually in plastic tubes with seal, 1 x instructions for use	83 19 910
SCSK – Saliva Confirmation Single Kit	Set consists of 10 sealed plastic bags, each containing 1 x DCD 5000 saliva sampler, 1 x tube of isopropanol, 1 x instructions for use included per set	83 26 300
SCK – Saliva Confirmation Kit	Set consists of 2 x DCD 5000 saliva samplers, 2 x tubes of isopropanol with packaging material: 1 x sealable transport bag for diagnostic material, 2 x safety seals, 1 x carbon copy form for laboratory use and instructions for use	83 23 880

Notes

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS
Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
www.draeger.com

REGION DACH
Dräger Safety AG & Co. KGaA
Revalstraße 1
23560 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE
Dräger Safety AG & Co. KGaA
Revalstraße 1
23560 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA
Dräger Safety AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA PACIFIC
Dräger Singapore Pte. Ltd.
25 International Business Park
#04-20/21 German Centre
Singapore 609916
Tel +65 6308 9400
Fax +65 6308 9401
asia.pacific@draeger.com

**REGION CENTRAL
AND SOUTH AMERICA**
Dräger Panama S. de R.L.
59 East Street, Nuevo Paitilla
House 30, San Francisco Town
Panama City, Panama
Tel +507 377 9100
Fax +507 377 9130
servicioalcliente@draeger.com

Locate your Regional Sales
Representative at:
www.draeger.com/contact

