

# CATALOG



**„Standing still means stepping backwards.  
Only by actively embracing change we will be successful.“**









- Model preparation
- Die preparation



## Die preparation



### SAKURA margin liner

Pencil for margin lines.

- reliable outline on all sorts of gypsum
- refills are free from graphite and wax, ideal for:
  - perfect casting results
  - long-term fixing with YETI CLEAR SPACER
- blue colour builds a high contrast to red cervical wax



SAKURA margin liner

1 pcs.

Art.-No. 986-0000

Refills, blue colour

free from graphite and wax

12 pcs., 0,5 mm

Art.-No. 986-1000

## SPACER

Our SPACER are available in 18 ml glass bottles with brush.

The bottles captivate by their elegant design, they adhere easily at Vario magnetic holder.

You constantly have tidy workplace. The SPACER can be applied more even due to improved brush quality.

18 ml




Magnetically  
adhering




„excellence line“  
Vario magnetic holder  
Art.-No. 312-9000


For more information order  
separate „excellence line“ catalogue  
or download it under [www.yeti-dental.com](http://www.yeti-dental.com)

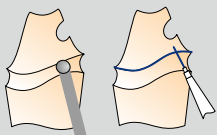


1. Prepare the margin line
2. SAKURA margin liner

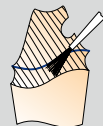

3. CLEAR SPACER

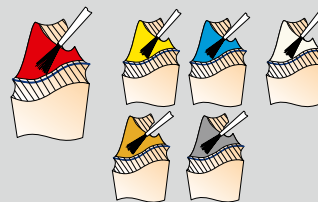

4. COLOR or DIE SPACER.  
Apply beginning 1 mm above  
the margin line

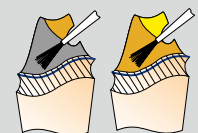

5. DIE SPACER.  
Apply 2 additional  
layers if necessary  
(use DIE SPACER  
only)



3 my

### CLEAR SPACER transparent 3 my

- Apply onto the complete plaster die including the area below the margin line.
- Fixation of margin line
  - A layer of 3 my on dry plaster
  - seals the die surface
  - base for a chemical bonding with further spacers



CLEAR SPACER transparent  
1 x 18 ml  
6 x 18 ml

Art.-No. 501-0018  
Art.-No. 501-6018

## Spacer

**Apply all Spacer beginning 1 mm above the margin line!**

7 my

### COLOR SPACER red 7 my

- Apply onto the die, beginning 1 mm above the margin line. A layer of 7 my on CLEAR SPACER
- the colour changes from orange to red during the drying process
  - best results for Inlay and Onlay Technique thanks to its fluency
  - transparent - any marking on the die remains visible
  - does not contain colour pigments - no shaking of the bottles necessary before use



COLOR SPACER red  
1 x 18 ml  
6 x 18 ml

Art.-No. 500-0018  
Art.-No. 500-6018

## DIE SPACER covering Spacer

Apply as first layer over CLEAR SPACER or if necessary as second additional layer (use only with DIE SPACER or CLEAR SPACER)

7 my

### DIE SPACER yellow 7 my

DIE SPACER yellow  
1 x 18 ml  
6 x 18 ml

Art.-No. 502-0018  
Art.-No. 502-6018



10 my

### DIE SPACER blue 10 my

DIE SPACER blue  
1 x 18 ml  
6 x 18 ml

Art.-No. 505-0018  
Art.-No. 505-6018





**10 my**

**DIE SPACER dentin 10 my**

- no colour irritations with full ceramic crowns  
(the colour shines through)




DIE SPACER dentin  
1 x 18 ml  
6 x 18 ml

Art.-No. 506-0018  
Art.-No. 506-6018

**13my**

**DIE SPACER gold 13 my**

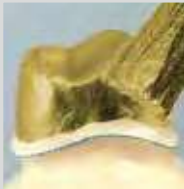

DIE SPACER gold  
1 x 18 ml  
6 x 18 ml

Art.-No. 504-0018  
Art.-No. 504-6018

**13my**

**DIE SPACER silver 13 my**



DIE SPACER silver  
1 x 18 ml  
6 x 18 ml

Art.-No. 503-0018  
Art.-No. 503-6018

**THINNER  
(not for light curing liquids and STEAM OFF)**

Universal thinner for all our spacers.



THINNER  
1 x 18 ml  
6 x 18 ml

Art.-No. 510-0018  
Art.-No. 510-6018

## LIGHT SPACER



3 my

Advantages of the light curing die spacer LIGHT SPACER:

- homogeneous curing
- creates a very hard, scratch-resistant and durable surface
- extremely smooth surface - no problem when lifting the coping
- excellent adhesion on plaster base
- unlimited application time cures only when put into the light curing unit

Curing in any available standard unit  
(We recommend PRECI NT SHUTTLE)

Curing time: approx. 2 - 5 min



10 my

Delivery: in black UVA protection bottles of 10 ml



**Die LIGHT SPACER can only be used with each other!**  
**Not to be mixed with other Spacers!**

LIGHT SPACER transparent	3 my	1 x 10 ml	Art.-No. 520-1000
LIGHT SPACER red	10 my	1 x 10 ml	Art.-No. 520-2000
LIGHT SPACER blue	10 my	1 x 10 ml	Art.-No. 520-3000



0 my


## DIE HARDENER 0 my

Completely absorbed by gypsum. Does not build up any dispersion layers.


- cure, seal and smooth the gypsum surface
- protects the plaster
- best results for the production of wax-copings
- stain and water repellent

**Not suitable as base because of insufficient bonding between spacer and hardener.**  
**Spacers can be affected by the solvents of the hardener.**

1. Prepare the margin line
2. SAKURA margin liner

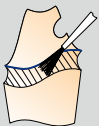

OR

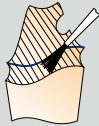

### 3. DIE HARDENER 0 my

- for curing the plaster surface
- only to be used below the margin line, if afterwards spacer (above the margin line) is used

### 3. DIE HARDENER 0 my

- curing of the plaster surface
- complete die, if no spacer is further used

DIE HARDENER			
Bottle with brush	1 x 18 ml		Art.-No. 540-0018
Bottle with brush	6 x 18 ml		Art.-No. 540-6018
Refill bottle	1 x 250 ml		Art.-No. 540-0250





**CUTEX BOND**  
block out wax



### YETI LUBE

YETI LUBE matches perfectly with all YETI waxes and spacers. It does not contain any alcohol and solvents, it is water soluble and does not react with the investment material. Even small quantities prove the great separating effect. YETI LUBE isolates wax against plaster, metal and wax. The separating effect remains for more than 24 hours.

**Do not use with STEAM OFF Spacer!**

#### YETI LUBE

Bottle with brush

1 x 18 ml

Art.-No. 550-0018

Bottle with brush

6 x 18 ml

Art.-No. 550-6018

Refill bottle

1 x 250 ml

Art.-No. 550-0250

### ISO QUICK

A handy separating pen filled with YETI LUBE

- thin and exact application without wasting material
- easy to refill
- economical and user friendly

1 Separating pen

Art.-No. 550-0000

1 x 250 ml Refill bottle

Art.-No. 550-0250

- 2 different tips:

1 x broad and flat

1 x pointed and fine

apply STEAM OFF Spacer:



remove applied material with a steamer:









## Waxes for crowns and bridges

quality made in germany







### DUO DIP Dipping Wax

DUO DIP is a ready for use wax mixture for everyday work. Smooth and rounded form of grains allows fluent melting of wax granulate. Dipping pot can be filled quickly, clean and fast with exact wax quantity.

- ashfree - burns without residue

Weight: 80 g



### MASTER DIP

MASTER DIP offers a new dimension in strength and stability. Rapid cooling allows for quick build-up of a stable and accurate coping. As dimensionally stable as a resin, Master Dip provides for a highly accurate casting.

A granulate wax available in natural transparent and orange. With an elastic memory, this wax is hard without being brittle.

- ashfree - burns without residue

Weight: 100 g



DUO DIP  
Dipping Wax orange (ashfree)

Art.-No. 741-0000

MASTER DIP  
elastic hard  
elastic hard

transparent  
orange

Art.-No. 741-0100  
Art.-No. 741-0110

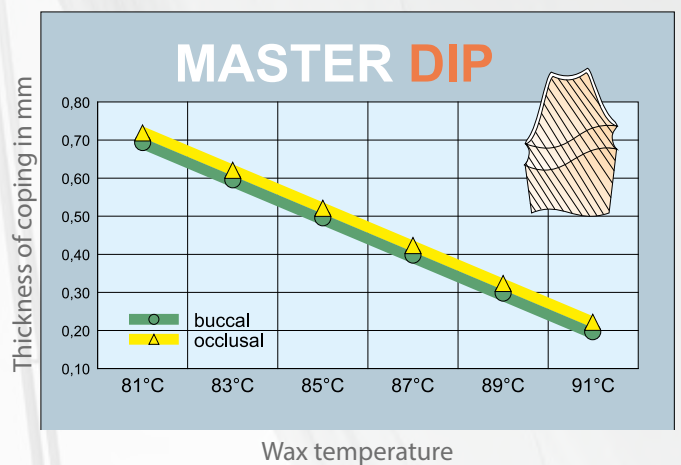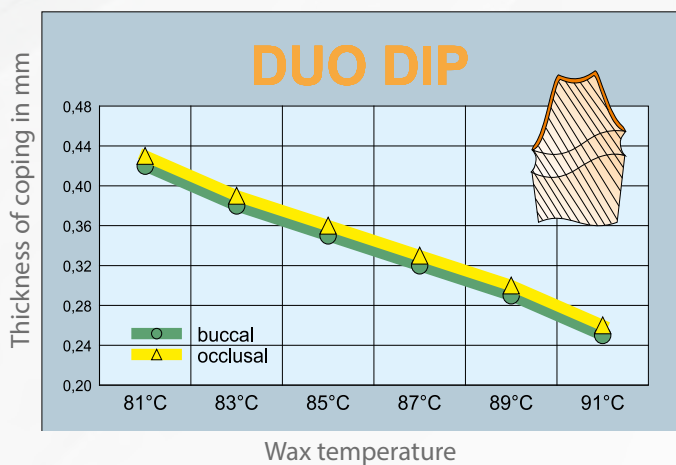



The modern modelling and cervical wax




- Double Line with integrated cervical wax inlay
- quick and time-saving
- ergonomical working height



- two foam rings for convenient removing



- firm base for stability
- transparent dust protection cover

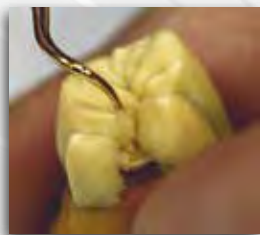

**Excellent wax quality - best modelling characteristics:**

- enables modelling of smallest fissures - natural looking occlusion
- allows drawing of the wax
- best stability in melted condition - easier building up of tooth shapes



The wax is:

- hard but not brittle - high internal strength
- good to scrape - clean filing
- tight fitting with few contractions





# WAX GIANT



WAX GIANT - **single line** -, 75 g  
 Cervical wax red (ash-free)  
 Modelling wax sand (opaque)  
 Modelling wax grey (opaque)  
 Modelling wax blue (ash-free)

Art.-No. 731-5000  
 Art.-No. 732-5001  
 Art.-No. 734-5001  
 Art.-No. 736-5000



WAX GIANT - **double line** -, 75g  
*(with integrated cervical wax inlay)*  
 Modelling wax sand (opaque)  
 Modelling wax grey (opaque)  
 Modelling wax blue (ash-free)

Art.-No. 733-5001  
 Art.-No. 735-5001  
 Art.-No. 737-5000



Jochen Peters shows in his special wax-up technique how it is possible to create high-class functional occlusal surfaces in a minimum of time following his efficient, systematic and easy-to-do working steps. Since in the concept by Jochen Peters build up and scrape down of the wax are combined, a wax which meets both criteria in a perfect way was required. Thanks to a good cooperation both the sculptural waxes "Ivory", "Sahara" and "Sunlight" as well as the suitable cervical wax "Wood" were created.



ivory



sahara



sunlight



wood

### WAX GIANT by Jochen Peters

#### The **sculpturing waxes "Ivory", "Sahara" and "Sunlight"**

are perfectly suitable for crowns and bridges, inlays, onlays, casting-, press- and milling techniques. Besides these, three universal waxes can be perfectly combined with brown cervical /underlining wax "Wood".

Modelling waxes captivate by the following characteristics:

- First-rate sculpturing properties. Waxes can be easily spread and have and still have first-class stability
- perfect wax hardness allows exact modelling as well as precise creation of cusps, cusp slopes and fine dental ridges
- after smoothing all wax surfaces are perfectly sealed without smearing
- excellent edge stability
- easy and clean to scrape down, no chipping
- exact fitting and a low contraction
- burns without residue

The colour for "Ivory", "Sahara" and "Sunlight" waxes are chosen to please the eyes and so that even smallest details can easily be seen

The **brown cervical / underlining wax "Wood"** can be used perfectly as cervical / underlining wax for crowns and bridges as well as sculpturing wax for inlays and onlays.

- fine wax consistency
- in a liquid state fills perfectly every small cavities
- soft but not sticky or smearing
- stable in form and utterly exact in fitting
- tension free
- burns without residue



WAX GIANT - **single line** -, 75 g

Cervical waxwood (ash-free)

Modelling wax sunlight (ash-free)

Modelling wax sahara (ash-free)

Modelling wax ivory (ash-free)

Art.-No. 738-5000

Art.-No. 738-5010

Art.-No. 738-5020

Art.-No. 738-5030



## Sculpturing wax



### THOWAX Sculpturing wax

THOWAX is a sculpturing wax suitable for the crown, bridge and inlay technique. It is a compound of natural waxes free from any synthetic materials.

Thanks to its great sculpturing features THOWAX has gained an excellent reputation worldwide.

- excellent sculpturing characteristics due to ingredients with extremely short molecules
- allows a drawing of the wax and the sculpturing of even smallest details
- quick cooling down facilitates rapid working
- excellent sculpturing characteristics without filling up the fissures
- opaque even when molten in order to have a better control during waxing-up
- high molecular density - easily carvable, leaves smooth surfaces
- matches perfectly with our special and dipping waxes
- free from synthetic materials
- available in 70 g tins, some colours also as wax chip



THOWAX Sculpturing wax	beige with brown	opaque	70 g	Art.-No. 700-0000
THOWAX Sculpturing wax	red	lightly opaque	70 g	Art.-No. 708-0000
THOWAX Sculpturing wax	blue	lightly opaque	70 g	Art.-No. 709-0000
THOWAX Sculpturing wax	beige	opaque	70 g	Art.-No. 710-0000
THOWAX Sculpturing wax	brown	lightly opaque	70 g	Art.-No. 711-0000
THOWAX Sculpturing wax	violet	lightly opaque	70 g	Art.-No. 712-0000
THOWAX Sculpturing wax	green	lightly opaque	70 g	Art.-No. 713-0000
THOWAX Sculpturing wax	grey	opaque	70 g	Art.-No. 714-0000

## THOWAX Special Waxes

Specialized waxes were developed to match perfectly in different application fields. They guarantee excellent quality in special working cases as well.

All waxes are in 70 g tins, some colours are available as wax chip.

**CERVICAL WAX, red**

- tension free, smooth-flowing wax for precise margins
- adapts perfectly to the margin line

CERVICAL WAX, red

Art.-No. 720-0000

**UNDERLINING WAX, bordeaux**

As first layer of waxing-ups for inlays, crowns, secondary parts and for sculpturing within the cervical area.

- high density enables the carving of sharp crown edges
- stable in form and exact fitting
- opaque for a better contrast to the die

UNDERLINING WAX, bordeaux

Art.-No. 722-0000

**MILLING WAX, grey opaque or green transparent**

- can be milled at a speed of 5000 rpm
- clean filing
- smooth surfaces
- good sculpturing features despite its hardness
- any unevenness of milled area is clearly visible

MILLING WAX, grey opaque

Art.-No. 727-0000

MILLING WAX, green transparent

Art.-No. 725-0000

**STICKY WAX, beige**

*not suitable for the casting, does not burn without residue*

- helpful for any repair of broken dentures or to solder metal parts
- very hard when cooled down and therefore immediately indicates any undesired slipping of the parts to be connected

STICKY WAX, beige

Art.-No. 724-0000

**STICK-ON-WAX, neon-pink**

To place casting channels, sprues and pontics

- perfect sticking, easy correction of position
- stays smooth and sticky
- burns without residue

STICK-ON-WAX, neon-pink

Art.-No. 726-0000





**VKS Waxes for all-ceramic systems**  
(for example for K2 press)

VKS waxes for all-ceramic systems are of THOWAX quality. They are filtered several times and are free from any opaque particles, which could possibly react with the investment material or casted objects. The excellent sculpturing features remain unchanged.



VKS Sculpturing Wax beige	70 g	Art.-No. 710-2000
VKS Sculpturing Wax grey	70 g	Art.-No. 714-2000
VKS Underlining Wax bordeaux	70 g	Art.-No. 722-2000



**SAMPLE SET - Liquids for die preparation and THOWAX**

A selection to introduce you to our products. Included are sculpturing and special waxes, spacer, separating liquid and high spot indicator. Detailed instructions are included and guarantee a correct, systematic and successful application of our products.

SAMPLE SET	Art.-No. 100-0000M
------------	--------------------

*Content: THOWAX Sculpturing wax grey, THOWAX Sculpturing wax beige, THOWAX Underlining wax bordeaux, THOWAX Cervical wax red, CLEAR SPACER transparent, YETI LUBE Separating liquid, BLUE MARKER High spot indicator*

## SPECIAL WAXES

### Cervical waxes

	Art.-No.	Colour	ash free/ opaque	Suitable for VKS	Solidification point	Page
IQ ASH-FREE Cervical wax	720-5000	black	ash free	X	60 °C = 140,0 °F	31
IQ COMPACT ASH-FREE Cervical wax	720-5100	black	ash free	X	60 °C = 140,0 °F	31
WAX GIANT Cervical wax single line	731-5000	red	ash free	X	63 °C = 145,4 °F	24
WAX GIANT Cervical wax single line	738-5000	wood	ash free	X	63 °C = 145,4 °F	25
THOWAX Cervical wax	720-0000	red	ash free	X	65 °C = 149,0 °F	27
Chips Cervical wax chip	720-1000	red	ash free	X	65 °C = 149,0 °F	29

### Underlining waxes

IQ ASH-FREE Underlining wax	722-5000	caramel	ash free	X	58 °C = 136,4 °F	31
IQ COMPACT ASH-FREE Underlining wax	722-5100	caramel	ash free	X	58 °C = 136,4 °F	31
THOWAX Underlining wax	722-0000	bordeaux	opaque		59 °C = 138,2 °F	27
Chips Underlining wax	722-1000	bordeaux	opaque		59 °C = 138,2 °F	29
THOWAX VKS Underlining wax	722-2000	bordeaux	ash free	X	58 °C = 136,4 °F	29

### Milling waxes

THOWAX Milling wax	725-0000	green	ash free	X	60 °C = 140,0 °F	27
THOWAX Milling wax	727-0000	green	opaque		61 °C = 141,8 °F	27
Chips Milling wax	725-1000	green	ash free	X	60 °C = 140,0 °F	29

### Stick on wax

THOWAX STICK-ON-WACHS	726-0000	neon-pink	ash free	X	68 °C = 154,4 °F	27
-----------------------	----------	-----------	----------	---	------------------	----

## Sculpturing wax



## Special wax





## Accessories for wax technique



### OCCLUSAL POWDER

Occlusal powder is a fine white powder which consists of wax only.

- easy to apply with a brush
- indicates frictions or undesired contact spots when using the articulator
- serves as a separating powder when applying hot wax
- burns without residue - washing not necessary

Occlusal powder white	17 g	Art.-No. 730-0000
-----------------------	------	-------------------



### WAX BRUSH

Wax brush to smooth and to even out wax surfaces

Wax brush	1 pcs.	Art.-No. 985-0000
-----------	--------	-------------------



### WAX DEBUBBLIZER

- reduces tension on all wax and silicone surfaces
- the investment material adapts precisely to the wax work

#### Spray:

- easy to spray on
- homogeneous
- colour: blue transparent

Wax debubblizer	100 ml / Spray	Art.-No. 142-0000
Wax debubblizer	1000 ml / Refill	Art.-No. 142-1000



### FINALE Wax Smoothing Liquid

- smoother surfaces of wax modellations, scratches and uneven surfaces disappear
- slim layer - approx. 2-3 my
- good adaption to the wax surface
- short drying time - depending on the layer thickness 1 - 3 minutes
- better casting results - smooth and homogeneous surfaces
- shorter finishing time


#### Fields of application:

- Smoothing of wax surfaces - occlusal areas
- reinforcement of contact points
- casting technique
- smoothing of wax profiles and contact points
- sealing of finest fissures between plaster and wax profiles - avoids the flowing in of investment material

FINALE Wax smoothing liquid	1 x 18ml	Art.-No. 507-0018
FINALE Wax smoothing liquid	6 x 18ml	Art.-No. 507-6018

## Preformed wax pattern

IQ Preformed patterns are helpful accessories in the dental laboratory, they save time and improve the quality of your casting.

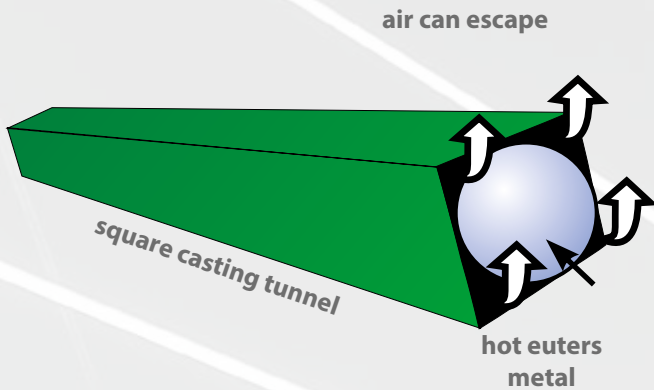

### IQ QUADRO STICKS



Why should a wax stick be square?  
After burn-out the wax stick leaves a square, hollow channel. Molten metal has round cross-section due to its high surface tension and simulates a rolling ball in motion. Displaced air escapes back through the unfilled corners when the molten metal enters the square channels, eliminating blow-back. The resulting castings are visibly smoother and imperfection free.

- 100% wax, burns without residue
- light green, hard
- Length: 110 mm
- Weight: 50g

IQ QUADRO STICK	diagonal $\varnothing$ 2,5 mm	Art.-No. 765-5000
IQ QUADRO STICK	diagonal $\varnothing$ 3,0 mm	Art.-No. 766-5000
IQ QUADRO STICK	diagonal $\varnothing$ 3,5 mm	Art.-No. 767-5000

### IQ STICKS round




More security when lifting your wax pattern.  
Stable connections for the whole casting channel system.

- 100% wax, burns without residue
- red, hard
- Length: 80 mm
- Weight: 250g



IQ STICKS round	$\varnothing$ 2,5 mm	Art.-No. 768-2500
IQ STICKS round	$\varnothing$ 3,0 mm	Art.-No. 768-3000
IQ STICKS round	$\varnothing$ 3,5 mm	Art.-No. 768-3500
IQ STICKS round	$\varnothing$ 4,0 mm	Art.-No. 768-4000
IQ STICKS round	$\varnothing$ 5,0 mm	Art.-No. 768-5000

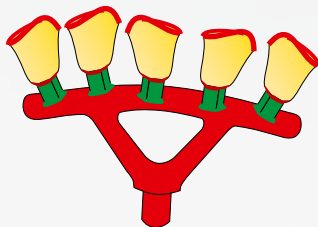

Casting object bonded to the IQ TRIANGEL with IQ QUADRO STICKS



### IQ TRIANGEL red casting bar, extra hard



- off-centered position of the casting bar due to the angle of inclination of the feeding channels
- the bent casting bar was adapted to the form of the muffle ( equal cooling of casting objects )
- saves time
- hardness: extra hard and stable
- burns without residue, pure wax, free from synthetic materials
- Package content: 50 pcs.

IQ TRIANGEL casting bar	$\varnothing$ 4,0 mm	Art.-No. 758-5000
IQ TRIANGEL casting bar	$\varnothing$ 5,0 mm	Art.-No. 759-5000





**DURON Wax wire green, hard**



Hard DURON wax quality is tension free and narrows while tearing off  
 - burns without residue  
 - time saving

DURON Wax wire green	Ø 2,0 mm	250 g	Art.-No. 746-1020
DURON Wax wire green	Ø 2,5 mm	250 g	Art.-No. 746-1025
DURON Wax wire green	Ø 3,0 mm	250 g	Art.-No. 746-1030
DURON Wax wire green	Ø 3,5 mm	250 g	Art.-No. 746-1035
DURON Wax wire green	Ø 4,0 mm	250 g	Art.-No. 746-1040
DURON Wax wire green	Ø 5,0 mm	250 g	Art.-No. 746-1050



**DETON Wax wire blue, medium**



DETON is a tension free wax wire without non-desirable rest forces  
 - middle hard wax wire with highest internal stability  
 - burns without residue

DETON Wax wire blue	Ø 2,5 mm	250 g	Art.-No. 741-2500
DETON Wax wire blue	Ø 3,0 mm	250 g	Art.-No. 742-3000
DETON Wax wire blue	Ø 3,5 mm	250 g	Art.-No. 743-3500
DETON Wax wire blue	Ø 4,0 mm	250 g	Art.-No. 744-4000
DETON Wax wire blue	Ø 5,0 mm	250 g	Art.-No. 745-5000




**CONSEQUENT Wax wire emerald green, soft**




- free of undesired tensions  
 - high inner stability  
 - good working properties also at low temperatures  
 - no splitt-off  
 - the diameter of the wire narrows while tearing off  
 - adheres to the wax model  
 - burns without residue

CONSEQUENT Wax wire	Ø 2,0 mm	250 g	Art.-No. 747-1020
CONSEQUENT Wax wire	Ø 2,5 mm	250 g	Art.-No. 747-1025
CONSEQUENT Wax wire	Ø 3,0 mm	250 g	Art.-No. 747-1030
CONSEQUENT Wax wire	Ø 3,5 mm	250 g	Art.-No. 747-1035
CONSEQUENT Wax wire	Ø 4,0 mm	250 g	Art.-No. 747-1040
CONSEQUENT Wax wire	Ø 5,0 mm	250 g	Art.-No. 747-1050

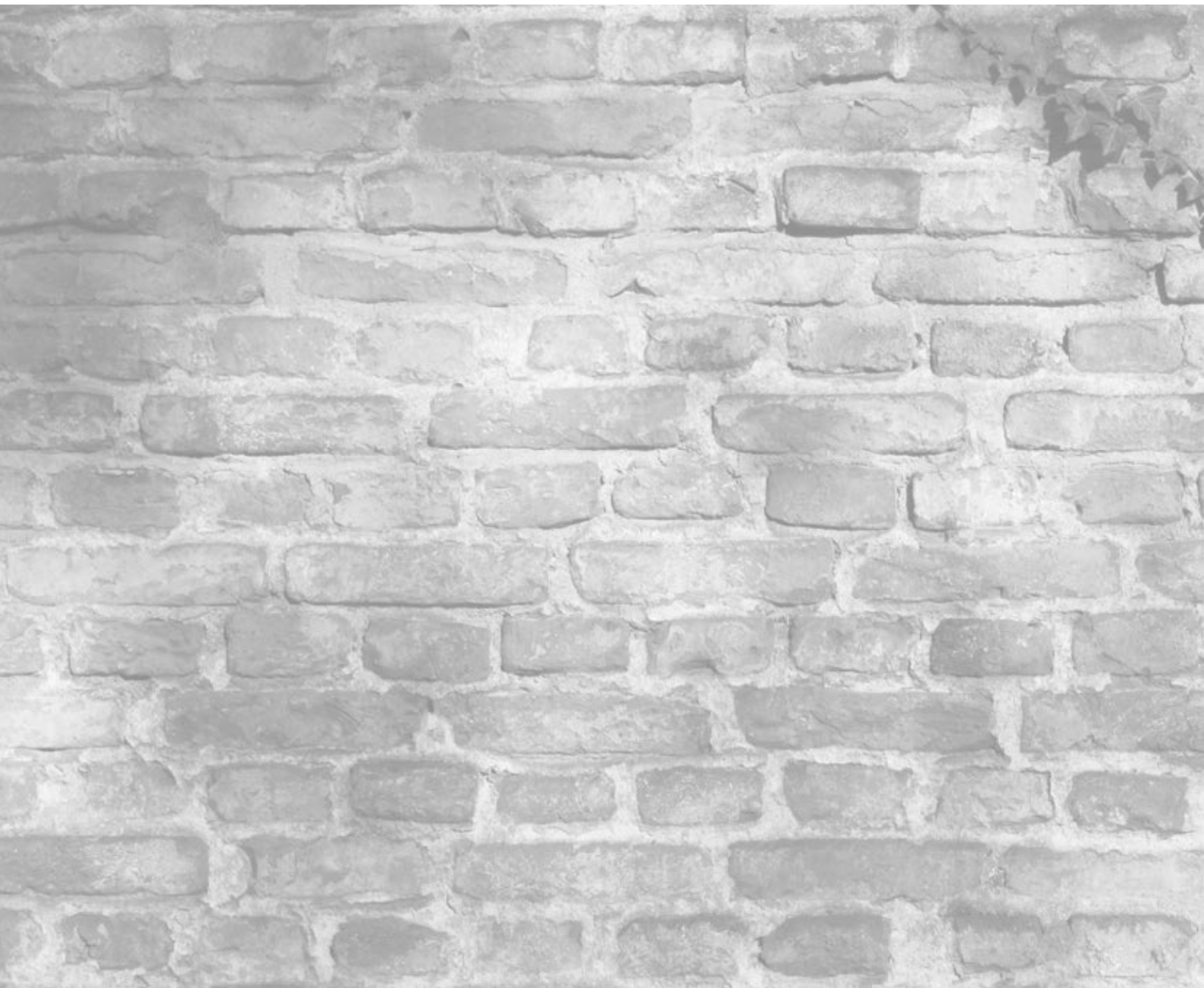




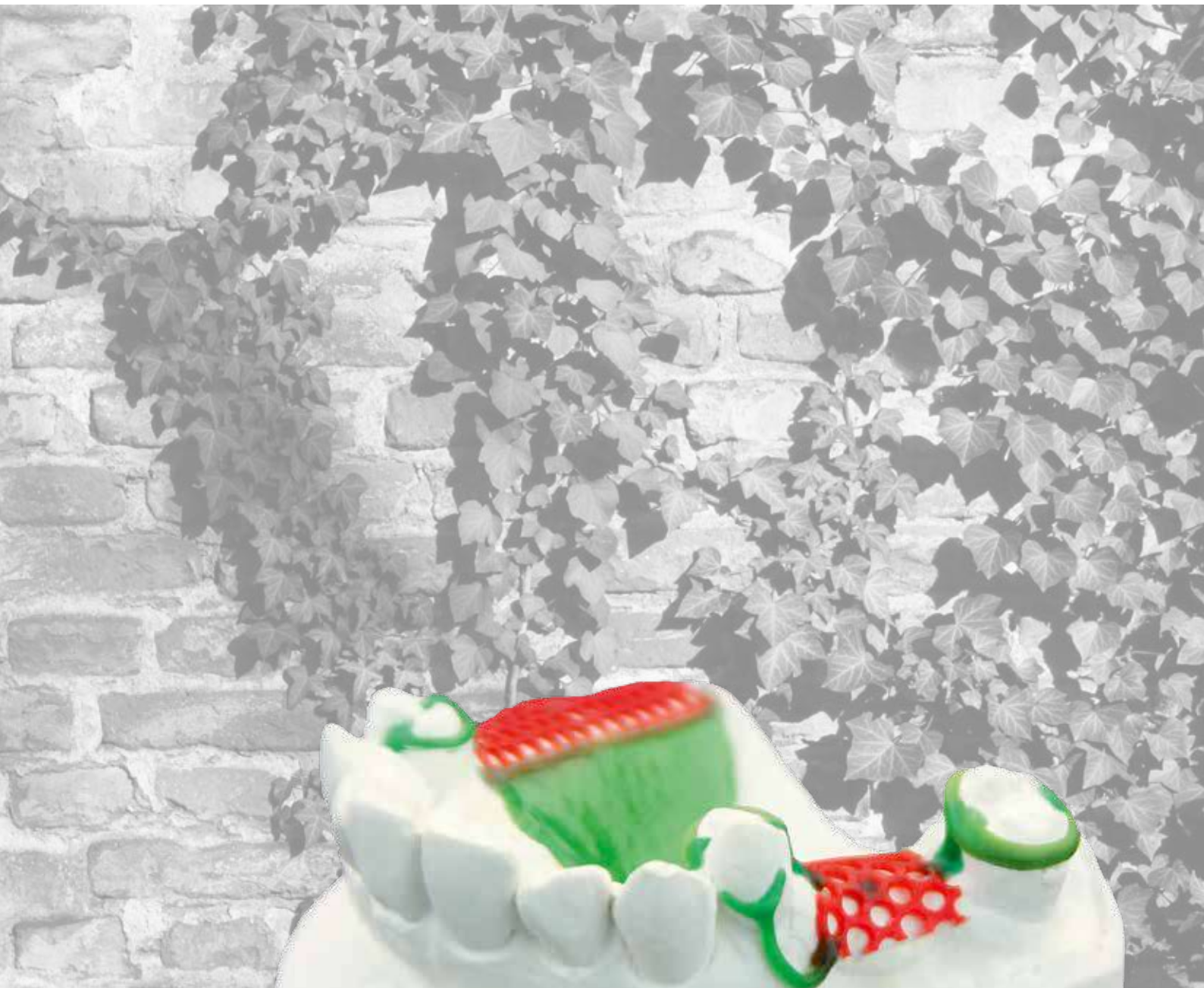

## Model casting

## Clasps profiles

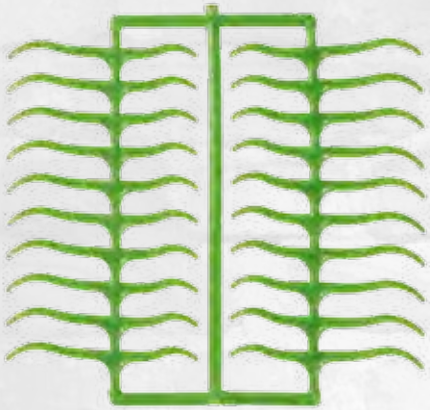


### CONSEQUENT multifunctional clasp

Wax from the CONSEQUENT product line can be easily processed even at low temperatures.

The wax is tension free with a simultaneous high inner solidness.

Stretching of the clasps allows making of more gracile clasp forms while keeping the same wax stability. This reduces the time of the refinishing processes and provides a remarkable saving of working time. In addition, by changing the form the clasps could be prolonged so that they can be used for both, molars and premolars.

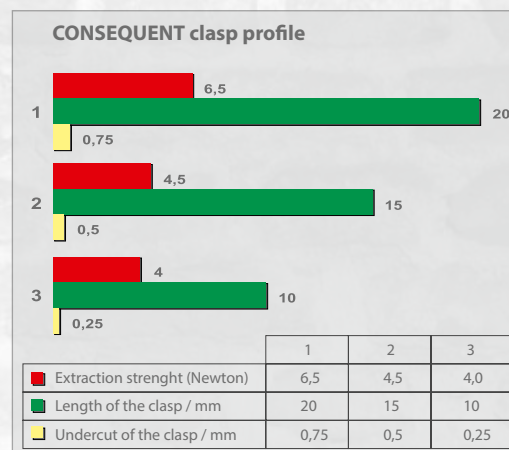

### CONSEQUENT multifunctional clasp

- the shifted ridge enables food residues to glide easily
- fine and filigree clasp form
- perfect anatomically shaped clasp
- high stability
- elongated clasp arms for molars of larger dimensions
- definitive saving of working time

Package content: 10 sheets = 200 clasps

Colour: emerald green


CONSEQUENT multifunctional clasp

Art.-No. 109-0000




## WAX CLASPS PROFILES / WAX TEMPLATES

- tension-free
- stable and self-adhesive
- package contents: 10 sheets = 200 pcs.
- colour: green

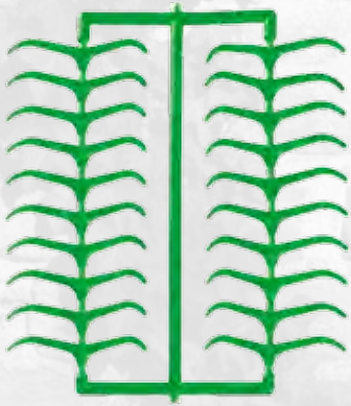


PK

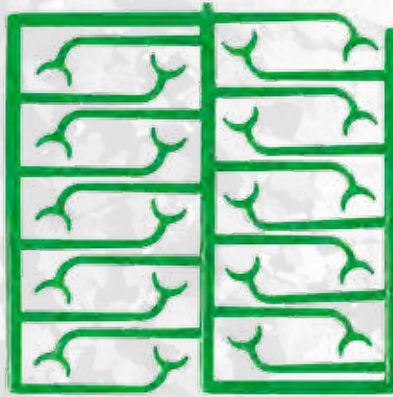


BH I




RK geb.

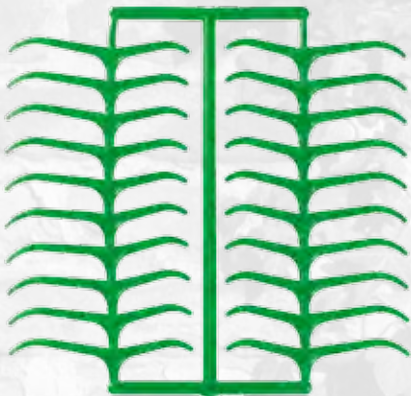

MK

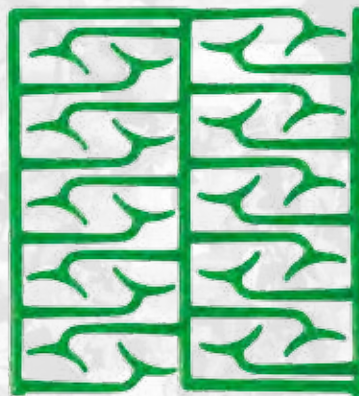

BH II

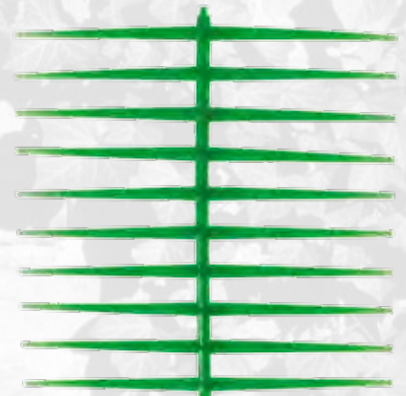

RK

PK	-	Premolar clasps	Art.-No. 110-0000
MK	-	Molar clasps	Art.-No. 111-0000
BH I	-	Bonyhard clasps small	Art.-No. 112-0000
BH II	-	Bonyhard clasps large	Art.-No. 113-0000
RK geb.	-	Ring clasps bent	Art.-No. 114-0000
RK	-	Ring clasps straight	Art.-No. 115-0000


## Wax profiles in bars

### CONSEQUENT Wax profiles emerald green

- **stability under pressure** = high inner stability, while adapting the profile, the diameter is not changing
- **tension-free** = profiles keep the desired position and do not bend up
- **easy to adapt** = easy to put on the model
- = good working properties also at low temperatures

### CONSEQUENT lower jaw bars

- CONSEQUENT wax quality
- when bending no folds and fissures
- tension-free, easy to adapt
- excellently suitable for titanium
- not necessary to reinforce
- ideal for tension-free connections
- length: 150 mm
- colour: emerald green
- hardness: medium

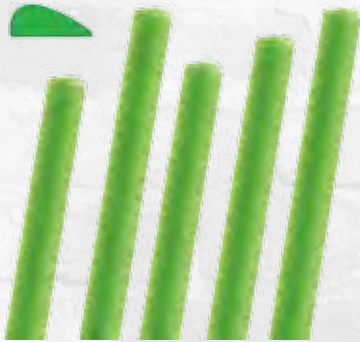

#### CONSEQUENT lower jaw bars

4,0 x 1,7 mm	65 g	Art.-No. 192-1700
4,0 x 2,0 mm	65 g	Art.-No. 192-2000
4,4 x 2,3 mm	65 g	Art.-No. 192-2300

### CONSEQUENT profiles in bars round

- in CONSEQUENT wax quality
- length: 150 mm
- colour: emerald green
- hardness: medium




#### CONSEQUENT lower jaw bars round

Ø 0,60 mm	12 g	Art.-No. 152-0600
Ø 0,80 mm	25 g	Art.-No. 152-0800
Ø 1,00 mm	38 g	Art.-No. 152-1000
Ø 1,20 mm	50 g	Art.-No. 152-1200
Ø 1,50 mm	50 g	Art.-No. 152-1500




### WAX PROFILES IN BARS


- smooth surface
- tension-free
- hardness: hard
- colour: blue
- length: approx. 150 mm

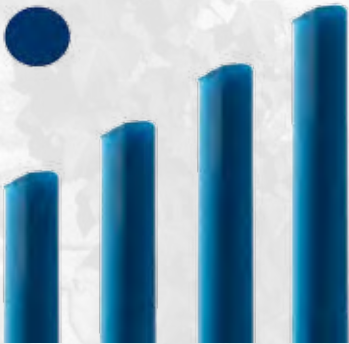

Wax wire round	Ø 0,60 mm	12 g	Art.-No. 151-0600
	Ø 0,80 mm	25 g	Art.-No. 151-0800
	Ø 1,00 mm	38 g	Art.-No. 151-1000
	Ø 1,20 mm	50 g	Art.-No. 151-1200
	Ø 1,50 mm	50 g	Art.-No. 151-1500

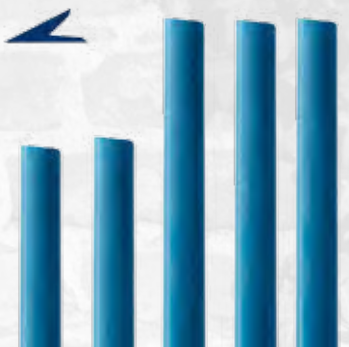

Finishing lines	2,0 x 1,0 mm	25 g	Art.-No. 171-1000
-----------------	--------------	------	-------------------

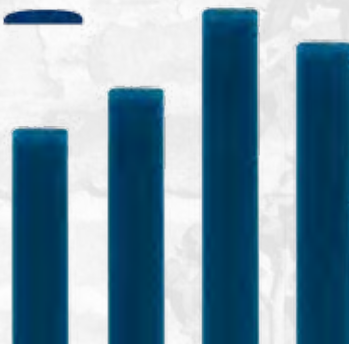

Upper jaw bar	5,0 x 1,3 mm	65 g	Art.-No. 181-1300
---------------	--------------	------	-------------------



Lower jaw bar	4,0 x 1,7 mm	65 g	Art.-No. 191-1700
	4,0 x 2,0 mm	65 g	Art.-No. 191-2000

## Wax retentions

### Wax retentions with border slat, red

- easy to adapt
- saving of time
- easy to process with rotary instruments
- available in two sizes
- 15 pcs.



#### Base plate retentions with border slat

Base plate retentions small

Art.-No. 116-0000

Base plate retentions large

Art.-No. 117-0000



#### Perforated retentions with border slat

7 x 127 mm

25 pcs.

Art.-No. 126-1000



#### Grid retentions with border slat

10 x 122 mm

25 pcs.

Art.-No. 125-1000



#### Oval retentions

130 x 8 mm

25 pcs.

Art.-No. 127-0000



#### Perforated retentions

135 x 6,5 mm


25 pcs.

Art.-No. 126-0000



### WAX RETENTIONS / MESHES red

- big assortment for individual applications
- tension-free
- easy to adapt

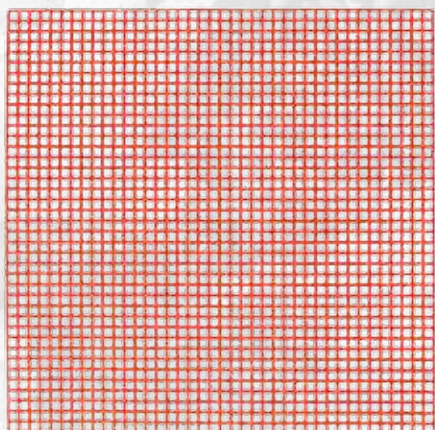


#### Wax meshes

100 x 100 mm

10 sheets

Art.-No. 120-0000

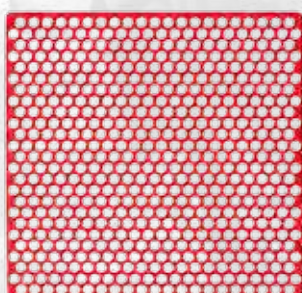


#### Grid meshes

71 x 70 mm

20 sheets

Art.-No. 121-0000

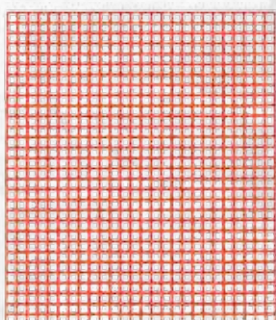


#### Wax meshes

75 x 65 mm

20 sheets

Art.-No. 122-0000

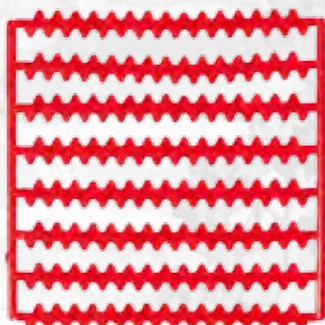


#### Bar retentions

75 x 65 mm

10 sheets

Art.-No. 123-0000

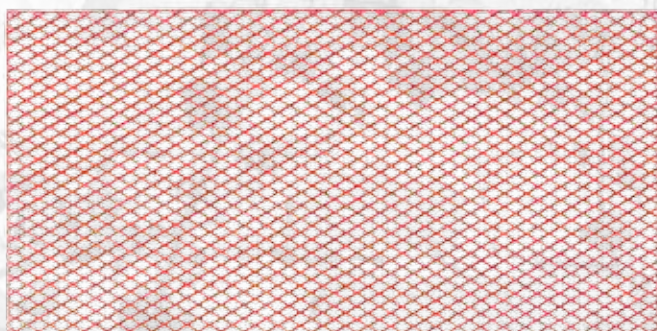

#### Wax grid diagonal

75 x 150 mm

10 sheets

Art.-No. 125-0000



## Model casting

### Casting wax green

- very malleable, transparent
- large sheet size: 17,5 x 8 cm
- Package content: 15 sheets



Casting wax green	0,30 mm	Art.-No. 130-3000
Casting wax green	0,35 mm	Art.-No. 130-3500
Casting wax green	0,40 mm	Art.-No. 130-4000
Casting wax green	0,50 mm	Art.-No. 130-5000
Casting wax green	0,60 mm	Art.-No. 130-6000
Casting wax green	0,70 mm	Art.-No. 130-7000

### Preparation wax pink

- well-proven large sheet size: 17,5 x 8 cm
- Package content: 15 sheets
- pink transparent



Preparation wax pink	0,50 mm	Art.-No. 135-5000
Preparation wax pink	0,60 mm	Art.-No. 135-6000
Preparation wax pink	0,70 mm	Art.-No. 135-7000



### Structured Casting wax

- easy to adapt
- tension-free
- adheres well to the model
- easily polishable textures
- sheet size: 17,5 x 8 cm

Package content: 15 sheets  
colour: green transparent



Sheets grained	0,40 mm	Art.-No. 133-4000
Sheets grained	0,50 mm	Art.-No. 133-5000
Sheets grained	0,60 mm	Art.-No. 133-6000



Sheets fine stippled	0,35 mm	Art.-No. 131-3500
Sheets fine stippled	0,40 mm	Art.-No. 131-4000
Sheets fine stippled	0,45 mm	Art.-No. 131-4500
Sheets fine stippled	0,50 mm	Art.-No. 131-5000
Sheets fine stippled	0,60 mm	Art.-No. 131-6000



Sheets coarse stippled	0,35 mm	Art.-No. 132-3500
Sheets coarse stippled	0,40 mm	Art.-No. 132-4000
Sheets coarse stippled	0,45 mm	Art.-No. 132-4500
Sheets coarse stippled	0,50 mm	Art.-No. 132-5000
Sheets coarse stippled	0,60 mm	Art.-No. 132-6000

## Model casting




### CONSEQUENT Block out wax

- does not remove from metal
- good visibility even when only a thin coat is applied
- good carving (with or without thermo-stick)
- does not react with monomer liquids and standard silicone
- colour: orange - is not absorbed by the plaster
- package size: 45 g tins

CONSEQUENT Block out wax

Art.-No.139-0000

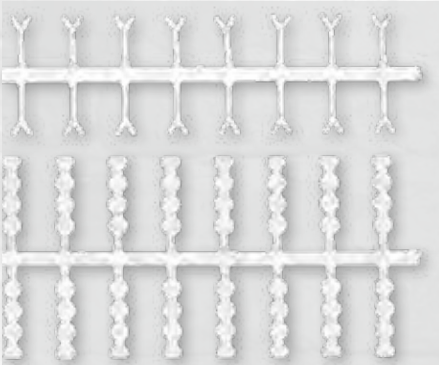

### PLASTIC RETENTIONS

- burn without residue
- different shapes - depending on available space
- made of plastic
- colour: white
- 500 pcs. per box

Plastic retentions  
Plastic retentions



Art.-No. 145-1000

Art.-No. 145-2000



### UNIVERSAL GLUE

- For fixing of retention pearls, preformed parts and plate waxes.
- burn without residue
  - colour: transparent

UNIVERSAL GLUE

18ml

Art.-No. 508-0018



### WAX RETENTION PEARLS blue

- good adhesion to the base
- burn without residue
- no swelling of the pearls

WAX RETENTION PEARLS blue

Ø 0,4mm

15g

Art.-No. 716-5010

WAX RETENTION PEARLS blue

Ø 0,2mm

15g

Art.-No. 717-5010



### CASTING FUNNEL

- made of plastic
- colour: white
- Package content: 50 pcs.

CASTING FUNNEL

Art.-No. 141-0000





### CONSEQUENT Universal wax

CONSEQUENT universal wax is mainly used within the casting technique but also suitable for the waxing-up. By using various colours good contrasts can be achieved in order to avoid too thin wax transitions.

- very good connecting wax
- melting point: 60°C
- softer wax quality
- cooling down slower
- 70 g tins

#### CONSEQUENT Universal wax

green  
red  
blue

Art.-No. 728-0010  
Art.-No. 728-0020  
Art.-No. 728-0030



### WAX DEBUBBLIZER

- reduces tension on all wax and silicone surfaces
- the investment material adapts precisely to the wax work

#### Spray:

- easy to spray on
- homogeneous

- colour: blue transparent

Wax debubbler  
Wax debubbler

100 ml / Spray  
1000 ml / Refill

Art.-No. 142-0000  
Art.-No. 142-1000



### FINALE Wax Smoothing Liquid

- smoother surfaces of wax modellations, scratches and uneven surfaces disappear
- slim layer - approx. 2-3 my
- good adaption to the wax surface
- short drying time - depending on the layer thickness 1 - 3 minutes
- better casting results - smooth and homogeneous surfaces
- shorter finishing time

#### Fields of application:

- Smoothing of wax surfaces - occlusal areas
- reinforcement of contact points
- casting technique
- smoothing of wax profiles and contact points
- sealing of finest fissures between plaster and wax profiles - avoids the flowing in of investment material

FINALE Wax smoothing liquid  
FINALE Wax smoothing liquid

1 x 18ml  
6 x 18ml

Art.-No. 507-0018  
Art.-No. 507-6018







## Duplication technique

## Dipping hardener

### YETIDUR

Bio-dipping hardener

Cold dipping hardener for the surfaces of the investment models.

- biological decomposition

Usage:

- dry the investment model for 25 - 30 min. at 180°C

- dip the warm model for 5 sec. in YETIDUR

- drain the model well and leave it to cool

- if needed put it into the oven to dry




YETIDUR Bio-dipping hardener

1000 ml

Art.-No. 943-1000



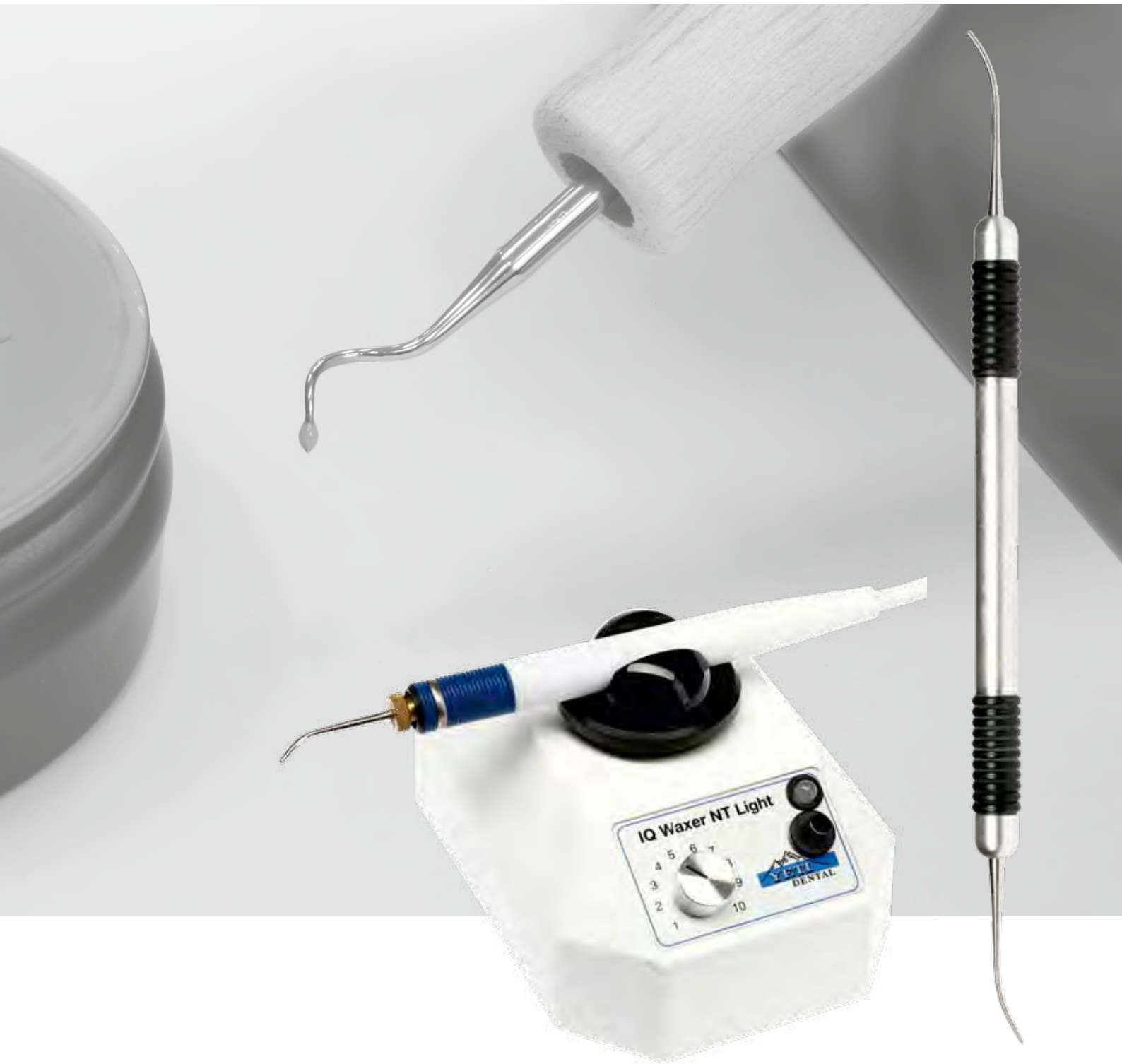

- Wax up instruments
- Units and accessories

# Wax up instruments



## ARTLINE Instruments

Modelling instruments are daily tools for dental technicians. Being aware of this importance our ARTLINE, ARTLINE plus and ARTLINE NAT instruments were created according to the latest knowledge of ergonomics, they are temperature isolated and non-slippery. The tips are made of high-quality steel, are hand-sharpened and polished. All ARTLINE instruments are " **MADE IN GERMANY** ".

ARTLINE Set in wooden case  
incl. (Artl.) Instruments

Art.-No. 850-0000

### ARTLINE

Carving sharp

(Artl.) ARTLINE

red

Art.-No. 850-1000



Carving round

Modelling tip medium

(Artl.) ARTLINE

silver

Art.-No. 850-2000




Modelling tip large

Modelling tip  
extra fine

(Artl.) ARTLINE

gold

Art.-No. 850-3000

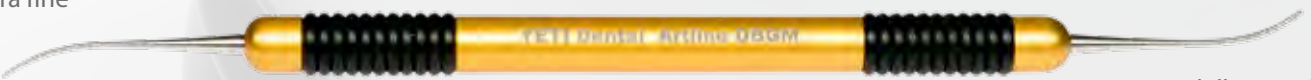

Modelling tip fine

### ARTLINE by Jochen Peters

ARTLINE Jochen Peters No.01, Lecron

Art.-No. 850-6001

my favorite



Moulding former

ARTLINE Jochen Peters No.02, PKT

Art.-No. 850-6002

Modelling tip fine



Modelling tip large



## Wax up instruments

### ARTLINE plus

Hollenback vertical

(AP) ARTLINE plus No.01 red Art.-No. 850-4001




Modelling top thin

(AP) ARTLINE plus No.02 green Art.-No. 850-4002



Universal moulder

(AP) ARTLINE plus No.03 yellow Art.-No. 850-4003

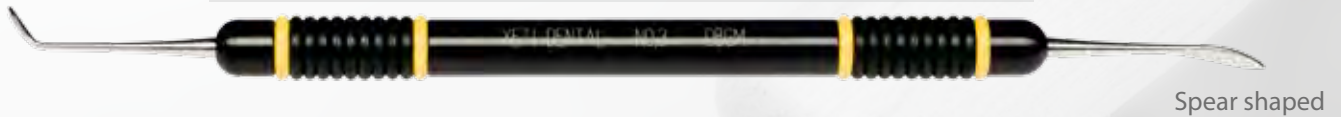

Wax knife large

(AP) ARTLINE plus No.04 blue Art.-No. 850-4004



### ARTLINE N-A-T

Universal knife

(AN) ARTLINE N-A-T No.01 red Art.-No. 850-5001



Cervical waxer

(AN) ARTLINE N-A-T No.02 green Art.-No. 850-5002



Moulding former

(AN) ARTLINE N-A-T No.03 yellow Art.-No. 850-5003



Base former

(AN) ARTLINE N-A-T No.04 blue Art.-No. 850-5004



## electric wax knife



### pro Waxer

Unit with analogue technology:

- very efficient:
  - adjustable 50°C - 225°C / 122°F - 437°F
- large display showing the set temperature
- the handpiece is:
  - extremely light, insensitive against temperature,
  - stays conveniently cool
- the connection cable to the handpiece is highly flexible and absolutely insensitive against temperature

pro WAXER	230 V	Art.-No. 412-0000
pro WAXER	115 V	Art.-No. 412-0115
<i>Delivery: incl. handpiece, handpiece clamp and modelling tip PKT 2 (411-0200)</i>		



### pro WAXER duo

Same as „pro Waxer“, with additional handpiece.

- every handpiece is separately adjustable
- no necessity of second unit

pro WAXER duo	230 V	Art.-No. 414-0000
pro WAXER duo	115 V	Art.-No. 414-0115
<i>Delivery incl.: 2 handpieces, 2 handpiece clamps and 2 modelling tips PKT 1 (411-0100) and PKT 2 (411-0200)</i>		

*Handpieces and modelling tips are not included in the warranty!*





**easy Waxer & easy Waxer duo**

- short heating time
- temperature range up to 190°C/ 374°F
- easy handling due to temperature control
- unbeatable price-performance ratio
- energy-saving on/off switch

easy Waxer	230 V	Art.-No. 413-0000
easy Waxer	115 V	Art.-No. 413-0115
<i>Delivery: incl.: handpiece, handpiece clamp and modelling tip PKT 2</i>		
easy Waxer duo	230 V	Art.-No. 415-0000
easy Waxer duo	115 V	Art.-No. 415-0115
<i>Delivery incl.: 2 handpieces, 2 handpiece clamps and 2 modelling tips PKT 1 and PKT 2</i>		

**duo Waxer NT**

**The first processor-controlled Waxer with a two-line LCD display offers unforeseen advantages which the technician only has known from the computer technology so far.**




- precise and stepless control 50°C - 225°C / 122°F - 437°F
- display shows set temperature of both handpieces at the same time
- **two-line LCD display** with alpha numerical presentation of both working temperatures
- possibility of intergrating personal data through CPU Software
- **possibility of memorising 2 temperatures per handpiece apart from the actual temperature - these temperatures can be changed and retrieved at any time**
- **automatical memorization of the latest data adjustments - also after power failure**
- the handpiece: extremely light, insensitive against temperature, stays conveniently cool
- the connection cable to the handpiece is highly flexible and absolutely insensitive against temperature
  
- **temperature control via high-quality contact sensors** - no use of wear and tear intensive parts, such as potentiometers or foil displays
- ergonomical design - **reduced in weight and dimensions** - no need of an additional transformer

duo Waxer NT	230V	Art.-No. 414-9000
duo Waxer NT	115V	Art.-No. 414-9115
<i>Delivery incl. 2 handpieces, 2 handpiece clamps and modelling tips PKT 1 (411-0100) and PKT 2 (411-0200)</i>		

*Handpieces and modelling tips are not included in the warranty!*

## electric wax knife accessories

### HANDPIECE

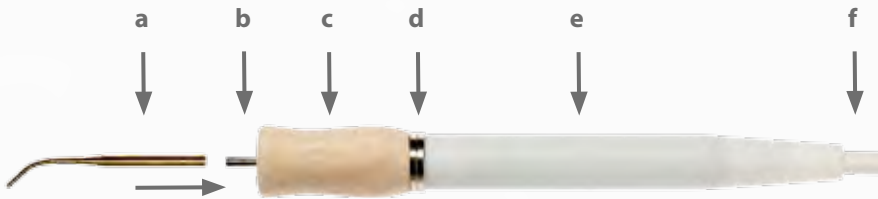

- the handpiece is extremely light and well insulated against heat, remaining cool even after prolonged use
- the handpiece has a highly flexible connection cable/lead



- a = modeling tip to be pinned
- b = heating pin
- c = ergonomically shaped wooden handle
- d = magnetic rings
- e = plastic holder
- f = cable/lead

Handpiece for all NT, Pro and Easy series with straight unipolar connector



Art.-No. 411-1001

Handpiece for Pro Waxer with angled 5-pin connector



Art.-No. 411-1000



### HANDPIECE CLAMP

- for safe fixing of the hanpiece
- just stick and use

Handpiece clamp black

Art.-No. 410-1003

### Modelling tips for Waxer NT and pro



(1) (2) (3) (4) (5) (6)

We offer an assortment of six different tips for all our Waxers series NT, pro and light. The tips can be individually modified by grinding.

The tips can be quickly changed thanks to the simple pin system. Their interior diameter is 2 mm.

- (1) PKT 1 thin Art.-No. 411-0100
- (2) PKT 2 thick Art.-No. 411-0200
- (3) Hollenback horizontal Art.-No. 411-0300
- (4) Knife shape Art.-No. 411-0400
- (5) Corkscrew thin Art.-No. 411-0500
- (6) Corkscrew thick Art.-No. 411-0600

Handpieces and modelling tips are not included in the warranty!





### IQ Waxer NT light

IQ Waxer NT Light with its additional „Touch Control“ function is a successful further development of electric wax knives. You can work as usual in the “continuous”-function, which means that the temperature of the modelling tip stays constant.

Furthermore, you have the possibility to use the extremely fast reacting “Touch Control” function. The modelling tip warms up and cools down in only few seconds. This is possible due to the micro-heating on the front point of the tip. The sensor contact can be controlled either with the contact at the handpiece or with the pedal (optional).

- adjustable 50°C - 230°C / 122°F - 446°F

The handpiece is extremely light, has a diameter of only 11 mm, is insensitive against temperature and stays conveniently cool. The connection cable is also highly flexible and insensitive against temperature. The distance from the touch sensitive ring to the modelling tip can be adapted individually by moving the ring.

IQ WAXER NT light	230 V	Art.-No. 410-9000
IQ WAXER NT light	115 V	Art.-No. 410-9115
<i>Delivery incl.: handpiece, handpiece clamp and modelling tip PKT 2 ( 420-0200 )</i>		
Pedal ( additional to the sensor contact on handpiece )		Art.-No. 409-9500

### MODELLING TIPS for IQ Waxer NT light

Thanks to micro heating elements and the operation of modern laser technology, the heating was placed at the point of the modelling tip. Therefore, the hottest point of maximal heat of the instrument lies at the point of the modelling tip, so that:

- the modelling tips heat up and cool down in fractions of a second
- there is no flow back of the wax drop

Besides, this new technology allows us to produce extraordinarily fine tips in new forms.



Assortment of modelling tips from left to right::

PKT 1 fine	Ø 1,0 mm	Art.-No. 420-0100
PKT 2	Ø 1,2 mm	Art.-No. 420-0200

*Handpieces and modelling tips are not included in the warranty!*

## dive monitoring device



### Dippy pro digital

- convenient device as alternative with analogue technology
- stepless temperature control 65°C - 105°C / 149°F - 221°F
- temperature precision  $\pm 0,10$  °C
- even temperature distribution thanks to direct heating around the tank
- the pot is integrated completely into the housing to guard against burns
- improved ventilation through ventilation slits
- electronic parts are treated with care - longer life of unit
- additional digital display showing the actual temperature in 3 numbers
- with measurement in °C or °F

DIPPY pro digital 230 V  
DIPPY pro digital 115 V

Art.-No. 405-1000  
Art.-No. 405-1115



Dipping wax



### Dippy easy



Small black Dipping Unit, of a superior quality, but economically priced.

Quick – you create accurate wax copings in seconds

Clean – easy to remove excess wax

Efficient – economic and time saving with minimal waste

- analogue technology
- temperature range 65°C - 105°C / 149°F - 221°F
- display accuracy  $\pm 1,0$  °C

DIPPY easy 230 V  
DIPPY easy 115 V

Art.-No. 402-1000  
Art.-No. 402-1115



### Dippy NT digital

**The first processor-controlled dipping pot with a two-line LCD display world-wide.**

With this unit you can read the set and the actual temperature at the same time.

- precise and stepless control 55°C - 105°C / 131°F - 221°F
- even temperature distribution thanks to direct heating system
- two-line LCD Display with alpha numeric indication, simultaneously indicating the set and current temperature
- temperature control via high-quality contact sensors
- no use of wear and tear intensive parts, such as potentiometers or foil displays
- automatical memorization of the latest data adjustments - also after power failure
- safety switch-off upon exceeding the highest admissible temperature
- no need of an additional transformator taking a lot of space
- ergonomical design - reduced in weight and size
- possibility of integrating personal data through CPU Software



Dippy NT digital 230 V  
Dippy NT digital 115 V

Art.-No. 405-9000  
Art.-No. 405-9115

## Light-curing equipment



### PRECI NT SHUTTLE II / IV UVA light curing unit ( with 2 or 4 light tubes )

UVA light curing unit for the polymerisation of light curing material and synthetic material.

The models can be safely positioned on the bottom of the drawer and be pushed into the interior of the unit.

The interior is big enough to take up to 4 models.

In case of closed drawer protection against UV rays by an anti-glare.

- already preset polymerisation time
- optical and acoustic signal after the polymerisation process is finished
- upon request with 2 or 4 light tubes with an emitter zone of 350 nm  
(PRECI NT SHUTTLE II = 2 light tubes,  
PRECI NT SHUTTLE IV = 4 light tubes)
- by changing the light tubes with another UVA spectral range, the unit can also be used for composites and LC block out material
- One-Touch switching
- no energy consumption during stand-by operation

Polymerisation time: about 5 min

Power:

PRECI NT SHUTTLE II approx. 20 Watt

PRECI NT SHUTTLE IV approx. 35 Watt

Size:

opening size 55 mm x 135 mm

drawer internal size 130 mm x 155 mm



PRECI NT SHUTTLE II	230V	Art.-No. 492-9000
PRECI NT SHUTTLE II	115V	Art.-No. 492-9115
PRECI NT SHUTTLE IV	230V	Art.-No. 494-9000
PRECI NT SHUTTLE IV	115V	Art.-No. 494-9115
Light tube	9W BL 350	Art.-No. 490-0001
Light tube for Composit	9W BL 450	Art.-No. 490-0010



## Light-curing equipment



### PRECI TRAY light curing material

- fits the upper and lower jaws
- easy to adapt to the model
- available in 3 different colours
- time-saving due to less reworking
- special filling materials were added to obtain an even layer without streaks
- sheet size: 70 x 100 mm
- thickness: 2,2 mm
- Package content: 50 pcs.

#### PRECI TRAY

rose  
blue  
white

Art.-No. 890-1000  
Art.-No. 890-1010  
Art.-No. 890-1020



### Build Up Resin LC

modelling light cure resin for casting technique

- easy to separate from primary construction
- high stability
- contouring accuracy
- very good casting results
- burns without residue
- outstanding burn out qualities - no break down of muffles
- easy to cut




#### Build Up Resin LC

2 x 3 g syringe (including 4 tips)

Art.-No.890-1300



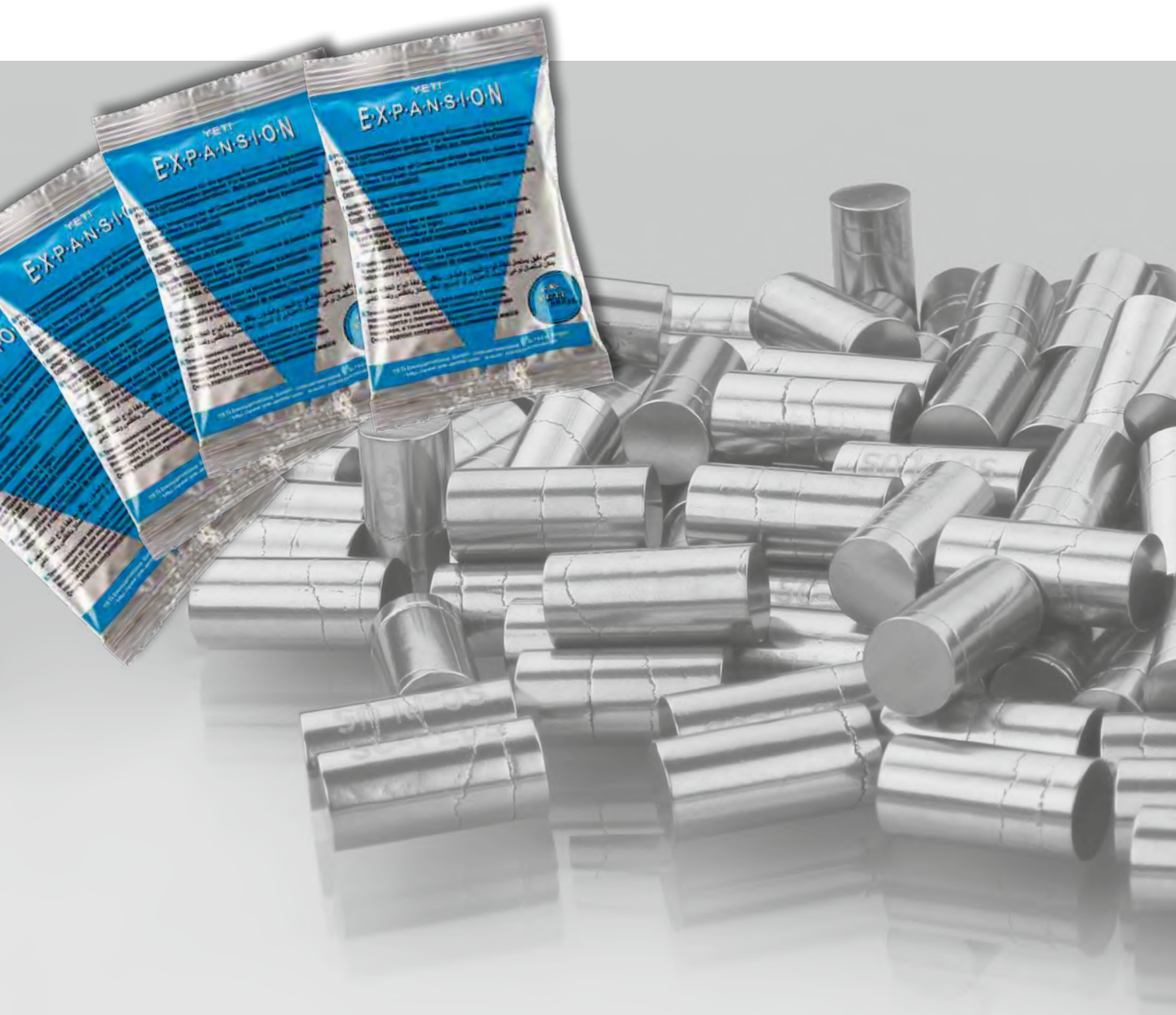





- Investment
- Alloys

## Investment



### YETIVEST

phosphate-bound precision investment for precious, metal-to-ceramic and non precious alloys

- fine grained consistency
- very smooth surfaces
- exact fitting

Casting: traditional or high-speed

YETIVEST Powder	4,0 kg (25 x 160 g)	Art.-No. 930-0160
YETIVEST Liquid	1000 ml bottle	Art.-No. 931-1000
YETIVEST Powder + Liquid	4,0 kg (25 x 160 g) + 1000 ml	Art.-No. 932-0000

**For all works where high expansion is required we advice YETI EXPANSION and YETI EXPANSION plus investment!**



### DUOVEST

Phosphate-bound precision investment for all dental alloys

- insensitive universal investment
- quality efficient alternative at a reasonable price

Casting: traditional or high-speed

DUOVEST Powder	4,0 kg (25 x 160 g)	Art.-No. 950-0160
DUOVEST Liquid	1000 ml bottle	Art.-No. 951-1000
DUOVEST Powder + Liquid	4,0 kg (25 x 160 g) + 1000 ml	Art.-No. 952-0000



### YETI EXPANSION

investment for precious and non precious alloys

- the expansion is regulated by the concentration of liquid and distilled water and delivers the desired results for both precious and non-precious alloys
- suitable for casting work, crowns and bridges technique, telescopic technique and implant technique
- high expansion values allow to make large non precious works
- exact control of expansion while making primary telescopic crowns
- cracking resistant, can be used also without muffle ring
- excellent price-performance ratio

Casting: traditional or high-speed

YETI EXPANSION Powder	4,5 kg (50 x 90 g)	Art.-No. 955-0090
YETI EXPANSION Liquid	1000 ml bottle	Art.-No. 955-1000
YETI EXPANSION Powder + Liquid	4,5 kg (50 x 90 g) + 1000 ml	Art.-No. 955-0000





### COBAVEST

phosphate-bound precision investment  
for partial denture technique

- suitable for both, silicone and gel duplication
- exact fitting
- excellent duplicating details
- fine grain consistency
- excellent fluidity

Casting: traditional or high-speed

COBAVEST Powder	20 kg (100 x 200 g)	Art.-No. 940-0200
COBAVEST Liquid	1000 ml bottle	Art.-No. 941-1000
COBAVEST Powder + Liquid	20 kg (100 x 200 g) + 3000 ml	Art.-No. 942-0000



### K2 Press-universal

Universal speed investment for all pressed ceramic

- also ideal for lithium disilicate
- low reaction layer
- exact fitting
- smooth and homogeneous press results
- fine flowing, creamy consistency
- easy devesting

Casting: high-speed

K2 Press-universal Powder	5 kg ( 50 x 100 g )	Art.-No. 360-0100
K2 Press-universal Liquid	1000 ml bottle	Art.-No. 360-1000
K2 Press-universal Powder + Liquid	5 kg ( 50 x 100 g ) + 1000 ml	Art.-No. 360-1005



### YETI EXPANSION plus

Universal investment for precious and non-precious alloys and for pressed ceramic

- expansion is regulated by the concentration of liquid and distilled water and delivering accurate results for precious, non-precious alloys and for pressed ceramic
- pressed ceramic restoration for inlays, veneers, press-to-metal and bonded techniques
- high expansion capability allows for large non-precious restorations to be cast accurately without the need for a separate high-expansion liquid
- economic but highly accurate
- easy to invest
- fine grain size and excellent flow characteristics result in extremely smooth casting surface
- 100g sachets

Casting: traditional or high-speed

YETI EXPANSION plus Powder	5 kg (50 x 100g)	Art.-No. 956-0100
YETI EXPANSION plus Liquid	1000 ml bottle	Art.-No. 956-1000
YETI EXPANSION plus Powder + Liquid	5 kg (50 x 100g) 1000 ml	Art.-No. 956-0000

### Silicone cylinder for 2g ingots, slightly conical, transparent



Conical cylinder form allows easy separation of investment.  
No additional base shaping for fitting into muffle is necessary!

SILICONE CYLINDER small  
Ø 37 mm, 100g +base complete

Art.Nr. 350-0701

SILICONE CYLINDER medium  
Ø 50 mm, 200g +base complete

Art.Nr. 350-0703



### Silicone ring for 5g ingots, tube form

SILICONE RING large, 300g  
Ø 50mm inner muffle base, Ø 60mm outer muffle base  
base, ring and muffle base shaper

SILICONE RING large

Art.Nr. 350-0706


### Plunger

Investment stamp for single use.  
Can be set into muffle without preheating.  
Cools down quickly,  
low possibility of fissuring.  
Available in 2g and 5g ingots.

Plunger (one-way) small (2g)  
Plunger (one-way) large (5g)

5 pcs.  
9 pcs.

Art.No. 350-0205  
Art.No. 350-0505

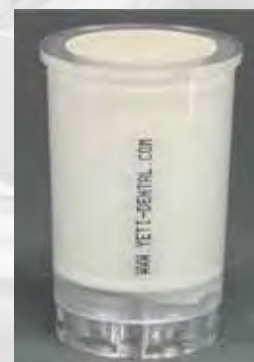



## BLUE MARKER high spot indicator

For indication and elimination of high spots when fitting crowns, inlays, etc.

- leaves a thin and even layer on the plaster die
- dries quickly and is well covering
- points of contact on the crown's inner sides are exactly indicated
- adheres to metal, ceramics, plaster and synthetics
- can be easily washed off from sealed plaster surfaces (to be washed off before porcelain firing or investing)
- 100% water soluble



### BLUE MARKER high spot indicator

Blue Marker	1 x 18 ml	Art.-No. 580-0018
Blue Marker	6 x 18 ml	Art.-No. 580-6018
Thinner for Blue Marker	1 x 18 ml	Art.-No. 580-2018
Thinner for Blue Marker	6 x 18 ml	Art.-No. 580-2618

## OCCLUSAL SPRAY

For the recognition of high spots

- targeted spraying
- does not blob
- 100% water soluble
- can be easily washed off from sealed plaster surfaces (to be washed off before porcelain firing or investing!)
- adheres to metal, ceramics, plaster and synthetics
- environmentally friendly: mixture of compressed gas, not harmful to the ozone layer
- Occlusive spray only 1 nozzle included

Content: 75 ml



OCCLUSAL SPRAY, green	Art.-No. 580-2007
OCCLUSAL SPRAY, red	Art.-No. 580-2008
OCCLUSAL SPRAY, blue	Art.-No. 580-2009
OCCLUSAL SPRAY, white	Art.-No. 580-2010

## SATI BRUSH

A matt finish of the occlusal surface can indicate possible friction points after a while inside the mouth. As soon as they are visible, they can be removed easily.

SATI BRUSH is used after polishing at slow speed on the occlusal surface.

SATI BRUSH can also be used to remove oxide layers on metal surfaces.



SATI BRUSH polishing wheels	22 mm	50 pcs.	Art.-No. 310-0000
-----------------------------	-------	---------	-------------------

## Alloys



### SOLIBOND N

metal-to-ceramic , Nickel-Chrome alloy

- beryllium free nickel-chrome alloy
- good manipulation when fluid, security is achieved by casting.
- fast polishing, shining surfaces as well as burning stability of the alloy
- exact fitting and easy manipulation with a hardness of 180 HV
- very thin oxide layer, esthetic construction of the ceramic
- optimum metal-ceramic bonding
- durable, high corrosion resistance, can be scanned with laser

Tip: best casting results in combination with the investments YETIVEST, DUOVEST, YETI EXPANSION and YETI EXPANSION plus YETIVEST, DUOVEST, YETI EXPANSION and YETI EXPANSION plus

Melting range	1.250 – 1.330°C
Casting temperature	1390°C
Vickers Hardness	180 HV
Ultimate elongation	25,0 %
CTE (20 - 600°C)	14,1 · 10 <sup>-6</sup> K <sup>-1</sup>
CTE (25 - 500°C)	13,9 · 10 <sup>-6</sup> K <sup>-1</sup>

Ni	Cr	Mo	Si	Fe
62,7%	24,5%	10,4%	1,4%	1,0%

SOLIBOND N	1000 g	Art.-No. 963-1000
SOLIBOND N	250 g	Art.-No. 963-0250

### LOT UNIVERSAL

universal lot for all non precious alloys

- durable in the patient's mouth, good colour resistance, easy to handle

LOT UNIVERSAL	5g	Art.-No. 960-0000
---------------	----	-------------------



**SOLIDUR CoCr**

partial denture alloy



- cobalt-chrome partial denture alloy, can be scanned with laser
- easy to handle due to the thin oxide layer
- excellent polishing and shine characteristics, exact fitting
- excellent fluidity and duplicating detail
- physical characteristics which permits the shaping of delicate structures
- high biocompatibility - durable in the patient's mouth

Tip: best casting results in combination with COBAVEST investment

Melting range	1360 -1410°C
Casting temperature	1470°C
Vickers Hardness	420 HV
Elongation	6,5%,
Ultimate elongation	8%
E-Modul	ca. 230.000 MPa

Co	Cr	Mo	Mn	C	Fe	Si
63%	29,4%	5,95%	0,6%	0,29%	0,06%	0,7%

SOLIDUR CoCr	1000 g	Art.-No. 960-1000
SOLIDUR CoCr	250 g	Art.-No. 960-0250

**SOLIBOND C plus**

Non precious metal-to-ceramic Cobalt-Chrome alloy



- nickel and beryllium free cobalt-chrome alloy
- very thin oxide layer, esthetic construction of the ceramic
- optimum metal-ceramic bonding
- good manipulation when fluid - security for casting results, fast polishing, shining surfaces as well as burning stability of the alloy
- exact fitting (time saving) - biocompatible - suitable to be scanned by laser
- improved flowing qualities
- improved casting results

Tip: best casting results in combination with the investments YETIVEST, DUOVEST, YETI EXPANSION and YETI EXPANSION plus


Melting range	1.310 – 1.370°C
Casting temperature	1430°C
Vickers Hardness	280 HV
Ultimate elongation	10,1 %
CTE (20 - 600°C)	14,0 · 10 <sup>-6</sup> K <sup>-1</sup>
CTE (25 - 500°C)	13,9 · 10 <sup>-6</sup> K <sup>-1</sup>

Co	Cr	Mo	W	Nb	Si
63%	24%	2,9%	8,1%	0,9%	1,1%

SOLIBOND C plus	1000 g	Art.-No. 969-1000
SOLIBOND C plus	250 g	Art.-No. 969-0250





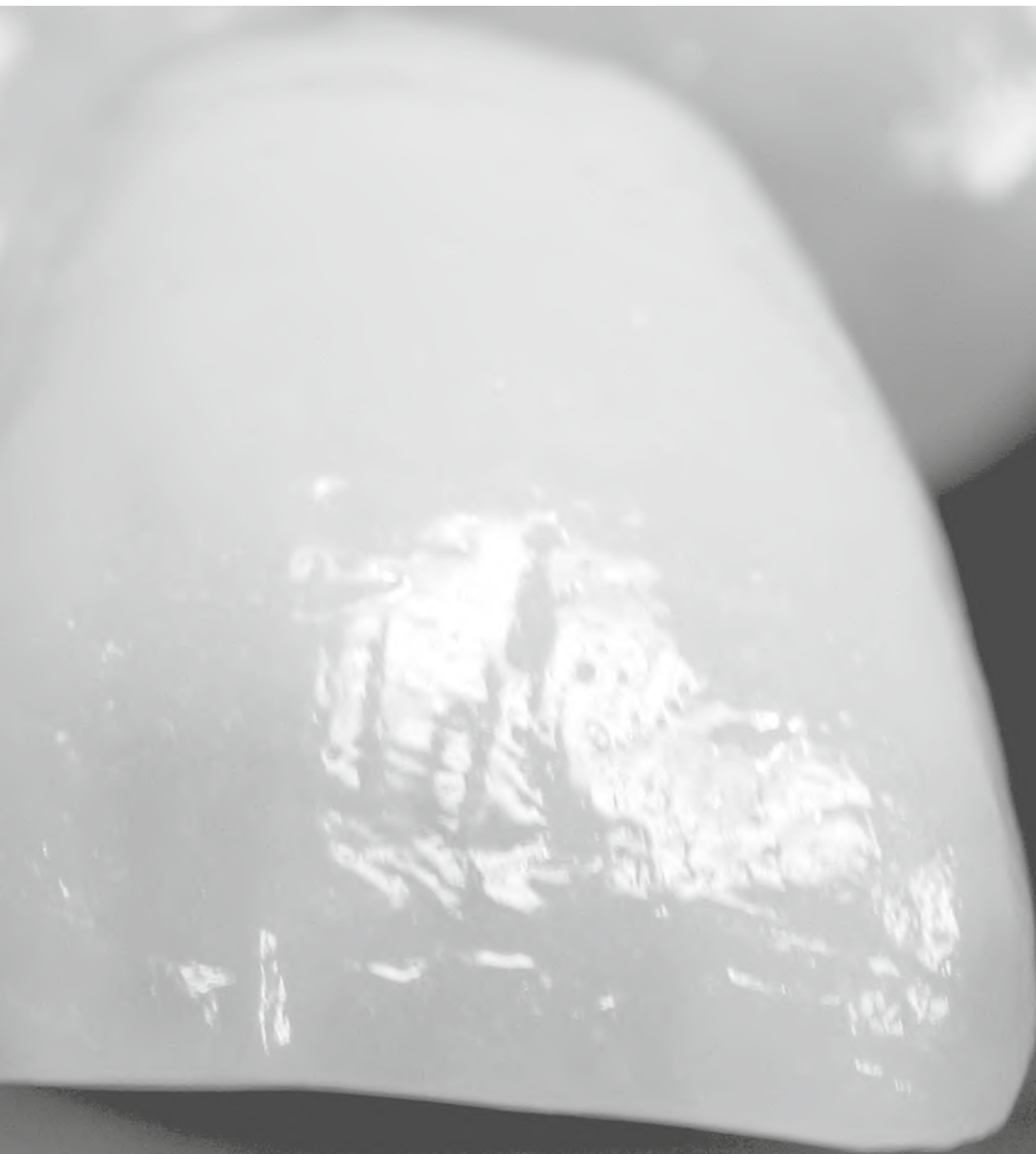

- Ceramic
- Ceramic accessories

Ceramic K2 LF







## K2 LF and K2 LF press

### natural, multi-purpose and easy to use

K2 LF is a highly condensed **2-phase leucite ceramic**.

It can be fired at 770°C and is suitable for all alloys

having a coefficient of thermal expansion **CTE between 13,8 and 14,9/25°-500°C**.

This ceramic system uses the nuances and reflection of natural teeth. It offers a vital result even in cases where there is insufficient space.

K2 is the first ceramic system developed as a comprehensive solution, including the Layering ceramic, Press ceramic and the Press over Metal (PoM) systems. The economical K2 ceramic system is revolutionary, the user needs only one ceramic layer for metal - as well as for full ceramics.

K2 assortments are available in the A1-D4 colours and include Opal, Fluor, Chroma and Intensive masses.

A higher shade measurement can be obtained by mixing or adding the various colours before firing. The K2 ceramic's press ingots may be used when working with full ceramic, inlays, onlays and PoM. There is a probability of 100% correspondence with the selected colour, even in the case of personalisation, as a result of the diversity of opalescence and fluorescence..



### Glaze paste for pressable ceramics

Only one glaze firing required

Content: 5 g

Glaze paste

Art.-No. 336-0010

# Ceramic K2 LF



K2 LF press: Table of contents for pressable ceramic ingots

V, I, O	K	PoM	Qty.	Fluor	Opazität	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4	Bleach D
X	XX	XX	4x2g	XX	75 %	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4	BD
Art.-No. 350-08..						..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16	..17
X	XX	X	4x2g	XX	60 %	TA1	TA2	TA3			TB1	TB2	TB3		TC1	TC2			TD2			
Art.-No. 350-08..						..18	..19	..20			..21	..22	..23		..24	..25			..26			
X	XX	X	3x5g	XX	60 %	TA1	TA2	TA3			TB1	TB2	TB3		TC1	TC2			TD2			
Art.-No. 350-15..						..01	..02	..03			..04	..05	..06		..07	..08			..09			

XX	X	X	4x2g	X	45 %	T1+	T2+	T3+	T4+
Art.-No. 350-08..						..27	..28	..29	..30
XX	XX	XX	4x2g	X	40 %	S57	S58	S59	S60
Art.-No. 350-08..						..35	..36	..37	..38
XX	X	X	4x2g	X	35 %	T1	T2	T3	T4
Art.-No. 350-08..						..31	..32	..33	..34

**V, I, O = Veneer, Inlay, Onlay**  
**K = Crowns**  
**PoM = Press over Metal**

## K2 LF Basis - Powder kit in cabinet

16 x 15 g Opaque powder ( A1 - D4 ), 16 x 15 g Opaque dentin ( A1 - D4 ), 16 x 20 g Dentin ( A1 - D4 ), 4 x 20 g Enamel ( E 57 - E 60 ), 4 x 20 g Enamel opal ( EO 57 - EO 60 ), 1 x 20 g Clear, 1 x 20 g Transpa normal, 1 x 10g Glaze, 1 x 5 g Bonder, 1 x 50 ml Opaque liquid, 1 x 50 ml Modelling liquid, 1 x 25 ml Shoulder liquid, 1 x 25 ml Glaze fluid

Art.-No. 320-0000

## K2 LF Basis - Paste kit in cabinet

16 x 15 g Paste opaque ( A1 - D4 ), 16 x 15 g Opaque dentin ( A1 - D4 ), 16 x 20 g Dentin ( A1 - D4 ), 4 x 20 g Enamel ( E 57 - E 60 ), 4 x 20 g Enamel opal ( EO 57 - EO 60 ), 1 x 20 g Clear, 1 x 20 g Transpa normal, 1 x 10g Glaze, 1 x 5 g Bonder, 1 x 50 ml Modelling liquid, 1 x 25 ml Shoulder liquid, 1 x 25 ml Glaze fluid

Art.-No. 320-0000/2

## K2 LF Confidential - A3 kit in paper box

1 x 15 g Opaque powder A3, 1 x 20 g Enamel E 59, 1 x 2 g Bonder, 1 x 20 g Dentin A3, 1 x 10 g Glaze, 1 x 20 g Transpa normal, 1 x 50 g Opaque liquid, 1 x 50 ml Modelling liquid, 1 x 25 ml Glaze fluid

Art.-No. 320-0001

## K2 LF Advanced - Powder kit in cabinet

1 x 15 g Bleach Opaque, 1 x 15 g Opaque gingiva powder, 7 x 15 g Shoulder powder ( SH1 - SH7 ), 1 x 15 g Shoulder gingiva, 1 x 20 g Bleach dentin, 1 x 15 g Gingiva, 4 x 15 g Chroma dentin ( CH A - CH D ), 1 x 20 g Transpa opal, 3 x 20 g Transpa ( T1 - T3 ), 3 x 15 g Mamelon ( MM1 - MM3 ), 4 x 20 g Enamel intensive ( T1, T2, O3, O4 ), 2 x 15 g K2 fluor ( PF 1, PF 2 )

Art.-No. 320-0002

## K2 press EF - Kit in paper box

Press ingots 4 x 2 g each ( A1 - A3, B1 - B3, C1, C2, D2, BD, TA1 - TA3, TB1 - TB3, TC1, TC2, TD2 T1+ - T3+, T4+, T1 - T4, S57 - S60 ), 10 x 2g Plunger (one-way)

Art.-No. 320-0004

## Enamel colours classification for K2 LF and K2 LF press

Color	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4
Dentine	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4
Enamel	E58	E58	E59	E59	E60	E57	E59	E59	E59	E60	E59	E59	E60	E60	E59	E59





Table of contents K2 LF

Color	Qty.	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4	Bleach	Gingiva		
Powder opaque O		<b>A1</b>	<b>A2</b>	<b>A3</b>	<b>A3,5</b>	<b>A4</b>	<b>B1</b>	<b>B2</b>	<b>B3</b>	<b>B4</b>	<b>C1</b>	<b>C2</b>	<b>C3</b>	<b>C4</b>	<b>D2</b>	<b>D3</b>	<b>D4</b>	<b>BO</b>	<b>OG</b>		
Art.-No. 320-15..	15g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16	322-1500	320-1517		
Art.-No. 320-50..	50g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16				
Paste opaque PO		<b>A1</b>	<b>A2</b>	<b>A3</b>	<b>A3,5</b>	<b>A4</b>	<b>B1</b>	<b>B2</b>	<b>B3</b>	<b>B4</b>	<b>C1</b>	<b>C2</b>	<b>C3</b>	<b>C4</b>	<b>D2</b>	<b>D3</b>	<b>D4</b>	<b>POBO</b>	<b>POG</b>		
Art.-No. 321-03..	3g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16	322-0300	321-0317		
Shoulder		<b>SH1</b>	<b>SH2</b>	<b>SH3</b>	<b>SH4</b>	<b>SH5</b>	<b>SH6</b>	<b>SH7</b>											<b>SH G</b>		
Art.-No. 324-15..	15g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16		..20		
Dentine D		<b>A1</b>	<b>A2</b>	<b>A3</b>	<b>A3,5</b>	<b>A4</b>	<b>B1</b>	<b>B2</b>	<b>B3</b>	<b>B4</b>	<b>C1</b>	<b>C2</b>	<b>C3</b>	<b>C4</b>	<b>D2</b>	<b>D3</b>	<b>D4</b>	<b>BD</b>	<b>G (15g)</b>		
Art.-No. 325-20..	20g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16	326-2000	338-1500		
Art.-No. 325-50..	50g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16	326-5000			
Chroma		<b>CH A</b>			<b>CH B</b>			<b>CH C</b>			<b>CH D</b>										
Art.-No. 327-15..	15g	..01			..02			..03			..04										
Art.-No. 327-50..	50g	..01			..02			..03			..04										
Opaque dentin OD		<b>A1</b>	<b>A2</b>	<b>A3</b>	<b>A3,5</b>	<b>A4</b>	<b>B1</b>	<b>B2</b>	<b>B3</b>	<b>B4</b>	<b>C1</b>	<b>C2</b>	<b>C3</b>	<b>C4</b>	<b>D2</b>	<b>D3</b>	<b>D4</b>				
Art.-No. 323-15..	15g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16				
Art.-No. 323-50..	50g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11	..12	..13	..14	..15	..16				
Enamel		<b>E 57</b>			<b>E 58</b>			<b>E 59</b>			<b>E 60</b>										
Art.-No. 328-20..	20g	..01			..02			..03			..04										
Art.-No. 328-50..	50g	..01			..02			..03			..04										
Enamel opal		EO 57			EO 58			EO 59			EO 60										
Art.-No. 330-20..	20g	..01			..02			..03			..04										
Art.-No. 330-50..	50g	..01			..02			..03			..04										
Enamel intensive		EI T1 (transpa)			EI T2 (transpa)			EI O3 (opaque)			EI O4 (opaque)										
Art.-No. 329-20..	20g	..01			..02			..03			..04										
Art.-No. 329-50..	50g	..01			..02			..03			..04										
Transparent		TN (normal)			TO (opal)			T1			T2			T3							
Art.-No. ....	20g	332-2000			333-2000			334-2001			334-2002			334-2003							
Art.-No. ....	50g	332-5000			333-5000			334-5001			334-5002			334-5003							
Clear	20g / 50g	<b>CL</b>			20g = Art.-No. 335-2000			50g = Art.-No. 335-5000													
Glaze	10g	<b>GL</b>			10g = Art.-No. 336-1000																
K2 Fluor		PF 1			PF 2																
Art.-No. 331-15..	15g	..01			..02																
Mamelon		<b>MM 1</b>			<b>MM 2</b>			<b>MM 3</b>													
Art.-No. 339-15..	15g	..01			..02			..03													
K2 Bonder	5g	<b>B</b>			5g = Art.-No. 337-0500																
K2 Liquid		<b>Opaque Liquid OL</b>			<b>Modelling Liquid ML</b>			<b>Shoulder Liquid SHL</b>			<b>Glaze Fluid GF UNIVERSAL</b>										
	25ml	Art.-No. 341-0050			Art.-No. 342-0050			Art.-No. 343-0025			Art.-No. 387-9150										
	50ml	Art.-No. 341-0250			Art.-No. 342-0050			Art.-No. 343-0050													
	250ml	Art.-No. 341-0250			Art.-No. 342-0250																
Stain UNIVERSAL		<b>ST-A</b>	<b>ST-B</b>	<b>ST-C</b>	<b>ST-D</b>	<b>ST-5 WHITE</b>	<b>ST-6 HONEY</b>	<b>ST-7 GREY</b>	<b>ST-8 BROWN</b>	<b>ST-9 BLUE</b>	<b>ST-10 PINK</b>	<b>ST-11 RED</b>									
Art.-No. 387-83..	2g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11									

# K2 *Li*





# K2 *Li*

K2 Li is the new Lithium disilicate ceramic in the K2 Ceramic assortment from Company Yeti Dental. It includes a Press system, as well as a layered ceramic, which can be also easily used for zirconium restorations.

This affects a cost-reducing for every laboratory.

The system impresses with:

- easy handling
- a natural fluorescence and translucency
- high stability
- natural aesthetics



## Ceramic K2 Li

full anatomic, Staining technique



Cut-back technique



Layering technique



transfer in K2 Li Press (see above) and finished work after layering and staining (see below)



### K2 Li Press ingots

The Press ingots ( 3g ) are offered in 3 different translucencies:

	Art.-No.	A1	A2	A3	B1	B2	C1	BL1	BL2	BL3	MO1
<b>HT</b> Highly translucent	387-80..	..01	..02	..03	..11	..12	..21	..31	..33	..32	
<b>LT</b> Low translucent	387-81..	..01	..02	..03	..11	..12	..21	..31	..32	..33	
<b>MO</b> Medium opaque	387-82..										..01

### indication

		HT Highly translucent	LT Low Translucent	MO Medium opaque
Processing technique	full anatomic / Staining technique	v	v	
	Cut-Back technique	v	v	
	Layering technique			v
indication	Thin veneer	v	v	
	Veneer	v	v	
	Inlay	v	v	
	Onlay	v	v	
	Partial Crown	v	v	
	Anterior Crown	v	v	v
	Posterior Crown	v	v	v



# K2 Li

K2 Li Enamel allocation table																	
Vita color	Bleach	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4
dentin	Bleach	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4
incisals	E-BL	E58	E58	E59	E59	E60	E57	E59	E59	E59	E60	E59	E59	E60	E60	E59	E59

Product overview K2 Li Layered ceramic																				
Color	Qty.	A1	A2	A3	A3,5	A4	B1	B2	B3	B4	C1	C2	C3	C4	D2	D3	D4	Bleach	Gingiva	
Dentin D- Art.-No. 387-10..	20g	A1 ..01	A2 ..02	A3 ..03	A3,5 ..04	A4 ..05	B1 ..06	B2 ..07	B3 ..08	B4 ..09	C1 ..10	C2 ..11	C3 ..12	C4 ..13	D2 ..14	D3 ..15	D4 ..16	BL1 ..17	BL3 ..19	G ..18
opaque dentine Art.-No. 387-20..	20g	OD-HO (HONEY) ..01									OD-WH (WHITE) ..02									
incisals Art.-No. 387-30..	20g	E-57 ..01			E-58 ..02			E-59 ..03			E-60 ..04			E-BL (Bleach) ..05						
intensive incisals Art.-No. 387-40..	20g	EI-WH (WHITE) ..01				EI-YE (YELLOW) ..02				EIO-NT (OPAL NEUTRAL) ..03				EIO-BL (OPAL BLUE) ..04						
Transparent materials Art.-No. 387-50..	20g	TN (NORMAL) ..01			TR-VL (VIOLET) ..02			T-BL (BLUE) ..03			T-YE (YELLOW) ..04			T-OR (ORANGE) ..05						
Clear	20g	CL Art.-No. 387-5006																		
Mamelon Art.-No. 387-60..	20g	MM-OR (ORANGE) ..01					MM-HO (HONEY) ..02					MM-IV (IVORY) ..03								
Modeling Liquid	50ml 250ml	ML UNIVERSAL Art.-No. 387-9050 ML UNIVERSAL Art.-No. 387-9060																		
Glaze fluid	50ml	GF UNIVERSAL Art.-No. 387-9150																		

Product overview Stains Universal												
Stain	Qty.	ST-A	ST-B	ST-C	ST-D	ST-5 WHITE	ST-6 HONEY	ST-7 GREY	ST-8 BROWN	ST-9 BLUE	ST-10 PINK	ST-11 RED
Art.-No. 387-83..	2g	..01	..02	..03	..04	..05	..06	..07	..08	..09	..10	..11

Product overview K2 Li Press-ingots										
Press-ingots HT High translucent	Qty.	HT A1	HT A2	HT A3	HT B1	HT B2	HT C1	HT BL1	HT BL2	HT BL3
Art.-No. 387-80..	5x3g	..01	..02	..03	..11	..12	..21	..31	..33	..32
Press-ingots LT Low translucent	5x3g	LT A1 ..01	LT A2 ..02	LT A3 ..03	LT B1 ..11	LT B2 ..12	LT C1 ..21	LT BL1 ..31	LT BL2 ..32	LT BL3 ..33
Press-ingots MO Medium opaque	5x3g	MO Art.-No. 387-8201								

## Mixing plates



### Ceramic Humidor

Ceramic mixing plate transparent

- constant consistency of ceramic
- **hydrophile PE inlay with defined width of pores, assures optimal absorption**
- **one inlay** without additional annoying membrane
- **easy to clean with a steamer or an ultrasonic unit**
- convenient height
- handy size (200 x 125 x 17 mm)
- transparent cover to protect against dust and quick drying
- Made in Germany

Ceramic Humidor  
Inlet grey

2 pcs.

Art.-No. 312-0000  
Art.-No. 312-0100



### Stain Humidor

Humidity plate for stains

With the further development of the Stain-Humidor you can treat the ceramic colours in the same way as the ceramic powder.

- unique capillary inlay with **15 hollows and 2 inserts**
- **automatically supply of the stains with glaze liquid**
- the stains can be used for months
- **mixed stains do not dry out**
- colours are provided with exact liquid quantity while mixing unlike on the usual mixing plates
- Stain Humidor works as Ceramic Humidor **without additional membrane**
- transparent cover protects humidity plate against dust
- **easy to clean with a steamer or an ultrasonic unit**

Stain Humidor

Art.-No. 312-0001





## Caramista

**Caramista - Ceramic brushes**  
fine brushes with wooden handle



no. 0 no. 1 no. O/B no. 4 no. 6 no. 8 no. 10



- **excellent water absorption and spreading** - even ceramic wetting while layering
- the natural elasticity of the hair allows the exact separation as well as precise application of ceramic
- **considerable prolongation of the lifetime of the brush hair**
- the brush **returns** to its original **point form** even after long time of use, rumped and spreaded brush points belongs in the past
- **exclusive brush hair** surpasses known standards
- convenient to use due to **ergonomically formed brush parts**
- **handmade in Germany**

Caramista Paint Brush no. 0	Art.-No. 312-3000
Caramista Paint Brush no. 1	Art.-No. 312-3001
Caramista Opaquer / Bonder Brush	Art.-No. 312-4000
Caramista Brush no. 4	Art.-No. 312-2004
Caramista Brush no. 6	Art.-No. 312-2006
Caramista Brush no. 8	Art.-No. 312-2008
Caramista Big Brush no. 10	Art.-No. 312-2010
Caramista Brush Set	Art.-No. 312-5000
(Content Caramista Brush Set: no. 1, no. O/B, no. 6, no. 8, no. 10)	

## Caramista N

**Caramista N – modelling brush for ceramic.**  
Natural hair brush from high-value Kolinsky sable hair.



no. 4 no. 6 no. 8 no. 10



- **Perfect tension**
- thick hair as water reservoir
- thin brush point for exact placement of ceramic
- in 4 different sizes
- **handmade in Germany**

Caramista N Brush no. 4	Art.-No. 313-2004
Caramista N Brush no. 6	Art.-No. 313-2006
Caramista N Brush no. 8	Art.-No. 313-2008
Caramista N Big Brush no. 10	Art.-No. 313-2010

## Ceramic accessories



### Magic-Glass-Stick for ceramic



- Glass-Stick for an easy and economic application of the opaquer powder onto framework
- with 2 different ball tips for exact quantity application
- even layering after firing
- safety glass with a special stability
- easy to work, saves time

Magic-Glass-Stick

Art.-No. 312-1000



### MAGIC DISC

- diamond-covered separating discs for ceramics
- ultra-thin and flexible
- break-resistant - glass fibre reinforced
- suitable for grinding of ceramics and zircon
- high efficiency while separation of pressed channels
- very high stability and durability




MAGIC DISC

10 pcs. / Ø 20 x 0,2 mm

Art.-No. 399-2002



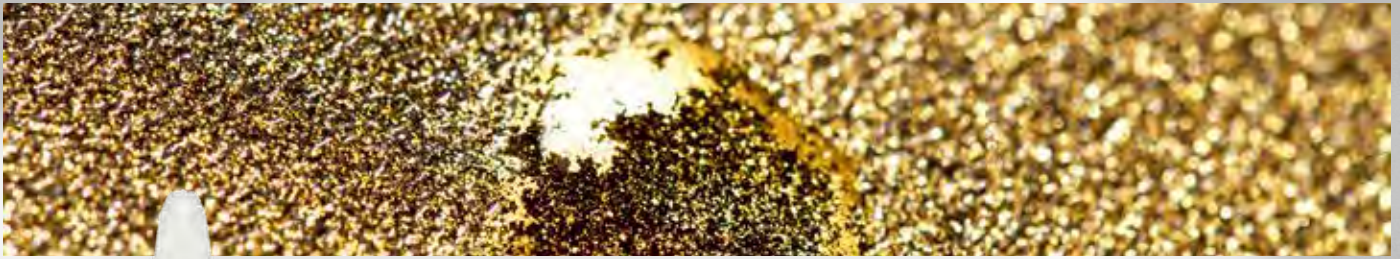

### SURFACE Marker

For the optical indication of the tooth surface structure. Works perfectly both with gypsum models and ceramics.

### Brilliant results ...with just one assistive product

Surface Marker - ideal product for sophisticated dental technicians and students, as well as for master examination

Please use the Liquid before applying the Surface Marker powder!  
This serves the purpose of providing easier cleaning of working model, better adhesion and equable distribution of powder.



- very cost-efficient
- very thin and smooth surfaces
- the best quality-price ratio



SURFACE Marker Set		Art.-No. 319-0000
Content: SURFACE Marker, SURFACE Liquid, Caramista Opaquer / Bonder Brush		
SURFACE Marker	5g	Art.-No. 319-0005
SURFACE Liquid	25ml	Art.-No. 319-0025
Caramista Opaquer / Bonder Brush		Art.-No. 312-4000

## Ceramic accessories



### Cabinet for ceramic (empty)

Cabinet for ceramic with 4 drawers grey/white (standard)  
Size:  
Total height 30,4 cm  
One drawer: L x W x H = 25,5 cm x 34,0 cm x 7,4 cm

Cabinet 4-drawers grey/white	Art.-No. 320-0010
1-drawer grey/white (as add on)	Art.-No. 320-0011



### Shade guide

Dentin	Art.-No. 320-0100
Shoulder/ Fluor / Intensive / Chroma	Art.-No. 320-0200
Transpa / Enamel / Mamelon	Art.-No. 320-0300
Press ingots	Art.-No. 320-0400



### K2 Instruction for use:

german	Art.-No. 310-0100
english	Art.-No. 310-0200
french	Art.-No. 310-0300
italian	Art.-No. 310-0400
spanish	Art.-No. 310-0500
russian	Art.-No. 310-0600

### K2 Lube ceramic separating liquid

This ceramic separating liquid ensures reliable separation of ceramic from gypsum.  
It burns without residue and causes no colour changes of ceramic.

K2 Lube	18ml	Art.-No. 345-0018
---------	------	-------------------





**DIA GLACE**

A polishing medium to obtain perfect surfaces of ceramic restorations. If ceramic is fired at too high temperature, in order to obtain a gleaming surface, it can appear artificial inside the mouth. Natural looks can only be achieved through polishing. We recommend the use of our lubricant-impregnated felt wheels with DIA GLACE.



Art.-No. 305-0100

- high diamond content
- high gloss within a few seconds at a speed of 7500 r/min
- economical in use thanks to its hard consistency, a very small quantity needed for one polishing procedure
- no wasted material, very clean while working

DIA GLACE can also be used on composites, veneer plastics and alloys. Use in combination with our non-impregnated felt wheels only on alloys.



Art.-No. 301-0100

DIA GLACE hard diamond paste 5 g	Art.-No. 300-0000
DIA GLACE Set, 5 g diamond paste + 100 lubricant-impregnated felt wheels	Art.-No. 302-1000
Felt wheels 100 pcs. lubricant-impregnated	Art.-No. 301-0100
Felt wheels 100 pcs. non-impregnated	Art.-No. 305-0100




**DIA QUICK**


Diamond-impregnated felt wheels

DIA QUICK are diamond-impregnated felt wheels for the quick polishing of ceramics, alloys and synthetic materials. When the wheels are in use, heat is generated by friction, releasing the impregnated diamond particles at the contact point. Consequently a permanent supply of diamond particles is created resulting in the perfect polish.  
- quick, clean, small and effective

DIA QUICK 10 pcs. diamond-impregnated felt wheels	Art.-No. 303-0000
---	-------------------

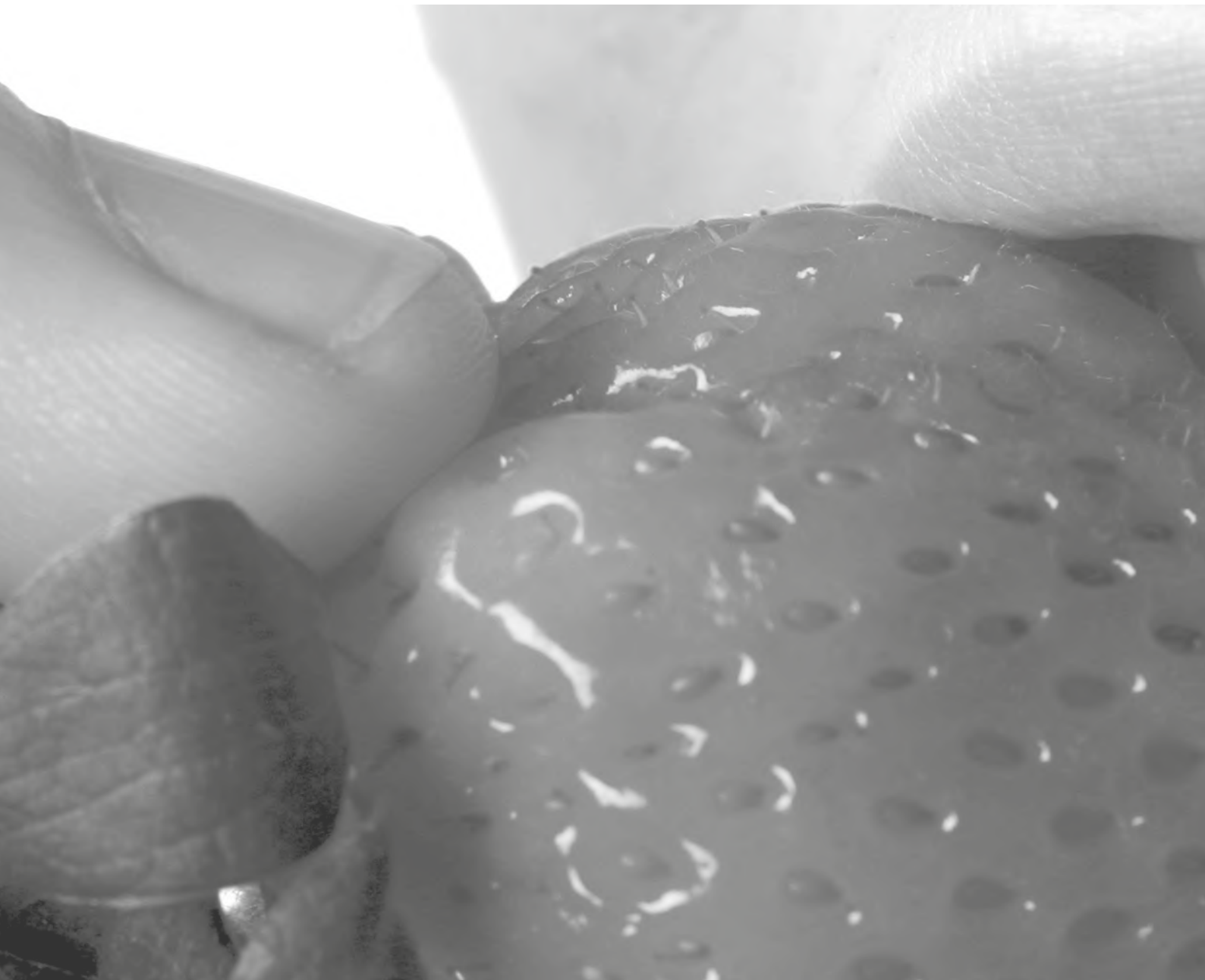



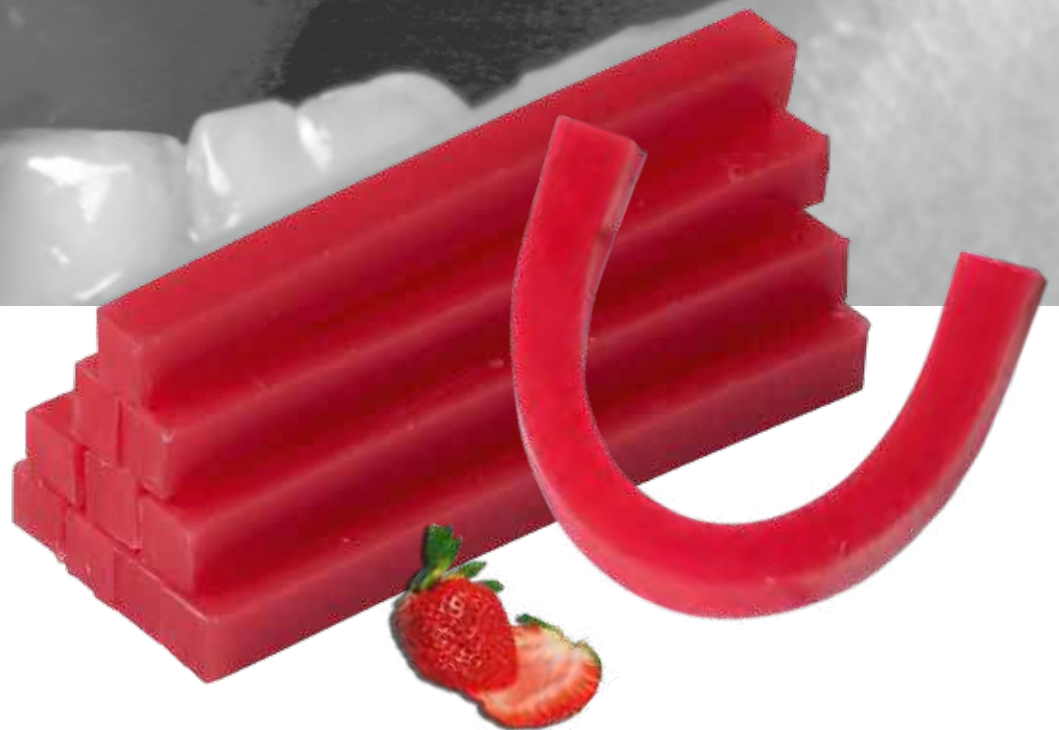

- **Modelling waxes (sheets)**
- **Bite registration waxes**

## Modelling waxes (sheets)

### BASE PLATE wax pink for full and partial denture

Our pink base plate waxes are available in four different qualities. You can choose between summer and winter quality 1,5 or 1,25 mm thickness, a higher or a lower solidification point, and between hard, soft or more elastic wax. All our base plate waxes are made of best raw materials to guarantee constant excellent quality.



#### SOLIDUS 80 & SOLIDUS 84

- both qualities have great stability which ensures they do not pull out of shape, even when fitting inside the mouth
  - they are extremely hard
  - very high solidification point
  - very good modelling features despite of their hardness
  - very precise
  - SOLIDUS is suitable also as bite wax
- A base plate wax for experts.

Differences:

- SOLIDUS 84 is harder - solidification point at 84 - 85°C
- SOLIDUS 80 is slightly softer - solidification point at 80 - 81°C

Size ca. 170 x 85 mm

SOLIDUS 80	1,5 mm	medium hard	500 g	Art.-No. 680-0000
SOLIDUS 84	1,5 mm	hard	500 g	Art.-No. 684-0000

Sculpturing  
SUPERPINK

Sculpturing wax  
SPEZIAL

IDEAL

SOLIDUS 80

SOLIDUS 84



softer



harder



### Sculpturing wax SPEZIAL



- pink base plate wax medium hard
- the "allrounder" among the base plate waxes
- summer quality with a solidification point of 56°C and winter quality with a solidification point of 54 °C
- size of the single plate: 175 x 80 mm
- sheet thickness: 1,5 mm or 1,25 mm
- package contents: 500 g or 2500 g
- our best-seller!

#### Sculpturing wax SPEZIAL

Summer	1,25 mm	500 g	Art.-No. 650-0500
Summer	1,25 mm	2500 g	Art.-No. 650-2500
Winter	1,25 mm	500 g	Art.-No. 652-0500
Winter	1,25 mm	2500 g	Art.-No. 652-2500
Summer	1,5 mm	500 g	Art.-No. 653-0500
Summer	1,5 mm	2500 g	Art.-No. 653-2500
Winter	1,5 mm	500 g	Art.-No. 655-0500
Winter	1,5 mm	2500 g	Art.-No. 655-2500

### Sculpturing wax SUPERPINK



- precise reproduction
- very elastic
- great stability
- solidification point: 57°C
- plate size: 175 x 80 mm

#### Sculpturing wax SUPERPINK

Sculpturing wax SUPERPINK	1,5 mm	500 g	Art.-No. 660-0500
Sculpturing wax SUPERPINK	1,5 mm	2500 g	Art.-No. 660-2500

### IDEAL



- natural pigmentation in gingival colour
- high breaking resistance of bending performance
- high accuracy and stability when fitting in the patient's mouth
- Size: 175 x 80 mm
- Solidification point: 58°C

#### IDEAL

IDEAL	1,5 mm	500 g	Art.-No. 651-0500
IDEAL	1,5 mm	2500 g	Art.-No. 651-2500

## Bite registration waxes




### **BITE REGISTRATION WAX flavoured**

Bite registration wax is very hard. It is very stable even at higher temperatures. The added orange flavour serves to relax the patient during the registration procedure.

Package content: 75 pcs.

Bite registration wax pink 2,0 mm	hard	Art.-No. 606-2000
--------------------------------------	------	-------------------

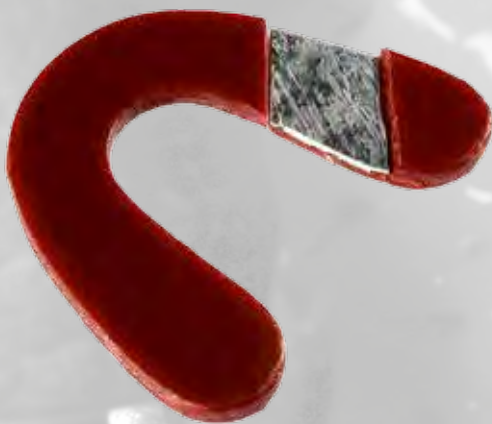

### **METAL FOIL WAX brown**

- for bite registration
- does not pull-out of shape as it is protected between two wax layers
- stable in shape

Package content: 60 pcs.

Metal foil wax brown 4,0 mm		Art.-No. 605-0060
--------------------------------	--	-------------------



### **BITE RIMS preformed non-flavoured**

Very useful bite rims for the warmer seasons of the year

- extrahard quality
- elevated solidification point


Package content: 100 pcs.

Bite rims preformed pink	extra hard	Art.-No. 614-0100
-----------------------------	------------	-------------------



## BITE WAXES flavoured sticks

Our bite waxes stay stable when warm and when cooled down, giving excellent bite registration. They guarantee the exact bite registration and are available in three different hardnesses. According to the clinical researches it is easier for a patient to accept a flavoured bite wax. This is the reason why we added popular flavours to our bite waxes such as strawberry, apple, orange and lemon. Your patient will feel comfortable and relaxed when using these bite waxes.

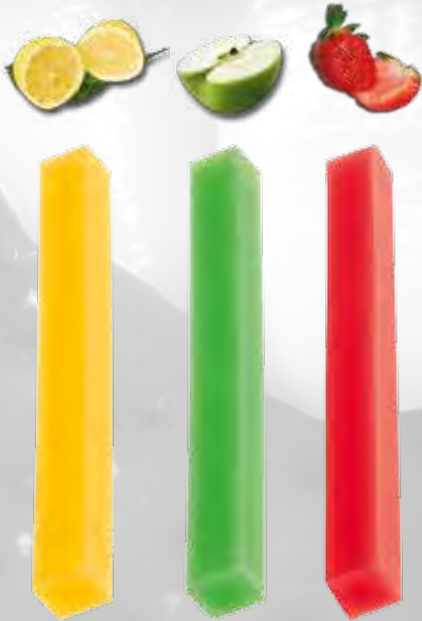

### BITE WAXES flavoured sticks

- size of one stick 108 x 10 x 10 mm
- colours correspond to their hardness:

yellow	=	lemon flavour	=	soft
green	=	apple flavour	=	medium
red	=	strawberry flavour	=	hard

Bite waxes				
yellow lemon flavoured	soft	490 g		Art.-No. 612-0490
yellow lemon flavoured	soft	1850 g		Art.-No. 612-1850
green apple flavoured	medium	490 g		Art.-No. 614-0490
green apple flavoured	medium	1850 g		Art.-No. 614-1850
pink strawberry flavoured	hard	490 g		Art.-No. 613-0490
pink strawberry flavoured	hard	1850 g		Art.-No. 613-1850



### BITE RIMS preformed flavoured

Use preformed bite rims to make your work easier. They are preformed in accordance to the jaw proportions. Available in three different hardnesses as the bite wax sticks.

Bite rims preformed				
yellow lemon flavoured	soft	96 pcs.		Art.-No. 610-0100
green apple flavoured	medium	96 pcs.		Art.-No. 612-0100
orange orange flavoured	medium	96 pcs.		Art.-No. 613-0100
red strawberry flavoured	hard	96 pcs.		Art.-No. 611-0100







- **digital technique  
and accessories**

## Scanner accessories



### Digiscan-Spray

Digiscan-Spray is an indispensable addition for the modern laboratory.

Ultra fine spray layer makes the correction after scanning much easier.

Digiscan-Spray can be applied onto silicone, diestone, wax or alloy.

Economic 300ml spray.

A distance of ca. 30 cm must be maintained between the spray and the object.

- eliminates reflections and reduce digital artifacts
- very fine spray
- one layer can be applied onto another
- smooth surface
- significantly better scan results
- easy to remove without leaving any residues

Digiscan-Spray	300 ml	Art.-No. 581-0300
----------------	--------	-------------------



### CAD / CAM Scan wax, beige

Laser-opaque and light-opaque sculpturing wax for CAD / CAM technique.

- guarantees excellent scanning results and exact fitting
- highly opaque wax for optimal scan data
- excellent as block out wax
- superb modeling qualities

CAD / CAM Scanwachs beige	45 g	Art.-No. 729-5000
------------------------------	------	-------------------




## exocad

Authorized distributor for leading dental CAD-Software (OEM)

Already the standard version of this software covers a wide variety of indications, making it an attractive choice both for expert users and newcomers in the dental CAD field.

The module system of exocad DentalCAD standard version is easily extensible and can be quickly adapted to individual needs. This software is the optimal solution for every user in accordance with personal requirements and budget.

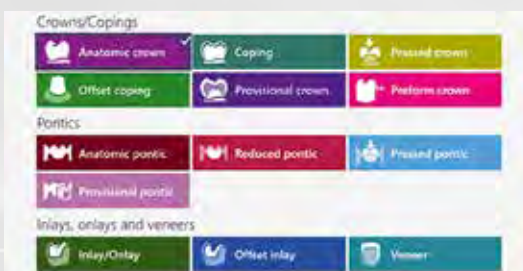

- easy to use even for newcomer
- modular system for highest flexibility
- variety of modules with additional functions
- purchase only what you need

Exocad basis software Further modules on request	Art.-No. 470-3200
---	-------------------

Monet Further modules on request





also as vinyl high resolution (HR) and vinyl open air available

SCAN BOX	Features	Vinyl 3D-Scanner
v	dental Scan	v
-	articulator scan	v
expandable	multiDie / multiCase	v
v	secondDie	v
-	Texture scan color	v
v	Texture scan monochrome	v
expandable	Triple Tray® Impression scan	v
-	Fully automated Z-axis	v
-	Touch display	v
v	Compatibility with exocad® DentalCAD	v

### Technical data

Dimension (WxDxH) 360 x 310 x 390 mm  
 Measurement field (X x Y x Z) 80 x 60 x 85 mm  
 Resolution 1,3 MP  
 Accuracy (ISO 12836) 6 µm  
 Weight 11 kg  
 Power Supply Voltage 100 - 240 VAC, 50/60 Hz  
 Interfaces 2x USB  
 Housing color white/silver  
 Scope of delivery: 1 SCAN BOX PRO, 1 object holder, 1x 3D calibration model, 2x USB cable, 1x mains cable, 1x data carrier

### Measurement Time

Complete Jaw  
 Scanning 18 sec  
 Postprocessing 14 sec  
 Total 32 sec  
 Single Tooth  
 Scanning 29 sec  
 Postprocessing 10 sec  
 Total 39 sec  
 3-unit bridge  
 Scanning 40 sec  
 Postprocessing 23 sec  
 Total 63 sec

### Technical data

Dimension (WxDxH) 455 x 430 x 420 mm  
 Measurement field (X x Y x Z) 80 x 60 x 85 mm  
 Resolution 1,3 MP  
 Accuracy (ISO 12836) 6 µm  
 Weight 23 kg  
 Power Supply Voltage 100 – 240 VAC, 50/60 Hz  
 Interfaces 1x USB  
 Housing color black / white  
 Scope of delivery: 1x Vinyl, 1x Object holder, 1x multiDie adapter, 1x Triple Tray® impression holder, 1x 3D calibrations model, 1x USB cable, 1x mains cable, 1x CD

### Measurement Time

Complete Jaw  
 Scanning 16 sec  
 Postprocessing 13 sec  
 Total 29 sec  
 Single Tooth  
 Scanning 25 sec  
 Postprocessing 9 sec  
 Total 34 sec  
 3-unit bridge  
 Scanning 33 sec  
 Postprocessing 22 sec  
 Total 55 sec

SCAN BOX

Art.-No. 470-8100

Vinyl 3D-Scanner

Art.-No. 470-8200

## Milling machines



also available as x and + variant  
**digi MILL<sup>5</sup> Dry milling machine**

- 5 simultaneously functioning axes
  - air pressure (min./max.): 7 bar/ 8 bar (recommended 7 bar)
  - automatic tool changer for 16 tools
  - user friendly CAM software
  - all open STL-files can be used
- Rotation range: 360° (A) | ± 35° (B)
- for blanks Ø 98,5 mm with a height of 10 - 40 mm
  - synchronous spindle rotational speed till 60.000 U/min
  - maximum power output (Pmax): 500 W
  - repetition accuracy ± 0,003 mm

- Working chamber illumination with color status indication

Plastics, wax, zirconium oxide, composites, CoCr, model plaster

### Technical Data

Dimensions (W/D/H): 450 x 530 x 630 mm

Weight: approx. 91 kg

digi MILL<sup>5</sup> Dry milling machine Art.-No. 481-5100



also available as + variant  
**digi MILL<sup>4</sup> wet - wet grinding machine**

- 4-axes wet grinding
  - Compressed air supply 4 bar · 35 l/min – 8 bar · 50 l/min
  - Automatic tool changer for 8 tools with 2 removable changer stations for easy equipment
  - User-friendly CAM-Software
  - All open STL-files readable
  - Rotation range: +190° to –10°
  - For 3 blocks with round shanks (total positioning range for workpieces: 60 x 40 x 20 mm)
  - Synchronous spindle with rotational speed range up to 60,000 RPM
  - Maximum power output (Pmax): 600 W
  - Repetition accuracy linear axes ± 0,003 mm
  - Machine base made of a solid cast aluminium body for maximum stability
  - 8 liquid nozzles on the spindle
  - Workspace lightning for status indication
  - No air extraction necessary
- Block materials such as glass ceramic, zirconium oxide, composites, titanium etc

### Technical Data


Dimensions (WxDxH) ca. 360 x 451 x 471 mm

Weight ca. 50 kg

digi MILL<sup>4</sup> wet - wet grinding machine Art.-No. 481-4600



## Milling machines

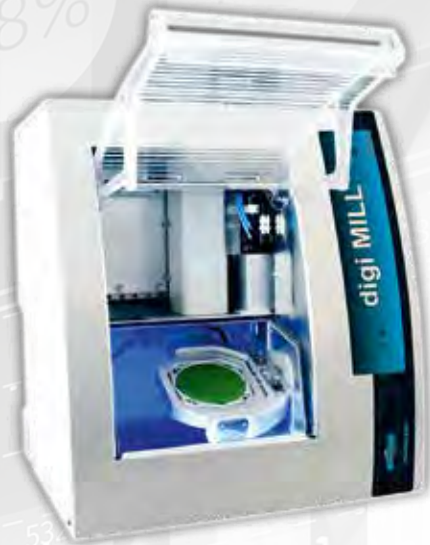

**digi MILL<sup>5</sup> milling machine**  
(without blank changer)

also available as x and + variant



**digi MILL<sup>5</sup> plus milling machine**  
(with blank changer for 8 blanks)

- 5 simultaneously working axes
- air pressure (min./max.): 7 bar/ 8 bar (recommended 7 bar)
- automatic tool changer for 16 tools
- user friendly CAM software
- all open STL-files can be used
- Rotation range: 360 ° (A) | ±30 ° (B)
- for blanks Ø 98,5 mm with a height of 10 - 30 mm
- synchronous spindle rotational speed till bis 60.000 U/min
- maximum power output (Pmax): 600 W
- Repetition accuracy ± 0,003 mm
- ready to be used with wet milling option
- reduced vibration due to stable construction

Plastic, wax, zirconium oxide, nano-composites, (glass ceramic\*), CoCr, (Titan\*)

digi MILL<sup>5</sup> milling machine  
without blank changer  
Dimensions (W/D/H): 490 x 445 x 540 mm  
Weight: ca. 75 kg

digi MILL<sup>5</sup> plus milling machine  
with blank changer for 8 blanks  
Dimensions (W/D/H): 692 x 445 x 540 mm  
Weight: ca. 95 kg

digi MILL<sup>5</sup> milling machine

Art.-No. 481-5000

digi MILL<sup>5</sup> plus milling machine

Art.-No. 481-5500

## Cutters digi MILL5



Universal Tools						
Ø Cutting edge	Form of cutting edge	Teeth	Length	Coating	Art.-No.	
0,30 mm	Radius	2	40 mm	no	481-0403	
0,50 mm	Flat	2	40 mm	no	481-0405	
0,60 mm	Radius	2	40 mm	yes	481-0506	
1,20 mm	Flat	2	40 mm	yes	481-0530	

Wax and Plastic (PMMA)						
Ø Cutting edge	Form of cutting edge	Teeth	Length	Coating	Art.-No.	
1,00 mm	Radius	1	40 mm	no	481-0611	
2,00 mm	Radius	1	40 mm	no	481-0612	
1,00 mm	Radius	2	40 mm	no	481-0610	
2,00 mm	Radius	2	40 mm	no	481-0620	
2,50 mm	Flat	1	40 mm	no	481-0613	

Zirkonium Oxide (ZrO <sup>2</sup> )						
Ø Cutting edge	Form of cutting edge	Teeth	Length	Coating	Art.-No.	
1,00 mm	Radius	2	40 mm	yes	481-0410	
2,00 mm	Radius	3	40 mm	yes	481-0420	

Zirkonoxid (ZrO <sup>2</sup> )						
Ø Cutting edge	Form of cutting edge	Teeth	Length	Coating	Art.-No.	
0,60 mm	Radius	2	40 mm	diamond	481-0306	
1,00 mm	Radius	2	40 mm	diamond	481-0310	
1,20 mm	Flat	2	40 mm	diamond	481-0412	
2,00 mm	Radius	3	40 mm	diamond	481-0320	

Composites						
Ø Cutting edge	Form of cutting edge	Teeth	Length	Coating	Art.-No.	
1,00 mm	Radius	2	40 mm	yes	481-0510	
2,00 mm	Radius	2	40 mm	yes	481-0520	
1,00 mm	Radius	1	40 mm	diamond	481-0509	
2,00 mm	Radius	1	40 mm	diamond	481-0511	



Non-precious alloys (CoCr) and Titanium						
Ø Cutting edge	Form of cutting edge	Teeth	Length	Coating	Art.-No.	
0,60 mm	Radius	2	35 mm	yes	481-0706	
1,00 mm	Radius	2	35 mm	yes	481-0710	
1,20 mm	Torus	2	35 mm	yes	481-0730	
2,00 mm	Radius	2	35 mm	yes	481-0720	
2,00 mm	Radius	4	35 mm	yes	481-0721	

Glass Ceramics (LiSi <sup>2</sup> )						
Ø Cutting edge	Form of cutting edge	Teeth	Length	Coating	Art.-No.	
0,60 mm	Radius	0	35 mm	no	481-1206	
0,60 mm	Torus	0	35 mm	no	481-1206-1	
1,00 mm	Radius	0	35 mm	no	481-1210	
1,20 mm	Torus	0	35 mm	no	481-1212	
2,40 mm	Radius	0	35 mm	no	481-1211	
2,60 mm	Torus	0	35 mm	no	481-1226	



	digi MILL <sup>4</sup> wet	digi MILL <sup>5</sup> Dry	digi MILL <sup>5</sup>	digi MILL <sup>5</sup> plus
		<b>also available as x and + variant</b>		
Machinable materials	Block materials such as glass ceramic, zirconium oxide, composites, titanium etc	Plastics, wax, zirconium oxide, nano-composites, (glass ceramic*), CoCr, (Titanium*) (*only digi MILL <sup>5</sup> and digi MILL <sup>5</sup> plus with wet milling option)		
Dimensions (W/D/H)	ca. 360 x 451 x 471 mm	ca. 450 x 530 x 630 mm	490 x 445 x 540 mm	692 x 445 x 540 mm
Weight	ca. 50 kg	ca. 91 kg	ca. 75 kg	ca. 95 kg
Rotary Axes	4-axes wet grinding	5 axes	5 axes	5 axes
	A axis: highest true running accuracy rotation range: +190° to -10°	· Harmonic Drive® gears free from backlash		
		· Rotary axis with highest true running accuracy		
		A Axis: 360 degree		
	B axis: rotation angle of ± 35 degrees	B axis: rotation angle of ± 30 degrees		
Axis measurement	Automatic axis measurement and compensation for absolutely exact results			
Repetition accuracy	linear axes ± 0,003 mm	± 0,003 mm	± 0,003 mm	
Fixing device	For 3 blocks with round shanks (total positioning range for workpieces: 60 x 40 x 20 mm)	Round universal blanks with a thickness of 10 to 40 mm (Can be processed: wax max 25 / 40 mm, alloy max. 15mm, zirconium max. 25 mm)		
		Diameter: 98,5 mm with step		
		No	For wet grinding option specific three-fold block holders available	
		No	· External wet grinding module (optional)	
Wet grinding option (*)	- 8 liquid nozzles on the spindle - removable drawer for cooling liquid and filter. Convenient to clean and fill - protection of mechanics, electronics and spindle by flexible rubber gaiter	No	3 liquid nozzles at the spindle for wet grinding	
		No	· External wet grinding module (optional)	
Exchange of workpieces	manual exchange	Manual exchange of blanks		Automatic blank changer for eight blanks
Collet chuck	pneumatic stainless steel collet chuck for tools with 3 mm shank diameter and max. 35 mm total length	For tools with 3 mm shank diameter and max. 40 mm total length		
Tool Change	automatic tool changer for 8 tools with 2 removable changer stations for easy equipment	Automatic tool change station for 16 tools		
	haptic tool length detection and tool breakage monitoring	· Tool length detection and tool breakage monitoring via pressure key		
	· Diamond-coated tools can be used			
	compressed air monitoring	· Compressed air monitoring for automatical tool changing		
Spindle	Synchronous spindle			
	· Nominal power under constant load (S1): 300 Watt			
	· Maximum power output (Pmax):			
	600 Watt	500 Watt	600 Watt	
	· 4-fold layering, · Hybrid ceramic ball bearings			
	· Rotational speed range up to 60,000 RPM			
	· Radial deviation at internal cone of the precision shaft < 1 µm			
· Blocking air prevents entering of foreign substances in the bearing area · Cone cleaning				
Air Extraction	no extraction necessary	· Opening at the housing side for air extraction		
		· Underpressure sensor for monitoring the air extraction		
		· 24 Volt output for connecting the switching unit		
Accessories		· Switching unit PSW 01-RSV for automatically switching a vacuum cleaner		
		· YETI digi MILL tool board without tools		
		No	· Wet grinding option	
		No	Block holders for glass ceramics and titanium blocks	
Art.-No.	481-4600	481-5100	481-5000	481-5500

## CAD / CAM Blanks

CAD / CAM Blanks							
Strength	Content	Wax					
		grey	beige	white	blue	green	
		normal hard (hardness 2)	hard (hardness 3)	hard (hardness 4)	hard (hardness 4)	elastic (hardness 1)	
		slightly elastic microwax					Model casting
Diameter in mm		Ø 98,5	Ø 95,0	Ø 98,5	Ø 98,5	Ø 98,5	Ø 98,5
Article number		791-....	791-....	792-....	793-....	794-....	795-....
14,0 mm	1 pcs.	....0014		....0014	....0014	....0014	
	12 pcs.				....0114	....0114	
16,0 mm	1 pcs.		....0116				
20,0 mm	1 pcs.	....0020	....0120	....0020	....0020	....0020	....0020
	12 pcs.				....0120	....0120	
25,0 mm	1 pcs.				....0025	....0025	
30,0 mm	1 pcs.				....0030	....0030	
							
Hardness 1 bis 4 = 1 elastic hard to 4 very hard							

CAD / CAM Blanks PMMA			
Strength	Content	PMMA	
		one colour	layered
		(colour indication follows article number)	
Diameter in mm		Ø 98,5	Ø 98,5
Article number		398-....	
18,0 mm	1 pcs.		....4118 (A1)
	1 pcs.		....4218 (A2)
	1 pcs.		....4318 (A3)
20,0 mm	1 pcs.	....0220 (A2)	
	1 pcs.	....0320 (A3)	
	1 pcs.	....1120 (B1)	
	1 pcs.	....8020 (transparent)	
25,0 mm	1 pcs.	....0225 (A2)	
	1 pcs.	....0325 (A3)	
	1 pcs.	....1125 (B1)	

CAD / CAM Blanks Metal		
Strength	Content	Metal
Diameter in mm		Ø 98,5
		Solibond C plus Blank
Article number		969-....
10,0 mm	1 pcs.	....9010
12,0 mm	1 pcs.	....9012
13,5 mm	1 pcs.	....9013
15,0 mm	1 pcs.	....9015
		CoCr sintered alloy blanks
Article number		969-....
10,0 mm	1 pcs.	....9110
12,0 mm	1 pcs.	....9112
14,0 mm	1 pcs.	....9114
16,0 mm	1 pcs.	....9116
18,0 mm	1 pcs.	....9118
20,0 mm	1 pcs.	....9120



## CAD / CAM Blanks



### CAD / CAM wax blanks

- Stable wax composition - allows the milling cutter to work on the narrowest space
- Homogenous surfaces – mirror-smooth milled surfaces
- Wax chips are not melting in the working area and the exhaust system remains efficient
- Produced from residue-free, temperature stable and investment compatible micro-waxes
- Quality - Made in Germany -



### CAD / CAM wax blanks

#### For partial dentures

- Stable in form and volumes
- Breakresistant
- Sharp contours



### CAD / CAM PMMA blanks In one color, transparent or layered

- For production of temporary crowns and bridges
- High break strength and good processing characteristics
- Very natural and aesthetic - layering is modeled on the natural characteristics (layered blanks)
- Available in different colors and thicknesses



### Solibond C plus blank alloy (non precious)

#### CoCr

- For ceramics
- Biocompatible
- Laser-capable
- Quick polishing and surface gloss
- Optimal bonding with ceramics

Vickers hardness: 330 HV  
Elongation: 6,5 %  
CTE (20 - 600°C) 14,4 · 10<sup>-6</sup>K<sup>-1</sup>

Co %	Cr %	Mo %	Si %
63	29	5,8	1,2



### CoCr sintered alloy blanks (non precious)

- longer working time of cutters
- minimum time for finishing
- perfect adhesion with ceramics
- can be furnace-fired, is biocompatible and laser-capable

Vickers hardness: 285 HV  
Elongation: 10,0 %  
WAK (20 - 600°C) 14,1 · 10<sup>-6</sup>K<sup>-1</sup>

Co %	Cr %	Mo %	Mn %	Fe %	C %	Ni %
Balance	27,5-29,5	5,5-6,5	<0,1	<0,1	<0,1	<0,1

# CAD / CAM Blanks



Ø 98,5 mm

Ø 95,0 mm

Pre-sintered and yttrium-stabilized zircon blanks for crown and bridge technique.

## K2 zircon blanks white classic

- opaque framework material
- perfect for covering of colored tooth stump as well as post & core

- suitable for large-span bridges till 14 units but with max. 2 pontics
- particular good biological compatibility
- homogenous density and very fine grain for smooth milling results
- high edge stability

Density	6.08 g/cm <sup>3</sup>
Flexural strength (biaxial)	>1100 MPa
Grain size	<0.4 µm
Breaking strength	5 MPavm
Vickers hardness HV10	1250
Translucency (1mm, D65)	35 %
CTE	10 10 <sup>-6</sup> /K
Thermal conductivity	2.5 W/mK
Dissolubility	≤50 µg/cm <sup>2</sup>

## K2 zircon blanks translucent

- natural-looking aesthetic results
- reduction of brightness level due to high translucency
- low light reflection
- natural adaptation to another teeth

Density	6.08 g/cm <sup>3</sup>
Flexural strength (biaxial)	>1100 MPa
Grain size	<0.4 µm
Breaking strength	5 MPavm
Vickers hardness HV10	1250
Translucency (1mm, D65)	41 %
CTE	10 10 <sup>-6</sup> /K
Thermal conductivity	2.5 W/mK
Dissolubility	≤50 µg/cm <sup>2</sup>

## K2 zircon blanks extreme translucent

- extremely highly translucent zirconium oxide
- perfectly suitable for monolithic restorations, particularly in the anterior zone
- the translucency is the same as of glass ceramic but features the high strength of zircon
- make highly aesthetic three-unit bridges possible in the anterior and posterior area
- excellent alternative to lithium disilicate, also from an economic point of view

Density	6.04 g/cm <sup>3</sup>
Flexural strength (biaxial)	>610 MPa
Grain size	<0.6 µm
Breaking strength	2.5 MPavm
Vickers hardness HV10	1250
Translucency (1mm, D65)	49 %
CTE	10 10 <sup>-6</sup> /K
Thermal conductivity	2.5 W/mK
Dissolubility	≤50 µg/cm <sup>2</sup>

Recommendation: DIP & BRUSH K2 zircon liquid for colorisation of milled works in suitable dental color

Strength	Content	Zircon					
		white classic		translucent		extreme translucent	
Diameter in mm		Ø 98,5	Ø 95,0	Ø 98,5	Ø 95,0	Ø 98,5	Ø 95,0
Article number		<b>396-....</b>	<b>396-....</b>	<b>397-....</b>	<b>397-....</b>	<b>397-....</b>	<b>397-....</b>
10,0 mm	1 St.	....0010	....1010	....0010	....1010	....2010	....3010
12,0 mm	1 St.	....0012	....1012	....0012	....1012	....2012	....3012
14,0 mm	1 St.	....0014	....1014	....0014	....1014	....2014	....3014
16,0 mm	1 St.	....0016	....1016	....0016	....1016	....2016	....3016
18,0 mm	1 St.	....0018	....1018	....0018	....1018	....2018	....3018
20,0 mm	1 St.	....0020	....1020	....0020	....1020	....2020	....3020
25,0 mm	1 St.	....0025	....1025	....0025	....1025	....2025	....3025





Ø 98,5 mm

Geo5 - K2 CAD / CAM Zirconia


- 5 colour gradient Layer
- 5 strength Layer
- 5 translucency Layer



before



finished glazed work

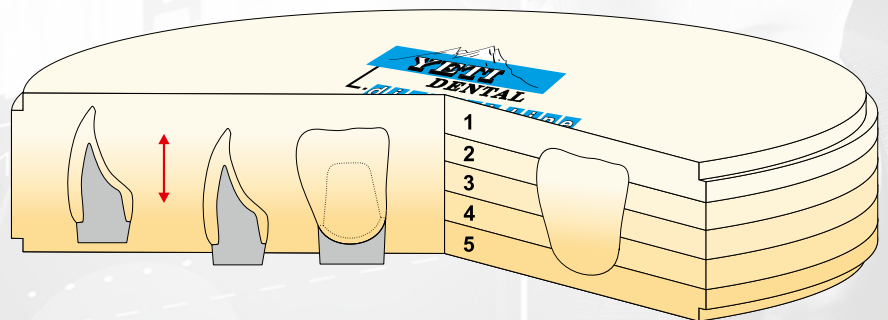

layer	height	translucency	strength
1	23%	47%	≥ 800
2	18%	46,5%	≥ 850
3	18%	46%	≥ 900
4	18%	45%	≥ 1000
5	23%	44%	≥ 1100

Geo5 K2 CAD / CAM Zirconia					
colour	Art-No.	height			
		14mm	16mm	18mm	20mm
A1	396-90..	..14	..16	..18	..20
A2	396-91..	..14	..16	..18	..20
A3	396-92..	..14	..16	..18	..20
A3,5	396-93..	..14	..16	..18	..20
B1	396-94..	..14	..16	..18	..20
B2	396-95..	..14	..16	..18	..20
C2	396-96..	..14	..16	..18	..20

# DIP & BRUSH



CARAMISTA-PURE  
excellence line  
Brush-SET (02/04) Art.-Nr. 314-5000





## DIP & BRUSH K2 CAD / CAM Zircon liquid


With K2 Dip & Brush it is possible to individualize parts of frameworks (Brush) or to colorize whole frameworks by dipping (Dip) before sinterisation.

Water based colours with 100 ml volume are available in 7 dentin and 2 enamel colours.





### Dipping technique:

Dentin colorization of frameworks by dipping them is applicable mainly in layering technique. It allows the reduction of zircon whiteness. Thus, less light is reflected and it leads to outstanding aesthetic results particularly in the thin areas.

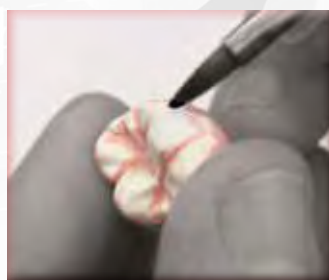

### Painting technique:

Painting technique (Brush) is suitable mainly for monolithic restorations which can be individualized with different dentin and enamel colours before sinterization. It makes the finalization with stains and glaze powder easier after sinterization.

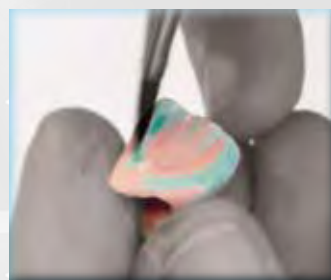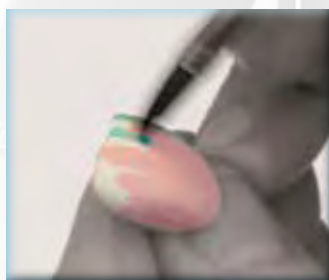

DIP & BRUSH  
K2 CAD / CAM Zircon liquid, 100 ml

A1 = Art.-No. 397-5001  
A2 = Art.-No. 397-5002  
A3 = Art.-No. 397-5003  
A3,5 = Art.-No. 397-5004  
B2 = Art.-No. 397-5007

C2 = Art.-No. 397-5011  
D2 = Art.-No. 397-5014  
BL1 = Art.-No. 397-5017  
GR3 = Art.-No. 397-5019

## Compressor / Sintering furnace

### YETI piston compressor

Fully automatic, air-cooled, quiet running ready to use compressor with electric motor which is installed on 20-liter pressure vessel.

#### Delivery contents

Basic configuration of dry-running compressors is a ready to operate system which allows installation without a special foundation for sound-insulated units.

Complete with safety equipment, pressure control, after cooler, condensate drain valve, manometer, individually adjustable pressure reducer, ventilation for pressureless start as well as motor safety switch.



Technical data	
Construction type	single-stage
Number of cylinders	2
Suction capacity	200 l/min
Effective air delivery	130 l/min
Maximum pressure	8 Bar
Compressor speed	1450 U/min
Volume of the container	20 l
Motor	1,5 kW
Voltage	230 V
Frequenz	50 Hz
Protection class	IP 44
Isolation classe	F
Sound level	55 dB(A)
Weight	110 kg

YETI piston compressor

Art.-No. 484-1000



### Sintering furnace Tabeo for Metal

30 single crowns or sintering discs

- Molybdenum-disilicide heating elements
- Sintering system, Ø 120mmx30mm; METAL (kit)
- Program control:
  - 7-Segment LED-
  - 4 stages
  - 9 program capacity
  - 4 fixed programs
- Shielding-gas supply with manual adjustment
- Timer function

Technical data	
Dimensions WxDxH	530 x 680 x 460mm
Heating-chamber height	92 mm
Heating-chamber capacity	1 tray Ø 120x 30 mm
Maximum temperature	1400°C
heating elements	4x MoSi2
Power max.	1600 W
Voltage range; frequency	220-240 V; 50/60 Hz
Maximum heating rate	40°C/min
Weight	80 kg



Sintering furnace Tabeo for Metal

Art.-No. 483-3000



## Sintering furnace

### Speed Sintering furnace for Zirconia

60 SINGLE CROWNS or SINTERING DISC



- High-quality molybdenum-disilicide heating elements PRO
- Program control
- LCD-Display
- 4 stages
- 30 program capacity
- Special functions
- Service programs
  - A-Temperature control (only in conjunction with test-kit)
  - C-Purge heating chamber
  - E- Regenerate heating elements
- Thermocouple PtRh-Pt, type S
- Temperature accuracy at 1500°C +/- 3°C
- Timer function

#### Technical information

Dimensions W/D/H)	390 x 790 x 500 mm
Max. heating-chamber capacity	2x Ø 120/30 mm
Maximum temperature	1650°C
Heating elements	4 x MoSi2
Power max.	3200 W
Voltage range / frequency	220 - 240 V, 50 / 60 Hz
heat-up rate up to	99 °C/min
Weight	60 kg

Speed Sintering furnace  
for Zirconia

Art.-No. 483-4000



### Sintering furnace Tabeo for Zirconium Oxide

- Molybdenum-disilicide heating elements
- Service programs
  - A-Temperature control\* (only in conjunction with test kit)
  - C-Purge heating chamber
  - E- Regenerate heating elements
- Program control
  - 7-Segment LED
  - 4 stages
  - 9 program capacity
  - Special functions
- Temperature accuracy at 1500°C +/- 3°C
- Timer function

#### Technical data

Dimensions WxDxH	400 x 600 x 400 mm	Heating-chamber height	42 mm
Heating-chamber capacity	1 tray Ø 100x 30 mm	Maximum temperature	1650°C
Maximum temperature	1650°C	Heating elements	4x MoSi2
Heating elements	4x MoSi2	Power max.	1500 W
Power max.	1500 W	Voltage range; frequency	220-240 V; 50/60 Hz
Voltage range; frequency	220-240 V; 50/60 Hz	Maximum heating rate	25°C/min
Maximum heating rate	25°C/min	Weight	55 kg
Weight	55 kg		



Sintering furnace Tabeo  
for Zirconium Oxide

Art.-No. 483-1000

# CAD/CAM Suction System




## CAD/CAM Suction System Z1 CAM

Single Suction for CAM milling machines

- Z1 CAM with guaranteed motor lifetime of 1500 hours
- Z1 CAM for CAM milling machines
- Integrated Bluetooth module and Zubler suction system app for easily understandable operation and individual setting options
- Designed with filter bag or permanent filter system according to the field of application
- Replacement of filter and motor without tools

Suction Power	1200W
Guaranteed motor lifetime (time is achievable by replacing the carbon brushes)	5000 h
Type of filter system (dust collection)	Filter bag
Interface for external control unit gate input	v
Interface for milling machines and devices with 24V output /PLC	v
Bluetooth interface for individual settings per Android-APP	v
Adapter for exhaust air to the outside (optional)	v
Exhaust air filter grade H14 for dust category H (optional)	v
Activated charcoal exhaust air filter (optional)	v
Quick exchange of the motor (no tools required)	v
Front color/material selectable (optional)	v



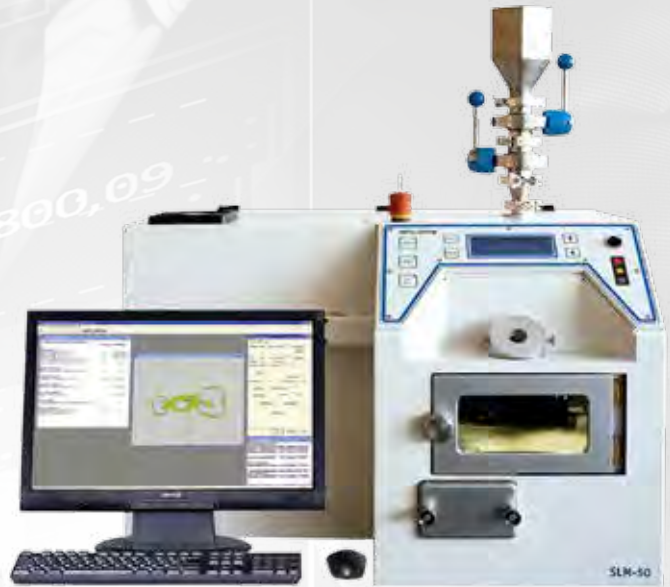

## SOLIBOND C plus Powder

Cobalt-Chrome alloy  
for laser based applications

An important component to SLM machines is the alloy powder, crucial for the quality of laser melted products. Solibond C plus powder with its spectrum of -10 / +63 my is produced in Germany with strict quality controlled production under inert gas. Each Lot is released to market, only after achieving the optimum test criteria. The popular alloy SOLIBOND C plus is registered worldwide and already enjoys a high level of customer confidence. SOLIBOND C plus is highly recommended for use with the SLM 50 machine and the powder stays constantly dry and free-flowing.

- exact fitting
- smooth surface of produced items
- fine and homogeneous alloy structure

Vickers hardness: 280 HV  
Elongation: 10,1 %  
CTE (20 - 600°C) 14,0 · 10<sup>-6</sup> K<sup>-1</sup>

Co %	Cr %	Mo %	W %	Nb %	Si %
63	24	2,9	8,1	0,9	1,1

## SLM 50

SLM 50 is the first SLM desktop machine for manufacturing components made of metal. The desktop device has been designed for the manufacturing of components with a diameter of up to 70 mm and a max. height of up to 80 mm. On the construction platform up to 30 individual components can be manufactured at the same time. In comparison to the conventional casting technique, the manufacturing with the SLM 50 becomes economical already starting at 100 units per month.


### Technical Data SLM 50

Dimensions	(H)500 x (W)800 x (D)700 mm
Construction volume	Ø70 mm, max. height 80 mm
Thickness of layer	20 µm - 50 µm
Laser type	Fibre laser 20W - 120W
Power connection	16 A, 230V
Power input	1,0 KW
Argon consumption	approx. 30 liter/hour.
Materials	CoCr, titanium, steel alloys etc.
Software	Realizer control software
Weight	ca. 100 kg

SOLIBOND C plus Powder 5000 g Art.-No. 969-5000

SLM 50

Art.-No. 460-1000

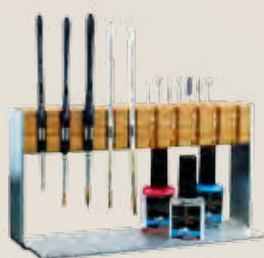

Article	Content	Article Number	
	Vario magnetic holder without instruments	1 pc. 312-9000	
	Vario S magnetic holder without instruments	1 pc. 312-9001	
	Zircon polishing set 1 HM Fissure bur 1 HM cutter 1 Smart Stone cone 1 Smart Stone flame 1 Arkansas 1 Abrasive paper cone (mounted on mandrell) 1 Rubber Flame	1 Set = 7 pcs. 440-0001	
Refill packages for zircon polishing set 	HM Fissure bur	5 pcs. 440-0300	
		HM cutter	1 pc. 440-0200
		Smart Stone cone	1 pc. 440-0400
		Smart Stone flame	1 pc. 440-0500
		Arkansas	5 pcs. 440-0600
		Abrasive paper cone (unmounted, without mandrell)	25 pcs. 440-0100
		Rubber Flame	5 pcs. 440-0700
	Aqua Spacer coral	1 x 18 ml 597-0001	
	Aqua Spacer ocean	1 x 18 ml 598-0001	
	Aqua Spacer stone	1 x 18 ml 599-0001	
	Aqua Sep Separating Liquid	1 x 18 ml 589-1001	
	Caramista pure Brush Set	2 brushes: 1 x no 2 1 x no 4 314-5000	
	K2 Lube	1 x 18 ml 345-0018	



## „excellence line“ - Product overview

Contact us to request the excellence line brochure!



# excellence line

	content / amount	item number
Contact base = base plate with wax scraper and dust cover	 <b>CONTACT freestyle</b> Wax dispenser Contact Basis	780-0100
	 <b>CONTACT traditional</b> Wax dispenser Contact Basis, sculpturing wax chip blue (ashfree) 25g, sculpturing wax chip grey (opaque) 25g, Cervical wax chip red (ashfree) 20g	780-0200
	 <b>CONTACT n. J. Peters</b> Wax dispenser Contact Basis, sculpturing wax chip ivory (ashfree) 25g, sculpturing wax chip sahara (ashfree) 25g, Cervical wax chip wood (ashfree) 20g	780-0300
	 <b>CONTACT esthetik</b> Wax dispenser Contact Basis, sculpturing wax chip blue transparent (ashfree) 20g, sculpturing wax chip orange transparent (ashfree) 20g, sculpturing wax chip sahara (ashfree) 25g	780-0400
	 <b>CONTACT NAT</b> Wax dispenser, 6 x 8g (ashfree) Contact Basis, sculpturing wax: blue, grey, neon-yellow, neon-green Underlining wax : caramel, Cervical wax: red	780-0500
	<b>Effect - Wachschiips CONTACT</b> (ashfree) magnetically adhesive blue transparent / orange transparent / incisal transparent each 20g	739-6320
	<b>Modellierwachschiips CONTACT sahara</b> (ashfree) magnetically adhesive by J. Peters, 3x25g	738-6325
	<b>Modellierwachschiips CONTACT ivory</b> (ashfree) magnetically adhesive by J. Peters, 3x25g	738-6335
	<b>Modellierwachschiips CONTACT beige</b> (opak) magnetically adhesive 3x25g	735-6325
	<b>Modellierwachschiips CONTACT grey</b> (opak) magnetically adhesive 3x25g	734-6325
	<b>Modellierwachschiips CONTACT emerald-green</b> (ashfree) magnetically adhesive 3x25g	737-6325
	<b>Modellierwachschiips CONTACT blue</b> (ashfree) magnetically adhesive 3x25g	736-6325
	<b>Modellierwachschiips CONTACT yellow</b> (ashfree) magnetically adhesive 3x25g	733-6325
	<b>Cervikalwachschiips CONTACT wood</b> (ashfree) magnetically adhesive by J. Peters, 3x20g	738-6320
	<b>Cervikalwachschiips CONTACT red</b> (ashfree) magnetically adhesive 3x20g	731-6320
	<b>IQ training set CONTACT</b> by G. Seubert (ashfree), 6 x 8g sculpturing wax: blue, grey, neon-yellow, neon-green Underlining wax : caramel, Cervical wax: red	102-0002M
	<b>phone holder CONTACT</b> magnetically adhesive 1 pc	780-0010
	<b>Wax scraper CONTACT</b> magnetically adhesive 6 pcs	780-0006

## Table of contents

<b>A</b>	Alloys	72 - 73	<b>C</b>	Cutex Bond	14
	Artline	56 - 57		Cutter	102
	Artline NAT	57	<b>D</b>	DETON wax wire	37
	Artline plus	57		DIA GLACE	89
<b>B</b>	Bite registration waxes	94 - 95		DIA QUICK	89
	Bite rims	95		Diagnostic wax	22
	Block out materials	14		Diamond paste	89
	Block out wax	14, 48		Diamond separating disc	86
	BLUE MARKER high spot indicator	71		Diamond-impregnated felt wheels	89
	Boxing wax, red	9		Die Bloc	14
	Brush for ceramic	85		Die Hardener	13
	Build Up Resin LC	65		Die Spacer	11 - 12
<b>C</b>	Cabinet for ceramic	88		digi MILL	100 - 103
	CAD / CAM Blanks	104 - 107		Digiscan-Spray	98
	CAD/CAM Suction System	112		DIP & BRUSH	108 - 109
	CAD / CAM scan wax	98		Dipping hardener	52
	CAD / CAM wax	104 - 105		Dipping wax	20 - 21
	CAD / CAM wax blanks	104 - 105		Dippy easy	63
	Caramista Ceramic brushes	85		Dippy NT digital	63
	Casting funnel	48		Dippy pro digital	62
	Casting sticks	36		dive monitoring device	62 - 63
	Casting wax sheets	46 - 47		DUO DIP	21
	Ceramic	76 - 83		duo Waxer NT	59
	Ceramic accessories	70, 84 - 89		DUOVEST investment	68
	Ceramic Humidor	84		Duplicating gel	53
	Ceramic mixing plate	84		Duplicating silicone	53
	Ceramic separating liquid	88		Duron wax wire	37
	Cervical wax	24, 25, 27, 29, 31	<b>E</b>	easy Waxer	59
	Clasp profiles	40, 41		easy Waxer duo	59
	Clear Spacer	10		ELASTO DIP	20
	COBAVEST investment	69		electric wax knife	58 - 61
	Color Spacer	10		electrical dive monitoring device	62 - 63
	Compressor	110		excellence line	114,115
	CONSEQUENT Block out wax	48		EXPANSION investment	68
	CONSEQUENT clasps	40		EXPANSION plus investment	69
	CONSEQUENT Universal wax	49	<b>F</b>	Felt wheels	89
	CONSEQUENT wax wire	37		Finale	34, 49
	CONTACT Wax dispenser	115	<b>G</b>	Glasstick	86
	Creation Set	22		Glaze paste	77
	Crown and bridge waxes	23 - 33, 49		Glue universal	48
	Cutex	14		Goldpowder for ceramic	87
				Gypsum hardener	13



<b>H</b>	Handpiece clamp	60	<b>M</b>	Milling wax	27, 29
	Handpiece for wax knife	60		Mixing plate ceramic	84
	Hardener	13		Model preparation	8-9
	High spot indicator	71		Model saw	8
	Humidity plate for stains	84		Modelling brush	85
	Humidor	84		Modelling instruments	56 - 57
<b>I</b>	Impression material light curing	65		Modelling light cure resin	65
	Investment	68, 69		Modelling tips	60, 61
	IQ Compact	30 - 31		Modelling wax sheets	91 - 93
	IQ Diagnostic wax dentin	22		Modelling waxes	23 - 33, 49
	IQ Quadro sticks	36		Monolith metal base & chip	29
	IQ Sculpturing wax	30 - 31		Monolith Thowax chip	29
	IQ Special wax	31		Multiflex saw blades	8
	IQ Sticks round	36		Multifunctional clasp	40
	IQ Triangel casting bars	36	<b>N</b>	Nawax compact	31
	IQ Waxer NT light	61	<b>O</b>	Occlusal powder	34
	ISO QUICK separating pen	15		OCCLUSAL SPRAY	71
	ISO STEAM	16 - 17	<b>P</b>	Periphery wax, red	9
<b>K</b>	K2 Ceramic	76 - 83		Plaster	8
	K2 Instruction for use	88		Plastic retentions	48
	K2 LF	76 - 79		Plunger	70
	K2 Lube ceramic separating liquid	88		PMMA blanks	104 - 105
	K2 Li / K2 universal (Liquid / stains)	80 - 83		Polishing paste	89
	K2 press	78, 81, 82		Polishing wheels	89
	K2 Zircon blanks	106, 107		Pontics	35
	K2 press universal	69		PRECI DIP	20
<b>L</b>	Lasermelting	112 - 113		PRECI NT SHUTTLE II / IV	64
	Light curing block out material	14		PRECI TRAY	65
	Light curing impression material	65		Precigel duplicating gel	53
	Light curing modelling resin	65		Precisil duplicating silicone	53
	Light curing spacer	13		Preformed pattern	35, 36
	Light curing unit	64		Preformed wax pattern	35, 36
	Light Spacer	13		Preparation wax	46
	LOT UNIVERSAL	72		Pressable ingots	78, 81, 82
<b>M</b>	MAGIC DISC	86		Pressable ingots LF	78
	Magic-Glass Stick	86		pro Waxer	58
	Margin liner	10		pro Waxer duo	58
	MASTER DIP	21	<b>R</b>	Retention pearl glue	48
	Metal	72,73, 112		Retention pearls	48
	Metal foil wax	94			
	Metal powder	112			
	Milling machines	100 - 103			

## Table of contents

<b>S</b>	Sakura margin liner	10	<b>U</b>	Underlining wax	27 - 29, 31
	Sample set	28		UVA light curing unit	64
	SATI BRUSH	71	<b>V</b>	Vario magnetic holder	10, 114
	Saw blades	8		Vinyl 3D-Scanner	99
	SCAN BOX	99		VKS wax for all-ceramic systems	28
	Scan wax	98	<b>W</b>	Wax brush	34
	Scanner	99		Wax chip	29
	Scanspray	98		Wax clasps	40, 41
	Separating disc	86		Wax Debubblizer	34, 49
	Separating materials	15, 17, 88		WAX GIANT	23 - 25
	Shade guide	88		Wax knife electrical	58 - 61
	Silicone cylinder	70		Wax meshes	45
	Silicone ring	70		Wax profiles	42 - 43
	Sintering furnace	110, 111		Wax profiles in bars	42 - 43
	SLM	113		Wax retention pearls	48
	SOLIBOND C plus	72		Wax retentions	44, 45
	Solibond C plus Blank	104 - 105		Wax sheets	46, 47, 92, 93
	SOLIBOND C plus Powder	113		Wax smoothing liquid	34, 49
	SOLIBOND N	73		Wax up instruments	56 - 57
	SOLIDUR CoCr	73		Wax wire	37
	SOLIDUS 80 & 84	92		Wire (wax)	37
	Spacer	10 - 13, 17	<b>Y</b>	YETI EXPANSION investment	68
	Special wax	24, 25, 27-29, 31		YETI EXPANSION plus investment	69
	Sprues	35, 36		Yeti Lube	15
	Square sticks	36		YETI piston compressor	110
	Stain Humidor	84		Yeti Rock	8
	STEAM OFF	16 - 17		Yeti Speed Sintering furnace for Zirconia	111
	Stick wax	27		YETIDUR Bio-dipping hardener	52
	Stick-On wax	27		YETIVEST investment	68
	Structured Casting wax	47	<b>Z</b>	Zircon blanks	106, 107
	SURFACE Marker	87		Zircon blanks GEO5	107
<b>T</b>	Tabeo Sintering furnace	110, 111		Zircon liquid	108 - 109
	Thinner	12		Zircon polishing set	114
	Thinner for Die Spacer	12			
	Thowax	26 - 29			
	Thowax sculping wax	26, 29			
	Thowax special wax	27, 29			
	Training set IQ	31			
	Transparent spacer	10			
	Triangel casting bar	36			

All information in this catalogue is subject to alteration without notice







Konzept & Design: Anke Biskupski, Engen

02.2020 V3 by YETI Dental



# YETI Dentalprodukte GmbH



# YETI Dentalprodukte GmbH

YETI Dentalprodukte GmbH · Industriestrasse 3 · D-78234 Engen · Germany  
<http://www.yeti-dental.com> · e-mail: [info@yeti-dental.com](mailto:info@yeti-dental.com)  
<http://www.dentalwax.com>