

Invacare Slings



OPTIMAL FIT AND COMFORT

PASSIVE LIFTERS | ACTIVE LIFTERS



**Making Life's
Experiences Possible**

Yes, you can.®

The Slings Range

Providing the optimal fit

Our new sling range has been designed to provide a safe, supportive and comfortable transfer for the client. The range provides options for various models, sizes and fabrics to ensure the sling fits perfectly to each person, whether they have some mobility impairment or require full body and head support.

Our team of experts and occupational therapists have spent time ensuring the slings are intuitive and easy to handle for the carer. This leads to reduced application times and improved confidence.

The range consists of eight models, each in up to six sizes and up to two fabric options to choose from.

This impressive selection ensures there is a sling suitable for everyone and for every use; daily transfer from bed to wheelchair/armchair, a toilet transfer or even for showering or bathing.



4

KEY FEATURES



◀ Optimal fit

- ▶ The right sling for the right client and transfer type
- ▶ Impressive range with 8 models, 6 sizes and 2 fabrics
- ▶ Well thought-out sling patterns supporting the right body parts
- ▶ All slings have loop attachment for a high spreader bar compatibility

Enhanced safety and confidence ▶

- ▶ Application videos easily accessible via QR code
- ▶ Clear sling labels with intuitive sling application drawings



◀ Optimal comfort

- ▶ Soft touch materials for increased comfort
- ▶ Padded back and leg sections for optimal weight distribution
- ▶ Cushioning under the tender armpit area

Easy to apply ▶

- ▶ Positioning pocket supports to slide between the client and backrest whilst sitting
- ▶ Two centre-line markers indicate the position over the spine
- ▶ Coloured loops facilitate the selection for the right positioning



INVACARE

Yes, you can.®

Selecting the right sling model

Always ensure the sling is specifically designed for the hoist you are using as failing to do so, may cause serious injury. There are different types and sizes of slings available. Selecting the most suitable sling for your client will ensure safety and improve comfort during the transfer. You should consider:

- Client's functional level, size, width and weight
- The type of transfer
- Sling shape matching body control e.g. U-shape or full-body support
- Any preferred position as the result of any medical condition e.g. stroke or amputation
- Material type and any padding
- Any sensitivity to pressure in any area



Some medical conditions such as stroke, orthopedic conditions, amputations or certain wounds may affect sling choice.

The main objective of the sling selection is to give the amount of body support that is needed. **Invacare Slings** are designed to give support in those areas where body control and muscle tone are reduced. The main body parts to assess are the legs, hip, trunk and head. The following chart gives an overview of which sling model is suitable depending on the client's physical ability:

Body Control assessment and Transfer Type:



Body Support Provided by Invacare Slings Model:

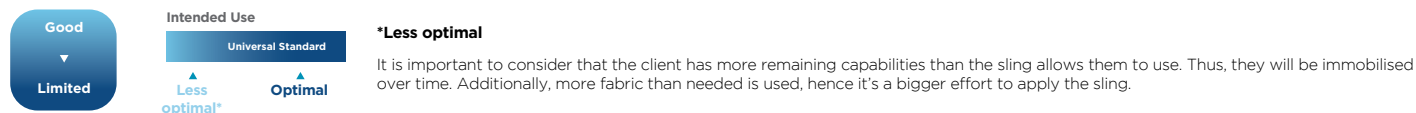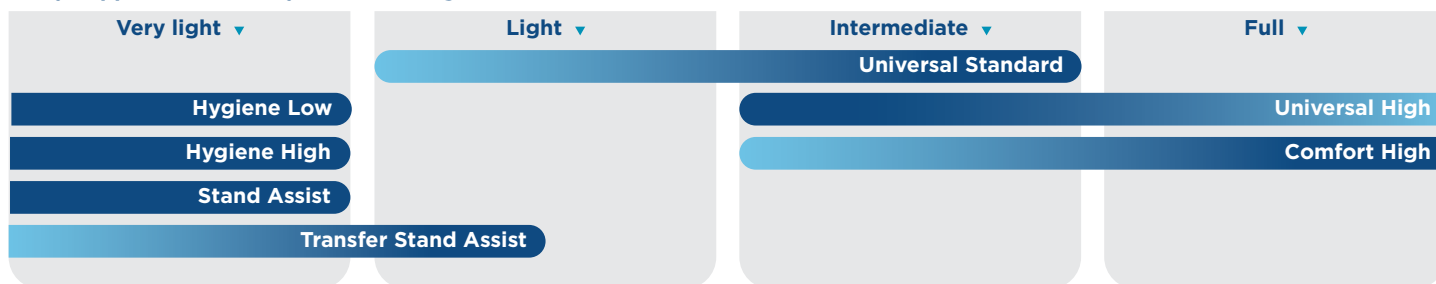


Selecting the right size

To select the correct sling size, the most important measurements are the hip width and back height in the sitting position. For slings with belts, the waist circumference is important, too.

The client's weight needs to be considered mainly for the sling's maximum safe working load. It's inferior to indicate the right sling size.

If no head support is needed, the correct back height is from the lower back until the shoulder/neck. With the need of head support the measurement needs to include the head, too.

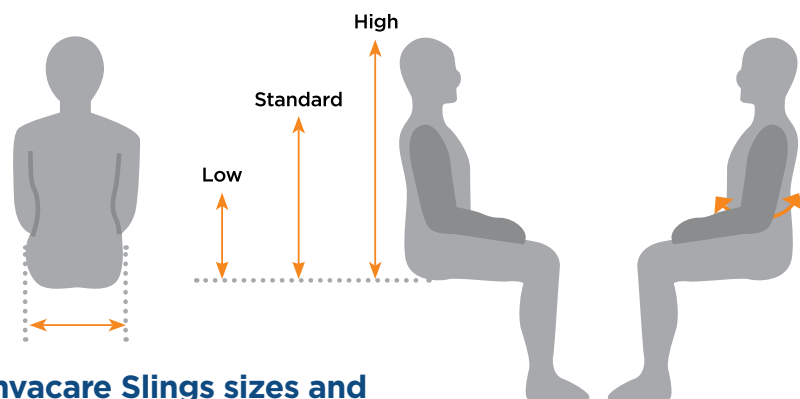

Invacare Slings sizes and safe working load:

Size	XS	S	M	L	XL	XXL
Safe Working Load	200 kg	200 kg	200 kg	200 kg	250 kg	300 kg

▶ If the sling is too big, the client may slip out. If the sling is too small, the client may still slip out, particularly with less body control. It could also exasperate the client's medical condition. When selecting a size, if your client is in-between sizes, having a slightly smaller sling, could keep your client more secure.

Sling features

Invacare slings have several features that enable you to position patients easily.



Invacare Slings for passive lifters


Universal Slings

All-round slings for everyday patient handling routines.



◀ Universal High

- Supports the whole body, including the head
- Client sits in slightly reclined position
- Suitable for those with limited control of the head, trunk and hips
- Can be used for transfers from a lying position, lifting from the floor and seated to seated transfers

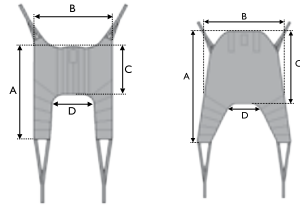

Universal family with regular opening


Universal High

Codes			A	B	C	D
Fabric	Mesh					
60124058	60124063	XS	1040	740	615	310
60124059	60124064	S	1080	845	770	350
60124060	60124065	M	1425	1000	900	375
60124061	60124066	L	1555	1120	1000	450
60124062	60124067	XL	1600	1190	1000	750
60125684		XXL	1650	1285	1000	600



◀ Universal Standard

- Supports the whole body, including the shoulders
- Suitable for those with good head control
- Can be used for transfers from a lying position, lifting from the floor and seated to seated transfers

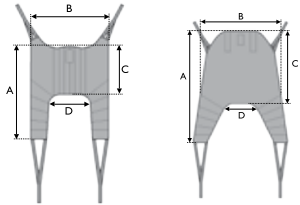

Universal family with regular opening

Universal Standard

Codes			A	B	C	D
Fabric	Mesh					
60124088	60124093	XS	965	800	405	365
60124089	60124094	S	1000	915	510	410
60124090	60124095	M	1085	1050	560	480
60124091	60124096	L	1180	1095	585	555
60124092	60124097	XL	1240	1240	590	720


Slings with a wider opening

Wide opening to facilitate toileting and easy to apply.



◀ Hygiene Low

- Suitable for those with hip, trunk and head control
- Belt offers increased support around the body and safe positioning
- Can be used for seated to seated transfers

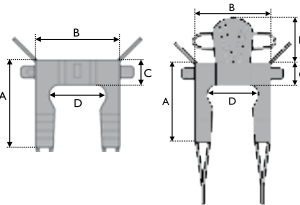

Wide opening


Hygiene Low

Codes			A	B	C	D
Fabric						
60123889		XS	840	620	160	380
60123890		S	870	710	200	430
60123891		M	950	800	230	495
60123892		L	1020	900	230	600
60123893		XL	1100	1080	270	760



◀ Hygiene High

- Same feature set as the Hygiene Low sling but comes with a head support for those with reduced head control

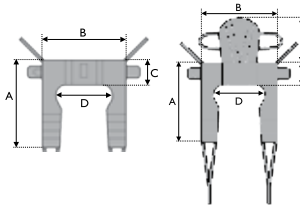

Wide opening

Hygiene High

Codes			A	B	C	D	E
Fabric							
60123894		XS	840	620	160	380	400
60123895		S	870	710	200	430	400
60123896		M	950	800	230	495	465
60123897		L	1020	900	230	600	465
60123898		XL	1100	1080	270	760	465


Comfort Slings

Hammock style slings, ideally suited to those with chronic pain. The sling is also useful for above-the-knee amputees.



◀ Comfort High

- Supports the whole body, including the head
- Suitable for those with limited control of the head, trunk and hips
- Including well positioned handles

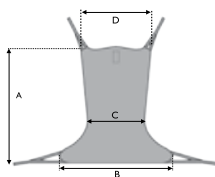

Comfort Slings

Comfort High


Codes			A	B	C	D
Fabric	Mesh					
60123856	60123861	S	1150	1050	550	1020
60123857	60123862	M	1395	1170	610	1160
60123858	60123863	L	1445	1290	675	1280
60123859	60123864	XL	1460	1330	745	1470

Invacare Slings for active lifters



◀ Stand Assist

- Suitable for clients with some weight bearing ability but cannot stand independently
- Client should be coherent, have good body and head control
- It's well-padded for extra comfort
- Comes with a non-slip material on the back

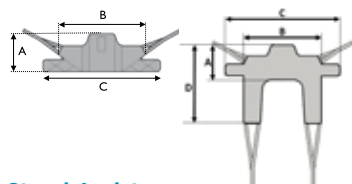


Stand Assist

Codes		A	B	C
Fabric				
60123941	S	260	840	950
60123942	M	330	895	1170
60123943	L	360	950	1270
60123944	XL	360	1020	1810



◀ Transfer Stand Assist

- Same feature set as the Stand Assist sling but comes with additional leg support

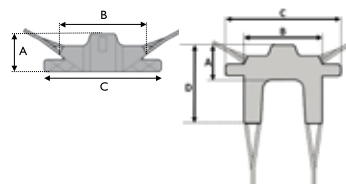

Transfer Stand Assist

Codes		A	B	C	D
Fabric					
60123945	S	260	840	950	965
60123946	M	330	895	1170	1060
60123947	L	360	950	1270	1160
60123948	XL	360	1020	1810	1260

Tech Specs



Max User Weight ▼

XS, S, M, L = 200 kg
XL = 250 kg
XXL = 300 kg



Max washing temperature ▼

90° C



Bleaching ▼

Do not bleach



Tumble dry ▼

Low temperature

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

All slings should be tested and adjusted to the individual client to ensure optimal comfort and safety during transfers.

A risk assessment must be carried out prior to use of one manufacturer's sling on another manufacturer's hoist.

Fabrics



Polyester - Solid

- ◀ A strong material that is hard-wearing and easy to apply. It slides smoothly into place and dries quickly after washing. This comes as standard on all models.



Polyester - Mesh

- ◀ Also known as net fabric, the mesh option is ideally suited for bathing or showering environments. It has quick drainage and drying properties.

Sling label

Each Invacare sling has an external label with information that supports you through application and patient positioning.



QR-Code

You can scan the QR (Quick Response) code using your Smartphone to access a sling application video and the user manual.


Microsite with all sling application videos

<https://marketing2.invacare.eu.com/slings>



Patient name

There is placeholder for the patient name on the label, allowing you to personalise the sling.

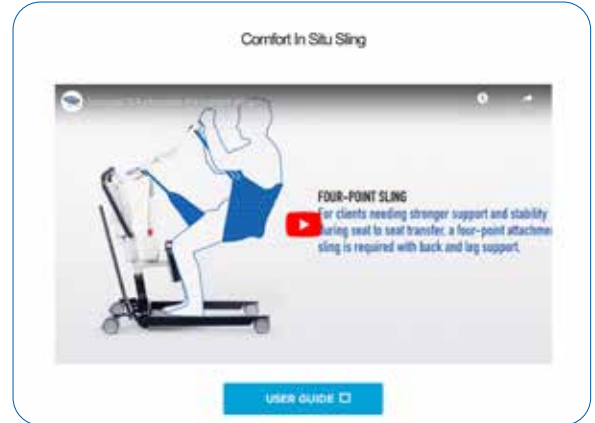

Various positioning options



Suspension strap

The suspension straps are equipped with colour coded loops which offer different lengths in order to place the patient in different positions.

Before applying the loops consider the users end position to choose the most appropriate loops. The pictograms on the sling's external label gives you quick support to remember the different positions. Additionally you can use loop markers to indicate the chosen loops.



For a transfer into a sitting position, place the patient as vertical as possible in the sling. Therefore, use shorter loops at the shoulders and longer loops at the legs.



For a semi-recumbent position in the sling use e.g. the blue loops on both the shoulders and on the legs.



For a transfer into a lying position use longer loops at the shoulders and shorter ones at the legs.

IMPORTANT!

For an even lift, always ensure the coloured loops selected on one side of the body match the other i.e. shoulder straps match and leg straps match accordingly.

The pivot sling range

Pivot clip release strap colour-coded by size

Upper body pivot clips

Care label with size and identification

Body section

Positioning handles

Padded leg section

Pivot clip release strap colour-coded by size

Leg section pivot clips



Selecting the correct size

Hygiene High Typically there is a degree of overlap between each sling size within a particular model. If the patient falls between two sizes, it is generally better to select the smaller size of sling, unless you are aware that the patient is likely to grow.

Measurements required for General Purpose with Head Support (Pivot and Yoke):

- Shoulder width
- Coccyx to top of head (General Purpose only)

Measurement required for Hygiene (Pivot and Yoke):

- Shoulder width

Measurements

One of the most important factors to perform a safe transfer is the match between hip/trochanter width and the measurement of the sling's opening between the leg support. If between sizes the smaller size may keep the client more secure.

Pivot General Purpose with Head Support

Codes		Shoulder Width (mm)	COCCYX TO TOP OF HEAD (mm)
SL146191	S	600 - 700	640 - 780
SL146291	M	700 - 810	750 - 900
SL146391	L	760 - 870	830 - 1000



Pivot Hygiene

Codes		Shoulder Width (mm)
SL247191	S	500 - 600
SL247291	M	500 - 600
SL247391	L	760 - 870



Tech Specs



Max User Weight ▼

S = 205 kg
M = 205 kg
L = 300 kg



Max washing temperature ▼

90° C



Bleaching ▼

Do not bleach



Tumble dry ▼

Low temperature

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.

All slings should be tested and adjusted to the individual client to ensure optimal comfort and safety during transfers.

A risk assessment must be carried out prior to use of one manufacturer's sling on another manufacturer's hoist.



Yes, you can.®

Cross reference guide

Invacare Item Code	New Description	Comments	Previous Item Code	Old Description	Comments
60123861	SLING COMFORT HIGH MESH S	No commode hole	1451097	SLING (GP-HYG-HS-FABRIC-M)	Hammock slings with commode hole
			1485968	AMPUTEE SLING FABRIC SMALL	
60123862	SLING COMFORT HIGH MESH M	No commode hole	1485969	SLING (COMFORT-MESH-M)	
			1451098	SLING (GP-HYG-HS-MESH-M)	Hammock slings with commode hole
			17336	SLING (HAMMOCK-MESH-MED)	
60123863	SLING COMFORT HIGH MESH L	No commode hole	1485970	SLING (COMFORT-MESH-L)	
			1451099	SLING (GP-HYG-HS-MESH-L)	Hammock slings with commode hole
			17129	SLING (HAMMOCK-MESH-LGE)	
60123864	SLING COMFORT HIGH MESH XL	No commode hole	1451100	SLING (GP-HYG-HS-MESH-XL)	Hammock slings with commode hole
60123889	SLING HYGIENE LOW FAB XS	NEW	-	-	-
60123890	SLING HYGIENE LOW FAB S		1478440	SLING(SUPREME-HYG-FABRIC-SM)	
			17530	SLING(SUPREME-HYG-FABRIC-SM)	
			THY-S	SLING HYGIENE FAB SM	
60123891	SLING HYGIENE LOW FAB M		1478441	SLING (SUPREME-HYG-FAB-MED)	
			17531	SLING (SUPREME-HYG-FAB-MED)	
			THY-M	SLING HYGIENE FAB M	
60123892	SLING HYGIENE LOW FAB L		1478442	SLING (SUPREME-HYG-FAB-LG)	
			17532	SLING (SUPREME-HYG-FAB-LG)	
			THY-L	SLING HYGIENE FAB L	
60123893	SLING HYGIENE LOW FAB XL	NEW	-	-	-
60123894	SLING HYGIENE HIGH FAB XS		17401	SLING GPHS HYGIENE FAB XSML	
60123895	SLING HYGIENE HIGH FAB S		1478444	SLING (DRESS HYG-HS-FAB-SM)	
			17402	SLING (DRESS HYG-HS-FAB-SM)	
60123896	SLING HYGIENE HIGH FAB M		1478445	SLING (DRESS HYG-HS-FAB-M)	
			17395	SLING (GP-HYG-H/S-FABRIC-MED)	
60123897	SLING HYGIENE HIGH FAB L		1478446	SLING (DRESS HYG-HS-FAB-L)	
			17396	SLING (GP-HYG-H/S-FABRIC-L)	
60123898	SLING HYGIENE HIGH FAB XL		1484780	SLING-DRESS TOILET-W/HD SP-XL	
			17403	SLING(HYGIENE-GP H/S-FAB-XLGE)	
60123941	SLING STAND ASSIST FAB S		1484685	SLING STAND ASSIST FAB S	
60123942	SLING STAND ASSIST FAB M		1484686	SLING STAND ASSIST FAB M	
60123943	SLING STAND ASSIST FAB L		1484687	SLING STAND ASSIST FAB L	
60123944	SLING STAND ASSIST FAB XL		1484688	SLING STAND ASSIST FAB XL	
60123945	SLING TRANS STAND ASST FAB S		1484689	SLING TRANSFER STAND ASSIST S	
60123946	SLING TRANS STAND ASST FAB M		1484690	SLING TRANSFER STAND ASSIST M	
60123947	SLING TRANS STAND ASST FAB L		1484691	SLING TRANSFER STAND ASSIST L	
60123948	SLING TRANS STAND ASST FAB XL		1485117	SLING TRANSFER STAND ASSIST XL	
60124058	SLING UNIVERSAL HIGH FAB XS		1483776	SLING (UNIVERSAL-HS-FAB-XS)	
			17209	SLING (GP H/S-FAB-XSML)	
60124059	SLING UNIVERSAL HIGH FAB S		1483777	SLING (UNIVERSAL-HS-FAB-SM)	
			17247	SLING (GP H/S-FAB-SML)	
60124060	SLING UNIVERSAL HIGH FAB M		1483778	SLING (UNIVERSAL-HS-FAB-M)	
			17067	SLING (GP H/S-FAB-MED)	
60124061	SLING UNIVERSAL HIGH FAB L		1483779	SLING (UNIVERSAL-HS-FAB-L)	
			17069	SLING (GP H/S-FAB-LGE)	
60124062	SLING UNIVERSAL HIGH FAB XL		1485949	SLING (UNIVERSAL-HS-FAB-XL)	
			17205	SLING (GP H/S-FAB-XL)	

Invacare Item Code	New Description	Comments	Previous Item Code	Old Description	Comments
60124063	SLING UNIVERSAL HIGH MESH XS	Dimensions similar to paed	1485950	SLING (UNIVERSAL HS MESH XS)	
			17260	SLING (PAED-GP-MESH-XSML)	
60124064	SLING UNIVERSAL HIGH MESH S	Dimensions similar to paed	1485951	SLING (UNIVERSAL-HS-MESH-S)	
			17126	SLING (GP H/S-MESH-SML)	
			17262	SLING (PAED-GP-MESH-SML)	
60124065	SLING UNIVERSAL HIGH MESH M	Dimensions similar to paed	1485952	SLING (UNIVERSAL-HS-MESH-M)	
			17127	SLING (GP H/S-MESH-MED)	
			17266	SLING (PAED-GP-MESH-MED)	
60124066	SLING UNIVERSAL HIGH MESH L		1485953	SLING (UNIVERSAL-HS-MESH-L)	
			17128	SLING (GP H/S-MESH-LGE)	
60124067	SLING UNIVERSAL HIGH MESH XL		1485954	SLING (UNIVERSAL-HS-MESH-XL)	
			17248	SLING (GP H/S-MESH-XLGE)	
60124088	SLING UNIVERSAL STD FAB XS	NEW	-	-	-
60124089	SLING UNIVERSAL STD FAB S		1485882	SLING (UNIVERSAL-STD-FAB-SM)	
			17105	SLING (GP-FABRIC-SML)	
60124090	SLING UNIVERSAL STD FAB M		1485883	SLING (UNIVERSAL-STD-FAB-M)	
			17041	SLING (GP-FABRIC-MED)	
60124091	SLING UNIVERSAL STD FAB L		1485884	SLING (UNIVERSAL-STD-FAB-L)	
			17043	SLING (GP-FABRIC-LGE)	
60124092	SLING UNIVERSAL STD FAB XL		1485885	SLING (UNIVERSAL-STD-FAB-XL)	
			17204	SLING (GP-FABRIC-XLGE)	
60124093	SLING UNIVERSAL STD MESH XS	NEW	-	-	-
60124094	SLING UNIVERSAL STD MESH S		1478411	SLING (UNIVERSAL-STD-MESH-SM)	
			17123	SLING (GP-MESH-SML)	
60124095	SLING UNIVERSAL STD MESH M		1478412	SLING (UNIVERSAL-STD-MESH-M)	
			17124	SLING (GP-MESH-MED)	
60124096	SLING UNIVERSAL STD MESH L		1478413	SLING (UNIVERSAL-STD-MESH-M)	
			17125	SLING (GP-MESH-LGE)	
60124097	SLING UNIVERSAL STD MESH XL		1478414	SLING UNIVERSAL STD MESH XL	
			17207	SLING (X-LARGE-MESH-GEN PURP)	
60125684	SLING UNIVERSAL HIGH FAB XXL	NEW	-	-	-
SL146191	SLING PIVOT GP HS MESH SM		17241	SLING (PIVOT-GP-MESH-SML)	
SL146291	SLING PIVOT GP HS MESH M		17242	SLING (PIVOT-GP-MESH-MED)	
SL146391	SLING PIVOT GP HS MESH L		17243	SLING (PIVOT-GP-MESH-LGE)	
SL247191	SLING PIVOT HYG MESH S	NEW	-	-	-
SL247291	SLING PIVOT HYG MESH M		17392	SLING (PIVOT-HYG-MESH-MED)	
SL247391	SLING PIVOT HYG MESH L		17393	SLING (PIVOT-HYG-MESH-LGE)	



FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

Asia: Phone +61 2 8839 5330, email anzexportsales@invacare.com or visit www.invacare.com

Thailand: Phone +66 0 2821 5515 or visit www.invacare.com.th



Yes, you can.[®]