


Ducted Fume Hood



Fume hood FH700: Negative pressure in work area protects operator and environment.

Air Flow Pattern:

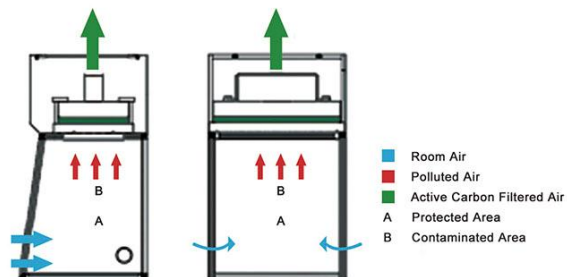

LED Display
Digital LED display: Airflow Level



Folding acrylic front window, down part with free stop function



Three Side Glass Windows
Back&Side glass windows maximize light and visibility, providing a bright and open working environment



Base Stand
Optional

ADVANTAGE

1. Folding acrylic front window, down part with free stop function.
 2. Microprocessor control system, LED display.
 3. With memory function in case of power-failure.
 4. Built-in centrifugal blower.
 5. Interlock function: UV lamp and blower, LED lamp. (Only for choosing UV lamp).
- Note:** Experiments with strong acids and alkali, flammable or explosive substances cannot be conducted with this type fume hood.

TECHNICAL PARAMETERS

Model	FH700	
External Size(W*D*H)	700*620*1150mm	
Internal Size(W*D*H)	640*550*700mm	
Max Opening	660mm	
Air Velocity	0.3~0.8m/s	
Noise	≤65dB	
LED Lamp	4W*2	
Blower	Built-in centrifugal blower; Speed adjustable	
Front Window	Folding acrylic window, down part with free stop function	
Material	Main Body	1.0mm cold-rolled steel with bacteria power coating
	Work Table	Chemical resistant phenolic resin
	Back&Side Windows	5mm toughened glass
Consumption	200W	
Power Supply	AC220V±10%, 50/60Hz; 110V±10%, 60Hz	
Standard Accessory	LED lamp, Active carbon filter, Blower, 4.3 meters PVC exhaust duct, Pipe strap	
Optional Accessory	15W UV lamp, HEPA filter, Base stand	
Gross Weight without Base Stand	100kg	
Gross Weight with Base Stand	115kg	
Package Size without Base Stand (W*D*H)	840*760*1400mm	
Package Size with Base Stand (W*D*H)	840*960*1400mm	