

# CONDYLAR HINGED KNEE (WRAP)

**EZYWRAP®** designed the Condylar Hinged Knee Wrap to be fully adjustable, exceedingly comfortable, and highly durable for active individuals. We use a breathable non-neoprene hypoallergenic material that keeps the patient cool during exercise and comfortable with anti-migration elements to minimize slippage during extended use.

Our Knee Orthosis has hidden, dual-axis hinges that are bendable and shapable to facilitate fitting and adjustment. The hinges include condyle pads for patient comfort and support. Additionally, the hinges are removable for ease of maintenance and patient cleaning. All while being designed with anterior closure for tailored compressive fit.

Approved **PDAC** **L1820**



## INDICATIONS FOR USE:

- Knee Ligament Injuries
- Meniscus Injuries
- Knee Sprains and Strains

## \*|FirstIce® FEATURES:

- Patented technology that contours to the body and lasts 5x longer than gel.
- Direct Stick™ global placement eliminates the need for pockets.
- Permanently bonded antimicrobial eliminates odor for the life of the product.




**P:** (800) 234-9004  
**F:** (800) 274-9006  
**E:** SALES@EZYWRAP.COM  
**W:** EZYWRAP.COM

FOR THE USE OF:  
CODE:

MADE IN THE USA     

# FEATURES & BENEFITS

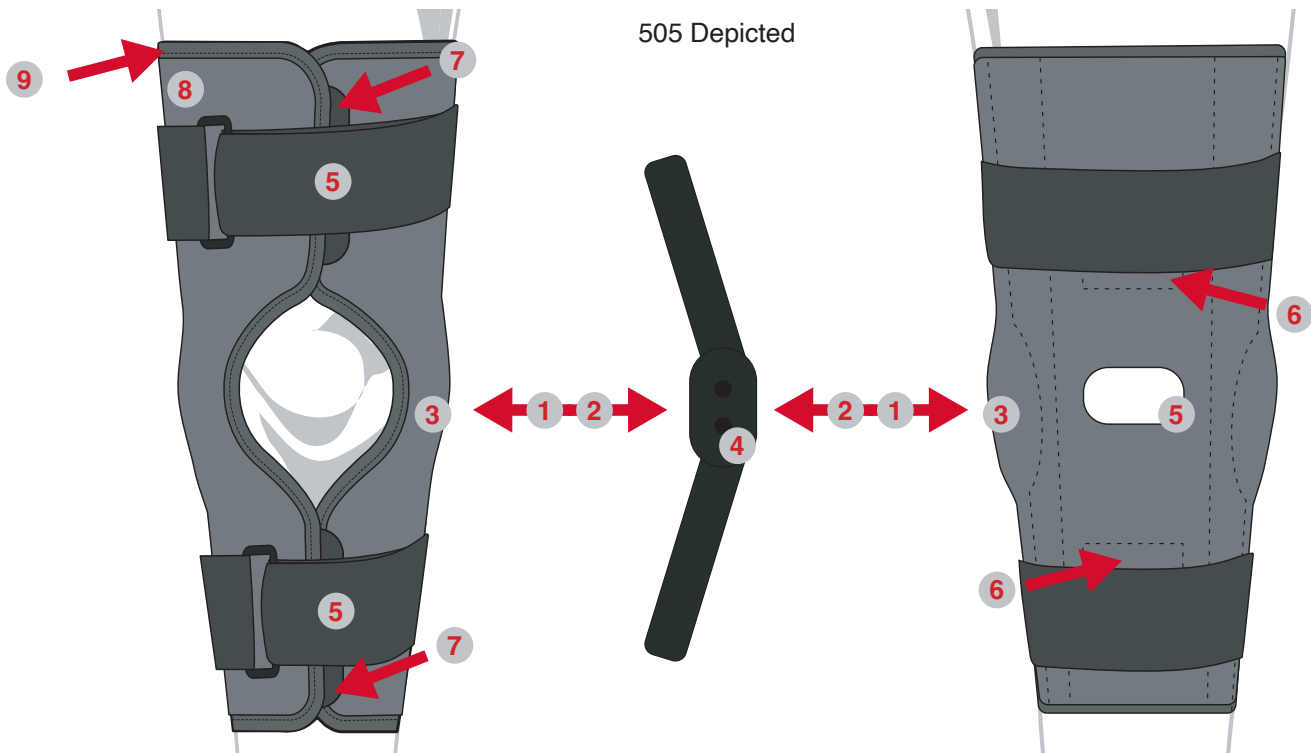

- 1. Shapable Dual Axis hinges for customization
- 2. Internal hinges for protection during activity
- 3. Condyle pads for support and comfort
- 4. Removable hinges for ease of maintenance and cleaning
- 5. Open popliteal and removable straps
- 6. Anti-migration elements to ensure effective support
- 7. Anterior closure for tailored compressive fit
- 8. Breathable non-neoprene hypoallergenic material
- 9. Soft edges for user comfort
- 10. FirstIcE® Direct Stick™ Cold Therapy application directly to the brace interior

## PRODUCT ORDERING

Product Name	Ref No.
Condylar Hinged Knee (Wrap) – Short	505-XX
Condylar Hinged Knee (Wrap) – Long	506-XX

## OPTIONAL COLD THERAPY ACCESSORIES

Product Name	Size	Ref No.
FirstIcE® Direct Stick™ Cold Therapy Pack	06 x 07"	101-0607
	06 x 10"	101-0610
	06 x 12"	101-0612

**Note:** Multiple FirstIcE® Direct Stick™ Packs may require increasing the brace size

## PRODUCT SIZING

Measure circumference of leg 6" above mid-patella with leg fully extended

Size	XX	Min	Max
XS	10	12" 31cm	16" 41cm
S	20	15" 39cm	18" 46cm
M	30	17" 44cm	20" 51cm
L	40	20" 51cm	23" 59cm
XL	50	23" 59cm	26" 66cm
XXL	60	26" 66cm	29" 74cm
XXXL	70	29" 74cm	32" 81cm

R&ID# \*285960\*  
REV A 6/22

FOR THE USE OF:  
CODE:

MADE IN THE USA     

**FOR ORDERING:**

**P:** (800) 234-9004

**F:** (800) 274-9006

**E:** SALES@EZYWRAP.COM

**W:** EZYWRAP.COM

*Scan for application  
and accessories.*

