

# EMSlim Neo Body Sculpture




**5000W STRONG POWER**  
**4 HANDLES WORK SIMULTANEOUSLY**

**+25% Muscle MORE**

**-30% Fat LESS**

**COMBINED 30-MIN. TREATMENT**




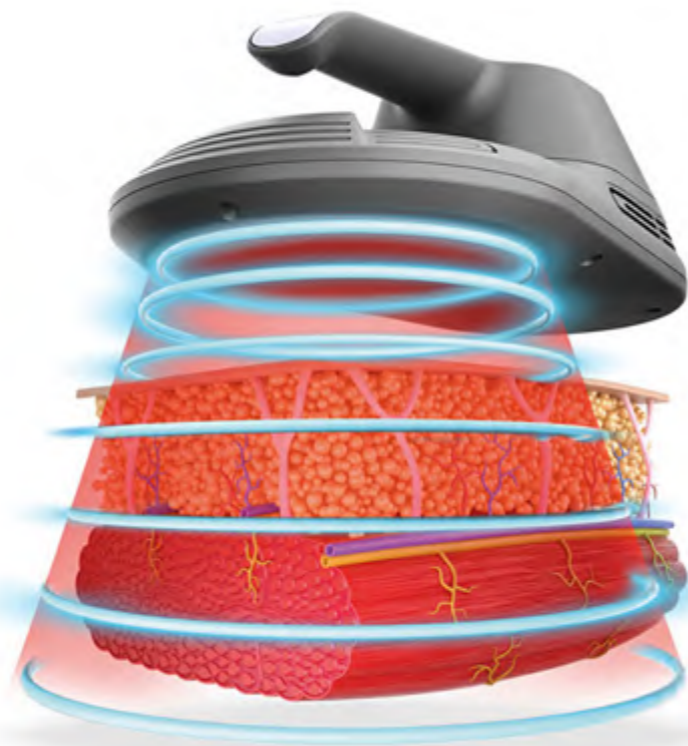

**RF**

**High Intensity Focused Electro-Magnetic**

**-30% Fat LESS**

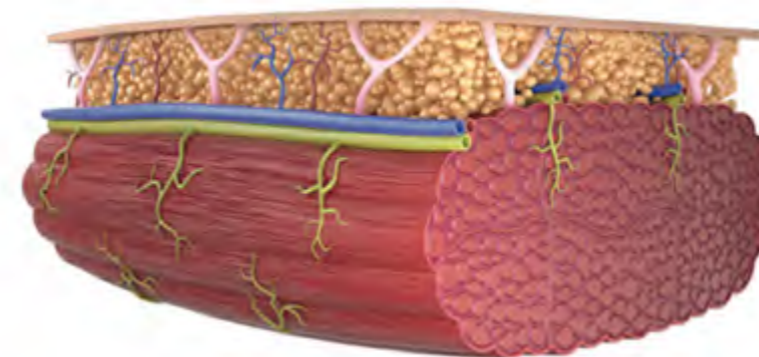

**+25% Muscle MORE**





MUSCULPTING NEO is based on an applicator simultaneously emitting synchronized RF and High Intensity Focused Electro-Magnetic energies.

Due to the radiofrequency heating, the muscle temperature quickly raises by several degrees. This prepares muscles for exposure to stress, similar to what a warm up activity does before any workout. In less than 4 minutes, the temperature in subcutaneous fat reaches levels that cause apoptosis, i.e. fat cells are permanently damaged and slowly removed from the body. Clinical studies showed on average a 30% reduction in subcutaneous fat.\*


Bypassing the brain limitations, High Intensity Focused Electro-Magnetic energy contracts the muscle fibers in the area at intensities that are not achievable during voluntary workout. Extreme stress forces the muscle to adapt, resulting in an increase in the number and growth of muscle fibers and cells. Clinical studies showed on average a 25% muscle volume increase.\*



# Mechanism of muscle contraction

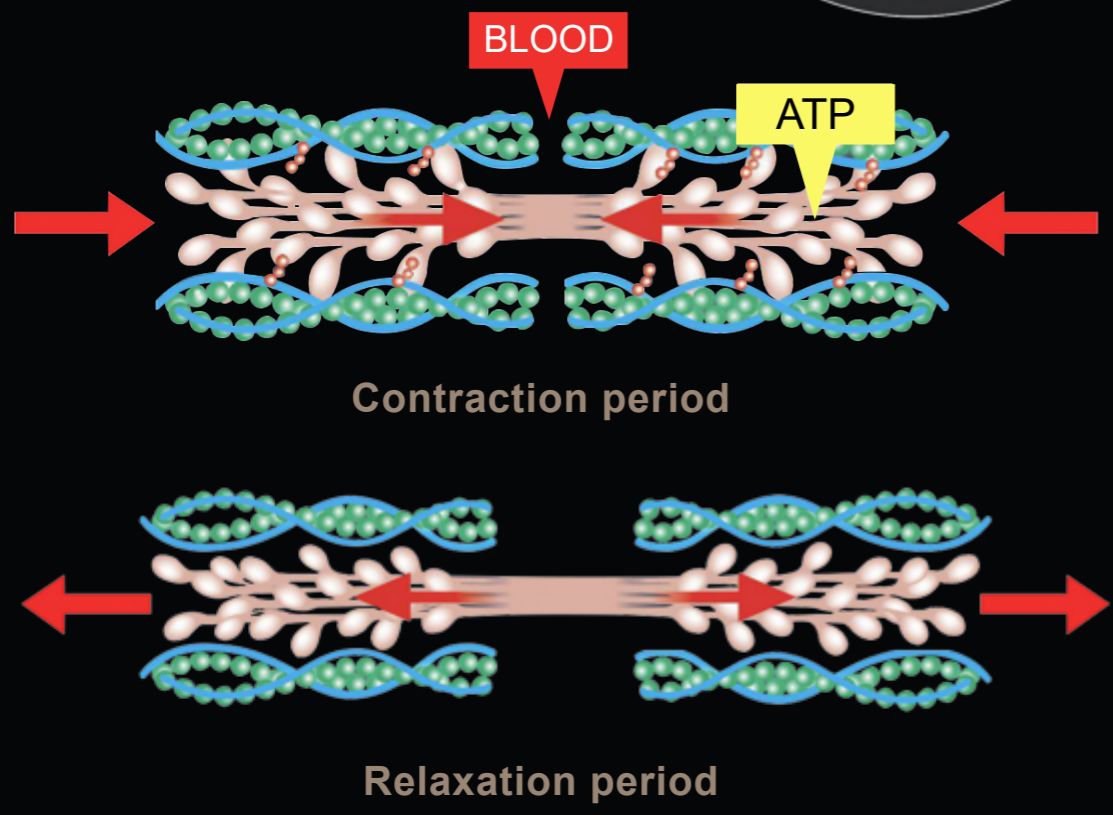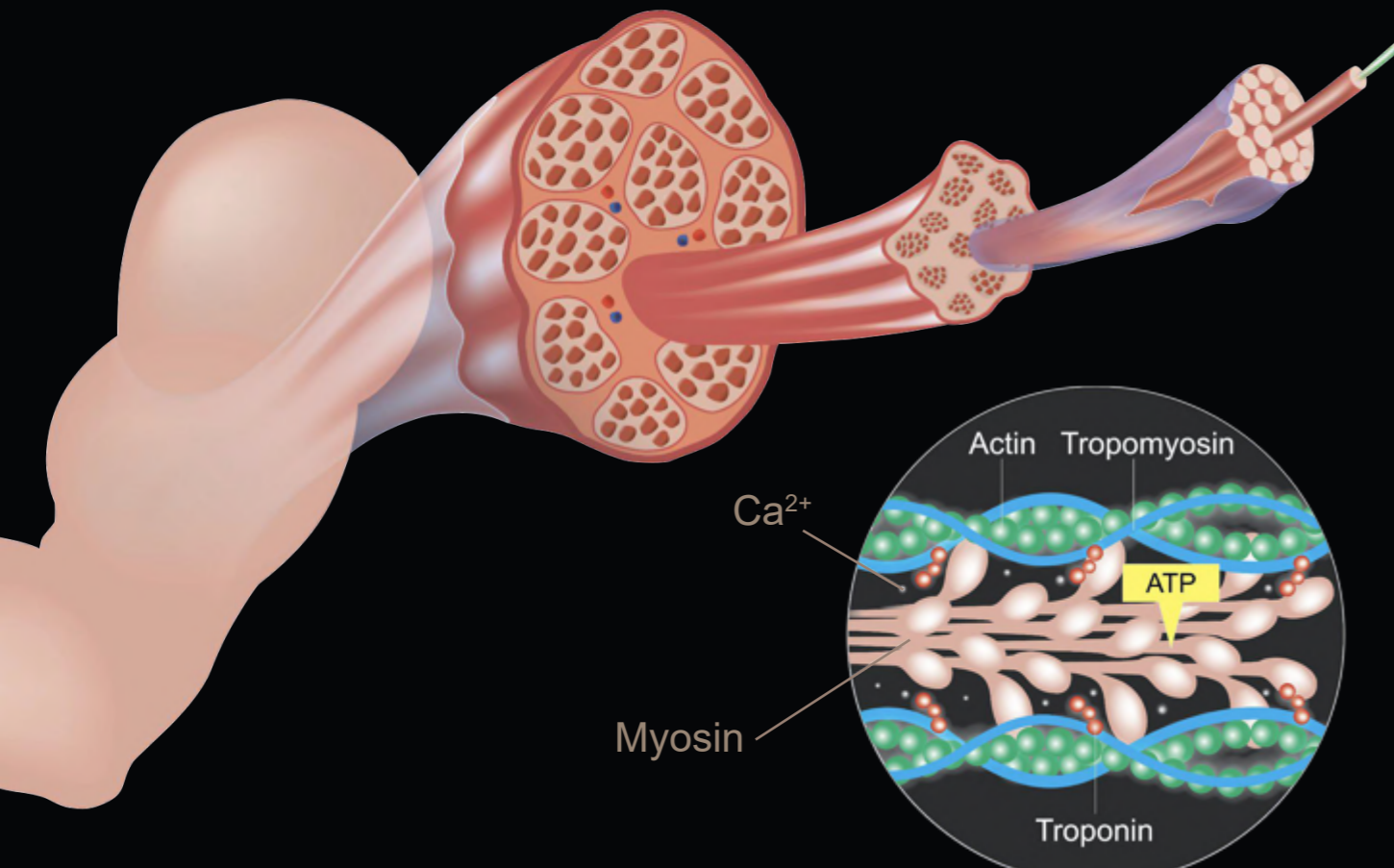

Nerve impulses set off a biochemical reaction that causes myosin to stick to actin

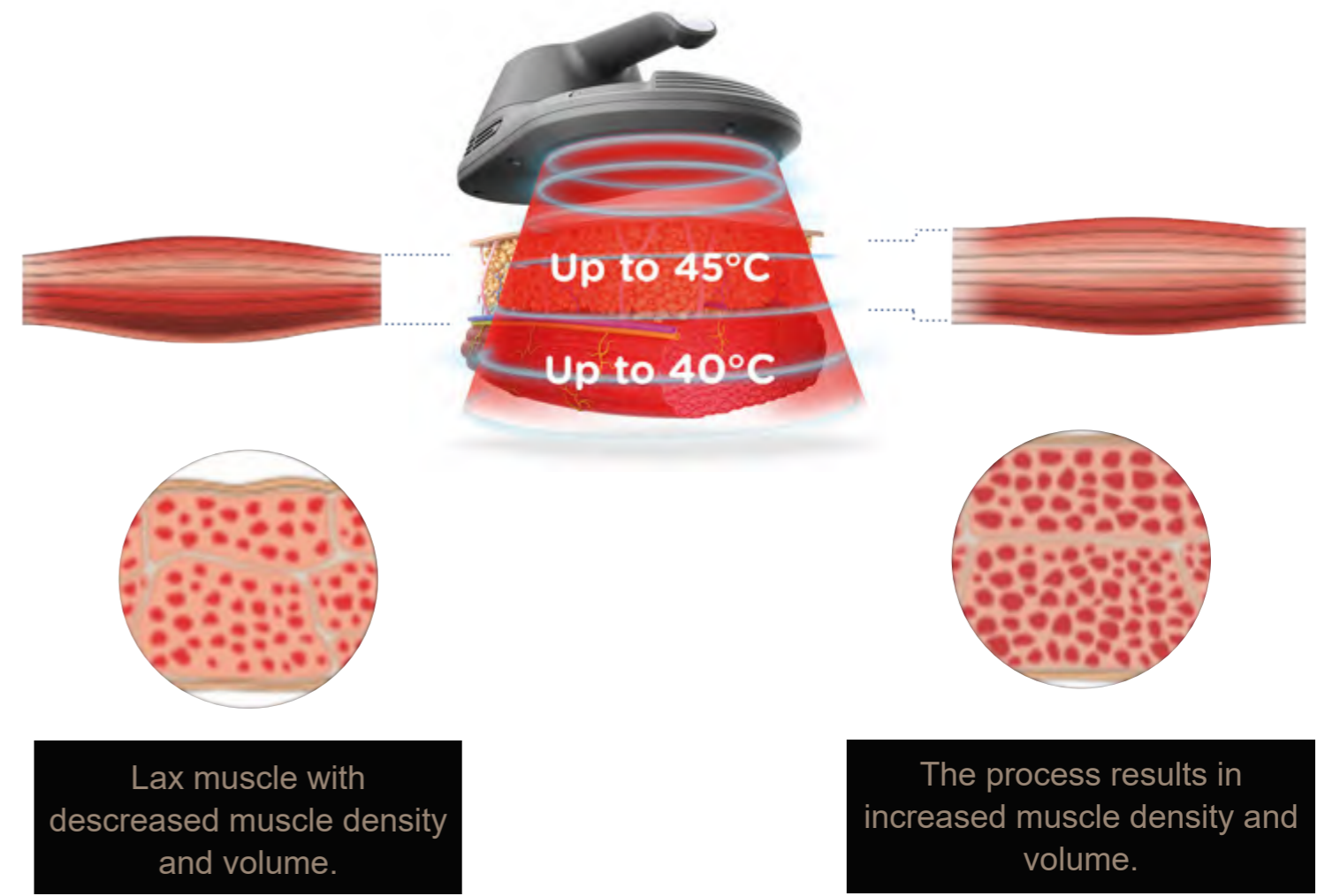

MUSCULPTING NEO induces very strong (supramaximal) muscle contractions, not achievable through workout. When exposed to these strong contractions, the muscle tissue is forced to adapt to such extreme condition. The muscle starts remodelling its inner structure by multiplication of fibers and their growth. The High Intensity Focused Electro-Magnetic uses a specific range of frequencies that does not allow muscle relaxation between two consecutive stimulations. The muscle is forced to remain in contracted state for multiple seconds. When repeatedly exposed to these high load conditions the muscle tissue is stressed and is forced to adapt. Recent studies reported that on average 15%- 16% increase in abdominal muscle thickness was observed in treated patients one to two months after High Intensity Focused Electro-Magnetic treatments.

## FLAT HANDLE



ABDOMEN

BUTTOCKS

THIGHS


BACK

## CURVED HANDLE



ARMS

CALVES

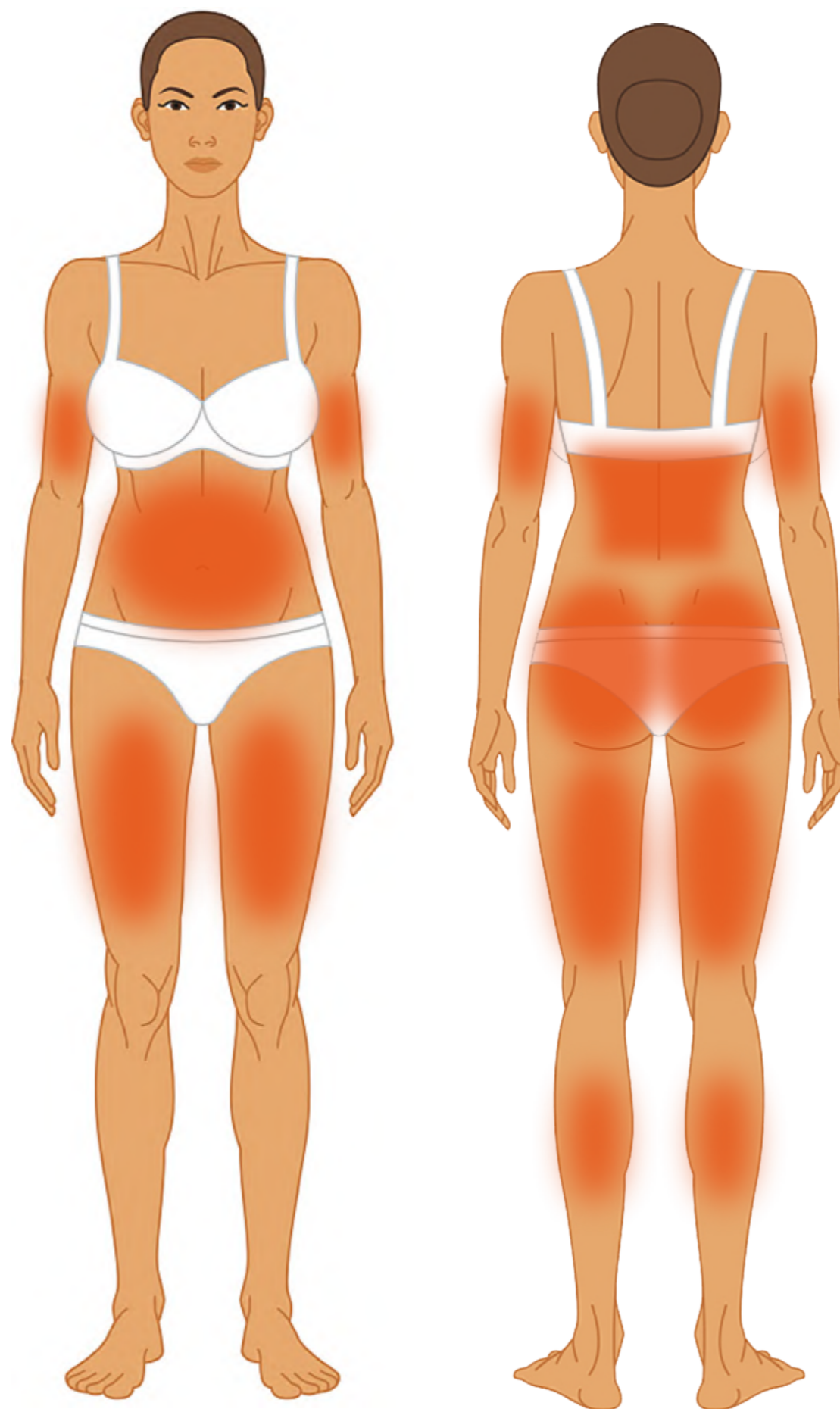


## PELVIC FLOOR HANDLE

Pelvic floor handle works to relieve incontinence by stimulating and strengthening the underlying muscles of the pelvic floor. During the treatment, you sit on a specially designed handle which delivers electromagnetic energy to the pelvic floor muscles, forcing them to contract. Over the course of your treatment, these contractions strengthen the weakened muscles and restore neuromuscular control of the bladder.





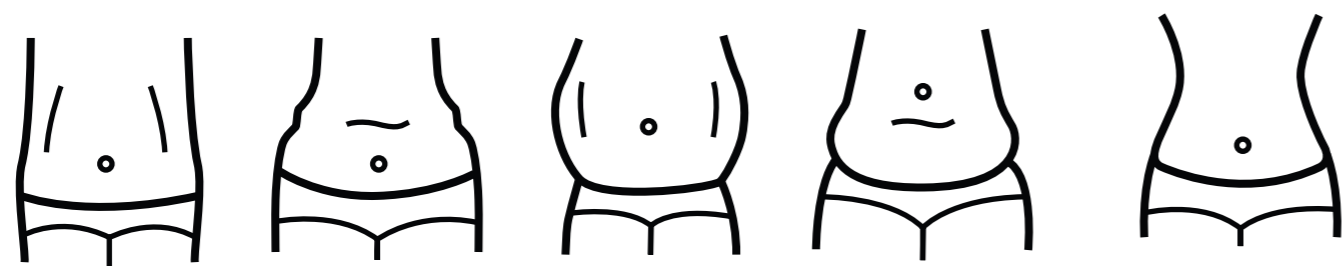


H Type

Folding Type

Fat Type

After Childbirth


Build abdominal muscles  
Perfect Firm abs  
and V line abs

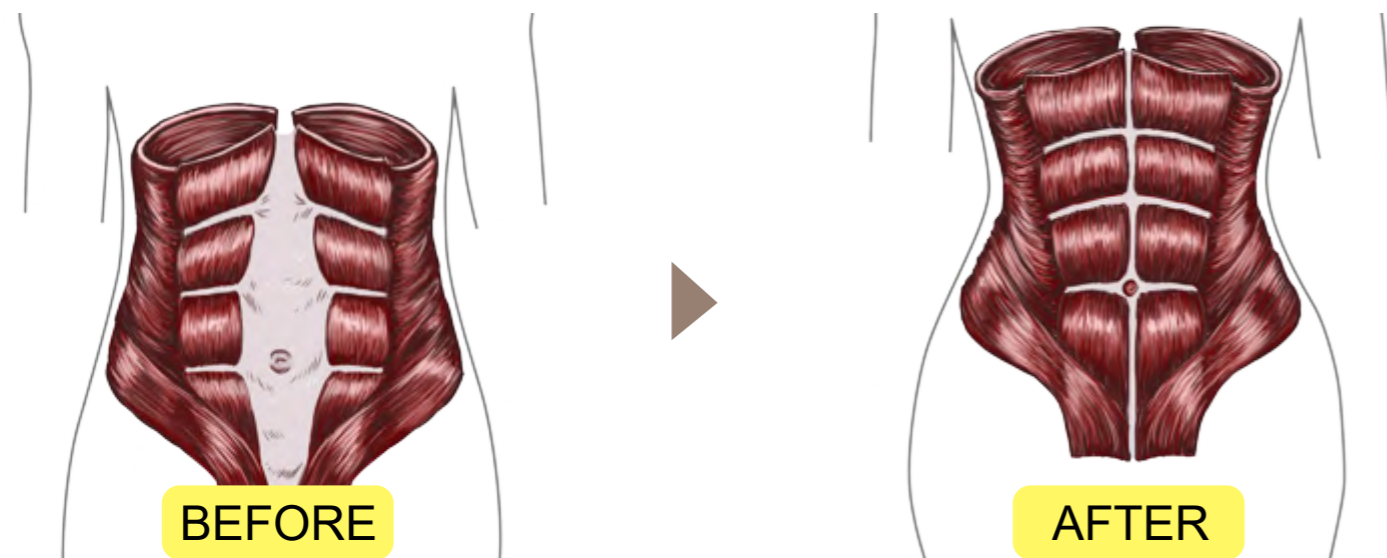

BEFORE



AFTER

The MUSCULPTING NEO works on the abdomen area contracting and tightening the muscles to both burn fat and sculpt the area. This treatment can help you improve the muscle definition that you are already seeing from your fitness routine. It is also ideal for addressing that stubborn inch of fat on your stomach that just won't go away.

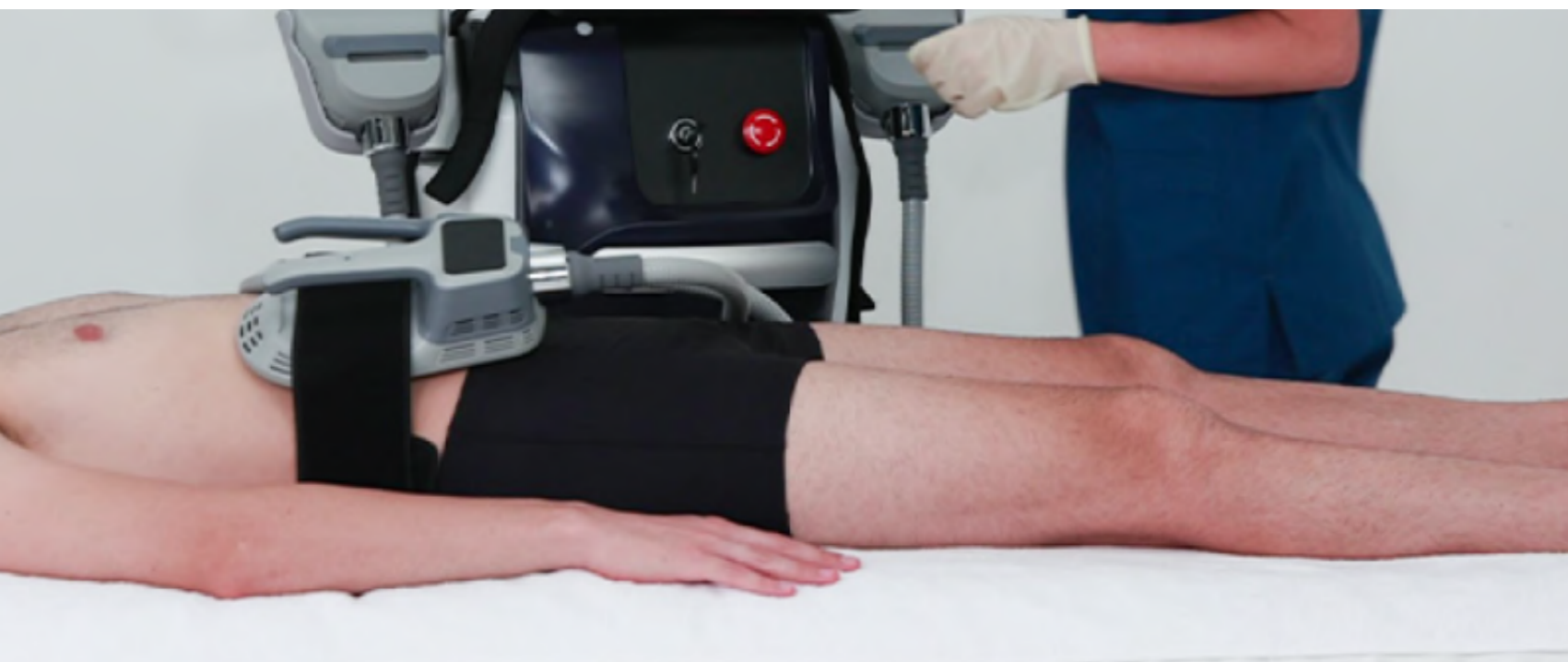

for the abdomen



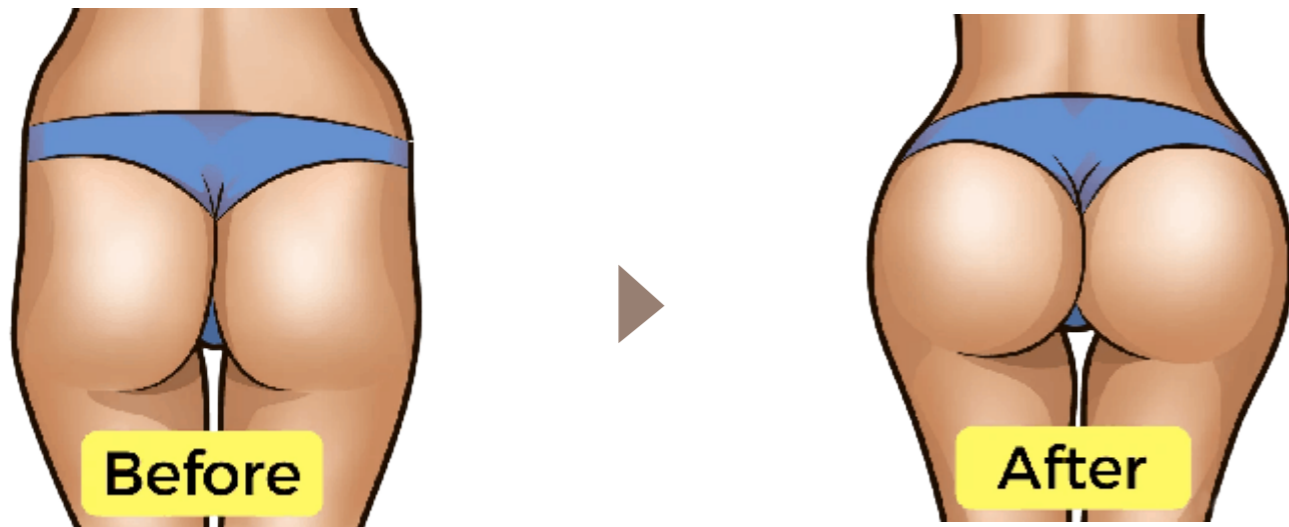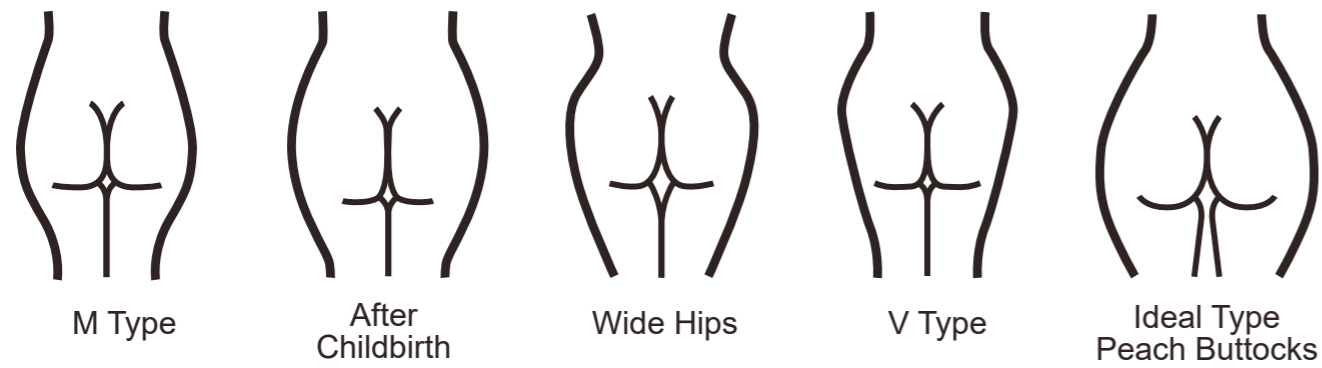

MUSCULPTING NEO works on the buttocks as well, using the same supramaximal muscle contractions as with the abdomen area. Many patients opt to try this treatment in place of an invasive Brazilian butt lift. High Intensity Focused Electro-Magnetic technology delivers the world's first noninvasive butt lift procedure. This treatment can give you the curvature, tone, and lift that you desire. It is an alternative to the Brazilian Booty Lift procedure and provides a defined shape instead of adding excess fat to the area.



for the buttock







## Dedicated multi-directional handle fixed arm

Multi-directional fixed support arm can be adjusted freely at any angle of 360°

The multi-directional fixed arm greatly facilitates the treatment of various parts of the body, and perfectly solves the problem that the handles of the shoulders, handles, shoulder blades and other parts are not easy to fix. Users can easily enjoy the treatment of muscle gain and fat reduction in various parts of the body.



# Intelligent Operating System

## GENDER SELECTION


Male | Female

## TREATMENT AREA SELECTION

Buttock | Abdomen | Thigh | Arm | Calf

## WORKING MODE SELECTION

Hiit | Strength | Muscle Building | Shape Up | Fat Burning

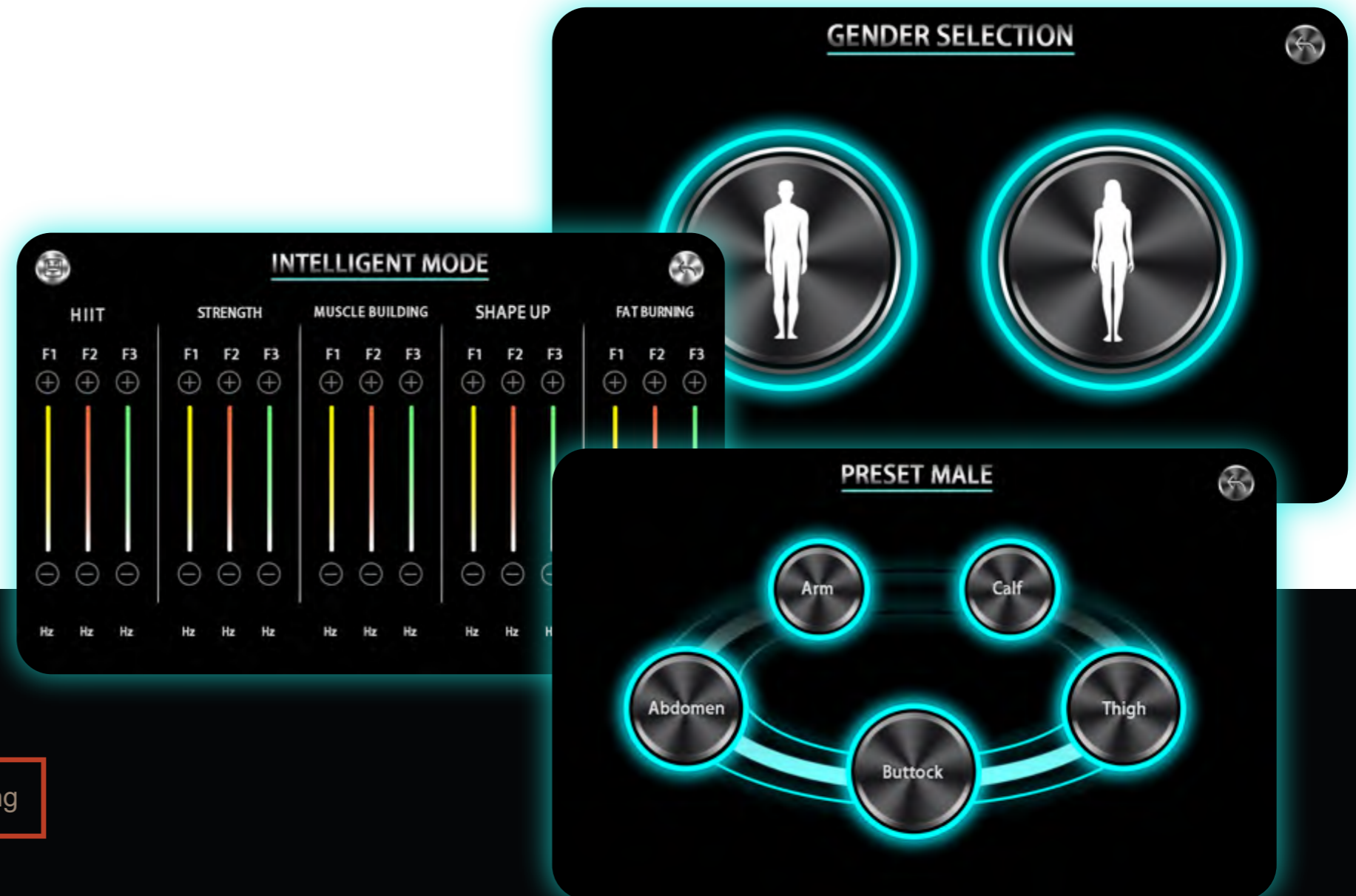

The new intelligent operating system can automatically match different treatment parameters according to different treatment parts; at the same time, it is equipped with five different working modes to directly match the corresponding treatment parameters according to the treatment needs of customers. Extremely easy to operate



# Comparison

## EMSlim Neo

## EMSlim

|                             |  |   |
|-----------------------------|--|---|
| Technology                  | High Intensity Focused Electro-Magnetic + RF | High Intensity Focused Electro-Magnetic |
| What can it do              | Build & tone muscle , Burn Fat, Tighten Skin | Build & Tone Muscle                     |
| What does it treat          | Muscle + Fat + Skin                          | Muscle                                  |
| How long is the treatment   | 30minutes                                    | 30 minutes                              |
| How many sessions           | 4 minimum                                    | 4 minimum                               |
| Contactless or not          | Needs to be on the skin tomake contact       | Can be on the skin or overclothes       |
| Discomfort                  | None   | None                                    |
| Recovery                    | Absolutely none                              | Absolutelly none                        |
| How often should one get it | Between 1-3 months, variable                 | Between 1-3 months, variable            |
| BMI requirement             | Up to BMI of 35                              | Up to BMI of 30                         |



Clinic

# ANYWHERE Shape Your Body ANYTIME



Home



Gym

MUSCULPTING NEO has a wide range of application scenarios, which can be used in beauty clinics, gyms and even offices or homes to shape your perfect body anytime, anywhere



Office

## SPECIFICATIONS

|              |   |
|--------------|---|
| Frequency:   | 1-150hz                                 |
| RF Power:    | 800W                                    |
| Power:       | 5000w                                   |
| Intensity:   | 7Tesla                                  |
| Screen:      | 12.4inch tough screen                   |
| Pulse width: | 310µs±10%                               |
| Cooling:     | Air Cooling                             |
| Voltage:     | AC220V±10% 10A 50HZ , 110v±10% 10A 60HZ |