

THREE TREATMENT MODES

SHR IPL ELIGHT three treatment modes to meet a variety of treatment needs








INTELLIGENT OPERATING SYSTEM

The operator only needs to select the part and skin color to be treated, and the system will automatically generate the recommended treatment parameters. Even an inexperienced operator can easily grasp the operation of the instrument.



SHR "In-motion" Hair Removal

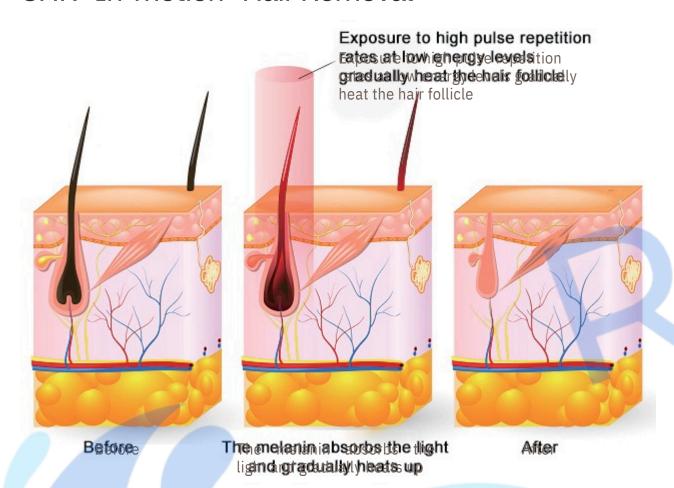

SHR stands for Super Hair Removal and it's the latest breakthrough innovation in IPL technology to launch onto the ever-advancing beauty and aesthetics market.













Faster

The SHR machine repetition rate is up to 10Hz. The handpiece could shoot 10 shots in 1 second.

SHR in-motion treatment method has much improved the treatment speed.




Unlike traditional IPL treatments, which use high energy levels that may cause discomfort and pain, SHR Skin Pulse treatments are mild and comfortable, SHR scans the treatment zone—gradually heating the skin with repetitive low fluence pulses, ensuring that you'll experience comfort even in the most sensitive of areas.



All Skin Types

The gradual heat distribution of SHR means treatments are totally comfortable, non-irritating and virtually pain free making it ideal for even sensitive dark skin types.


IPL ELIGHT SKIN TREATMENT

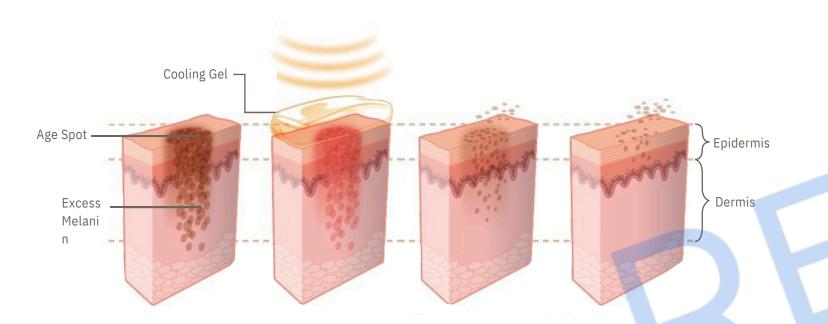

IPL technology is a type of Intense Pulsed Light that generates intense bursts of red or blue light.

These heat the dermal tissue without harming or removing the outer epidermal layers of the skin.

The body's natural processes then work to break down and remove these marks. The appearance of any skin imperfection is then visibly reduced.

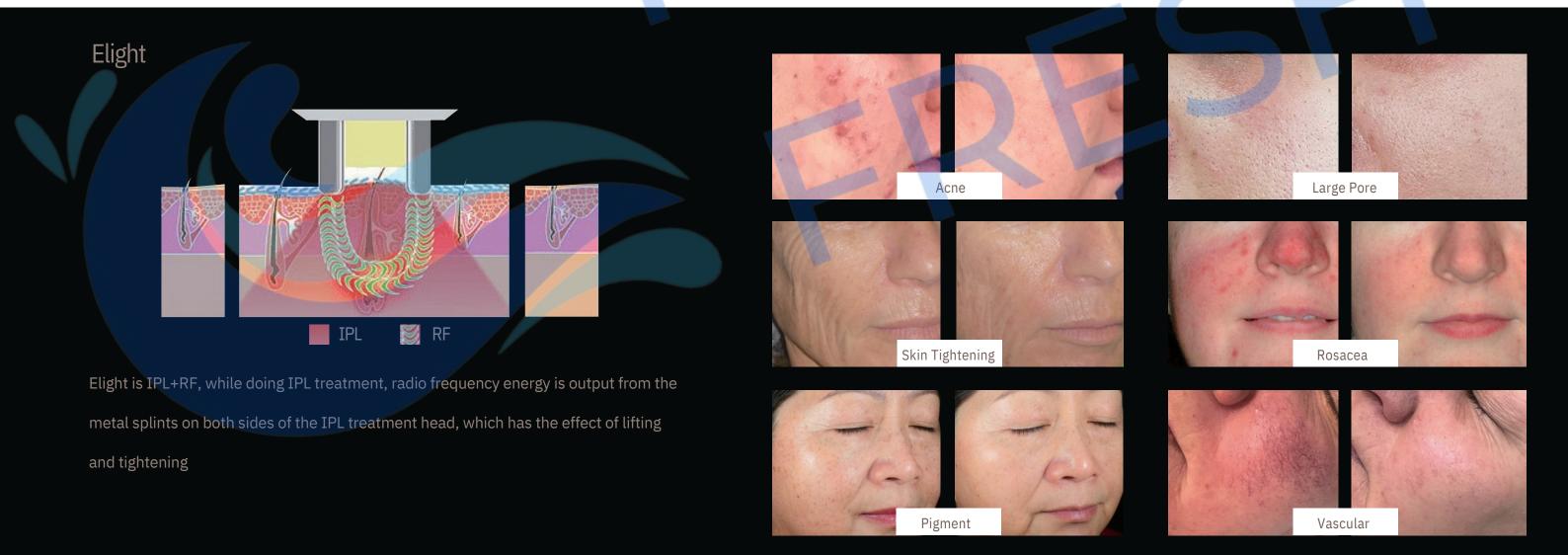





evel detection

15*50mm larger spot size



REAL-TIME INTELLIGENT MONITORING SYSTEM

20 intelligent detection functions, including water flow, water temperature, handle temperature, handle connection, water level

etc. Ensuring the stability of the machine and find out the possible problems of the machine in time




Ensuring handle surface is always cool and the treatment always comfortable.

At the same time, making sure that the machine works continuously for 24 hours

without damage by the high temperature.

Six different filters

One handle is equipped with 6 filters with different wavelengths, which can meet various treatment needs

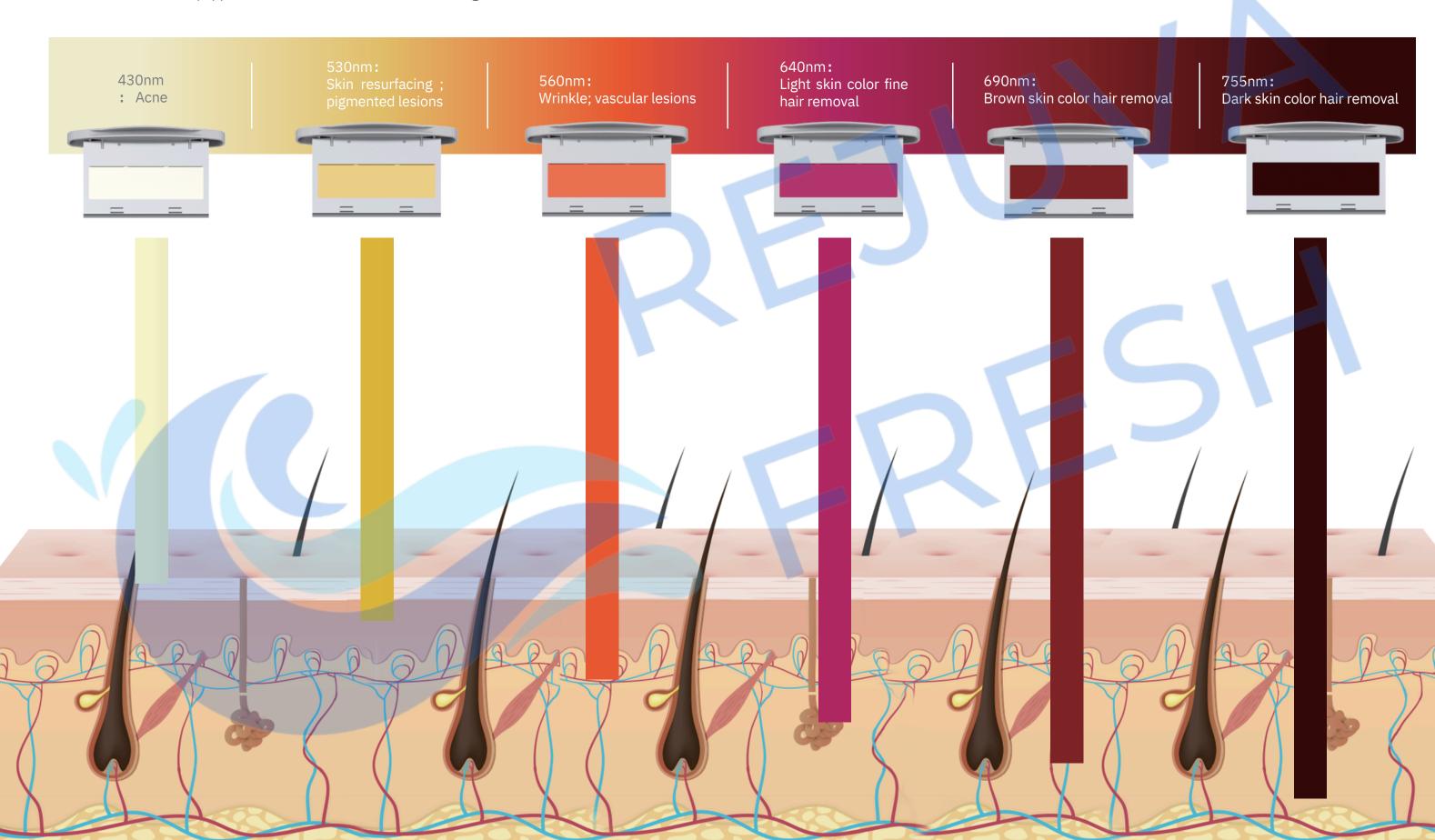



Handle Cooling

The handle system adopts 108W high-power refrigeration film, which is twice as much as other machines on the market, ensuring comfortable and safe treatment



Silver Reflector

The reflector is made of silver, which is not easy to be oxidized, has higher energy and longer life.



Sapphire Crystal

The handle crystal adopts sapphire crystal, which has better conduction of cold and can ensure the comfort of treatment







The UK xenon lamp, diameter 9MM, high energy and long life. Warranty 1 million shots.



Power	1500W
Wavelength	430nm,530nm,560nm,640nm,690nm,755nm
working modes	IPL, SHR, E-light
IPL energy	10-60j/cm2
Spot size: SR	15*50mm
SHR energy	3-15j/cm2
SHR frequency	1-10Hz
SHR pulse width	1-10ms
RF power	160W
RF energy	1-50j
Screen size	12.1 inches
Voltage	AC220V±10% 10A 50HZ, 110v±10% 10A 60HZ

