

KAHUNA  CHAIR®

LEG & FOOT MASSAGER

PAU

KAHUNA **MANUAL**





KAHUNA MASSAGE CHAIR®

BEST 3D ZERO GRAVITY MASSAGE CHAIR®









CONTENTS

Please read and follow all safety and operating instructions listed in this user manual prior to use. Basic precautions should always be used when operating an electrical appliance.

Safety Instructions	2
Operation Guide	6
Manual control	8
Names of parts	9
Maintenance	10
Troubleshooting	11
Product Specifications	12

SAFETY INSTRUCTION





Thank you for purchasing this PAU portable leg & foot Massager. Please read the operation manual thoroughly for the correct methods of usage and take special notice of the safety section before you start. This manual should be kept on hand for easy reference on safe and efficient operation of the massager.

 Warning	This sign indicates warning of causing possible death or hurt to body.
 Caution	This sign indicates caution of causing possible pain or damage to body.
	This sign indicates forbidden behavior.
	This sign indicates must follow instruction to operate.

- The following should be kept for any further references.

SAFETY INSTRUCTION

Thank you for purchasing this PAU portable leg & foot Massager. Please read the operation manual thoroughly for the correct methods of usage and take special notice of the safety section before you start. This manual should be kept on hand for easy reference on safe and efficient operation of the massager.

 Warning	This sign indicates warning of causing possible death or hurt to body.
 Caution	This sign indicates caution of causing possible pain or damage to body.
	This sign indicates forbidden behavior.
	This sign indicates must follow instruction to operate.

- The following should be kept for any further references.



Warning



- To ensure safe and correct use of the massager, do not operate it without carefully reading this operation manual.
 - Close supervision is necessary when the product is used by any person who is an invalid, handicapped, or with any form of physical, sensory or mental incapability and with limited experience and knowledge of the product (including children).
 - This product is not intended for children, ensure that they do not play with the product.
 - There are heat generation components inside, the surface of the product could be hot, so the people who is not sensitive to heat should be more careful when using this product.
 - If the wire is damaged, to avoid danger, one should ask manufacturer or authorized maintenance institution/team to change the wire.
 - The massager should be kept away from sharp items. Do not use the massager when the fabric cover is torn in order to prevent toes be drawn into the machine causing injuries. One should ask manufacturer or authorized maintenance institution/team to change or repair before use again.
 - Please unplug the machine before cleaning or after finish using.
 - Do not cover the massager during operation. This may result in a rise in motor temperature and subsequent malfunction or damage.
 - Stop using immediately should you experience any discomfort.
 - When the machine malfunctions, please send it back to manufacturer or authorized maintenance department to repair before you use it again.
 - Please do not use the spare parts unless it's recommended by original manufacturer.
 - Never drop or insert any object into any opening of the product.
 - Do not operate the massager with wet hands.
 - Do not expose the control switch to water or any liquids.
 - To disconnect, switch off all controls and remove plug from electrical outlet.
 - This product is not intended for self-diagnosis or self-treatment for conditions that should be managed by a qualified health care provider.
-
- Patients that are under treatment or post-surgery patients.
 - Patients who suffer from high blood pressure, cardiac disease, cerebral arteriosclerosis or cerebral thrombosis.
 - Patients who suffer from malignant tumor
 - Pregnant women or women during menstruation.
 - Patients who suffer from osteoporosis.
 - Person with cardiac pacemakers or have implantable medical device.
 - Patients who suffer from acute diseases.

SAFETY INSTRUCTION



Warning



- Do not use the massager in an excessively dusty or humid environment such as a bathroom since this may result in malfunction, electric shock or damage to the product.
- Do not use the massager in a room with temperature of 40°C (104°F) or higher. Do not expose the massager to heaters, stoves or direct sunlight.
- Do not use the product when the ambient temperature is below 5°C (41°F) for the thickening of the lubricant could increase abrasive wear of the gears, which will shorten the operational life span of the product.
- Recommended not use the product outdoor or place it under direct sunlight, it might cause damage to the paintcoat of the product.
- When moving the machine from low temp environment to relative high temp to use, the inside metal parts might malfunction or even damaged because the water vaporization. We suggest placing the machine under room temperature for more than one hour before power it on.



Caution



- Use this massager 1-2 times a day, 5-10 minutes each time would effectively eliminate muscle fatigue, improve blood circulation. Use this massager for 5- 15minutes before sleep, would effectively improve sleep quality.
- Do not fall asleep during using.
- Please do not use this product when your are inebriated or when you are not feeling well.
- Please do not massage over 30mins, warm foot do not over 2 hours.
- If the cloth cover is broken, please don't use it.
- Do not use this product one hour before or after meal, it might cause opposite effect.
- When moving the machine from low temp environment to relative high temp to use, the inside metal parts might malfunction or even damaged because the water vaporization. We suggest placing the machine under room temperature for more than one hour before power it on.



Caution



- Place the massager on a clean floor and check if the machine's foot cloth cover are loose or broken
- First plug the end of the adapter into the massager outlet, the other end connected to the power cord, and then plug the power cord into the mains socket, the massager power indicator light starts flashing, the massager into the standby state, people sitting in the rear of the massager, put the feet into the product foot cloth cover (note not stand to use).
- According to the needs of massage, you can choose to flip part fold or unfold, flip folding can massage the foot, flip unfolding can massage the foot and calf at the same time.
- The upper operation panel is invalid when the flip is folding, the lower operation panel is invalid when the flip is unfolding.
If you want to change the flip state in the using, you need to press the
- ON/OFF button, until the product stops working and then lift the flip.
- In the working state press the power button or the end of the timer time, the massager will stop working, the power indicator flashes, enter the standby mode.




Warning : The massager is finished using in the clamshell state(or stop in the middle), should not be flip folding immediately, the massager should be allowed to stand for one minute, until the gas is completely released and then folded.

After use
Tip

- After use, a cup of warm water is better.
- To prevent children accidentally start the massager caused an accident or affected parts, please unplug the power cord from the electrical outlet when you don't use the massager, and proper storing.

| Features |

- This product is developed by AJX,LLC product R&D team. The design is elegant and with assorted colors. Classy lacquer surface radiates fashionable qualities. It is also light-weighted and portable.
- Ergonomic design, enable you to enjoy professional leg & foot massage at any time. The ergonomic design, unique styling and acupoint distribution all fit to the human leg & foot structure. Just like bringing a professional masseuse at home, you will enjoy ultimate relaxed massage on your leg & foot while watching TV program, or before sleep.
- Double safety protection makes you feel secure when use this product. When the product is overheating or the current is overload, the circuit will break off automatically. This double protection is safe and reliable, ensure you to enjoy the massage joyfully and carefree.

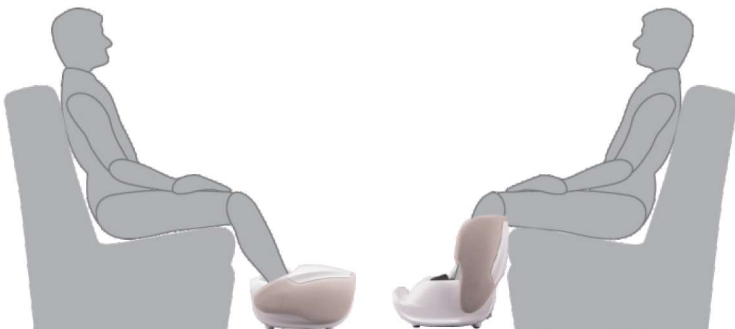

- Combined various massage ways of kneading, scrapping, shiatsu, airbag press, warm foot etc.
- Unique calf massage and care function.
- Nano removable washable foot cloth cover, keep you clean at all time.
- Equipped with double sensitive touch panel.





ON/OFF

Press the ON/OFF button on the front panel of the massager to start the massager and enter the automatic working state.



HEATING

Press the HEATING button to turn off or turn on the foot heating function.



STRENGTH

Press the STRENGTH button, you can choose the appropriate air massage strength, the button indicator light according to the weak, medium and strong, synchronized display for the blue, blue and red, red; long touch the button for one second, you can turn off the air massage function.



SCRAPING

Press the SCRAPING button, you can turn off or turn on the foot massage function; long touch the button for one second, you can turn off the foot massage function and switch to a soft foot air care experience.



VIBRATION


Press the VIBRATION button, you can turn off or turn on the foot vibration function and select the vibration education experience for different strength modes.




LEGS

Press the LEGS button, you can turn off or turn on the leg foot massage function. Long touch the button for one second, you can turn off the foot massage function and switch to soft leg foot air care experience(This function is not available when the clamshell is folded).

NAMES OF PARTS

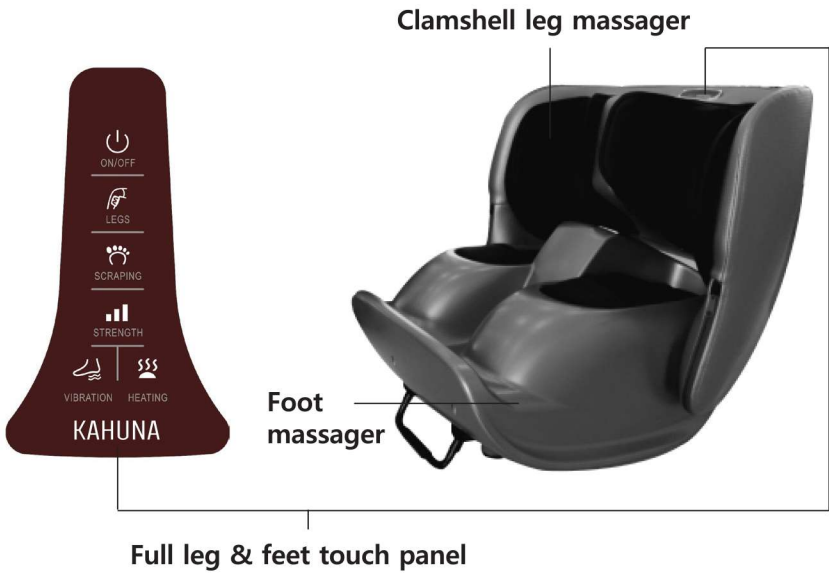


MAINTENANCE

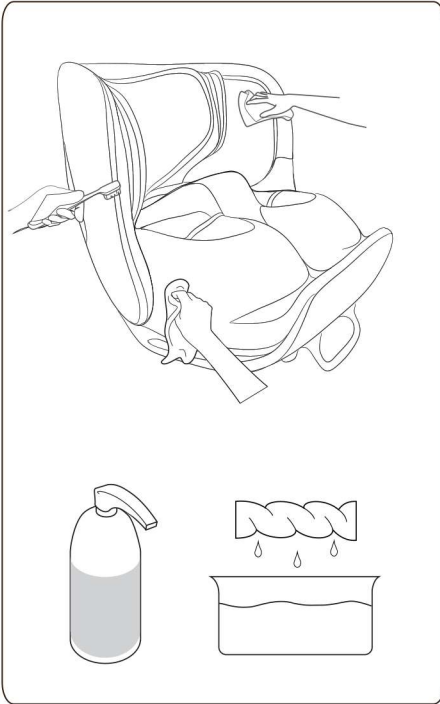

| Use dry cloth to clean |

- Before cleaning, make sure the power switch is off and unplug the power cord.
- Keep the massager clean and away from high temperature and humidity place. If there is a long time not to use the massager, please use the dust cover trap to prevent dust from falling into.

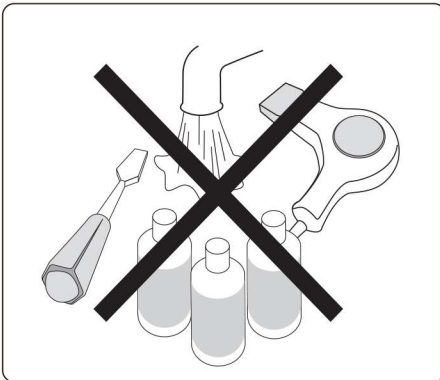

- Please use a dry cloth to clean the surface of the leg & foot massager. Do not use corrosive cleaners, such as benzene or thinner to clean the massager.



This leg & foot massager has a temperature protection device, it will automatically stop working when the temperature exceeds the safe temperature, please unplug the power to allow its natural cooling and then use.

TROUBLE SHOOTING

If you have problems using this leg & foot massager, please refer to the following tips. If this massager is still in a failed state. Please contact our service department.

855.646.9500 / 1.714.388.2625, service@kahunachair.com

Problem	Solution
<ul style="list-style-type: none">• Not power on• The function button is not turned on	<ul style="list-style-type: none">• Connect the power supply or check whether the power supply is connected properly• Turn on the ON/OFF button
<ul style="list-style-type: none">• Power off• Massage time more than 15 minutes• Continuous using time is too long, massager overheating protection• Force too much, massage ball stalled	<ul style="list-style-type: none">• Check if the power connector is connected• Turn on the ON/OFF button again• Turn on the ON/OFF button again after cooling the massager• Relax the feet, reduce the pressure load on the massager



KAHUNA CHAIR PRODUCT WARRANTY REGISTRATION
<https://www.kahunachair.com/services/>

PRODUCT SPECIFICATIONS

Description	Specifications
Product Name	Kahuna Leg & Foot Massager
Model No.	PAU
Rated Voltage	DC24V
Rated Power Consumption	55W
Timing time	15 minutes
Dimension (LxWxH) Inches	19.7 x 18.9 x 13 inches
Weight (LBS)	Net Weight: 24.25 lbs



Product Leg & Foot Massager

Model PAU


Distributed by AJX Distribution LLC

11125 Knott Ave, Suite D

Cypress, CA 90630

Tel 1-855-646-9500

Email INFO@KAHUNACHAIR.COM

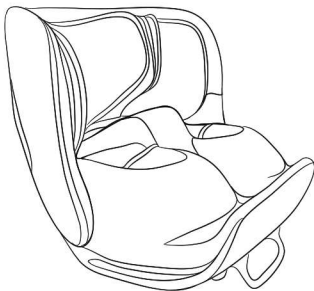

You must register warranty at www.kahunachair.com/services/

INSTRUCTION MANUAL
MODEL : PAU

Plaease read this user manual thoroughly
prior to operating the leg & foot massager.
Store this booklet for future reference.