


# EB03 High Frequency Surgical Unit

## Features

- Multiple working modes, including 4 monopolar cut modes (pure cut, blend cut 1, blend cut 2, and blend cut 3) and 2 monopolar coagulation modes (soft coagulation and point coagulation) and 1 bipolar coagulation mode
- High power cut mode up to 350W , applicable for cutting thick body fat or vaporized cutting operation
- Automatic detection for high frequency current leakage, avoiding burns on patients
- Real-time monitoring the stick area between return electrode and skin, automatically cutting off output in case stick area drops to dangerous level
- Tissue responding technology, keeping consistency of cutting quality regardless difference in tissues, electrodes or cutting range
- Memory of setting mode and output power, automatically restored for the next operation
- Audio and visual alarms for errors of connection, operation, and the parts of main units, with error codes display
- Subdivided monopolar coagulation modes - “soft coagulation” for precise control of coagulation and “point coagulation” for quick hemostasia
- Current flowing just between two electrodes of forceps in bipolar coagulation mode, without return electrode, applicable in operations of sensitive tissues like brain and oviduct
- Widely applicable in general operations, like gynecological, facial, orthopedic, urological, and dermatological surgeries, etc.

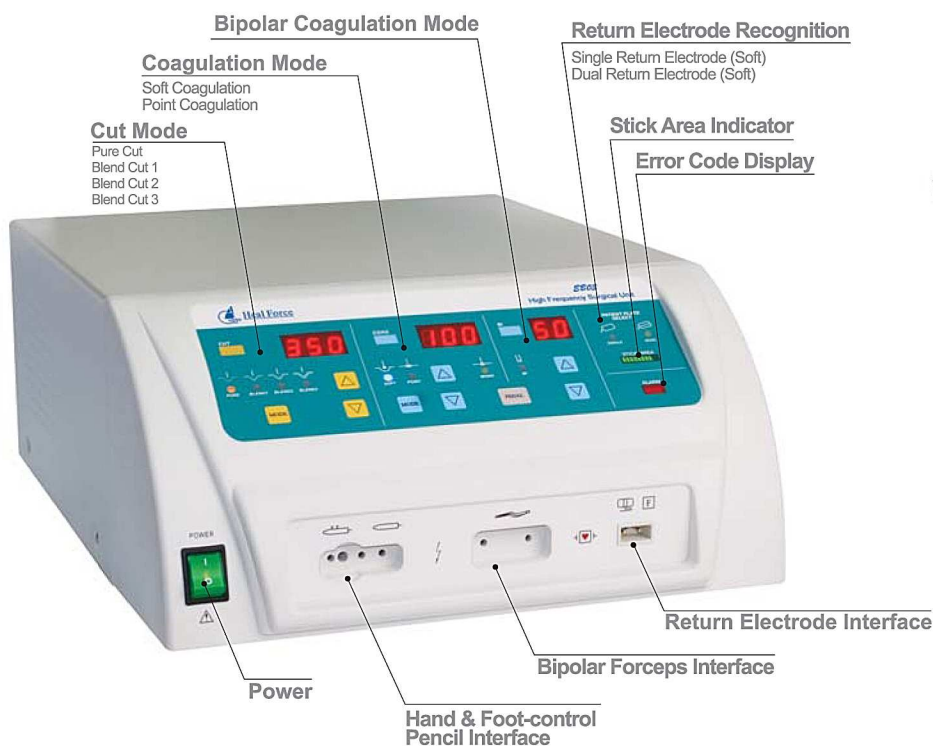

EB03 (Cart for Option Only)

# EB03 High Frequency Surgical Unit

## EB03 Technical Specifications

Working Modes	Output Power	Max. Open Voltage (Vp-p)	Crest Factor	Rated Impedance
Pure Cut	350W	1800V	1.5	500Ω
Blend Cut 1	300W	2700V	1.9	500Ω
Blend Cut 2	200W	2400V	2.2	500Ω
Blend Cut 3	150W	2900V	3.2	500Ω
Soft Coagulation	100W	3600V	4.0	500Ω
Point Coagulation	100W	4600V	3.9	500Ω
Bipolar Coagulation	50W	240V	—	100Ω
<b>General Parameters</b>				
Power Supply	220V±10% 50/60Hz, optional for 110V/60Hz			
Dimensions (W×D×H)	395mm×510mm×205mm			
Net Weight	18Kg			

## Standard Accessories\*



Hand-Control Pencil



Foot Switch



Dual Return Electrode (soft)



Single Return Electrode (soft)



Single Return Electrode with Cable (hard)



Foot-Control Pencil



Bipolar Coagulation Forceps (for EB03)




Connecting Cable of Bipolar Coagulation Forceps (for EB03)



Straight Knife Electrode



Standard Knife

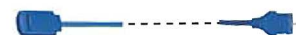

Connecting Cable of Soft Return Electrode

\*Please refer to ordering information for complete accessory list.



Heal Force Bio-Meditech Holdings Group  
Nison Instrument (Shanghai) Limited

International Sales Office: Floor 15, No.2 Hua Shan Road, Shanghai 200040, P.R. China

Tel: +86 21 62728646

E-mail: [export@healforce.com](mailto:export@healforce.com)

Fax: +86 21 62710529

Website: [www.healforce.com](http://www.healforce.com)

Information is subject to change and/or updating without notice.  
Copyright © 2010 Heal Force. All Rights Reserved.