



Zetuvit® E highly absorbent wound dressing pad

For treatment of highly exuding wounds



Hypoallergenic fibres permeable to air and water

- Prevents irritation of surrounding skin
- Enables skin's natural functions to continue

Filling of fluff pulp and fluid dispersing tissue

- Provide high level of absorbency

Blue, hydrophobic backing

- Provides a barrier against contamination



PAUL HARTMANN LTD
 Heywood Distribution Park
 Pilsworth Road
 Heywood
 Lancashire
 OL10 2TT
 Tel: 01706 363200
 Fax: 01706 363201


Web: www.hartmann.co.uk
E-mail: info@uk.hartmann.info

Zetuvit® E highly absorbent wound dressing pad

For treatment of highly exuding wounds

Zetuvit® E is a non-irritant absorbent dressing pad. The covering layer consists of hydrophobic, non woven polypropelene material. The dressing itself is made up of a highly absorbent fluff filling, fluid dispersing tissue and a hydrophobic tissue backing which acts as a barrier to protect the wound from contamination.

The dressing should be worn with the blue backing upper most, facing away from the wound and is available sterile and non sterile.

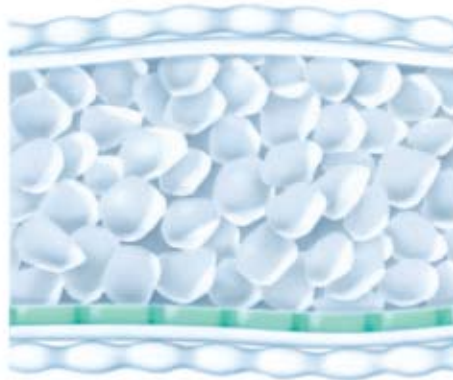

Indication

Zetuvit® E can be used to treat heavily exuding wounds in all areas of medicine and also for cushioning and protecting wounds.



Product	Size	Pack	Hartmann Code	NHS Code	PIP Code
Zetuvit® E sterile	10cm x 10cm	25 pcs	413770	EJA025	322 7618
	10cm x 20cm	25 pcs	413771	EJA026	322 7584
	20cm x 20cm	15 pcs	413774	EJA027	322 7592
	20cm x 40cm	10 pcs	413776	EJA028	322 7600
Zetuvit® E non sterile	10cm x 10cm	50 pcs	413860	EJA021	328 6085
	10cm x 20cm	50 pcs	413861	EJA022	328 6093
	20cm x 20cm	50 pcs	413864	EJA023	328 6101
	20cm x 40cm	30 pcs	413866	EJA024	328 6119