


# CASCADE 3171 WALL MOUNT HANDPIECE DELIVERY SYSTEM

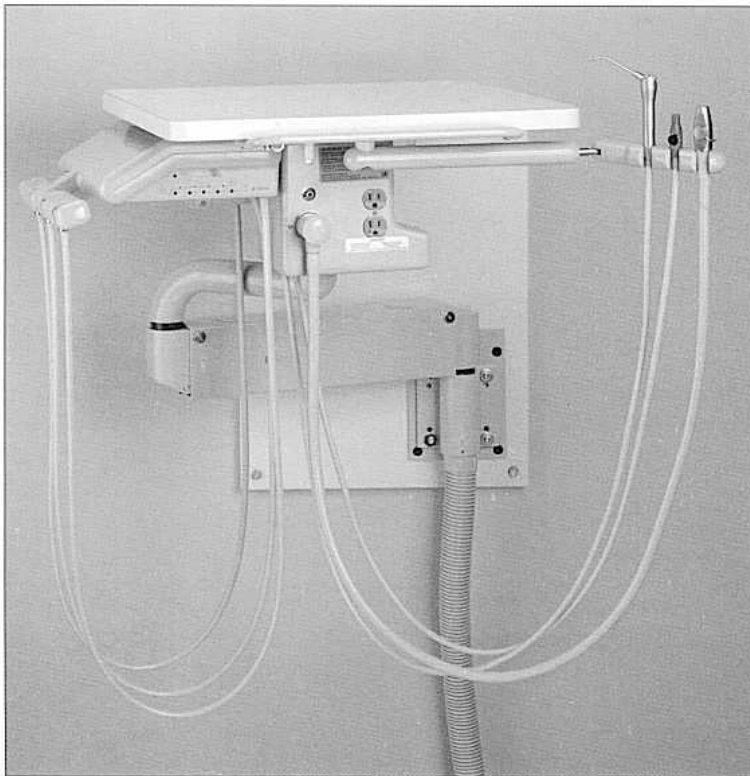

Cascade 3171 Wall Mount


## STANDARD PRODUCT FEATURES

- **Cascade Handpiece Delivery System** — Control for three handpieces with the proven reliable Century Plus control blocks and unitized holders.
- **Cascade Assistant's Instrumentation** — Solids collector, Autoclavable QD HVE, an Autoclavable syringe and a Saliva Ejector.
- **8 Conductor Cable Prewiring** — Allows economical addition of electrical accessories.
- **Foot Control** — Wet/dry disc foot control with 7' (2134 mm) of tubing.
- **Air Actuated Brake** — Toggle-activated brake holds the system where you want it.
- **Work Surface with Laminate Finish** — 19" by 22" (483 x 559 mm) work surface for instrument tray; laminated for easy cleaning.
- **Left/right Convertible** — Handpiece controls and assistant's instruments can be pivoted for left/right conversion.
- **Flex Arm Assembly** — 8.5" (216 mm) vertical travel; 23" (584 mm) reach.
- **Medium Stainless Steel Floor Box** — Includes air and water filters, pressure regulators, and shutoff valves.
- **Umbilical** — 7' (2134 mm) for wall mount and 1' (305 mm) for Preference mount.
- **Backboard** — 18" x 23.5" (457 x 597 mm) for mounting on wall (Wall mount only).
- **Duplex** — 120 volt duplex outlet for electrical accessories.
- **NOTE:** Chair touch pad and handle not applicable to this unit.




**72**

Wall Mount 15-0450-00  
with Crystal Laminate  
Backboard only

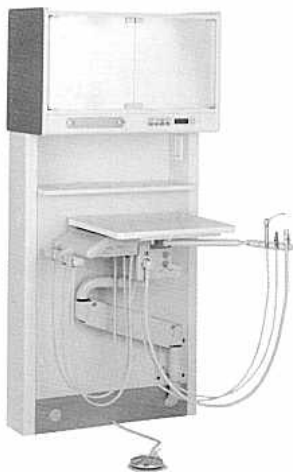

Cascade 3171 shown on  
Preference Model 5561

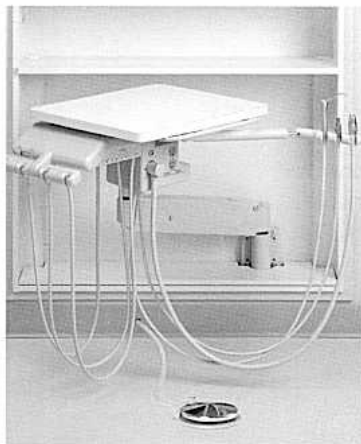

Cascade 3171 shown on  
Preference Model 5661 42"

