

W E T - R I N G  
V A C U U M S



# C L A S S I C

## Midmark ClassicSeries® Wet-Ring Dental Vacuums

Midmark ClassicSeries® Wet-Ring Vacuums bring together the best features and most proven technologies into one uniform offering. These vacuums satisfy those who have grown to appreciate the high standards created by Apollo and Matrix brand wet-ring vacuums, offering:

- High output, proven performance
- Durable aluminum bronze construction
- Exclusive internal water recycling
- Water recycling by-pass
- Sound reducing mufflers standard
- NFPA 99C level 3 compliant
- UL/CUL classified
- 100% long term burn-in testing

When it comes to selecting a wet-ring dental vacuum, Midmark is the natural choice. ClassicSeries® vacuums evolved from the Apollo and Matrix brand wet-ring vacuums, both of which have provided years of proven vacuum performance to thousands of satisfied dental professionals.


ClassicSeries® Wet-Ring Dental Vacuums come standard with many of the same value-driven features as their predecessors. Standard internal water recycling is an exclusive design that contributes to minimizing fresh water consumption – saving over three times the amount of fresh water consumed vs. non-water recycling vacuums. In addition, a valuable water recycler by-pass valve is utilized to help protect against the potentially harmful effects of water mineral deposits.

Aluminum bronze, rated harder than stainless steel (Brinell Hardness Scale), is used exclusively to ensure ClassicSeries® wet-ring vacuums provide many years of durable operation. Stainless steel motor shafts and heavy-duty 1HP/2HP motors assure extra strength and reliability. Dual voltage (1.25 HP vacuums) allows for versatility and ease of installation in all mechanical rooms.


Along with the best features brought forward from Apollo and Matrix, several new features elevate the ClassicSeries® above all competitors. The component layout and design has been

simplified and re-organized to allow for extremely easy, convenient installation and service. From the strategically-designed left-hand/right-hand twin vacuum installation layout to the convenient all-access front electrical panels, ClassicSeries® Wet-Ring Dental Vacuums definitely raise the bar in fit, finish and serviceability.

You can depend on Midmark ClassicSeries® Wet-Ring Dental Vacuums. Proven. Reliable. Trusted.

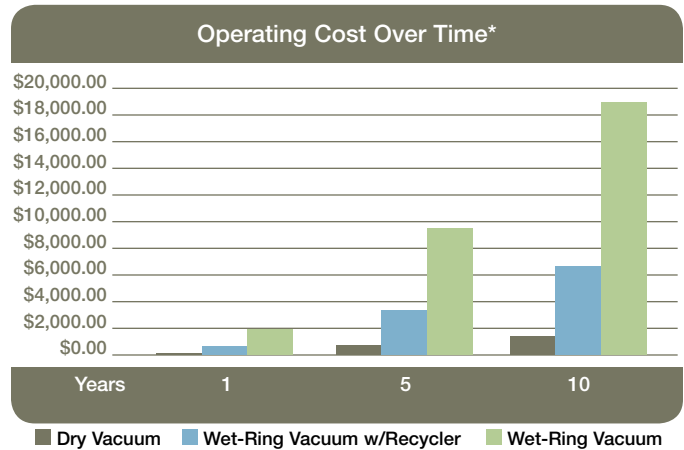

\* 2,000 hours/yr; 57,600 gal/yr (Wet-Ring); 17,280 gal/yr (Wet-Ring w/Recycler)

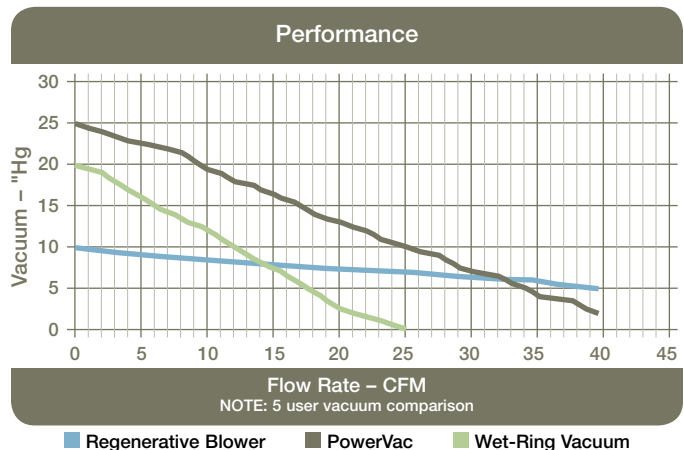



Model CV5R

Model CV3R

Model CV3



Redesigned intake manifold allows for smooth air flow and greater efficiency.



Separate twin model on/off switches and hour meters provide operational versatility.



Quick connect hose connections allow for tool-free serviceability. Water recycler by-pass valve provides extra protection against water mineral deposits.

## Specifications

Model	Max Users*	No. Vacuums	Total HP	Vacuum ("Hg)	Voltage	Amps (Per Vac)	Water Recycling	Overall Dimensions H x W x D	Shipping Weight (lbs.)
CV3	3	1	1.25	up to 12"	115/208-230	15/7.5	No	15" x 12" x 13"	70
CV3R	3	1	1.25	up to 12"	115/208-230	15/7.5	Yes	15" x 12" x 13"	70
CV5	5	1	2	up to 12"	208-230	7.5	No	17" x 12" x 13"	76
CV5R	5	1	2	up to 12"	208-230	7.5	Yes	17" x 12" x 13"	76
CV6	6	2	2.5	up to 12"	115/208-230	15/7.5	No	18.5" x 25.5" x 20"	188
CV6R	6	2	2.5	up to 12"	115/208-230	15/7.5	Yes	18.5" x 25.5" x 20"	188
CV10	10	2	4	up to 12"	208-230	7.5	No	18.5" x 25.5" x 20"	208
CV10R	10	2	4	up to 12"	208-230	7.5	Yes	18.5" x 25.5" x 20"	208

\* HVE = 1 vacuum user; SE = .5 vacuum user; N<sub>2</sub>O Scavenger = .5 vacuum user

## Warranty

- 5 years/10,000-hour product warranty

Complies with NFPA 99C level 3



Precision Clense Plus  
Vacuum System Cleaner



Keep your vacuum thoroughly clean, naturally disinfected and deodorized

Dispos-A-Bowl



The safe, smart way to dispose of used vacuum filters and waste.

Air Water Separators



Air Water Separator/  
Muffler for single and twin vacuums.

Midmark is an ISO 13485 and ISO 9001 Certified Company. For more information or a demonstration, contact your Midmark dealer or call: 1-800-MIDMARK Fax: 1-877-725-6495 Outside the U.S.A. call: 1-937-526-3662 Fax: 1-937-526-8214 or visit [midmark.com](http://midmark.com)

© 2007 Midmark Corporation  
60 Vista Drive P.O. Box 286 Versailles, Ohio 45380-0286  
Products subject to improvement changes without notice  
Litho in U.S.A. 007-0652-00 Rev. C1 (7/10)



Because we care.