


# Ambu® aScope™ 4 Cysto

The single-use flexible cystoscope  
Item no. 600001000

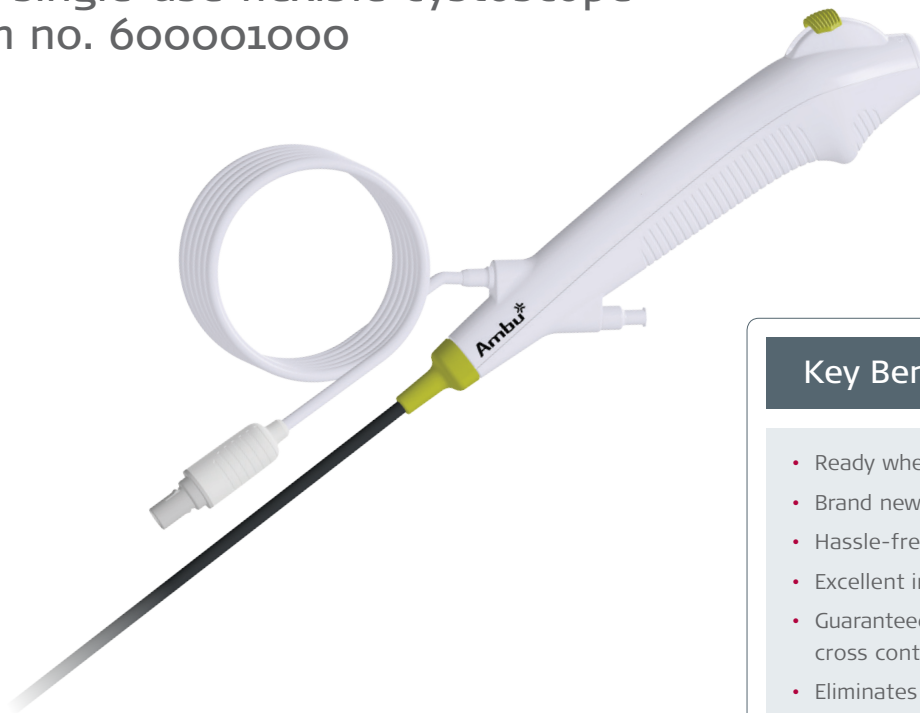

## Key Benefits

- Ready when you are
- Brand new every time
- Hassle-free portable solution
- Excellent imaging & manoeuvrability
- Guaranteed sterility with no risk of cross contamination
- Eliminates reprocessing and costly repairs

## Ambu® aScope™ 4 Cysto

**aScope 4 Cysto is a single-use flexible endoscope solution that allows you to take control of your schedule and be more productive - without compromising on the quality of your work.**

It is always available and portable making it easy to manage your schedule and deal with inhouse consult procedures.

With the aScope 4 Cysto you can count on clear, sharp images, making it easy to identify anatomical structures. Bending angles of 210°/120° enable you to manoeuvre and navigate smoothly and allow for complete retroflexion to inspect the bladder neck, with or without forceps inserted. It offers consistent quality without any deterioration of image or bending quality because you get a brand-new cystoscope for every procedure.

The sterile aScope 4 Cysto eliminates need for reprocessing, costly repairs, and the risk of cross-contamination. As a result, aScope 4 Cysto simplifies workflow, frees up resources, and enables you to treat more patients and you can assure your patients that you will be performing the procedure with a single-use sterile cystoscope opened straight from the package just for them.

The aScope 4 Cysto makes it easy to plan and manage your day. From the outpatient clinic to inpatient consult procedures, you can take the lightweight single-use cystoscope and portable monitor with you under your arm. And when you're done, you simply dispose of it; there's no more cleaning hassle.

# Specifications

Optical system	
Field of View	120°
Depth of Field	3 - 100 mm
Illumination method	LED
Insertion portion	
Bending angle*	210° up, 120° down
Insertion cord diameter	16.2 Fr / 5.4 mm (0.21")
Distal end diameter	16.2 Fr / 5.4 mm (0.21")
Maximum diameter of insertion portion	18 Fr / 6.0 mm (0.24")
Outer diameter (mm)	Max. 6.0 mm
Inner diameter (mm)	Min. 2.2 mm
Working length	390 mm (15.4")
Working channel	
Minimum instrument channel width**	6.6 Fr / 2.2 mm (0.086")
Operating environment, storage & transportation	
Temperature - transportation	10 ~ 40° C
Temperature - operating environment	10 ~ 40° C
Temperature - storage***	10 ~ 25° C
Relative humidity	30 ~ 85%
Atmospheric pressure	80 ~ 109 kPa
Altitude	≤2000 m
Sterilisation	
Method of sterilisation	ETO




aScope 4 Cysto offers consistent quality every time



Bending angles of 210°/120° enable complete retroflexion.

## Order specifications

Product	Item no.	Quantity
Ambu® aScope™ 4 Cysto	600001000	5 pcs

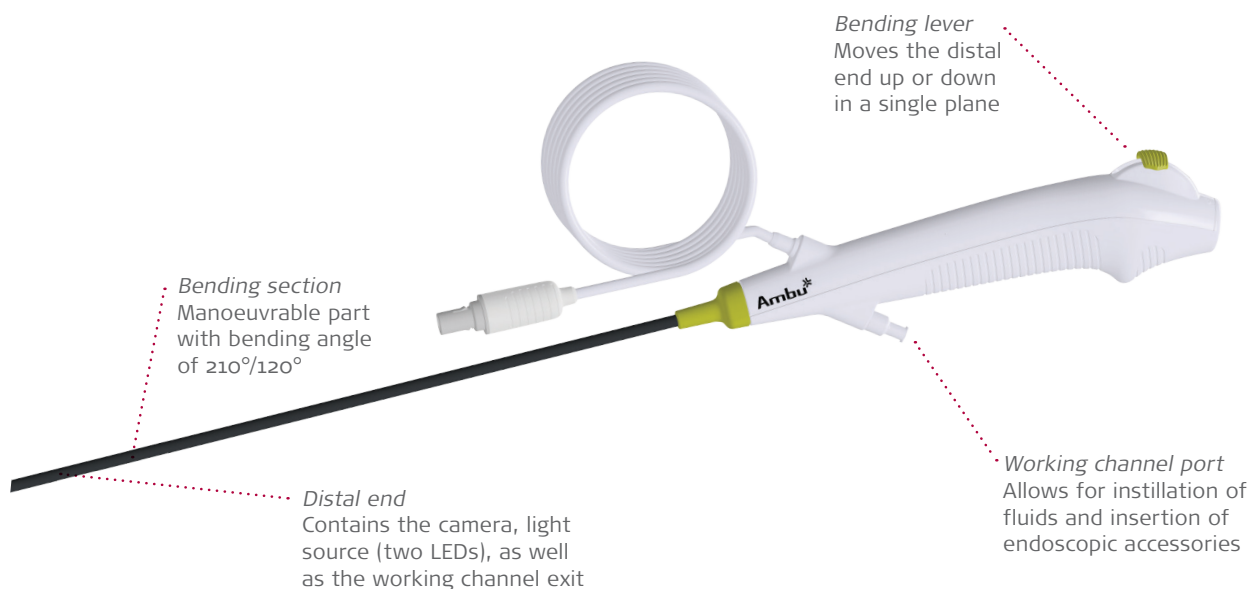

**Ambu Australia Pty Ltd**  
Unit 5, 2 Daydream St  
Warriewood, NSW 2102  
Australia  
Tel. 1300 233 118  
aussales@ambu.com  
www.ambuaustralia.com.au

\* Please be aware that the bending angle can be affected if the insertion cord is not kept straight or have inserted endoscopic accessories.  
\*\* There is no guarantee that endoscopic accessories selected solely using this minimum working channel width will be compatible in combination.  
\*\*\* Storage under higher temperatures may impact shelf life.