

COCOON

Convective patient warming system



THE COCOON CONVECTIVE PATIENT WARMING SYSTEM IS A MICRO PROCESSOR CONTROLLED COCOON CWS 4000 WARMING MACHINE WITH A FULL RANGE OF COCOON WARMING BLANKETS TO SUIT VARIOUS APPLICATIONS IN PATIENT CARE FOR MAINTAINING PATIENT TEMPERATURE CONTROL (NORMOTHERMIA).

The development of the micro pore blanket system for the delivery of air gives a constant air flow, thereby providing a uniform distribution of warm filtered air across the patient. The universal inlet port ties ensure a secure hose connection to the blanket, enabling our blanket's to deliver a warm soothing cocooning effect to your patients.

The Cocoon CWS 4000 is light and compact and has been developed to be user friendly for both Health Clinicians and Biomedical Engineers. We take care to exceed essential regulatory standards.



TEMPERATURE MANAGEMENT SYSTEM

Convective Warming Machine

Microprocessor controlled, filter status indicator, digital temperature sensor system, washable protective hose cover, built-in safety monitoring system, retractable bed hook brackets, external testing port.



Electrical

Supply Voltage	230-240 Vac
Supply Frequency	50 Hz

Classification

Applied Part Classification/Device Classification	Type BF
TGA - Medical devices & MDD 93/42/EEC	Class II b device

Performance

Temperature settings measure the average air temperature at the end of the delivery hose.	Ambient, 34°C, 37°C, 40°C, 43°C or 46°C
Temperature accuracy of delivery air, except for the ambient temperature setting.	+2°C with delivery hose cover


Physical

Dimensions	29cm x 22cm x 40cm
Weight	6.5 kg
Filter Type/Life	Bacteriological HEPA filter/1000 hours 0.2 micron



Care Essentials Pty Ltd
25 Slevin Street, North Geelong, Victoria, Australia, 3216
Phone: +61 3 5277 1455 Fax: +61 3 5277 1411
Email: queries@careessentials.com.au

CE 0805

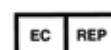






QPF 058 - Version: 1 - Supersedes: 0

Disposable Patient Warming Blankets





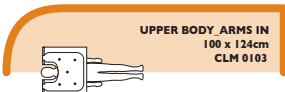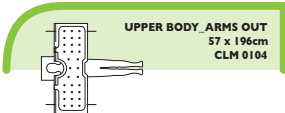

INTRAOPERATIVE BLANKETS

Cocoon patient warming provides effective temperature management with our full range of warming blankets. Our unique microscopic perforations of over 26 holes/cm² provide for even air distribution without the risk of residue.

<p>Model - CLM 0101 Full Body Size - 100 x 195cm Hose ports - 1</p>	<ul style="list-style-type: none"> • Ample dimensions to accommodate all patients. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	 <p>FULL BODY 100 x 195cm CLM 0101</p>
<p>Model - CLM 0102 Lower Body Size - 100 x 120cm Hose ports - 1</p>	<ul style="list-style-type: none"> • Generous amount of tape to secure blanket over abdomen. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	 <p>LOWER BODY 100 x 120cm CLM 0102</p>
<p>Model - CLM 0103 Upper Body - Arms In Size - 100 x 124cm Hose ports - 2</p>	<ul style="list-style-type: none"> • Generous amount of tape to secure blanket over abdomen. • Two hose inlets for operator convenience. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	 <p>UPPER BODY - ARMS IN 100 x 124cm CLM 0103</p>
<p>Model - CLM 0104 Upper Body - Arms Out Size - 57 x 196cm Hose ports - 2</p>	<ul style="list-style-type: none"> • Two hose inlets for operator convenience. • Arm board ties to secure blanket once in position. • Clear face shield to retain warmth around the head whilst providing a clear view. • Universal inlet port ties ensure the hose remains firmly in the blanket. • Ample tape (latex free) to secure blanket to torso. 	 <p>UPPER BODY - ARMS OUT 57 x 196cm CLM 0104</p>
<p>Model - CLM 0122 Adult Tube Size - 100 x 217cm Hose ports - 1</p>	<ul style="list-style-type: none"> • Full surgical access as patient lies between two tubular sections. • Warm air is directed both across and down the length of the patient. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	 <p>ADULT TUBE 100 x 217cm CLM 0122</p>

PAEDIATRIC BLANKETS


Cocoon patient warming provides effective temperature management for the infants and children in any medical environment.

<p>Model - CLM 0105 Paediatric - Neonatal Size - 87 x 89cm Hose ports - 2</p>	<ul style="list-style-type: none"> • Full surgical access as patient lies on top of quilted area. • Two hose inlets for operator convenience. • Universal inlet port ties ensure the hose remains firmly in the blanket. • 2 sheets of clear plastic (for additional warmth, if required). 	 <p>PAEDIATRIC NEONATAL 87 x 89cm CLM 0105</p>
<p>Model - CLM 0106 Paediatric - Tube Size - 117 x 78cm Hose ports - 1</p>	<ul style="list-style-type: none"> • Full surgical access as patient lies between two tubular sections. • Warm air is directed both across and down the length of the patient. • 1 sheet of clear plastic (for additional warmth, if required). • Universal inlet port ties ensure the hose remains firmly in the blanket. 	 <p>PAEDIATRIC TUBE 117 x 78cm CLM 0106</p>
<p>Model - CLM 0107 Paediatric - Child Size - 130 x 92cm Hose ports - 2</p>	<ul style="list-style-type: none"> • Full surgical access as patient lies on top of quilted area. • Two hose inlets for operator convenience. • Universal inlet port ties ensure the hose remains firmly in the blanket. • 2 sheets of clear plastic (for additional warmth, if required). 	 <p>PAEDIATRIC CHILD 130 x 92cm CLM 0107</p>
<p>Model - CLM 0112 Multi Purpose Air Sock Size - 100 x 12cm Hose ports - 1</p>	<ul style="list-style-type: none"> • Full surgical access as patient lies next to the tubular air Sock. • By placing a covering sheet over the Air Sock warm air is retained and directed both across and down the length of the patient. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	 <p>MULTIPURPOSE HOSE SOCK 100 x 12cm CLM 0112</p>



Care Essentials Pty Ltd
25 Slevin Street, North Geelong, Victoria, Australia, 3216
Phone: +61 3 5277 1455 Fax: +61 3 5277 1411
Email: queries@careessentials.com.au

CE 0805

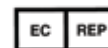

QPF 058 - Version: 1 - Supersedes: 0

Disposable Patient Warming Blankets



SPECIALITY & CARDIAC BLANKETS

Cocoon provides speciality blankets to suit every procedure.

<p>Model - CLM 0108 Cardiac Access Size - 100 x 195cm Hose port - 1</p>	<ul style="list-style-type: none"> • 40cm slit from the neck area down the length of the blanket over the chest and abdomen allowing for easy access. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	<p>CARDIAC ACCESS 100 x 195cm CLM 0108</p>
<p>Model - CLM 0118 Abdominal Access Size - 190 x 180cm Hose ports - 3</p>	<ul style="list-style-type: none"> • Cut out section over the chest and abdomen allowing for easy access. • Ample dimensions to accommodate all patients and positions. • Three hose inlets for operator convenience. • Arm board ties to secure blanket once in position. • Ample tape (latex free) around fenestration to secure blanket to torso. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	<p>ABDOMINAL ACCESS 190 x 180cm CLM 0118</p>
<p>Model - CLM 0121 Multi Access Size - 100 x 195cm Hose port - 1</p>	<ul style="list-style-type: none"> • Flexible design to allow surgical access by folding back sections of the blanket. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	<p>MULTI ACCESS 100 x 195cm CLM 0121</p>
<p>Model - CLM 0124 (New) Surgical Access - Arms In Size - 100 x 195cm Hose port - 1</p>	<ul style="list-style-type: none"> • Larger cut out section over the chest and abdomen allowing for easy access. • Ample dimensions to accommodate all patients and positions. • Ample tape (latex free) around fenestration to secure blanket to torso. • Universal inlet port ties ensure the hose remains firmly in the blanket. 	<p>SURGICAL ACCESS 100 x 195cm CLM 0124</p>

COCOON - TROLLEY STAND

Model C 0005 A	670 mm (height)
Dimensions pole	25mm Dia
Castors	5, anti-static, latex free
Weight	2.25 kg



COCOON - TROLLEY CART

Model C 0005 B	800 mm (height)
Dimensions pole	25mm Dia
Castors	4, anti-static casters with brakes
Weight	3 kg



Care Essentials Pty Ltd
25 Slevin Street, North Geelong, Victoria, Australia, 3216
Phone: +61 3 5277 1455 Fax: +61 3 5277 1411
Email: queries@careessentials.com.au

