



Bio-Therapeutic

TECHNOLOGY ACCELERATED SKIN CARE™


bt-micro® FUSION EXFOLIATION SYSTEM

QUICK START GUIDE



BT-MICRO FUSION
MODEL MF-69

Read operation manual prior to using the equipment.
Do not use the unit if any contraindications are present.

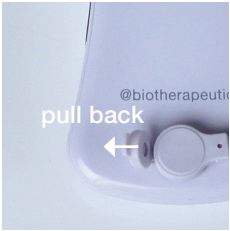

CHARGING THE UNIT

Pull back charging port cover and plug round end of charging cord into the back of the unit. You will feel the magnet engage when cord is secure. Plug the charging cord into AC adapter. Plug AC adapter into standard wall outlet.



The LED lights cycle on and off from left to right while the unit is charging. All 3 LED's will turn on solid when charging is complete.



1. Turn the unit on or off by pressing and holding the **Power Button** until the unit powers on or off.

NOTE: To conserve battery power the unit is designed with an auto shut off feature. If the unit turns off during a service simply turn it back on and resume your service.



2. Press **MODE** to select mode, EXFOL is active by default.




3. Press +/- to start mode in **LOW**. Press +/- to cycle between **HIGH** and **LOW**.



MODE: EXFOL

Apply exfoliating product of choice. Hold the bt-micro fusion at a 45 degree angle with the buttons facing outward. Gently smooth with upward motions. Periodically wipe away any debris from the applicator.

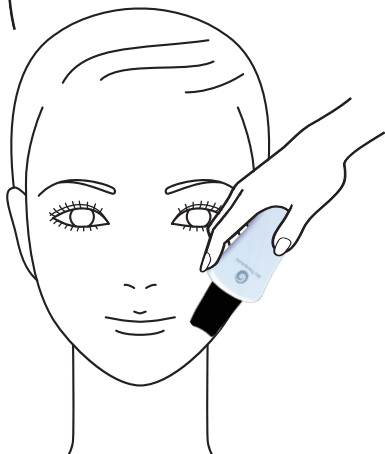

MODE: CLEAR

Apply deep cleansing product of choice. Hold the bt-micro fusion with buttons facing inward and the applicator flat against skin. Gently smooth with upward motions.

MODE: APPLY

Apply serum, cream, or ampoule of choice. Hold the bt-micro fusion with buttons facing inward and the applicator flat against skin. Gently smooth with upward motions.

bt-micro® FU



Bio-Therapeutic

TECHNOLOGY ACCELERATED SKIN CARE™



bio-therapeutic.com

@biotherapeutic

1-206-938-5800

watch the video: bio-therapeutic.com/fusion

Models and items shown are not available in all countries. The bt-micro fusion does not penetrate the living layers of the skin, and is not intended in the diagnosis of disease or other conditions, or in the cure, mitigation, treatment, or prevention of disease, and is not intended to affect the structure of the human body.



Distributed by: Micro Current Technology, Inc.
2244 1st Ave South Seattle, WA 98134

