





SUPER POLARIZED.

With the world's first **variable polarization system** in a hybrid dermatoscope, DermLite® DL5 lets you uncover a whole new layer in handheld dermoscopy. Not only can you toggle between polarized and non-polarized modes, now you're able to visualize extremely superficial structures under parallel polarization and, as you scroll through the entire polarization range, optically explore a lesion's visual depths. And on that note, you'll love DL5's optics, which are stunning up-close and outstanding at a distance: Its **32 mm lens system** with *true* 10x magnification reveals the supertiniest of details — better than ever before.

CROSS-POLARIZED



NON-POLARIZED



PARALLEL-POLARIZED



Simulated images






SUPER CLINICAL.

DermLite DL5 supports your clinical workflow unlike any other dermatoscope. It's the world's first hybrid dermatoscope with **365nm UV**, allowing you to high-light fluorescent structures — at full 10x magnification.

And, when screening your patients at a distance, DL5's **high-powered torch LED** provides phenomenal illumination to larger areas or hard-to-reach spots.

Need to measure a lesion? DL5 comes with an on-board, magnetically attached 100 mm ruler. Simply slide it off and its circular markings make it easier than ever to assess a lesion or choose the proper biopsy punch.

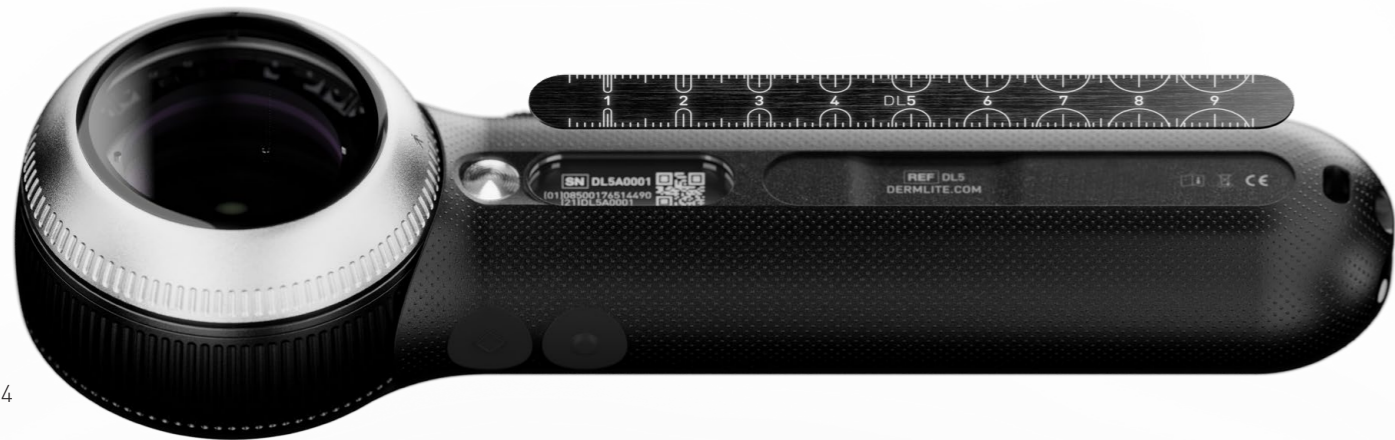



SUPER CHARGED.

DL5 lets you work substantially longer on a single charge, thanks to a battery that's more than **twice the capacity of DL4**. Between uses, keep it charged up on the included desktop charging base — which doubles as home base for your IceCap® supply. Or, use the supplied USB-C cable while on the go; regardless, DL5's 4-level charge indicator keeps you up to date.



D DL5

POLARIZATION Variable polarization (incl. polarized/non-polarized toggle)

OPTICS ø32mm, 10x

ILLUMINATION 16x white, 4x UV (365nm), 4x PigmentBoost, 1 torch LED
(each mode is brightness-adjustable)

POWER 3000mAh 3.7V lithium-ion, 4-level charge indicator,
charging base included

SMARTPHONE USE Via included MCC™ smartphone adapter

CONTACT PLATE Removable via bayonet (incl. 10mm reticle)

INFECTION CONTROL IceCap® disposable covers (5 included)

INCLUDED ACCESSORIES Magnetic eyepiece, MCC smartphone adapter, 100 mm
steel ruler, leather belt pouch, protective holster, 5 IceCap
disposable faceplate covers, 2m USB charging cable,
microfiber cleaning cloth

SIZE & MASS 182x58x36mm, 285g



SN | DL5A0000
10110850017651449
21DL5A0001