


ASEPSIS SELF-CONTAINED STANDARD SINGLE WATER SYSTEMS #8181 & #8183



Installation Instructions

INSTALLATION

- A. Remove the water system from packaging and inspect.
- B. For Wall Mounting, select a location, keeping in mind that the system must be accessible for easy bottle replacement and servicing.
- C. For Post Mounting use a 2" post mount DCI #8137 or #8162.
- D. Connect the Gray 1/4" line to an air source. Pressure should be between 40 PSI and 120 PSI.
- E. Connect the Red 1/4" line to the Dental Unit water supply line.

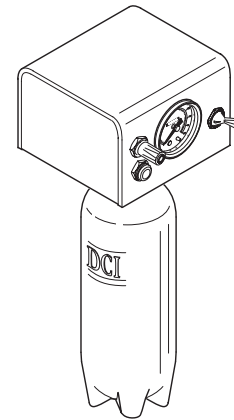

NORMAL OPERATION

- A. Install the bottle with water or solution you choose, in the selected bottle position.
- B. Place Master Toggle located on the bottom right rear of chassis in On position (to front).
- C. Check pressure gauge on front panel to assure that pressure reads 40 PSI, (Factory set at 40 PSI). If pressure needs to be adjusted, use regulator located on the bottom left rear of the chassis. Clockwise to increase pressure and counterclockwise to decrease pressure.

****CAUTION: Pressure should not exceed 60 PSI****

- D. The Quick Disconnect located on the lower left corner of the front panel will provide solution to auxiliary equipment. The Quick Disconnect source is from the bottle only. The flow control valve above it controls the flow through the quick disconnect. Clockwise reduces flow and counter clockwise increases flow.

AIR PURGING OPERATION

To air purge water lines in the dental unit use the procedure below.

- A. Place Master Toggle located on the bottom right rear of chassis in Off position (to rear).

- B. Activate all Dental Unit instruments that use self-contained water, i.e. air water syringe, handpieces, cup filler, and any other water using attachments. To flush the QD connection, you must insert a #0014 male QD. Attach open-ended poly-flow tubing and drain into proper receptacle.
- C. Activate Purge Toggle located on upper right corner of front panel. To air purge, place toggle in Purge position (down).
- D. Flush water lines until all liquid is gone.
- E. Return Purge Toggle to Off position (up).
- F. Follow Normal Operation instructions to continue use.

CHANGING BOTTLE

- A. Place Master Toggle located on the bottom right rear of chassis in Off position (to rear).
- B. Remove empty bottle.
- C. Replace with a filled bottle.
- D. Place Master Toggle located on the bottom right rear of chassis in On position (to front).

SURFACE CLEANING

- A. The water system should be wiped down with a mild detergent solution.
- B. Harsh detergents or chemicals should not be used.

STANDARD SINGLE WATER SYSTEMS SCHEMATIC 91609 FOR #8181 & #8183

