


ICare ActiveX Overlay

Designed to provide added pressure relief and comfort on an existing mattress

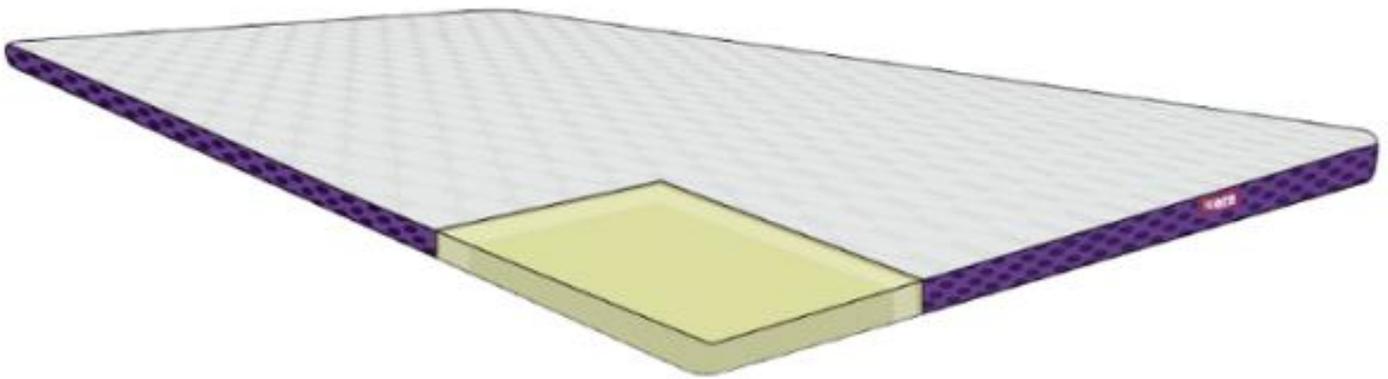

This overlay is heat and pressure sensitive causing them to soften around pressure points and inflamed areas. A mattress overlay is an affordable way to provide relief for the short-term; especially for short-term illness, recovery from surgery or to provide added comfort to an existing mattress.

Pressure Redistribution

The overlay features the ActiveX 84 material, which is a high grade of elastic type foam. ActiveX is responsive to temperature around inflamed areas and will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.

Ventilated Side Fabric

Vent fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. This is on the side wall of this overlay cover to greatly increase air movement.

Medical Covers Available

Zerotec fabric is a medical barrier fabric that has high elasticity in every direction, which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties.

- Safe to use on an alternating mattress.
- High Risk Category (rating based on average peak pressure of <18mmHg).

We care about
your comfort!

www.capeabilities.com.au

 26 Strickland Street, Bunbury WA 6230
 08 9717 6622
 bunbury@capeabilities.com.au

 Unit 1/64 West Street, Busselton WA 6280
 08 9751 1446
 sales@capeabilities.com.au

Product Code	Size	Mattress Dimensions	Maximum Weight Capacity	Risk Category	Feel Rating	Warranty
ICTS	Single	1900 x 900 x 50mm	215kg	High	Soft	5 years
ICTLS	Long Single	2030 x 900 x 50mm				
ICTKS	King Single	2030 x 1070 x 50mm				
ICTD	Double	1900 x 1350 x 50mm				
ICTQ	Queen	2030 x 1520 x 50mm				
ICTSQ	Split Queen	2030 x 750 x 50mm				

We care about
your comfort!

www.capeabilities.com.au

 26 Strickland Street, Bunbury WA 6230
 08 9717 6622
 bunbury@capeabilities.com.au

 Unit 1/64 West Street, Busselton WA 6280
 08 9751 1446
 sales@capeabilities.com.au