

Functional characteristics

1. 24 Led lights:

12 lights on one side, 24 lights on both sides available. The front is red & infrared light ,and the back is blue light.

2. Integrated a variety of wavelengths:

Double-sided available wavelengths of 450nm, 660nm, 850nm, 950nm a total of four.

3. Efficient light therapy:

The most effective light output power for human body is 200mw/cm².The optical output power of this product both exceeds 200mw/cm² at 0 distance and 1 inch (2.5CM).

4. Portable design:

Small size lamp, weight only 243g (0.53lb).

5. Built-in booster, flashlight use time up to 3 hours.

6. Excellent detail design:

Metal appearance, beautiful and fashionable, faster heat dissipation, won't be hot.

How to use

Please read the manual in detail

Quality and Technology Certification

1. Product batteries are supplied by factories that have passed the international MSDS battery certification.

2. Product is entrusted to SHENZHEN YIHONG LIGHTING CO., LTD by Scienlodic Ltd. for production.

3. FDA owner number is 10082798.

4. The product has passed the European patent and the US patent is under application. Please do not steal and imitate without authorization.

Warning

1. Products are not waterproof, pay attention to waterproof when using.

2. The higher the light power generated heat will also be higher, pay attention to the use of distance and use time, to avoid skin burns.

3. Please use the matching charging cable, it is recommended to use a 5A1V plug for charging.

4. Battery power is low when using the product will affect the performance of the work, pay attention to timely charging.

5. If you have any questions during use, please consult with the seller in a timely manner.

After Sales

If you have any questions in using, please ask the seller or professional guidance for help.

All Scienlodic products with any quality problems can get free reissue or replacement within 2 year. Please check and use the product in time to avoid missing the after-sales period.




Usage Demonstration



Brightness adjustment
Turn on/off & Switch

Long press the second button to turn on/off the power. After turning it on, you can press this button to switch between red and blue light.

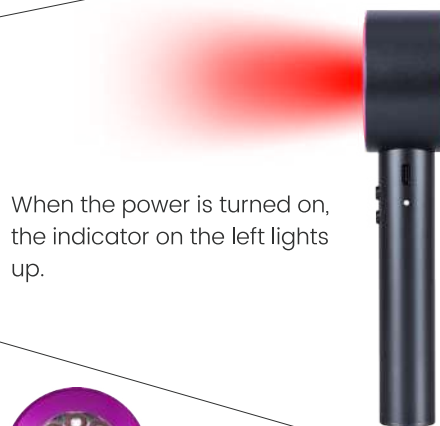

When the power is turned on, the indicator on the left lights up.



After use, please charge it in time. The indicator light will be white when charging, and the indicator light will be red when fully charged.

Product name	Double-sided multi-wavelength light therapy lamp	
Illuminating surface	Front	Back
Wavelength	660nm&850nm	450nm&950nm
Number of led light	660nm*9+850nm*3	450nm*9+950nm*3
Integrating sphere test result	When the voltage is 2.1V and the current is 400mA the photosynthetic effective radiation flux PAR WATT: 562.05mW	
Product size	Handle: 13*3*3cm (5.1*1.2*1.2inch) Lamp head: 6.5*6.5*4cm (2.5*2.5*1.5inch) Overall: 19*6.5*4cm (7.4*2.5*1.5)	
Product power	6-12W	
Battery information	3.7V 5000mAh 18.5Wh	

Product in small size has no high voltage and no EMF