


BRUSHES

CLEANING and STERILIZATION

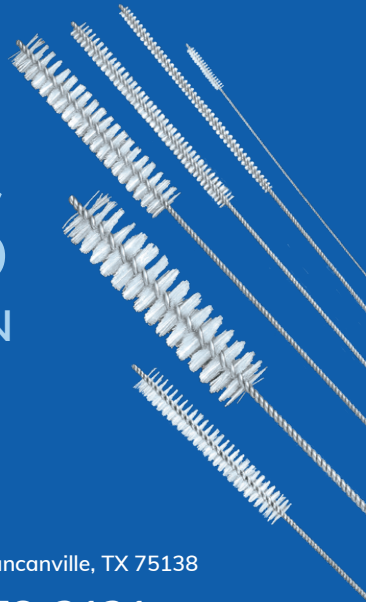

P.O. Box 382120 • Duncanville, TX 75138

1(800)858-2121

www.qualityaspirators.com

CLEANING & STERILIZATION

INSTRUCTIONS

All Quality Aspirators cleaning brushes can be sterilized by dry heat, autoclave, and chemiclave. They can also be cleaned using any standard method of cleaning for medical devices.

Typical steam sterilization methods of 121° C for 20 minutes at 1.036 bar (15.03 psi) and 134° C for 3-4 minutes at 2.026 bar (29.41 psi) are acceptable for all Quality Aspirators unwrapped suction tips.

All Quality Aspirators cleaning brushes can be sterilized without causing damage.