



The Blom-Singer® FISTULA PROSTHESIS INSERTION & REMOVAL (I&R) CATHETER

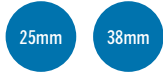


Insertion & Removal (I&R) Catheter
for use with 25mm and 38mm prosthesis

- Material > Semi-rigid thermoplastic
- For insertion and removal with these sizes only, 25mm and 38mm prosthesis
- 8 Fr. catheter diameter gradually increasing to 11 Fr. in the keyhole end
- Designed for non-surgical insertion of the fistula prostheses in the OFFICE or O.R.
- Sterile, 3 year shelf life
- Single use, disposable

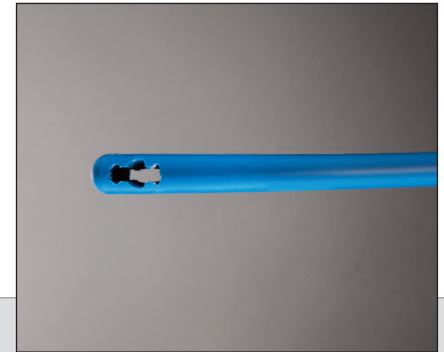

Retrograde Insertion

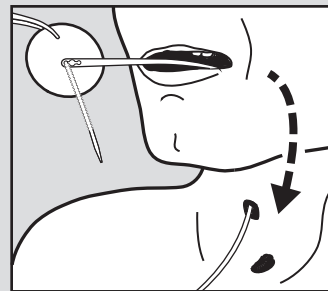

Anterograde Removal



ORDERING INFORMATION


Call InHealth Technologies Customer Service, 800.477.5969, and provide the REF number(s) below.

SIZE	REF NUMBER	FLANGE DIAMETER	OVERLAP of FLANGE to FISTULA DIAMETER
SMALL I&R Catheter provided	FM1003	25mm	8-10mm
MEDIUM I&R Catheter provided	FM1004	38mm	8-10mm
LARGE	FM1005	50mm	12.5mm
CATHETER	FM1010	Blom-Singer® Fistula Prosthesis Insertion and Removal Catheter *	

[* Fistula Prosthesis Insertion and Removal Catheter available separately]

Introducing the Blom-Singer® ADJUSTABLE BI-FLANGED FISTULA PROSTHESIS

For the management of hypopharyngeal fistulae

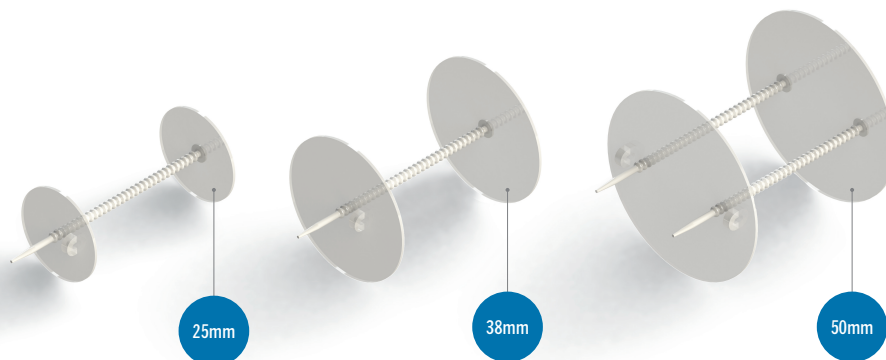

Designed by Eric D. Blom Ph.D.

MANAGEMENT OF THE HYPOPHARYNGEAL FISTULA – QUALITY OF LIFE

- Enhances maintenance of a dry fistula tract depending on size of fistula and condition of tissue
- May reduce leakage of saliva, food and drink
- NPO Instructions may be alleviated
- Possible avoidance/reduction of enternal/parenteral nutrition
- Potential for Post-laryngectomy speech with TEP

PROMOTES COST SAVINGS

- Insertion & Removal Catheter (I&R) provided
- Can be placed in OFFICE or O.R. (Choice in Anesthesia protocol/necessity)
- Anesthesia and facility cost may be reduced

WHAT IS INCLUDED IN THE OFFERING?

- Insertion & Removal Catheter (I&R) provided with the 25mm and 38mm Prosthesis
- The 50mm Adjustable Bi-Flanged Fistula Prosthesis does not require the I&R Catheter

IS IT STERILE?

Yes.


SHELF LIFE?

3 years from date of manufacture.

CAN THIS DEVICE BE REPROCESSED FOR MULTI-PATIENT USE?

No, Single patient use only.

FISTULA PROSTHESIS PLACEMENT & FUNCTION

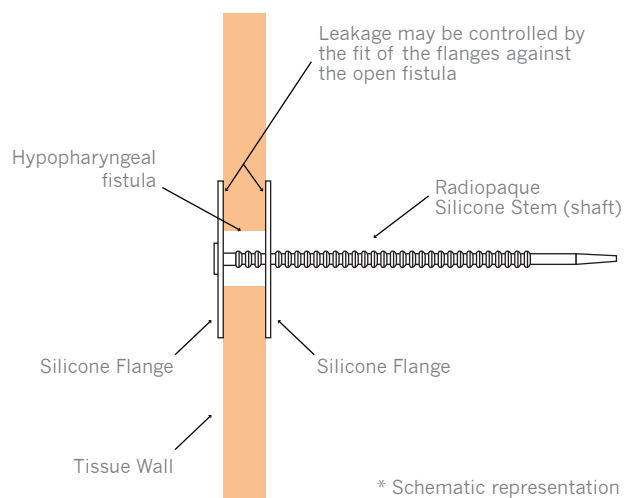

MATERIALS

FLANGES

- Soft, flexible
- Medical grade silicone
- Translucent to allow Fistula tract visualization

BEADED STEM

- Medical grade silicone, radiopaque
- Elastomeric (stretchable)
- Locking loop prevents flange movement