


# Protekt Akra-SE™ Semi-Electric Bed

The **Akra-SE™** Homecare bed has a heavy-duty bed frame constructed from reinforced steel for maximum stability and durability. The head and foot position feature motorized electric adjustment for customizable patient comfort and support. In the event of a power outage, a 9V battery backup can be used to allow movement to the lowest position. The snap on or off low-voltage DC motor box ensure quick, efficient set-ups and easy routine maintenance.



Half Rails

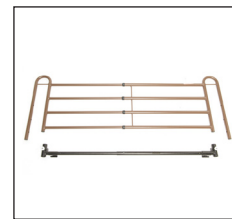

Full Rails



Pendant

## FEATURES

- The entire bed, including the rails comes complete in only (2) boxes
- Quick set-up & easy maintenance with motor & electronics that can be removed for service or cleaning
- Easy-to-use control adjusts upper body, knees, and height of bed
- Head and foot sections can be raised independently or simultaneously
- Large button, simple to use hand pendant
- Manual crank at the foot adjusts the frame height
- Scratch-resistant bed ends keep the bed looking new
- 9V battery back-up lets user go to the lowest bed position during power failures (battery not included)
- Headspring and foot assembly are recessed to help reduce pinching and injuries
- Beds include 2 locking and 2 non-locking casters
- Quiet, efficient ETL-approved DC low voltage single motor box
- Optional full-electric conversion kit

## SPECIFICATIONS

Surface Height:	15" min. to 23" max. (without mattress)
Sleep Surface:	80"L x 36"W
Overall Bed:	86"L x 36"W
Weight Capacity:	450 lbs.
Warranty:	Limited Lifetime Frame 2 Years on Motor

## SHIPPING DIMENSIONS

Component	Carton Size	Shipping Weight
Foot/Head Springs, Motors, Pendant	47"x35.5"x9"	101 lbs. w/HR 106 lbs. w/FR
Head/Foot Bed Ends, Casters	38"x26"x5"	57 lbs.

Item #	Description	HCPCS Code
PMSEBHR	Semi-Electric Bed w/Half Rails (No Support Surface)	E0261
PMSEBFR	Semi-Electric Bed w/Full Rails (No Support Surface)	E0261
PMSEBHRM	Semi-Electric Bed w/Half Rails w/Protekt 100 Foam Mattress (Group 1)	E0260
PMSEBFRM	Semi-Electric Bed w/Full Rails w/Protekt 100 Foam Mattress (Group 1)	E0260
PMSEB	Semi-Electric Bed w/OUT Rails (No Support Surface)	E0295
BP-HR	Half Length Side Rails (pr)	E0305
BP-FR	Full Length Side Rails (pr)	E0310
BP-FEKIT	Full-Electric Conversion Kit (Incl. Motor, Drive Shaft & Pendant)	N/A