

SPAN PressureGuard
Easy Air[®] XL



Bariatric Microclimate
Management and Alternating
Pressure Mattress System

Industry-best
microclimate
management
performance

Uniquely equipped to
meet the extreme skin
integrity challenges of
bariatric care

Safety Edge[™]
for added fall
prevention
and stability

science, safety, simplicity


 **savaria**[®]

SPAN PressureGuard Easy Air® XL


Packed with the design innovations that have made Span a recognized pioneer in high-performance microclimate management and pressure management systems.

Treatment Flexibility

Select one of two therapy modes:

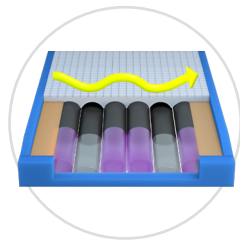

Alternating Pressure
Gently changes loading across the surface (22-minute cycle).

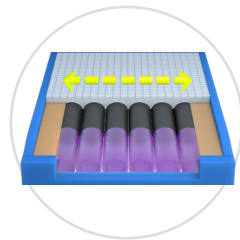

"Float" Mode
Provides active, powered flotation therapy.


PressureGuard Design

Integrated air system/engineered foam shell with "Safety Edge" bolster system delivers superior edge-of-bed safety, stability and comfort.



Indications

- Patients weighing 350–1,000 lbs.
- Patients with excessive skin moisture due to sweating, decreased mobility, inability or refusal to be repositioned
- Patients with macerated skin due to any of the above
- Patients with increased skin or body temperature due to infection, sepsis, or other conditions
- Patients with a combined score ≤ 5 on Braden Moisture, Activity and Mobility scales
- Treatment of stage 1–4 pressure injuries
- Prevention of pressure injuries for high-risk patients
- Some patients with stabilized flaps/grafts (see span algorithm for parameters)

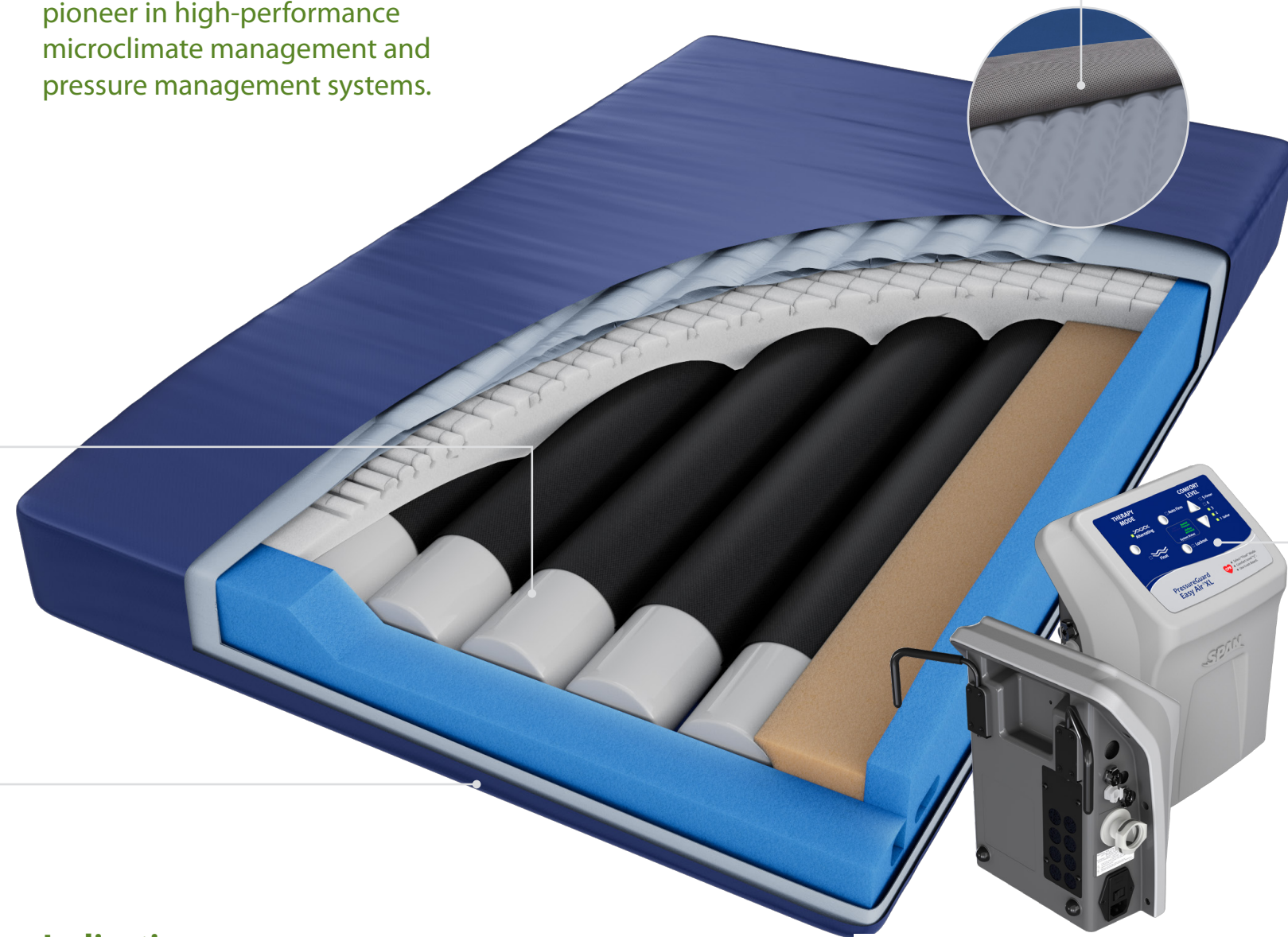


Proprietary Air Diffusion Matrix™

Three-dimensional liner helps ensure uninterrupted air flow in both cover and coverlet.

Control Unit Features

- Does not deflate during power interruption and transport
- Timed (30-minute) Auto Firm mode for added stability during transfers and ADL
- Simple comfort adjustment, including ultra-soft "Max Immersion" setting for users under 120 lbs
- Lockout function eliminates accidental or unintentional changes
- Normal/Low Pressure indicator
- Streamlined, lightweight design

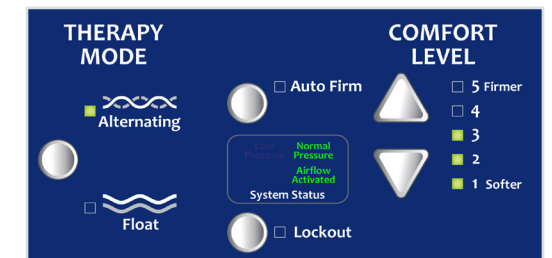


- Efficient: Uses less than 15 watts
- Generates less heat, noise than traditional roll-up surfaces
- Recessed, shielded connectors for protection from impact, damage
- Detachable power cord for easy replacement

Outer Coverlet

- Highly vapor-permeable
- Superior air flow
- Wipes clean in place with standard hospital-grade cleaners
- Zippered, machine-laundryable (air dry)

Inner Air Delivery Cover

- Dedicated air supply
- Directs constant stream of air
- Sweeps away moisture vapor before it can re-form as liquid
- Reduces maceration

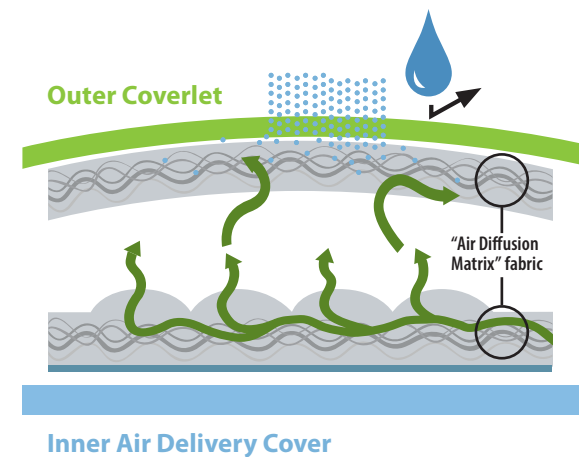

WEIGHT RANGE

Appropriate for users from 350–1,000 lbs.

Note: May not provide adequate pressure redistribution or comfort for patients weighing less than 350 lbs.

FLAMMABILITY

Conforms to NFPA 101 (Life Safety Code); Cal. TB 117; ASTM E1590; and 16 CFR Parts 1632 and 1633.

COMPLETE SYSTEM

Includes Mattress, Control Unit, and Coverlet

To order coverlet only, add "CLT-" to part number (example: CLT-L8042XL). Mattress weight approximately 32 lbs. Two-year warranty (all components).

CONTROL UNIT

Dimensions	12" x 9.5" x 7"
Weight	6 lbs.
Leakage Current	< 250 micro-amps
Voltage	120 AC
Max Current	1.0 Amp
Frequency	60 Hertz
Conforms to:	ANSI/AAMI ES60601-1:2005 + C1:2009 + A2:2010 + A1:2012 CAN/CSA-C22.2 NO. 60601-1:2014; EN 60601-1:2006 + A1:2013; ANSI/AAMI HA60601-1-11:2015; CAN/CSA-C22.2 NO. 60601-1-11:2015; IEC 6060101-2:2014 (4th ed.); EN 60601-1-2:2014 (4th ed.)

ITEM CODE

DESCRIPTION

L8042XL-29	Easy Air XL System 80" L x 42" W x 7" H
L8048XL-29	Easy Air XL System 80" L x 48" W x 7" H
L8053XL-29	Easy Air XL System 80" L x 53" W x 7" H
L8442XL-29	Easy Air XL System 84" L x 42" W x 7" H
L8448XL-29	Easy Air XL System 84" L x 48" W x 7" H
L8453XL-29	Easy Air XL System 84" L x 53" W x 7" H
62507	Easy Air XL System 86" L x 42" W x 7" H
8310	Easy Air XL Control Unit



Get more information

800.888.6752 (Span)

savariapatientcare.com



Span-America

| 800.888.6752

| PO Box 5231 Greenville SC 29606

Span-America is a division of Savaria Corporation, a global leader in accessibility. Learn more about Savaria at savaria.com.