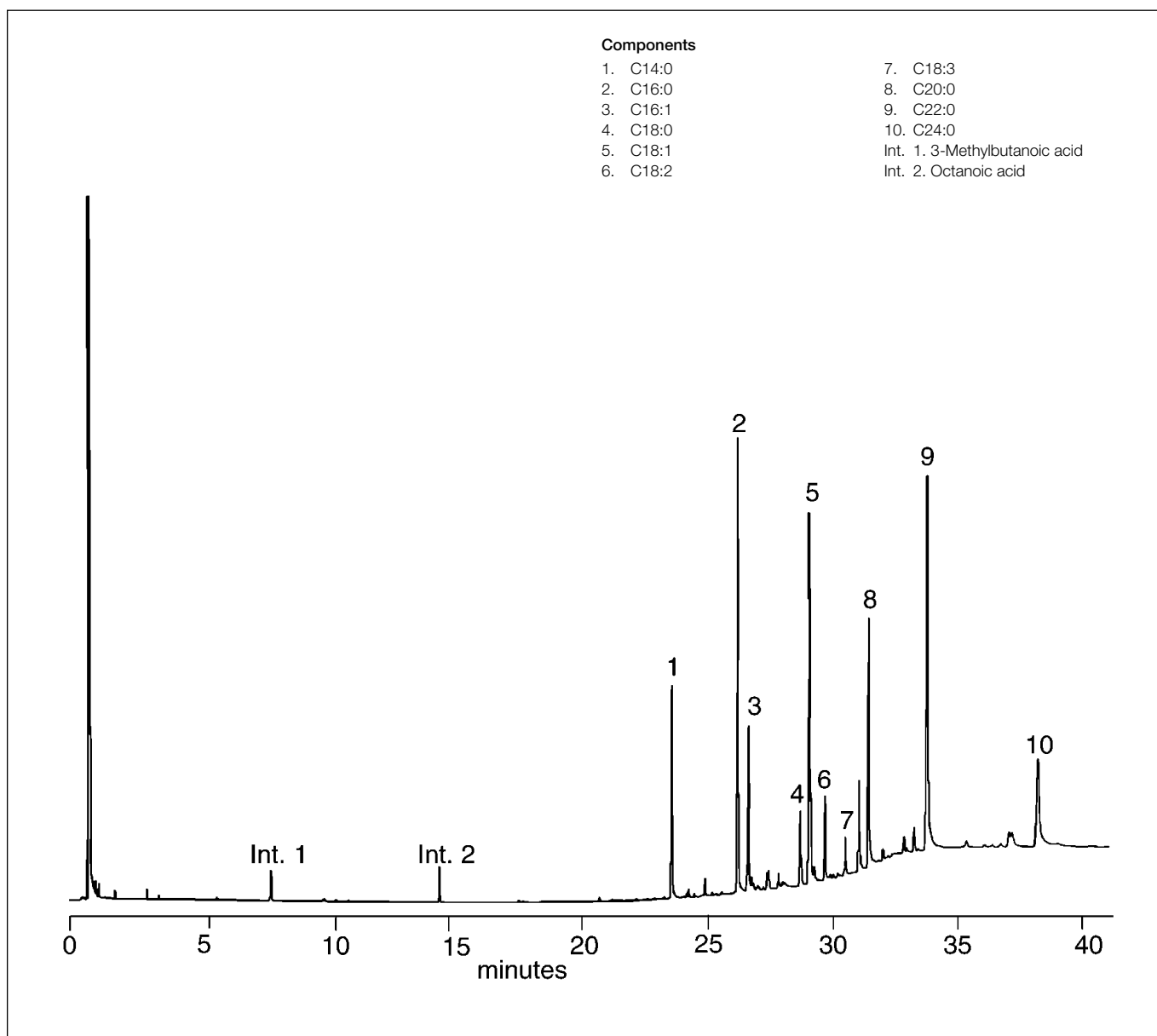


Analysis of cod oil fatty acids (underivatised)

BP21 (FFAP)

Column part number	054468		
Phase	BP21 (FFAP)	Detector	FID, 280°C
Column	25 m x 0.32 mm x 0.25 µm	Injection mode	Split (50:1)
Initial temperature	70°C, 0.5 min	Carrier gas	He, 15 psi
Rate	6°C/min	Injection volume	1 µL
Final temperature	260°C, 10 min	Injection temperature	280°C


 For more information visit www.trajanscimed.com or contact techsupport@trajanscimed.com