

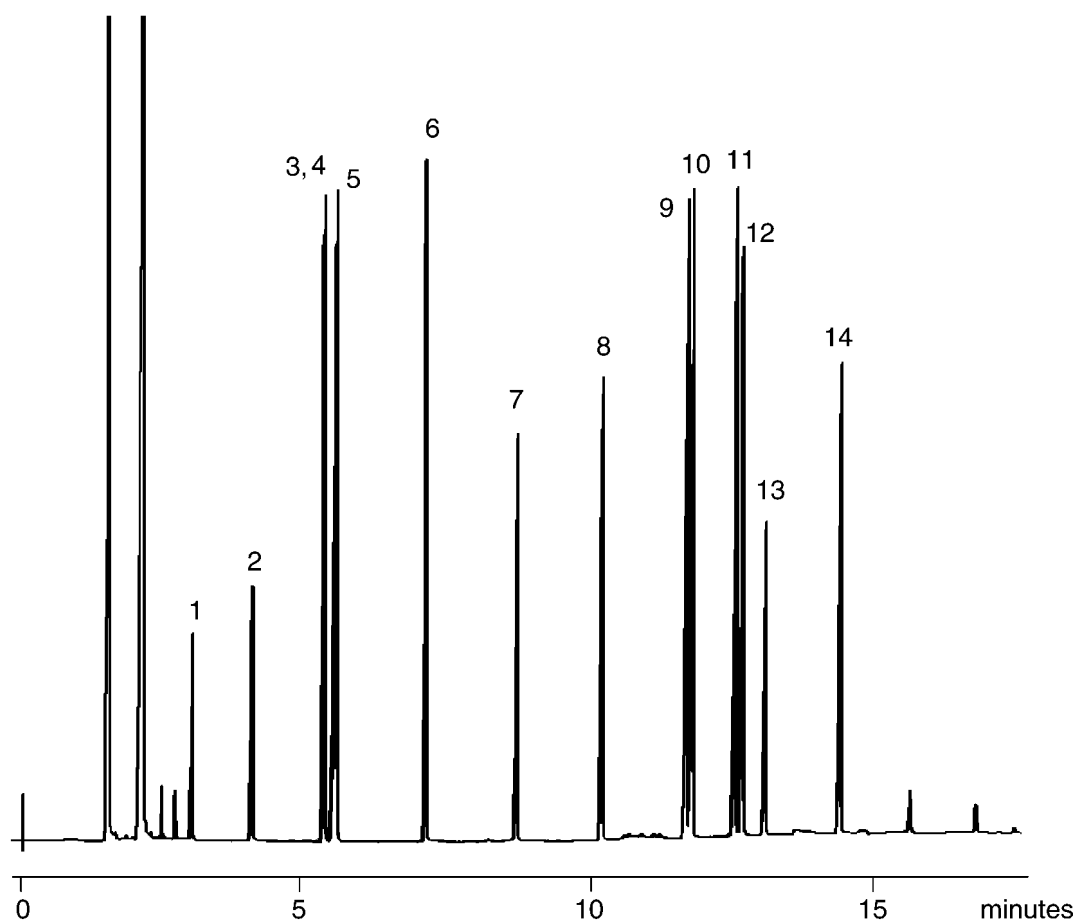
Analysis of dimethyl esters of dicarboxylic acids

BP1

| | | | |
|---------------------|--------------------------|-------------------|--------------|
| Column part number | 054047 | | |
| Phase | BP1 | Final temperature | 240°C, 0 min |
| Column | 25 m x 0.22 mm x 0.25 µm | Detector | FID |
| Initial temperature | 60°C, 0 min | Injection mode | Split |
| Rate | 10°C/min | | |

Components

- | | | | |
|-----------------------|-----------------------|----------------------------|-----------------------|
| 1. Dimethyl oxalate | 5. Dimethyl maleate | 9. Dimethyl suberate | 13. Dimethyl azelate |
| 2. Dimethyl malonate | 6. Dimethyl glutarate | 10. Dimethyl phthalate | 14. Dimethyl sebacate |
| 3. Dimethyl fumarate | 7. Dimethyl adipate | 11. Dimethyl terephthalate | |
| 4. Dimethyl succinate | 8. Dimethyl pimelate | 12. Dimethyl iso-phthalate | |


 For more information visit www.trajanscimed.com or contact techsupport@trajanscimed.com