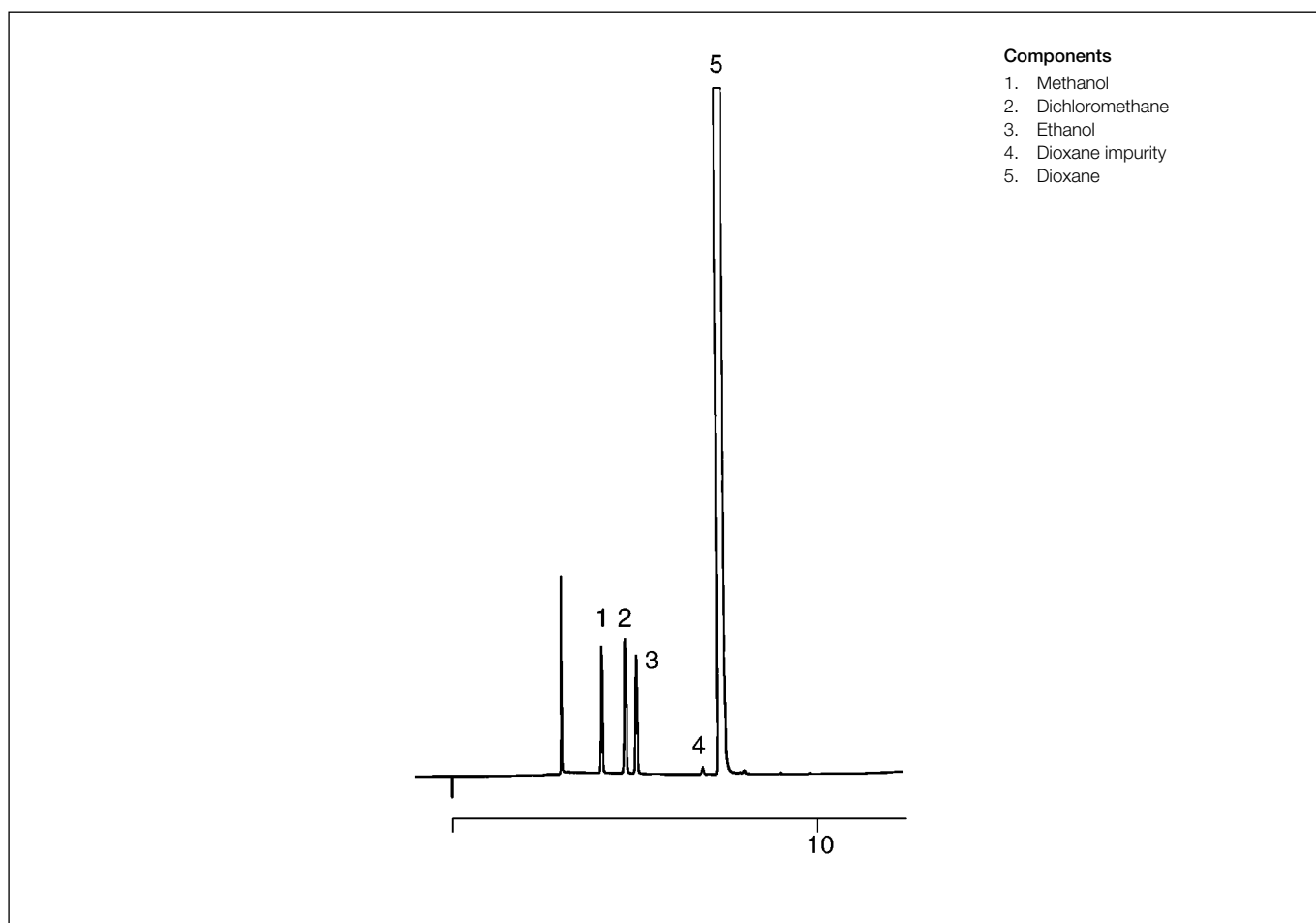


Analysis of solvent impurities in dioxane (10 ppm)

BP20 (WAX)

Column part number	054448		
Phase	BP20 (WAX)	Detector	FID, 280°C
Column	25 m x 0.53 mm x 1.0 µm	Carrier gas	H ₂ , 2 psi
Initial temperature	40°C, 2 min	Injection mode	Split (30:1)
Rate	10°C/min	Injection volume	0.2 µL
Final temperature	120°C		


 For more information visit www.trajanscimed.com or contact techsupport@trajanscimed.com