

Geuder®

UNO Colorline Mach 4 Vitreous Cutter

Four Cuts In a Row

Low Traction - High Flow

Up to
24,000 cpm

Quad-point
precision

Maximum
duty cycle



Geuder®
Germany

Optimized Vitrectomy

Maximum Flow


Reduced viscosity of aspirated vitreous body through four times higher cutting rate

Optimal aspiration by up to 24.000 cpm in combination with Geuder megaTRON S4 HPS system

Low Traction

Permanently open aspiration port allows for fast core vitrectomy

Large aspiration port and optimized cutting properties provide constant flow and precise vitreous shaving

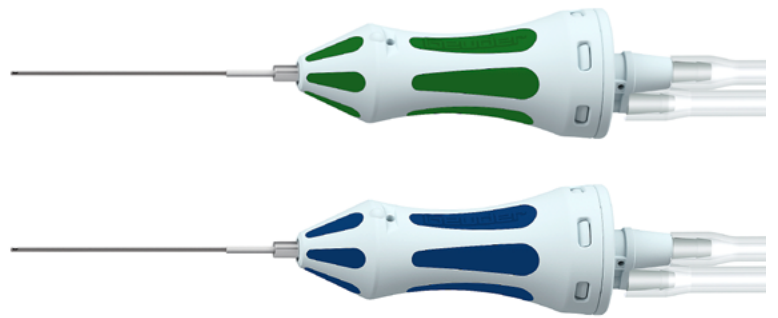

Einmal-High-Speed-Vitrektor UNO Colorline

Vierfachschnide, pneumatisch, mit Geuder Sicherheitskonnektor, steril
VPE 6 Stück

☒ -46304

Vitrektor 23 G Mach 4
23 Gauge/ 0,64 mm

☒ -46504

Vitrektor 25 G Mach 4
25 Gauge/ 0,5 mm

Single-Use High Speed Vitreous Cutter UNO Colorline

Quadruple-blade, pneumatic, with Geuder safety-connector, sterile
6 pcs. per box

Vitreous Cutter 23 G Mach 4
23 gauge/ 0.64 mm

Vitreous Cutter 25 G Mach 4
25 gauge/ 0.5 mm



Vierfachschnide
Quadruple-Blade



Geuder® UNO Colorline Mach 4

Technical Data

Compatible with Geuder megaTRON, megaTRON S3, megaTRON S4 und megaTRON S4 HPS

24,000 cpm (in combination with megaTRON S4 HPS)
One blade with four cutting points in two ports
incl. safety connector