

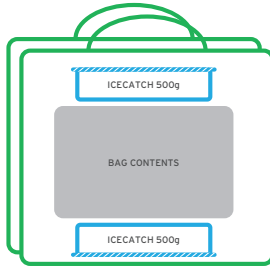




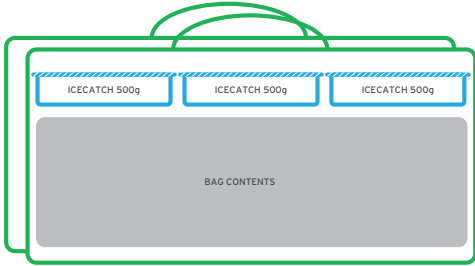




Blood Carriers User Guide

Tested to maintain temperature of 2-8 degrees

ICECATCH must be pre-conditioned to -18 degrees before use

Product Name/Code	Carrier	ICECATCH	Placement	Ambient Temperature	*Temperature maintained for
Small Blood in Transit Medical Carrier (BLTD1)		 500g x2		 21°C 30°C	 21hrs 29min 18hrs 12min
Large Blood in Transit Medical Carrier (BLTD2)		 500g x3		 21°C 30°C	 28hrs 43min 17hrs 20min

*Study by TuV GmbH. ICECATCH preconditioned to -18°C for 24 hours. Ensure that placement of ICECATCH is correct with non-print insulated side in contact with contents. For more information please contact us.

