

FRED PA-1

Simplicity as a priority



FRED PA-1

Dedicated to public access, SCHILLER's latest defibrillator, the FRED PA-1, has been designed so that even untrained users can save lives. The FRED PA-1 guides the rescuer step-bystep during the entire resuscitation process.






IT CONCERNS ALL OF US

More than 7 million people worldwide fall victim to a sudden cardiac arrest every year.

A sudden cardiac arrest can hit anyone, anywhere and at any time.

Defibrillation consists in delivering electrical energy to the heart in order to restore a proper heart rhythm. A defibrillator in combination with well-performed chest compressions is the only therapy that can save a patient's life.

The FRED PA-1 is a well-designed lifesaving equipment, which is highly visible and clearly recognizable in all environments.

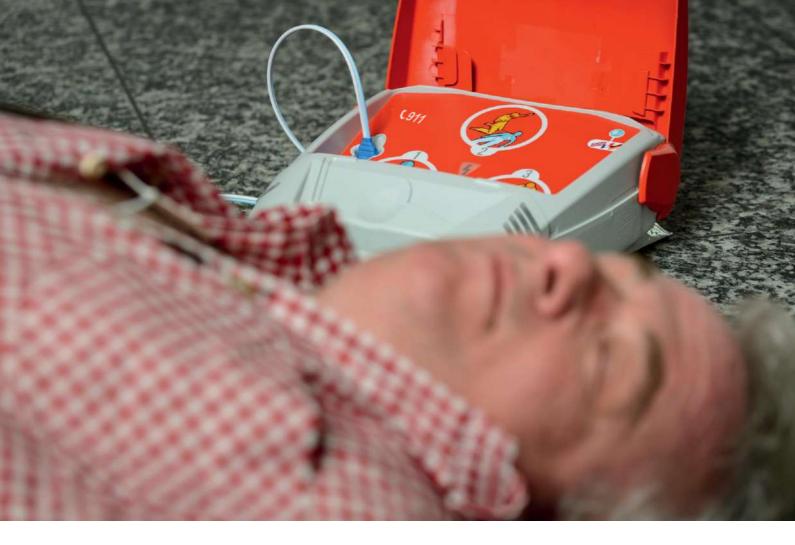

HIGHEST USER-FRIENDLINESS

In emergency situations, procedures that would normally be clear and easy to follow can suddenly become confusing. Every detail counts and can make a difference in saving lives. That's why a defibrillator needs to be as simple as possible.

1-2-3 STEPS

Not only is the user interface streamlined with 1-2-3 steps, but spoken instructions and clear pictures also guide the userthrough the process.

NO DELAY

When you lift up the device cover, the FRED PA-1 starts up immediately. The electrodes are pre-connected: apply them to the patient's chest and let the FRED PA-1 guide you.

SHOCK DELIVERY POSSIBLE AS SOON AS A SHOCK IS RECOMMENDED

Resuscitation Guidelines emphasise the importance of minimizing the time-to-shock after Cardio Pulmonary Resuscitation (CPR). In order to reduce the hands-off time (time while no CPR can be performed), the FRED PA-1 features a fast analysis system that systematically pre-charges the correct energy so that the shock can be delivered as soon as it is recommended.



EMERGENCY NUMBER

The FRED PA-1 first advises you to call the rescue services and indicates the appropriate phone number (customisable).

91



FAST SET-UP

Pre-connected electrodes allow faster application of the electrodes to the patient's chest; analysis starts automatically.

ALWAYS OPERATIONAL

The FRED PA-1 automatically performs a regular self-test and indicates its ready-to-use status.

EASY TO CARRY

Grab the device with one hand to run to the scene of the incident.



NO DELAY

The FRED PA-1 starts up immediately and automatically when the cover is lifted.

ELECTRODE SLOT

Electrodes are immediately accessible and easily recognizable.

"REPEAT" FUNCTION

Should you have missed information given by the device, don't worry: just push the "Repeat" button to listen again to the last instruction.

CLEAR USERINTERFACE

1-2-3 steps

CLEAR OPERABILITY INDICATIONS

When a check is necessary, the FRED PA-1 provides warnings: acoustic signals as well as LEDs indicate the state of the device and the actions that have to be taken (replacing the battery, the electrodes...).



USED ANYWHERE

Sudden cardiac arrest can hit anyone, anywhere. The FRED PA-1 is designed for all locations, wherever you might need it.

With the ability to choose one of 3 languages¹, the FRED PA-1 can be installed in airports, hotels, railway stations, in any international environment.

With its high IP rating (IP55-protected against dust and water projection), the FRED PA-1 withstands even unfavourable conditions.

AUTOMATIC OR SEMI-AUTOMATIC

The FRED PA-1 analyses the patient's heart rhythm and determines whether a shock is necessary or not. Rescuers can fully rely on the result of the analysis, which is based on proven algorithms and decades of experience.

The FRED PA-1 is available either as semi-automatic or fully automatic device.

Delivering a shock to a victim of a sudden cardiac arrest should be done as soon as possible. In its automatic version, the FRED PA-1 delivers the shock without any action of the rescuer. In its semi-automatic version, the FRED PA-1 prompts the rescuer to deliver the shock by pressing the shock button.

DATA MANAGEMENT

The FRED PA-1 can store up to 2 hours of intervention data, including 500 events. Using a simple SD card, the retrieval of the intervention data as well as the software update are very easy to perform.

¹ The three available languages can be customised at the time of order.

IT'S THE DETAILS THAT MAKE THE DIFFERENCE

EXTENDED WARRANTY PERIOD, UP TO 10 YEARS¹

SCHILLER provides peace of mind with up to 10 years of warranty in the case of a warranty, the device is replaced with a new one.

READY WHEN YOU NEED IT

The FRED PA-1 has a long battery life (6 years) in order to be ready the day you need it. Battery level is checked during the self-test and a LED warns of low battery level.

ELECTRODE SHELF LIFE CHECK

Thanks to a highly innovative feature, the device is able to determine when the shelf life has expired and prompts the user to replace the electrodes.

1 Contact your local sales representative for more information

FOR ALL TYPES OF PATIENTS

Used with adult or children electrodes, the FRED PA-1 automatically recognises the electrode type and adjusts the defibrillation energy accordingly.

At the back of the device, a set of children electrodes can be safely stored.

CPR FEEDBACK

Defibrillation is not enough. In order to keep bringing oxygen to the organs, the blood flow must be kept. This is the purpose of chest compressions, and the reason why it's important to do it right. The FRED PA-1 helps you to achieve the highest efficiency in chest compressions by:

- indicating the right positioning of the hands
- setting a regular rhythm (metronome)
- providing a live feedback on the frequency of the chest compressions based on impedancemetry, without any additional sensor (FreeCPR).







Phone +603 6272 3033 Fax +603 6272 2030 sales@schiller.com.my www.schiller-asia.com

Phone +36 1 383 4780 / +36 1 460 9491 Fax +36 1 383 4778 info@schillerhungary.hu www.schillerhungary.hu

India

SCHILLER Healthcare India Pvt. Ltd. Mumbai - 400 059, India Phone +91 22 6152 3333 / 2920 9141 Fax +91 22 2920 9142 sales@schillerindia.com www.schillerindia.com

Poland SCHILLER Poland Sp. zo.o. PL-02-729 Warszawa Phone +4822 843 20 89 / +48 22 647 35 90 Fax +48 22 843 20 89 schiller@schiller.pl www.schiller.pl

Russia & C.I.S. AO SCHILLER.RU 125124 Moscow, Russia Phone +7 (495) 970 11 33 Fax +7 (495) 956 29 10 mail@schiller.ru www.schiller.ru www.schiller-cis.com

Serbia SCHILLER d.o.o. 11010 Beograd Phone +381 11 3979 508 Fax +381 11 3979 518 info@schiller.rs

www.schiller.rs

Slovenia SCHILLER d.o.o.



Spain SCHILLER ESPAÑA, S.A. ES-28232 Las Rozas/Madrid Phone +34 91 713 01 76 Fax +34 91 355 79 33 schiller@schiller.es www.schiller.es



Turkey SCHILLER TÜRKIYE Okmeydani-Sisli - Istanbul Phone +90 212 210 8681 (pbx) Fax +90 212 210 8684 info@schiller.com.tr www.schiller-turkiye.com



C*

SCHILLER UK Bellshill, ML4 3PR Phone +441698744505 Fax +44 1698 744 474 sales@schilleruk.com www.schilleruk.com





Schiller medicinski instrumenti d.o.o.

10000 Zagreb Phone +385 1 309 66 59 Fax +385 1 309 66 60 info@schillerzg.hr www.schiller.hr

Austria SCHILLER Handelsgesellschaft m.b.H.

AlfredSchiller (Beijing)Medical Equipment

France SCHILLER Médical S.A.S. F-67162 Wissembourg Phone +33 3 88 63 36 00 Fax +33 3 88 63 36 49 info@schiller.fr www.schiller-medical.com

100015 Beijing, China Phone +86 010 52007020 Fax +86 010 52007020-8016 info@schillerchina.com www.schillermedical.cn Croatia

A-4040 Linz Phone +43 732 70 99

Fax +43 732 757 000

sales@schiller.at

www.schiller.at

China

Co..Ltd.

