



**SAM**  
MEDICAL

# SAM<sup>®</sup> Soft Shell

## SHAPE. STRENGTHEN. SUPPORT.

### FLEXIBLE COMFORT


SAM<sup>®</sup> Soft Shell Splint is engineered with the same revolutionary C-Curve™ Principle as the original SAM<sup>®</sup> Splint, building strength from curves. With unrivaled flexibility, SAM<sup>®</sup> Soft Shell Splint bends into any simple arc, becoming exponentially stronger and more supportive, immobilizing sprains, strains and minimally displaced fractures.



### FEATURED TECHNOLOGY

#### ADAPTABLE C-CURVE™ PRINCIPLE

SAM<sup>®</sup> Soft Shell Splint leverages the revolutionary C-Curve™ Principle to mold its foam and aluminum layers into structural arcs, strong enough for pre- or post-cast care, while remaining lighter and more transportable than traditional splints that rely on heavy, rigid materials. Bend it back into its original form and it's ready to be remolded.


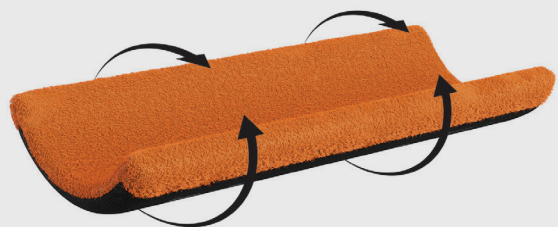




## ENGINEERED FOR VERSATILE SUPPORT

- 1 MAXIMUM-COMFORT SURFACE**  
Absorbent-brushed terry cloth fabric covers skin-contact surface, providing enhanced comfort.
- 2 MOLDABLE ALUMINUM**  
Flat 0-temper aluminum inside the splint's soft outer layer molds easily, allowing for a wide spectrum of applications.
- 3 RADIOLUCENT MATERIAL**  
Medical professionals are able to take X-rays and CT-scans without removing the splint from the fractured or injured limb.

- 4 SIZE THAT FITS**  
Full product line covers most sizes needed, including XL and sugar tong.
- 5 RESISTANCE CONTROL EDGE**  
Firm and extra-firm ends allow different levels of resistance to motion, according to needs. (Not available on the Sugar Tongs.)

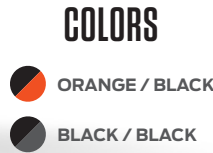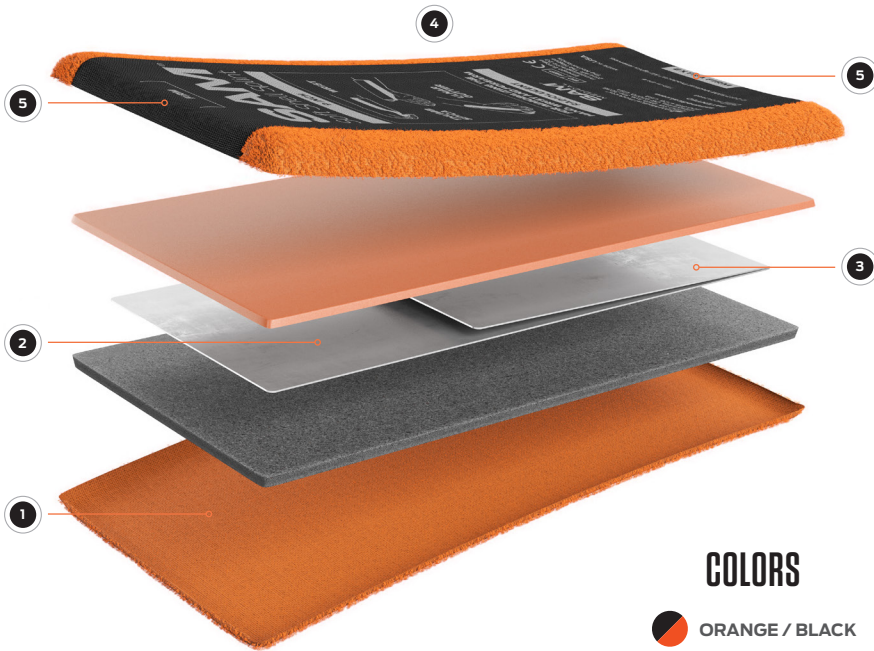

## MODELS



**STANDARD**  
4.25 in wide



**XL**  
5.5 in wide



**SUGAR TONG**  
Pediatric, Adult, Adult XL

## PRODUCT SPECIFICATIONS

### MODELS AVAILABLE\*

9" (SC042668)	15" (SC121802)
9" XL (SC042668XL)	15" XL (SC121802XL)
12" (SC122069)	
12" XL (SC122069XL)	

### SUGAR TONG MODELS

Pediatric (ST052106P)
Adult (ST052106)
Adult XL (ST052106XL)

### GLOBAL SALES

U.S. Toll Free: 800.818.4726  
 U.S. Local: 503.639.5474  
 U.S. Fax: 503.639.5425  
[sammedical.com](http://sammedical.com)

Emergo Europe, Prinsessegracht 20, 2514 AP, The Hague, The Netherlands	 SAM Medical Products® 27350 SW 95th Ave, Ste 3038 Wilsonville, OR 97070 USA	

\* Standard width 4.25 in wide, XL 5.5" wide.  
 SAM® is a registered trademark in the United States and other countries.  
 ©2017 SAM Medical Products. All rights reserved.