


T-Quick 2540MKA Autoclave

The Only Manual Autoclave with Closed Door Drying



Post Covid-19: Chairside sterilization is no longer recommended due to aerosols!
The T-Quick autoclave provides:

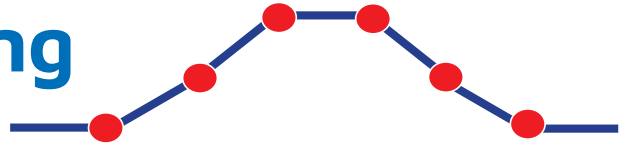




- Faster cycles that deliver dry instruments and no wicking
- Closed door drying
- More partial loads when needed vs. full loads



Features:

- T-Quick allows you to easily track each phase of the sterilization process
- Automatic safety shutoff prevents overheating
- Water detector with automatic shutoff
- HEPA air filter provides sterile, bacteria-free air for drying.
- Aligns with AAMI guidelines regarding sterilization cycles.

Sterilization Programs

Icon	Cycle Program	Load Type	Cycle Time (min.)*
	Unwrapped 134	Unwrapped instruments	11 Minutes
	Wrapped 134	Wrapped instruments	15 Minutes
	Unwrapped 121	Unwrapped delicate instruments	25 Minutes
	Wrapped 121	Wrapped delicate instruments	25 Minutes
	Unwrapped	Single Instruments	9 Minutes

* Cycle times are based on: 1 kg standard loads (not including drying), Warm chamber (and subject to load size)

Panel Features



- ① Cycle Status Indicator
- ② Pressure Indicator
- ③ Power Switch
- ④ Sterilization Timer
- ⑤ Water Pump Switch
- ⑥ Drying and Sterilization Timer
- ⑦ Night Mode Switch
- ⑧ Temperature Selector

Technical Specifications:

- Chamber Dimensions W x L, 10" x 19"
- Chamber Volume 23L / 6.1gal
- Standard Cassettes Capacity 3 Full/3 Half
- No. of Trays (4)
- Tray Dimensions WxHxD, 16.3" x 6.7" x 0.8"
- Voltage(V) / Freq.(Hz) 230V / (50/60 Hz)
- Current(A) / Power(W) 9.6A / 2200W
- Overall Dimensions WxHxD, 21.5" x 20" x 14.4"



Tuttnauer
Innovation · Legacy · Partnership

Learn from our experts.
Visit our blog: tuttnauerusa.com/blog

Tuttnauer USA Co. Ltd.
25 Power Drive, Hauppauge, NY 11788
Tel: +800 624 5836, +631 737 4850, Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com, www.tuttnauerUSA.com

